

Maestro 7.5

Tutorial

Copyright © 2006 Schrödinger, LLC. All rights reserved. CombiGlide, Epik, Glide, Impact, Jaguar, Liaison, LigPrep, Maestro, Phase, Prime, QikProp, QikFit, QikSim, QSite, SiteMap, and Strike are trademarks of Schrödinger, LLC.

Schrödinger and MacroModel are registered trademarks of Schrödinger, LLC.

The Visualization Toolkit (VTK) is a copyrighted work (1993-2002) of Ken Martin, Will Schroeder, Bill Lorensen. All rights reserved.

Python is a copyrighted work of the Python Software Foundation. All rights reserved.

The C and C++ libraries for parsing PDB records are a copyrighted work (1989) of the Regents of the University of California. All rights reserved.

The NCSA HDF5 (Hierarchical Data Format 5) Software Library and Utilities is a copyrighted work (1998-2004) of the Board of Trustees of the University of Illinois. All rights reserved.

See the [Copyright Notices](#) for full copyright details.

To the maximum extent permitted by applicable law, this publication is provided “as is” without warranty of any kind. This publication may contain trademarks of other companies.

Please note that any third party programs (“Third Party Programs”) or third party Web sites (“Linked Sites”) referred to in this document may be subject to third party license agreements and fees. Schrödinger, LLC and its affiliates have no responsibility or liability, directly or indirectly, for the Third Party Programs or for the Linked Sites or for any damage or loss alleged to be caused by or in connection with use of or reliance thereon. Any warranties that we make regarding our own products and services do not apply to the Third Party Programs or Linked Sites, or to the interaction between, or interoperability of, our products and services and the Third Party Programs. Referrals and links to Third Party Programs and Linked Sites do not constitute an endorsement of such Third Party Programs or Linked Sites.

Revision A, April 2006

Contents

Document Conventions	vii
Chapter 1: About This Tutorial.....	1
Chapter 2: Starting Maestro and Viewing Molecules	3
2.1 Copying the Files You Need To Run this Tutorial	3
2.2 Starting Maestro	5
2.3 Features of the Maestro Main Window	5
2.4 Importing a Maestro File	8
2.5 Viewing a Structure in the Workspace	10
2.6 Including and Excluding Entries in the Project Table.....	12
2.7 Saving the Maestro Project.....	14
Chapter 3: Building Molecules Using Maestro	17
3.1 Opening the Build Panel	18
3.2 Clearing the Workspace and Building with Fragments	19
3.3 Adjusting Angles	21
3.4 Deleting and Labeling Atoms and Adjusting Formal Charge	23
3.5 Saving Your Structure to the Project Table	24
3.6 Copying and Renaming a Project Entry	25
3.7 Deleting Bonds and Reattaching Fragments	26
3.8 Adding Hydrogen Atoms to a Structure.....	28
3.9 Substituting Functional Groups in Place Mode.....	29
3.10 Replacing Functional Groups in Grow Mode.....	30
3.11 Adjusting a Dihedral Angle	31
3.12 Duplicating a Project Entry	32
3.13 Changing Bond Orders	33

3.14	Changing Elements.....	34
3.15	Adding Atoms by Freehand Drawing	35
3.16	Moving Atoms.....	38
3.17	Cleaning Up Structures.....	40
3.18	Changing the Chirality	41
3.19	Exporting Structures to a File.....	42
Chapter 4: Changing the Appearance of Structures.....		45
4.1	Getting Ready for the Exercises in this Chapter	45
4.2	Displaying and Undisplaying Atoms	46
4.3	Changing the Atom Coloring Scheme.....	51
4.4	Changing the Appearance of Atoms and Bonds.....	53
4.5	Displaying the Protein Backbone as a Ribbon	55
4.6	Creating and Displaying a Molecular Surface.....	57
Chapter 5: Measuring, Analyzing, and Superimposing Structures....		61
5.1	Getting Ready for the Exercises in this Chapter	61
5.2	Measuring Distances	62
5.3	Measuring Angles	64
5.4	Measuring Dihedrals.....	65
5.5	Superimposing Multiple Structures	67
5.6	Displaying Hydrogen Bonds	69
Chapter 6: Beyond the Basics		73
6.1	Advanced Maestro Features.....	73
6.2	Programs Accessible from Maestro.....	74

Chapter 7: Getting Help	77
Copyright Notices	79

Document Conventions

In addition to the use of italics for names of documents, the font conventions that are used in this document are summarized in the table below.

Font	Example	Use
Sans serif	Project Table	Names of GUI features, such as panels, menus, menu items, buttons, and labels
Monospace	<code>\$SCHRODINGER/maestro</code>	File names, directory names, commands, environment variables, and screen output
Italic	<i>filename</i>	Text that the user must replace with a value
Sans serif uppercase	CTRL+H	Keyboard keys

In descriptions of command syntax, the following UNIX conventions are used: braces { } enclose a choice of required items, square brackets [] enclose optional items, and the bar symbol | separates items in a list from which one item must be chosen. Lines of command syntax that wrap should be interpreted as a single command.

In this document, to *type* text means to type the required text in the specified location, and to *enter* text means to type the required text, then press the ENTER key.

About This Tutorial

This document is designed to show you how to use the basic capabilities of Maestro, in easy steps. In the early parts of the book, we will explain everything in great detail—what to do with the mouse, where to click, where buttons are, and so on—so you aren’t missing any part of the procedure. At the same time, we’ll provide a more compact notation for explaining what to do. Later in the book, we’ll use the compact notation, with references back to the explanation, and still later, we’ll just use the compact notation. By that time, you’ll have learnt what it means and how to perform the actions, and won’t need the detailed explanation. But it’ll always be there to refer back to.

This tutorial is designed for you to work through from the beginning, but you can start at any of the chapters if you so desire. However, the explanations in the later chapter are more compact, so you might have to go back to the earlier chapters to find out how to perform an action.

If you are reading this book online, you may notice that there are words that are marked in an indigo color, like this: [Contents](#). These are links to other parts of the book. If you click on the word in your PDF reader, it takes you to the other part of the book. To go back to the place you came from, click the Back button in your PDF reader. Try it with the word [Contents](#), above. It will take you to the table of contents for this book.

There are also indigo-colored links to various sections of the *[Maestro User Manual](#)*. These links are provided so you can read more about the feature you are working with, if you are interested. Clicking on these links takes you to the section referred to. Try clicking the words *Maestro User Manual* in the first sentence of this paragraph, then click the Back button.

Starting Maestro and Viewing Molecules

Maestro™ is a freely available, full-featured molecular visualization environment that also serves as the interface to all of Schrödinger's computational chemistry software. When coupled with Schrödinger software such as Glide™, Prime™, or Phase™, Maestro is a powerful tool for interpreting, managing, and sharing the results of computational experiments. As a standalone program, Maestro is an easy-to-use tool for building, visualizing, and sharing 3-dimensional chemical models. In this tutorial, you'll learn how to perform these simple tasks using Maestro's intuitively designed interface. Working through the exercises in this chapter, you'll start Maestro, import structure files, and learn how to view them in the Maestro Workspace.

Maestro is designed to run on Linux, Windows, and SGI machines, but can be displayed over a network to many commonly used platforms. With a little effort, Windows users can also install one of many Linux distributions on their PCs. Not only do these distributions include easy-to-operate graphical user interfaces, but they can be installed on a computer without removing the Windows operating system.

2.1 Copying the Files You Need To Run this Tutorial

To perform the exercises in this tutorial, you'll need to make sure that you have the structures easily accessible. These structures can be copied from the Maestro tutorial subdirectory in the installation to your own directory.

To copy the files under Windows:

1. Click Start, then click My Computer.

An explorer window opens.

2. Move to the Schrödinger software installation.

This installation is not in Programs and Settings, because the path to the software cannot have spaces in it. It's likely to be on your C drive, so start looking there if you don't know exactly where it is.

3. When you get to the installation, open the `maestro-vversion` folder.

Here, version is the 5-digit Maestro version number, for example, 75111. This folder should contain a folder labeled `tutorial`.

4. Open another explorer window and move to the place you want to keep the tutorial files.
For example, click Start, then click My Documents. You can create a folder to keep the files in if you want to.
5. Drag the tutorial folder from the first window to the second.

To copy the files under UNIX:

1. Open a terminal window so that you have a command-line prompt.

The steps to open a terminal window vary slightly between operating systems, but it is generally possible to right-click in your desktop and select an option like New Terminal from the menu that appears.

2. Next, set the SCHRODINGER environment variable.

This environment variable contains the path to your installation of Maestro. How you set the SCHRODINGER environment variable depends on the “shell” used in the terminal window. If you are using either the bash or ksh shell, enter the following command:

```
export SCHRODINGER=installation-directory
```

If you are using either the csh or the tcsh shell, enter the following command:

```
setenv SCHRODINGER installation-directory
```

In the commands above, *installation-directory* should be replaced with the path to the directory where Maestro is installed. If you’re not sure which shell you’re using you can type the command `which $SHELL`, which will return the location and name of the shell that you’re using.

3. Change to a directory for which you have write permission.

```
cd my-directory
```

4. If you like, create a subdirectory for the tutorial, and change to that subdirectory.

```
mkdir -p my-subdirectory  
cd my-subdirectory
```

5. Copy the tutorial files to your current working directory with the following command:

```
cp $SCHRODINGER/maestro-vversion/tutorial/*.* .
```

Here, version is the 5-digit Maestro version number, for example, 75111.

Many UNIX systems have a graphical file manager that operates much like a Windows system. You can use this file manager to copy the files instead of typing the commands.

2.2 Starting Maestro

After you've installed Maestro, starting the program is a simple matter.

To start Maestro under Windows:

- Simply double-click the Maestro icon on your desktop.

The Maestro interface should appear after a short delay.

To start Maestro under UNIX:

1. Make sure that you've set the `SCHRODINGER` environment variable and moved to the directory you want to work in

If you haven't done this, follow the instructions under [To copy the files under UNIX:](#) in the previous section.

2. Start Maestro by entering the following command:

```
$SCHRODINGER/maestro &
```

The Maestro interface should appear after a short delay.

The “&” at the end of the line means that Maestro is run in the “background”, and you can continue to use the terminal window. Without it, you would have to wait until Maestro has finished running to type any other commands. (This has no effect on how well Maestro runs.)

Here, *version* is the 5-digit Maestro version number, for example 70113. From here on, the terminal and the command line are no longer used—all tasks can be performed within the Maestro interface.

2.3 Features of the Maestro Main Window

Once you've started Maestro, you'll want to take note of various features in the Maestro main window, as shown in [Figure 2.1](#).

- At the center of the interface is the Maestro Workspace, which remains empty until you import structures, build molecules, or open a “project”.
- At the top of the Maestro interface is the Maestro menu bar, which gives you access to all kinds of options, features, and programs.
- On the left-hand side of the interface is the main Maestro toolbar—all of the frequently used Maestro tasks can be performed using the toolbar.

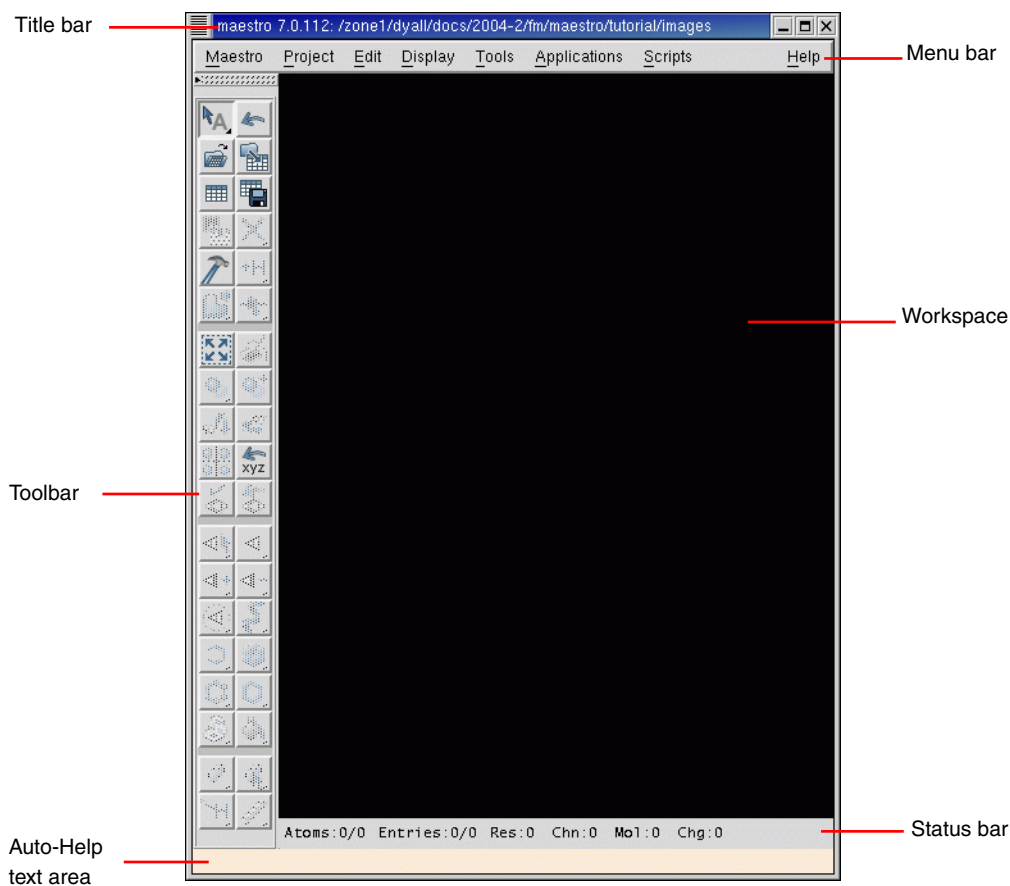


Figure 2.1. The Maestro main window.

- The status bar near the bottom of the window gives information about the number of atoms and molecules in the Workspace, and also displays more specific structural information when you pause the pointer over an atom in the Workspace.
- At the bottom of the window is the Auto-Help text area. This area automatically displays a hint about the task you're doing.

There are other features that can be added to the main window, that aren't displayed by default. You can display these features using the Display menu. For more information, see [Section 2.2](#) of the *Maestro User Manual*.



In this tutorial, you'll be using the toolbar a lot, so it's useful to know a bit about it beforehand. Notice that in [Figure 2.1](#), most of the buttons are dimmed, which means they can't be used. That's because there's nothing in the Workspace. The toolbar is shown to the left with all the buttons active.

Many of the toolbar buttons have a little triangle in the lower right corner—we'll discuss those shortly. The ones that don't are buttons that perform a simple action, like opening a panel or changing the view of what's in the Workspace. All you need to do is click the button to perform the action.

The buttons that have little triangles in the lower right corner are *menu buttons*. If you click and hold on any of these buttons, a menu appears to the lower right. This menu is called a *button menu*. You can choose an item on the menu by continuing to hold the mouse button down, moving the mouse until the item you want is highlighted, then releasing the mouse button. The menu disappears and the item you highlighted is selected. We'll call this action "choosing an item from a button menu". It's just like choosing an item from any other menu. If you accidentally click on the wrong button, just move the pointer out of the menu and release the mouse button.

If you chose an item that represents a group of atoms or a structural feature—Atoms, Bonds, Molecules, Residues, and so on—the toolbar button is indented after you release the mouse button. This indicates that the button's action will be applied to the group of atoms you pick. When you move the pointer into the Workspace, it turns into a small square with a letter beside it. The letter tells you what you chose on the menu. The button indentation and the appearance of the pointer tells you that Maestro expects you to click on an atom ("pick" an atom) or bond in the Workspace to select some atoms. The action represented by the button is then performed on these atoms. To get out of this "picking" mode, just click on the button again. The button is no longer indented, and the pointer has its normal appearance in the Workspace. When you display the button menu again, you'll see that the item you chose has a red diamond beside it.

Other items on the button menus can be groups of atoms to which the action is applied immediately, options that are set, or properties that can be applied to structures. Some items open dialog boxes in which you can perform an action or select atoms to apply the button's action to. Some examples of button menus are shown in [Figure 2.2](#).

One of the useful features of the Maestro interface is the availability of help features. In addition to the Auto-Help text area, in the upper right-hand corner of the menu bar, you'll see a Help menu that can be used to find information about Maestro features. You may have also noticed that when you pause the pointer over one of the buttons in the toolbar, a pop-up

“tooltip” appears, with text that describes the function of the button. We’ll use this text to refer to the buttons. (When you move the mouse again, it disappears.) All Maestro windows, called *panels*, have a Help button at the bottom that can be used to display relevant information. If you can’t find the information you need, you can also email help@schrodinger.com to get an answer to your questions.

2.4 Importing a Maestro File

The most common operation within Maestro is the simple manipulation of structures in the Workspace—they can be rotated, translated, and made to appear larger or smaller. But before you can experiment with these manipulations, you’ll need to import some structures. To do so:

1. Open the Import panel by clicking on the Import Structures button on the toolbar. (You can also open the panel with the keyboard shortcut CTRL+I.)



Take a moment to study the layout of the Import panel. In the top is the Path text box, which displays the location of the directory that you’ll be importing structures from. By default, the path is first set to your current working directory, but you can change this setting by navigating through the Directories list. To the right is the Files list, which displays all of the files of the selected Format that exist in the selected directory. If you correctly copied the tutorial structure files to your current working directory, you should see the file `small_molecule.mae` in the Files list.

2. To import the file `small_molecule.mae`, either double-click on the file name, or click the file name, and click Import.

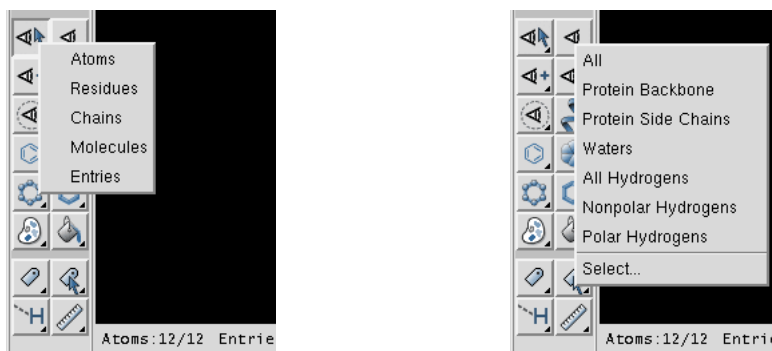


Figure 2.2. Display only selected atoms *and* Display only *button menus*.

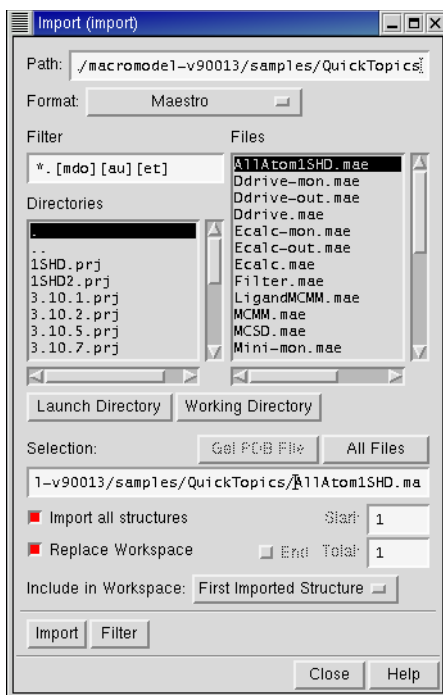


Figure 2.3. The Import panel.

You'll see that a structure appears in the Workspace. Take a moment to observe the color scheme used. Carbon atoms are colored gray, nitrogen atoms are colored blue, oxygen atoms are colored red, and hydrogen atoms are colored white. This commonly used color scheme makes it easy to identify the elements in the Workspace.

3. Close the Import panel by clicking the Close button in the lower right corner of the panel.

If you're curious about some of the options in the Import panel that weren't used, you can either click on the Help button in the Import panel, or see [Section 3.1](#) of the *Maestro User Manual*.

2.5 Viewing a Structure in the Workspace

After importing a molecule into Maestro, you're ready to manipulate the structure in the Workspace. In this exercise, you'll learn how to translate, rotate, and zoom in or out on a molecule. You'll also want to note that these manipulations require the use of a three-button mouse. If you don't have a three-button mouse (or a two-button mouse capable of three-button emulation), you may want to consider obtaining one.

But first, a word about the axis system used by Maestro. The x axis is the horizontal axis, and the y axis is the vertical axis. The z axis is the axis coming out of the screen towards you (this is the positive z direction).

Now, let's try some manipulations. And, by the way, when you do these manipulations, you aren't really changing the position of the molecule. It's as if you were operating a video camera, and you are changing the camera angle to get a different view of the molecule.

Rotating the molecule:

- To rotate the molecule left or right (around the y axis), place the pointer anywhere in the Workspace, hold down the middle mouse button, and move the mouse to the left or right.

The pointer changes to an arrow going around in a circle, with the letters xy above it, to show that you are rotating around the x (horizontal) or y (vertical) axis.

- To rotate the molecule up or down, (around the x axis) repeat the above step, but move the mouse up or down.

You can rotate in both directions simultaneously by moving the mouse diagonally. If you want to make sure you rotate only in the x direction or the y direction, hold down the SHIFT key along with the middle mouse button.

- To rotate the molecule in the plane of the screen (around the z axis), place the pointer anywhere in the Workspace, hold down the CTRL key and the middle mouse button, and move the mouse to the left or right.

The pointer changes to an arrow going around in a circle, with the letter z above it, to show that you are rotating around the z axis (the one perpendicular to the screen).

Translating the molecule:

- To translate the molecule in the Workspace (change its position in the Workspace), place the pointer anywhere in the Workspace, hold down the right mouse button, and move the mouse left, right, up, or down.

The pointer changes to a pair of double-headed arrows perpendicular to each other, to show that you are translating the Workspace structure.

Zooming in or out:

- To zoom in on the structure, place the pointer anywhere in the Workspace, hold down the middle and right mouse buttons, and move the mouse to the right.

The pointer changes to a magnifying glass, to indicate that you are zooming.

- To zoom out on the structure, repeat the above step, but move the mouse to the left instead of to the right.
- It's also possible to automatically "fit" the structure to the Workspace, so that it fills the screen. To do so, simply click the Fit to screen button on the toolbar.



Saving the view:

You can also save the orientation of the Workspace, and return to it later. Here's an example.

1. First, click the Save view button on the toolbar.



2. Next, rotate the structure toward you by clicking on the Rotate around X axis by 90 degrees button.



3. Rotate the structure 90 degrees to the right by clicking on the Rotate around Y axis by 90 degrees button.



4. Click on the Restore view button to return to the view that you saved in [Step 1](#).



Congratulations, you're now manipulating structures in the Workspace. Feel free to practise these operations until you are familiar with them. You can always restore the original view of the molecule by clicking the Reset Workspace button:



2.6 Including and Excluding Entries in the Project Table

Now that you know how to manipulate structures in the Workspace, you're ready to learn how to include and exclude structures from the Workspace using Maestro's Project Table. The Project Table is where Maestro stores structures that you've imported. It can also be used to manage the results of computations, and store structure-associated data. For now, though, we'll just use the Project Table to control which molecules are displayed in the Workspace. For more information about the Project Table, see [Chapter 8](#) of the *Maestro User Manual*.

1. Open the Import panel by clicking the Import structures button on the toolbar.



2. From the Format menu, choose SD.

This action causes all files in SD format (with the extension `.sdf`) to appear in the Files list. There are many file formats used to store molecular structures, and Maestro can read most of the common ones.

3. Double-click on `multi-structure_file.sdf` to import it.
4. Click the Close button in the bottom right-hand corner of the Import panel to close it.
5. Click the Open/Close Project table button on the toolbar to open the Project Table panel.



In the Project Table, you'll notice that several "entries", represented by various rows with text in them, are highlighted in yellow. These highlighted entries are the structures that you just imported. You'll also notice a tan row with a "-" in the In column. This is the label row for an *entry group*. When you imported the structures, they were added as a group, named after the file they came from. You can open and close the group by clicking its In column.

You can include the structures in the Workspace, or exclude them from the Workspace, by clicking on the diamonds to the left of the title text (in the In column). There is no limit to the number of entries that can be included in the Workspace. An "X" through the diamond indicates that the structure is currently included in the Workspace. By default, only the first new entry (in this case `sd_molecule_01`) is included in the Workspace when you import multiple structures.

6. Include `sd_molecule_02` in the Workspace by clicking the diamond next to its title.
(Later on, we'll just say "Include *entry* in the Workspace".)

When you click, the molecule that was originally in the Workspace, `sd_molecule_01`, is excluded from the Workspace, and `sd_molecule_02` is included in the Workspace. This is the normal behavior: including one entry in the Workspace excludes all others.

7. To include both `sd_molecule_02` and `sd_molecule_03` in the Workspace, control-click the diamond next to the entry for `sd_molecule_03`.

Holding down the control key when including an entry just affects that entry: other included entries remain included.

8. To clear the Workspace, control-click the diamond next to `sd_molecule_02`, and then control-click the diamond next to `sd_molecule_03`.

Control-clicking the diamond for an entry that is already included in the Workspace excludes it from the Workspace. This only affects whether or not the structure appears in the Workspace—no atoms are deleted when excluding an entry.

9. To include `sd_molecule_01`, `sd_molecule_02`, and `sd_molecule_03` in the Workspace, click the diamond next to `sd_molecule_01`, then shift-click the diamond next to `sd_molecule_03`.

Shift-clicking includes all entries between the two most recently clicked entries in the Workspace.

10. Clear the Workspace by clicking the Clear Workspace button on the main toolbar.



This automatically excludes all entries from the Workspace. (Later on, we'll just say "Clear the Workspace".)

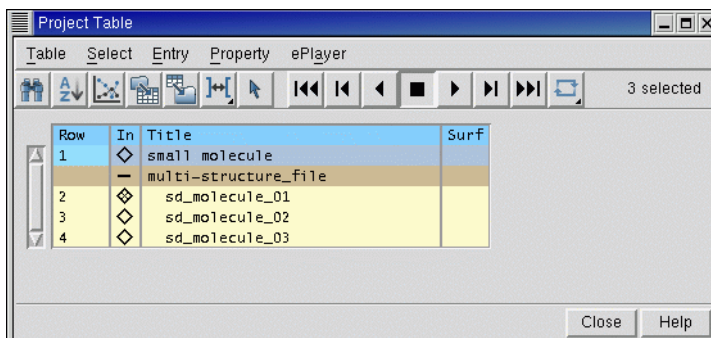


Figure 2.4. The Project Table *panel*.

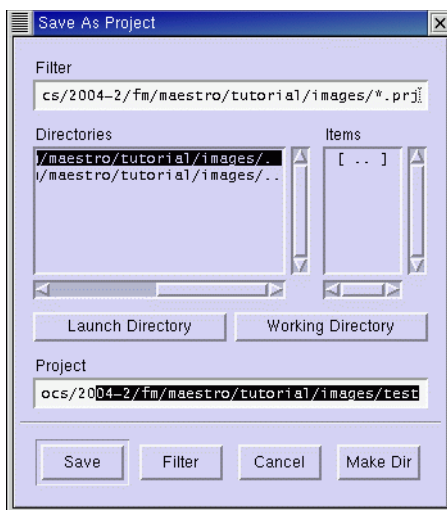


Figure 2.5. The Save As Project *dialog box*.

2.7 Saving the Maestro Project

Maestro “projects” are collections of structures and data that you can save in a designated location and then return to later on. Because structures that you work with in Maestro are copies that are saved to disk in a designated directory, structural changes that you make in Maestro will not affect the actual files that you imported. If you’ve saved your Maestro project, these copies are saved in the project directory that you specified. Otherwise, the copies are saved in a temporary directory. In this exercise, you’ll save your project so that you can return to it later.

1. Click the Save as button in the main toolbar.



The Save As Project dialog box appears.

2. In the Project text box, type a name for your project, such as tutorial.
3. Click the Save button.

Changes you make in the project are automatically saved, so there’s no need to save the project again later.

If you quit Maestro and want to re-open this project later, you can use the Open project button on the Maestro toolbar.



This button opens a file chooser in which you can navigate to the project and select it.

Building Molecules Using Maestro

Once you've mastered the process of starting Maestro and importing structures, you're ready to begin building structures of your own. Maestro's Build panel supports a wide variety of 3D building operations that can be used to easily create any reasonably sized chemical structure that you're interested in. You can even build polypeptides and DNA sequences, although experimentally resolved structures of proteins and other macromolecules are freely available online at <http://rcsb.pdb.org/>. (For protein sequences without experimentally resolved structures, Schrödinger's Prime program can be used to predict the naturally occurring low-energy protein structure.)

In this chapter, you'll work through the process of creating three familiar analgesic drugs—aspirin, tylenol, and ibuprofen—using Maestro's Build panel. For a full description of the Build panel and all its various features, see [Chapter 4](#) of the *Maestro User Manual*. If you'd like to learn about using Schrödinger software to prepare experimentally resolved protein structures for calculations, see [Chapter 4](#) of the *Glide User Manual*.

Building 3D structures is a multi-step process, and for any given structure there are many different ways to create it. This tutorial suggests just one way of building these structures, and is designed to introduce most of the useful facets of the Build panel rather than demonstrate the quickest possible way to build this particular molecule. After you complete the exercises in this chapter, feel free to rebuild these molecules without the aid of the tutorial, and see how the process can differ.

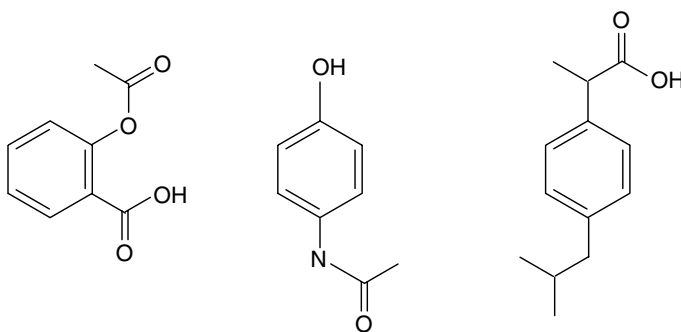


Figure 3.1. Aspirin, tylenol, and ibuprofen.

If you make an error while working through the exercises in this chapter, keep in mind that you can always undo your last action by clicking on the Undo button in Maestro's toolbar:



3.1 Opening the Build Panel

To begin the exercises in this chapter, we'll simply open the Build panel and survey its contents. To get started:

- Either click on the Open/Close Build Panel toolbar button, or use the CTRL+B keyboard shortcut to open the panel.



You'll see the Build Panel open—along with the main toolbar, it contains all the tools you need to perform the structure-building exercises in this chapter.

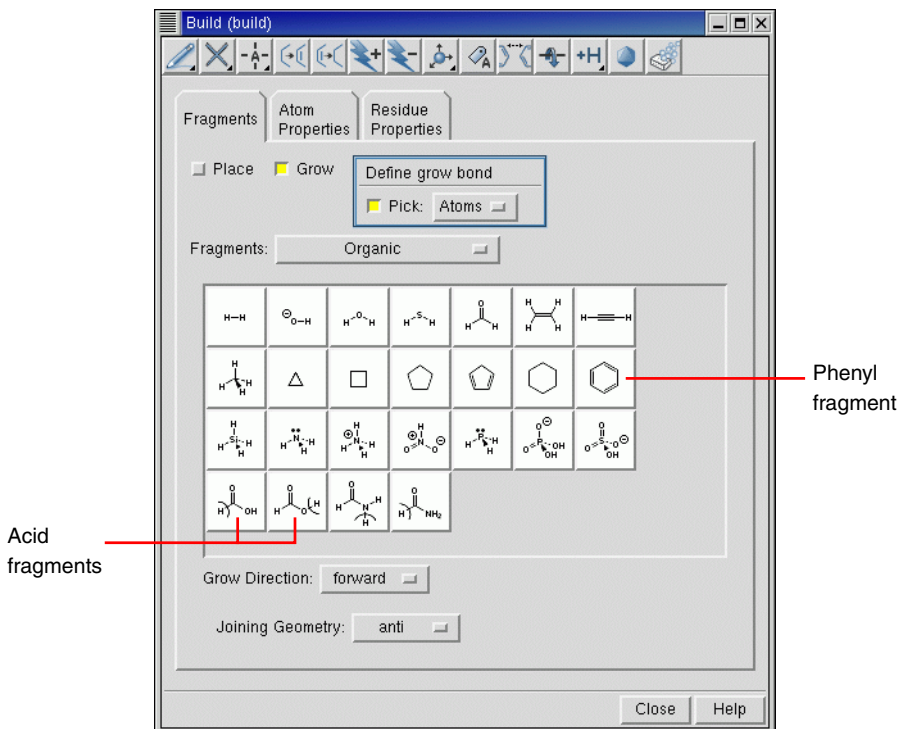


Figure 3.2. The Build panel.

At the top of the Build panel is a toolbar that allows you to quickly perform many common build operations. Below the toolbar are three folders, Fragments, Atom Properties, and Residue Properties. The Fragments folder allows you to simply build common, predefined substructures in the Workspace with one click—anything from simple organic fragments to more complicated heterocycles. The Atom Properties folder allows you to set element types, formal charges, and so on, although much of this behavior can be performed using the toolbar. The Residue Properties folder can be useful when building polypeptides or refining protein structures, but is generally not needed when building small molecules.

3.2 Clearing the Workspace and Building with Fragments


Now that you've opened the build panel, you can begin building the structure using Maestro's fragment libraries. These pre-defined groups of substructures and functional groups can be placed in the Workspace or attached to existing structures, eliminating the need to manually add every individual atom in a structure.

1. Begin by clicking the Clear workspace button on the main toolbar.



This simply “excludes” any structures that may currently be in the Workspace—no atoms or molecules in existing structures are deleted. This button is dimmed if there's nothing in the Workspace, so you don't need to click it.

2. Click the Fragments tab to display a selection of fragment libraries.
3. Choose Organic from the menu labeled Fragments.
4. Click on the phenyl fragment in the library in the center of the panel (see [Figure 3.2](#)).

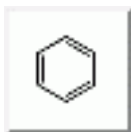


5. To place the fragment, click near the center of the Workspace.

You'll see the complete phenyl ring appear in “wire frame” representation. Carbon atoms are colored gray, and hydrogen atoms are colored white.

Next, we'll add an acid group to the ring. There are two acid fragments in the Organic fragment library, one for adding at the C–H bond, and the other for adding at the O–H bond. The bond that is replaced is indicated by a curved line.

6. Click on the carboxylic acid fragment whose C–H bond is divided by a curved line.

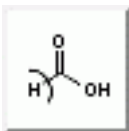


7. To place this fragment, click on any of the white hydrogen atoms in the Workspace.

When you move the pointer into the Workspace, it turns into a yellow square with an “A” beside it. The square is the “target area”—when an atom is inside the square and you click, that atom is selected for the operation you are doing.

Once you click, the hydrogen atom is replaced with a carboxylic acid group. Now you have a benzoic acid molecule in the Workspace.

To complete the structure of aspirin, we need to add an acetate group. We can do this by adding the other carboxylic acid group and a methyl group, in two steps.


8. Click on the other carboxylic acid fragment, whose O–H bond is divided by a curved line.

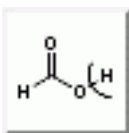



9. Click on the ring hydrogen atom nearest to the OH of the acid group, as indicated in the figure to the right.

The acid is added to form an ester, but it is a formate ester, not an acetate ester. The last step is to add a methyl group in place of the formate hydrogen.

10. Click on the methyl fragment.

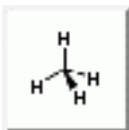

11. Now click on the hydrogen of the formate group, as indicated in the figure to the right.
12. The formate group is converted into acetate. You should now have a structure that looks like the one in [Figure 3.3](#).

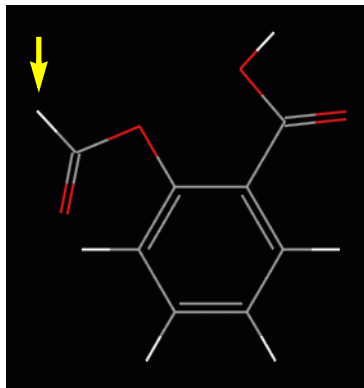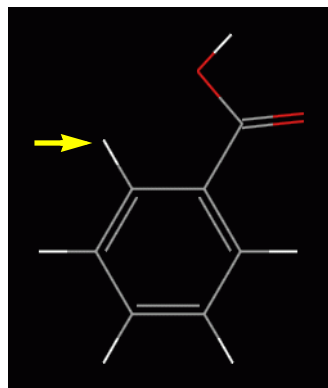

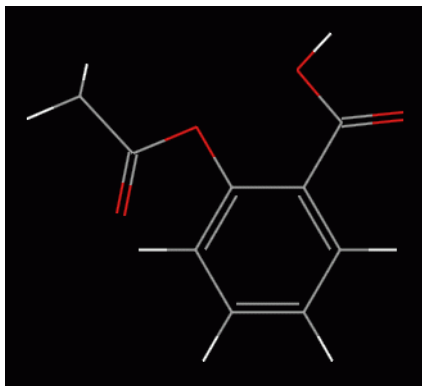


Figure 3.3. The completed aspirin structure.

Congratulations! You’ve just successfully build aspirin (acetylsalicylic acid) in a few easy steps using Maestro’s fragment library. In the following exercises, you’ll see how to adjust the positions of the atoms in a molecule and transform this structure into other molecules.

3.3 Adjusting Angles

Some of the angles in the structure you’ve just built don’t look quite right. In this exercise, we’ll adjust the C–O–C angle of the ester to make the structure look a bit better.

1. Click and hold on the **Adjust** button on the main Maestro toolbar (*not* the Build panel toolbar), and choose **Angle** from the menu that appears.



(If you don’t remember how to do this, see the explanation on [page 7](#).)

Notice that the pointer changes to a yellow square with an “A” beside it, to indicate that you are picking atoms, in this case to define the angle.

2. Click on the carbon atom in the ring, then the oxygen atom, then the carbon atom of the carboxylate (see [Figure 3.4](#)).

At each click, a purple box appears around the atom. When you have clicked the third atom, the boxes disappear and the angle is marked with a green line over the bonds and an orange arrow, with the value of the angle next to it (see [Figure 3.4](#)).

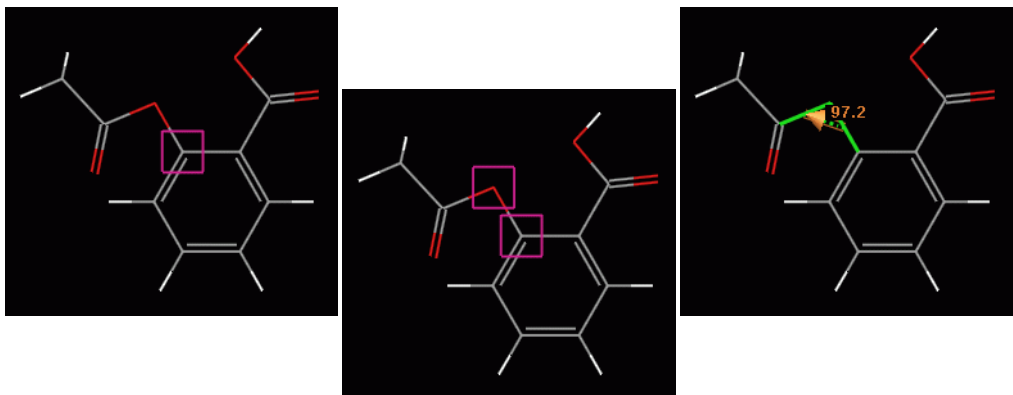


Figure 3.4. Selecting an angle.

3. Hold down the middle mouse button and move the mouse horizontally until the angle is about 120° (this is called “dragging horizontally with the middle mouse button”).

You’ll see the angle and its value change as you move the mouse. The atoms that move are the ones attached to the last atom you picked in [Step 2](#).

4. Click and hold on the Adjust button on the main toolbar, and choose Delete Adjustments from the button’s menu.



The lines and the arrow that marked the angles disappear from the Workspace.

5. Drag horizontally with the middle mouse button until the angle is about 0°.
6. Choose Delete Adjustments from the Adjust button menu on the main toolbar. (Click and hold on the Adjust button, and choose Delete Adjustments from the button’s menu).

The lines and the arrow that marked the angles disappear from the Workspace.

7. Finally, click the Adjust button (click and release immediately, not click and hold) to leave adjustment mode.

3.4 Deleting and Labeling Atoms and Adjusting Formal Charge

Acetylsalicylic acid is partially ionized in solution. In this exercise, you'll delete the acidic hydrogen, label atoms with structural information, and adjust the formal charge of an atom.

1. Choose Atoms from the Delete button menu on the main toolbar.



(If you don't remember how to do this, see the explanation on [page 7](#).) When you move the pointer into the Workspace, it turns into a square with an "A" beside it, to show that you are picking atoms, which will be deleted.


2. Click on the acidic hydrogen atom in the Workspace.

The hydrogen atom is removed. If you make a mistake when you delete an atom, you can always click the Undo button to restore it.

3. Next, choose Formal charge from the Label atoms button menu on the main toolbar.



You'll notice that no new information appears in the Workspace—that's because the formal charge of all the atoms is still zero. The oxygen atom of the acidic group needs to have a formal charge of -1 for a carboxylate anion.

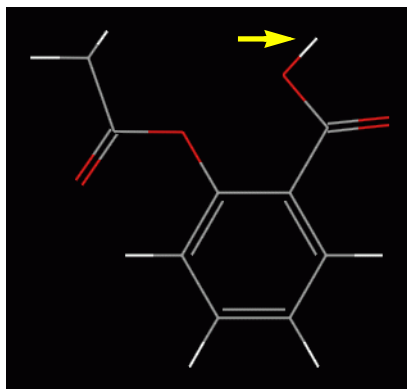

4. Click the Decrement formal charge button on the Build panel toolbar.



The button is indented. When you move the pointer into the Workspace, it turns into a square with an "A" beside it, to show that you are picking atoms, which will have their formal charge decremented.

5. Click on the oxygen atom that the hydrogen atom was attached to.

You'll see that the atom is then labeled with formal charge information (see [Figure 3.5](#)).

6. To undisplay the atom labels, choose Delete Labels from the Label atoms button menu.

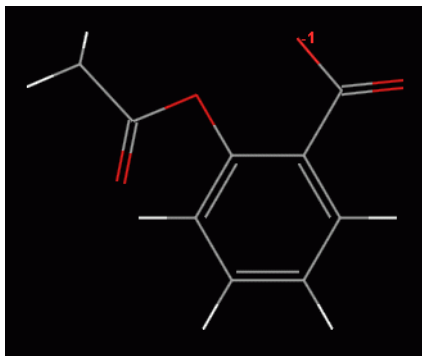


Figure 3.5. The ionized aspirin structure showing the formal charge.

Congratulations, you've just adjusted the geometry and formal charge of a structure. In the next exercise, you'll place the structure in the Maestro Project Table.

3.5 Saving Your Structure to the Project Table

Once you've created your structure, it's time to place it in the Project Table. Placing your structure in the Project Table allows you to include and exclude the structure from the Workspace, save it to disk, copy it, and add properties to it. First, we'll open the Project Table (if it's not already open from the last chapter), then add the structure to it.

1. If the Project Table panel isn't open, click the Open/Close project table button on the main toolbar.



2. Click the Create entry from Workspace button on the main toolbar.



The Create entry from Workspace dialog box opens.

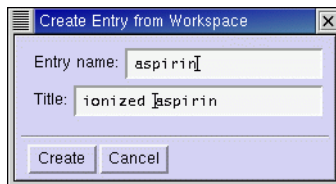


Figure 3.6. The Create Entry from Workspace dialog box.

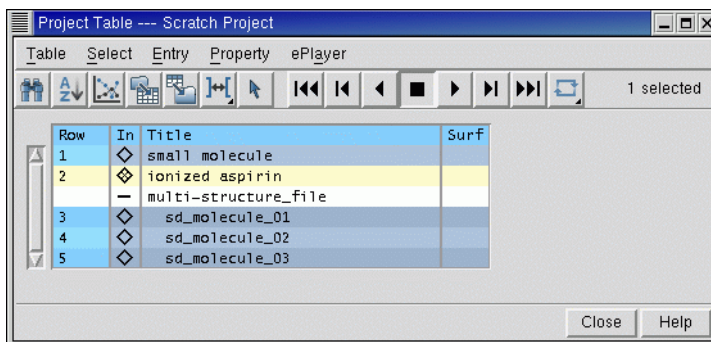


Figure 3.7. The Project Table panel with the aspirin entry.

3. Set the Entry name to “aspirin” and the Title to “ionized aspirin”.
4. Click Create to add the structure to the Project Table.

The structure is added to the Project Table as an “entry” in the last row.

3.6 Copying and Renaming a Project Entry

For the next few exercises, you’ll need a copy of the aspirin structure that you just created. We’ll create a copy by duplicating the aspirin entry and renaming it.

1. Make sure the aspirin entry is selected by clicking on its row number.

The entry should be highlighted in yellow, as in [Figure 3.7](#), to show that it is selected.

2. From the Entry menu on the Project Table menu bar, choose Duplicate.

Another row appears in the Project Table, with the same title. This isn’t a problem for Maestro, because the entry also has an “entry name” which is unique, but which is normally hidden. (If you want to see it, choose Show Entry Name from the Property menu.)

Now we’ll change the title, so we can tell which entry is which. (There’s no need to change the entry name.)

3. Click in the Title column of the second row of the Project Table.

The table cell turns a slate blue, and a cursor appears at the end of the text, to indicate that you can edit the contents of the table cell.

4. Change the title to “tylenol” and press ENTER.

The table cell turns back to yellow, and the cursor disappears. The title is now changed.

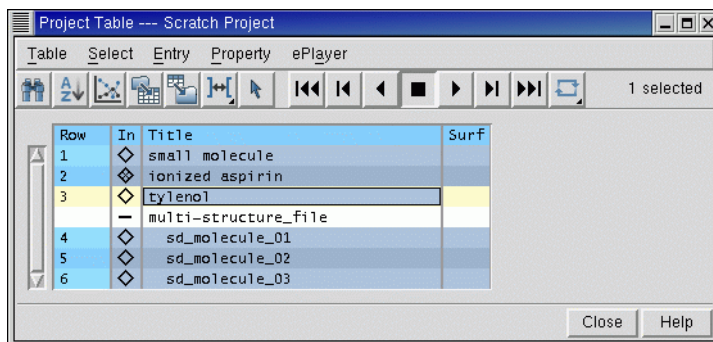


Figure 3.8. The Project Table panel with the title highlighted in the duplicated entry.

3.7 Deleting Bonds and Reattaching Fragments

Now that you have duplicated and renamed the aspirin entry, you can modify the new entry to create another molecule. This process can be useful when you want to make a series of molecules that have a similar core structure but with different groups attached.

In this exercise, you'll move the acetate group from the ortho position to the para position. First, we want the new entry in the Workspace.

1. Click the diamond for the tylenol entry in the Project Table.

An “X” appears in the diamond, but no change is visible in the Workspace. That's because we haven't made any change in the structure yet.

2. Choose Bonds from the Delete button menu on the main toolbar.



(If you don't remember how to do this, see the explanation on [page 7](#).) When you move the pointer into the Workspace, it turns into a square with an “B” beside it, to show that you are picking bonds, which will be deleted.

3. Click on the bond between the ring and the acetate oxygen (see [Figure 3.9](#)).

The bond disappears, and there are two fragments in the Workspace.

4. Click the Delete toolbar button again to exit Delete mode.

Next, we want to move the acetate fragment. Normally, when you move something in the Workspace, everything moves. To move this fragment and leave everything else where it is, we must perform a “local transformation”.

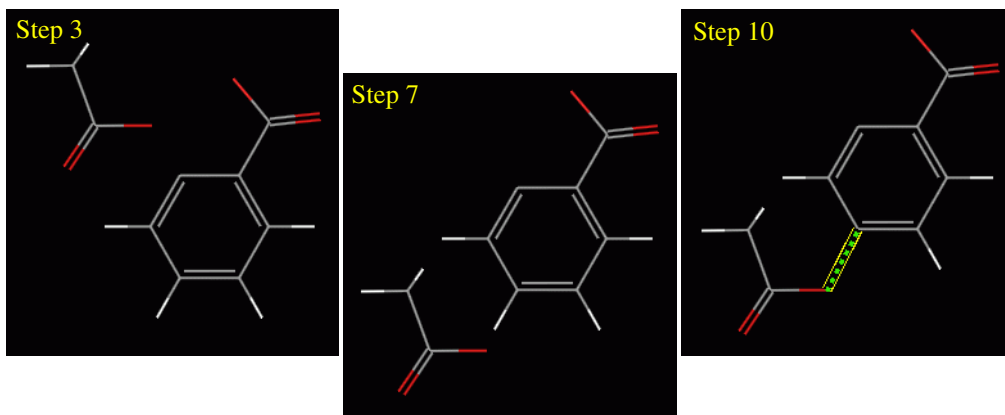


Figure 3.9. Detaching, moving and reattaching the acetate fragment.

5. Choose Molecules from the Local transformation button menu on the main toolbar.



6. Click an atom in the acetate fragment.

This selects the acetate molecule for a local transformation. You can now move it without moving the rest of the structure.

7. Drag the fragment with the right mouse button so that the oxygen atom that was attached to the ring is near the para hydrogen of the ring (see [Figure 3.9](#)).

If you see a menu appear in the Workspace, it means that you were a bit too slow in dragging the fragment. Just release the right mouse button and press it again to drag. If you pressed it while the pointer was over an atom, a yellow dot appears on the atom. To get rid of the dot, click the Workspace selection toolbar button, then click in the Workspace; if you didn't finish moving the fragment, click again on the Local transformation button. These menus are shortcuts to various commands that you can apply to the Workspace or to the selected atoms. We won't be using them, but feel free to explore them.

Now we'll reconnect the acetate fragment to the ring.

8. Click the Display Connect & Fuse panel button on the Build panel toolbar.



The Connect & Fuse panel opens.

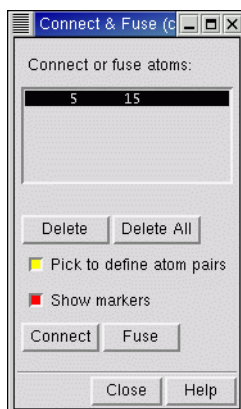


Figure 3.10. The Connect & Fuse panel.

9. Make sure that Pick to define atom pairs and Show markers are both selected.
10. Click on the ring carbon atom in the para (or 4) position, then on the singly-bonded oxygen atom of the acetate fragment (see [Figure 3.9](#)).

A purple box appears around the carbon atom on the first click, then on the second click, a green dotted line appears between the two atoms, with yellow full lines on either side (see [Figure 3.9](#)). In the Connect or fuse atoms list in the Connect & Fuse panel, the atom numbers of the two atoms are displayed.

11. Click Connect.

A bond appears between the acetate oxygen and the ring carbon, and the hydrogen that was attached to the ring is deleted. The methyl group of the acetate is very close to the ring. We will adjust the structure later.

If you want to connect two molecules with more than one bond, you can select more atom pairs. When you click Connect, a bond is created for each atom pair.

3.8 Adding Hydrogen Atoms to a Structure

Now that the acetate group is joined back on, we need to add a hydrogen atom to the place where it was removed. Although you can add hydrogen atoms one-by-one to a structure, this is unnecessarily tedious. Maestro allows you to automatically add hydrogen atoms wherever they're needed with just a few clicks of the mouse. In this exercise, we'll add the missing hydrogen atom to the ring, where we cut off the acetate group.

Missing hydrogen atoms can be added to any structure in the Workspace by double-clicking the Add hydrogens button on either the Build panel toolbar or the main Maestro toolbar.

- Double-click the Add hydrogens button to build the necessary hydrogen atoms.




The ring carbon atom now has a hydrogen atom attached to it. For more information on the Add Hydrogens feature, see [Section 4.8](#) of the *Maestro User Manual*.

3.9 Substituting Functional Groups in Place Mode

When you built the aspirin structure, you added functional groups by clicking on hydrogen atoms. The new group replaced the hydrogen atom. You can also substitute functional groups in the middle of a molecule, provided they have similar sorts of connections. For example, you can convert a hydrocarbon into an ether by replacing a $\text{-CH}_2\text{-}$ group with an -O- group. In this exercise, you'll convert the acetate group into an acetamide group.

1. In the Build panel, click on the amine fragment.

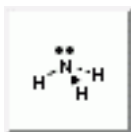

2. Click on the oxygen atom that connects the acetate group to the ring.

The oxygen atom is replaced with an NH group: the acetate group is now an acetamide group.



Figure 3.11. The structure after conversion of acetate to acetamide.

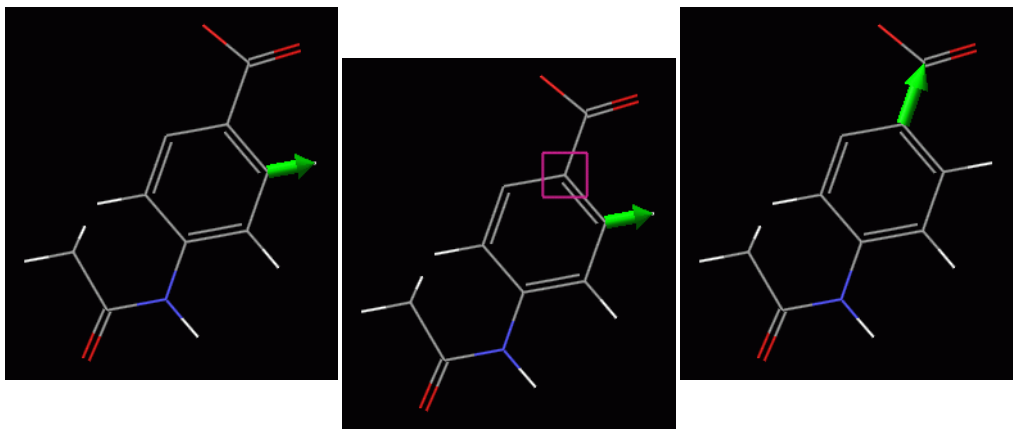


Figure 3.12. *Moving the grow bond.*

3.10 Replacing Functional Groups in Grow Mode

Maestro has two modes for adding or replacing fragments in a molecule, “Place” mode and “Grow” mode. So far, you have been using Place mode, in which you click to place a fragment in the desired location. In Grow mode, the location of the fragment is determined by the “grow bond”, which is displayed in the Workspace as a green arrow. This mode is very useful for rapidly adding fragments to a molecule. All you need to do is to click on the fragment you want to add next, and it’s added to the structure. You can also use Grow mode to replace fragments, which is what you’ll do in this exercise.


1. In the Build panel, select Grow.

A green arrow appears on the C–H bond next to the carboxylate group (see [Figure 3.12](#)). This isn’t the place we want to add a fragment, so we’ll move the grow bond. (It’s possible that the green arrow might be somewhere else, but it doesn’t matter if it is.)

2. Make sure that Pick is selected in the Define grow bond section, and that Atoms is selected on the Pick menu.
3. Click on the ring carbon that’s attached to the carboxylate group, then on the carboxylate carbon (see [Figure 3.12](#)).

The green arrow moves to the bond you just picked, and is pointing towards the carboxylate. If you had picked the atoms in the other order, the bond would be pointing in the opposite direction. The atom that the arrow head points to is the one that is replaced.

- Click on the hydroxyl fragment in the Build panel.

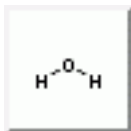

The carboxylate is replaced with a hydroxyl group.

- Deselect Grow in the Build panel.

3.11 Adjusting a Dihedral Angle

Finally, the acetyl group needs to be rotated around the C–N bond because the methyl group is too close to the ring hydrogen. We'll change a dihedral angle to do this.

- Choose Dihedral from the Adjust button menu on the main toolbar.



- Click on the carbon atom in the ring that the nitrogen atom is attached to, then the nitrogen atom, then the carbonyl carbon atom, and finally, the oxygen atom (see [Figure 3.13](#)).

At each click, a purple box appears around the atom. When you have clicked the fourth atom, the boxes disappear and the angle is marked with a red line over the bonds and a turquoise arrow, with the value of the angle next to it. The angle should be about 180° (it might be plus or minus).

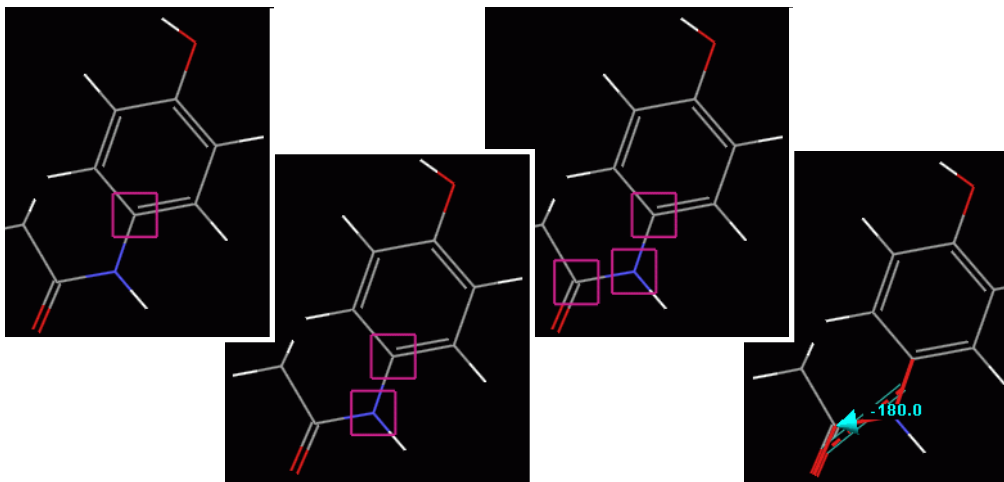


Figure 3.13. Selecting a dihedral angle.

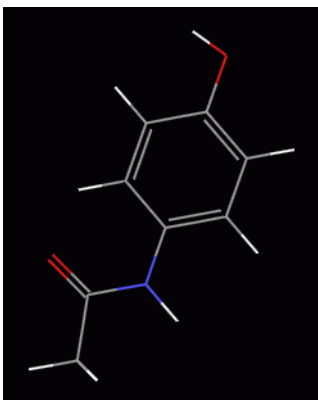


Figure 3.14. The completed tylenol structure.

3. Hold down the middle mouse button and move the mouse horizontally until the angle is about 0°.

You should see the acetyl group rotate around the C–N bond as you move the mouse.

4. Choose Delete Adjustments from the Adjust button menu on the main toolbar.



The lines and the arrow that marked the angles disappear from the Workspace.

5. Finally, click the button again to leave adjustment mode.

Congratulations! You have just turned aspirin into tylenol. You've also learned how to use Grow mode, replace functional groups in both Grow and Place modes, and detach and reattach fragments.

In the remaining exercises, you'll convert tylenol into ibuprofen.

3.12 Duplicating a Project Entry

Before we go ahead with the structural changes, we need to duplicate the tylenol entry and rename it.

1. Make sure the tylenol entry is selected in the Project Table by clicking on its row number.

The entry should be highlighted in yellow.

2. From the Entry menu on the Project Table menu bar, choose Duplicate.

Another row appears in the Project Table, with the same Title.

3. Now change the title to “ibuprofen”. If you aren’t sure how to do this, see [Section 3.6 on page 25](#).
4. Finally, click the In column of the new entry to place it in the Workspace.

We now have a new entry, and we will convert the structure in it to ibuprofen.

3.13 Changing Bond Orders

So far, we haven’t had to change any of the bond orders. The functional groups we added already had the correct bond orders. In this exercise, the first thing to do is to change some bond orders. Maestro has a very easy way of changing bond orders. We’re going to experiment with adjusting bond orders while we prepare the structure for the next step.

1. In the Build panel, click the Decrement bond order toolbar button.



2. In the Workspace, click on the carbonyl C=O bond.

The double bond turns into a single bond.

3. In the Build panel, click the Increment bond order toolbar button.



4. In the Workspace, click on the C–N bond between the carbonyl carbon and the nitrogen.

You’ll see the bond order increase, and the number of hydrogen atoms is automatically adjusted.

5. Click again on the carbon-nitrogen bond.

The bond is now a triple bond, but the bonds to the other atoms aren’t changed. Maestro lets you violate the valences of atoms when adjusting bonds, because it’s sometimes more convenient to increase the bond order of one bond before decreasing the bond order of other bonds. Before you finish changing bond orders, don’t forget to make sure that the valences are correct. Now we’ll set the bond orders the way they need to be.

6. In the Build panel, click the Decrement bond order toolbar button.



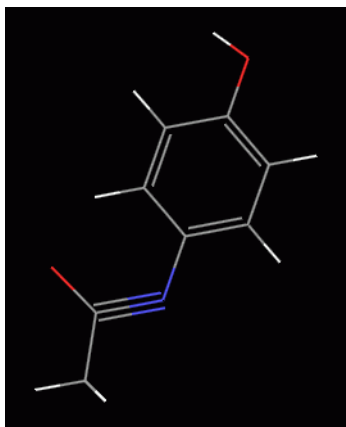


Figure 3.15. The incremented C–N bond showing the valence violation.

7. In the Workspace, click twice on the carbon-nitrogen bond.

The triple bond turns back into a single bond. The hydrogen atoms aren't added back again, so the nitrogen is missing a hydrogen atom. We'll add hydrogens later, when we have finished changing the rest of the structure.

3.14 Changing Elements

Once you've built a molecule, it's possible to change the atoms in the molecule to atoms of a different element without deleting and then rebuilding a substructure. In this exercise, you'll see how to change elements using the Build panel. This is the next step in converting tylenol to ibuprofen.

1. Choose C from the Set element button menu on the Build panel toolbar.



2. Click on the nitrogen atom, then on the two oxygen atoms.

The blue nitrogen atom and the red oxygen atoms are replaced with gray carbon atoms, changing the side chains into hydrocarbon chains. No hydrogen atoms were added to the carbons—we'll add these in a later exercise.

For less commonly modeled elements that aren't available from the toolbar, you can click the Atom Properties tab, choose Element from the Property option menu, and then select an element in the periodic table by clicking on it. Clicking on an atom in the Workspace then changes it to the selected element.

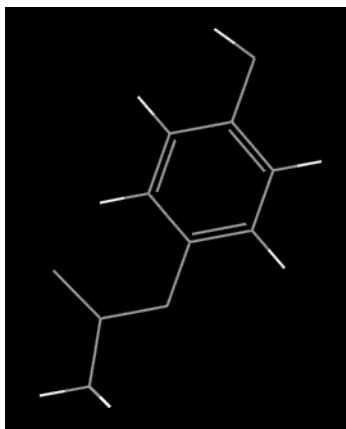


Figure 3.16. The structure after conversion of N and O to C.

3.15 Adding Atoms by Freehand Drawing

Although the Build panel contains a variety of common fragments that can be combined to form complex molecules, it is sometimes useful to build parts of your structure manually. In freehand drawing mode, Maestro allows you to do just that. This exercise will introduce you to the freehand drawing tool as you continue to convert tylenol to ibuprofen. For more information on freehand drawing in Maestro, see [Section 4.3](#) of the *Maestro User Manual*.

1. Click and hold the Draw structures button on the Build panel toolbar, and choose C from the menu that appears.



This action sets the element type to carbon. The menu allows you to choose from a variety of common elements, including hydrogen, oxygen, nitrogen, sulfur, and phosphorous.

2. Click on the carbon atom attached to the ring that only has a hydrogen attached to it (see [Figure 3.17](#)).

A purple box appears around the atom, indicating that you've selected this atom (the "active" atom) as the starting point for the next drawing operation.

3. Click near the hydrogen atom to place a new carbon atom (see [Figure 3.17](#)).

You'll see that a new carbon atom is placed and bonded to the first carbon atom. It's only a carbon atom—it doesn't have any hydrogens attached to it. We'll add those later.

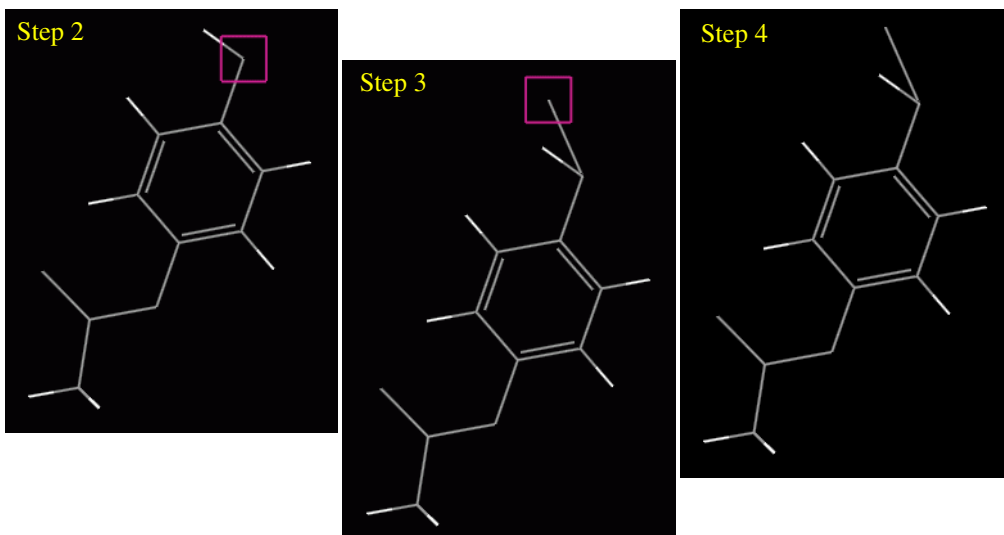


Figure 3.17. Adding a carbon atom by drawing.

4. Click again on the new carbon atom.

The purple box disappears, to show that you've finished drawing at this location.

We need to add a carboxylic acid group to the carbon atom. We could do this with the fragments, but this time we'll do it by drawing.

5. Click on the same carbon atom as you clicked on in [Step 2](#) (see [Figure 3.18](#)).

The purple box around the atom reappears.

6. Click near this carbon atom, but away from the carbon and hydrogen atoms (see [Figure 3.18](#)).

Another new carbon atom is placed and bonded to the first carbon atom.

Now, we'll change the element we're drawing with to oxygen, so we can draw the acid group.

7. Choose O from the Draw structures button menu on the Build panel toolbar.



Notice that the purple box is still on the carbon atom you placed. This means that drawing is still active, and you can click in the Workspace to place a new atom that is bonded to this atom, but this time it will be an oxygen atom.

8. Click near the active carbon atom (see [Figure 3.18](#)).

A red oxygen atom appears, bonded to the carbon atom, but the bond is only a single bond. We could use the toolbar buttons to increment the bond, but this time we'll draw in the second bond.

9. Click back on the carbon atom that's bonded to the oxygen (see [Figure 3.18](#)).

The bond changes from a single bond to a double bond, and the carbon atom becomes the active atom. Notice that, even though we are drawing with oxygen, the carbon atom didn't turn into an oxygen atom. If you click on an atom that's already there, it doesn't change to a different element. That's very useful for drawing in bonds—you don't have to worry about which element you have selected.

10. Click again near the active carbon atom, away from the other atoms (see [Figure 3.18](#)).

Another oxygen atom is placed. New atoms are always atoms of the selected element; old atoms keep their elemental identity.

You can also join one atom to another in Draw mode. We'll do this, and then undo the operation.

11. Click on the carbon atom you placed first.

A bond is created, to form a 4-membered ring.

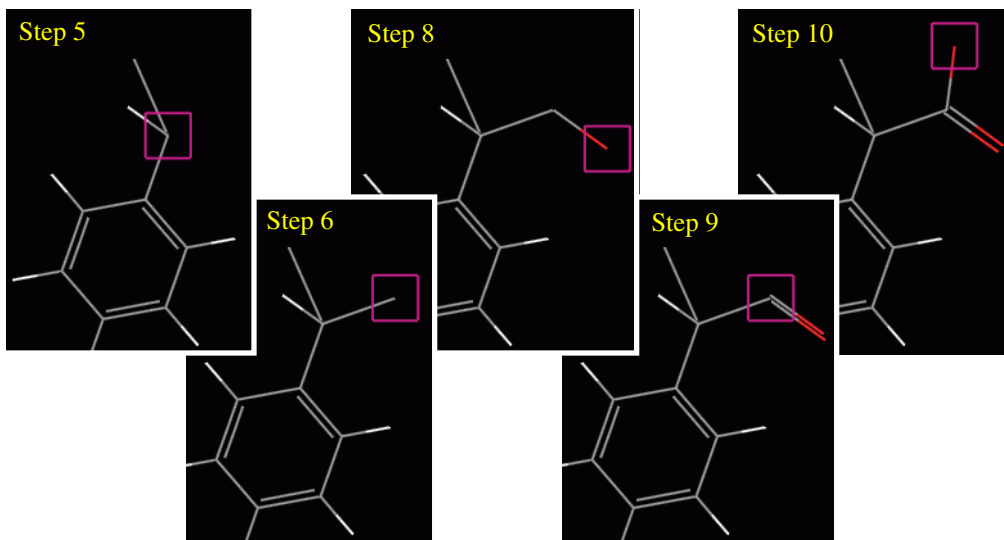


Figure 3.18. Completing the carboxylate group by drawing.

- Click the Undo button on the main toolbar.



The bond disappears, and the purple cube is now on the oxygen atom.

- Click the active oxygen atom to finish drawing.

Finally, we need to add hydrogen atoms. We could draw in the hydrogen atoms on this structure (and feel free to do so if you like), but as we saw in [Section 3.8](#), Maestro provides a convenient way of automatically adding hydrogen atoms wherever they're needed.

- Double-click the Add hydrogens button to add the necessary hydrogen atoms.



Hydrogen atoms are added to all the carbon atoms that are missing them, and to the singly bonded oxygen in the carboxylate group.

3.16 Moving Atoms

When you've been drawing atoms, it's usually the case that they aren't in the optimum location—the bond lengths and angles can be far from what they are in the real structure. Also, when you are drawing atoms, they are always placed in the same plane. If you want to place atoms in a different plane, you can rotate the structure and start placing them again, but if they are already in place, you might want to move them. The Build panel provides a tool for moving atoms that have already been placed.

In this exercise, you'll move some of the atoms so that they are in a more realistic location. In the next exercise, you'll use another tool to get much better structures.

First, let's move one of the CH hydrogens out of the plane.

- Choose $-Z$ from the Move button menu on the Build panel toolbar.



- Click on the CH hydrogen atom in the $-\text{CH}(\text{CH}_3)\text{COOH}$ group (see [Figure 3.19](#)).

Although it looks like nothing happened, the hydrogen moved 0.5 \AA away from you (in the $-z$ direction) each click. You can rotate the structure to see what happened. Click the Reset Workspace button when you have finished looking at the rotated structure.

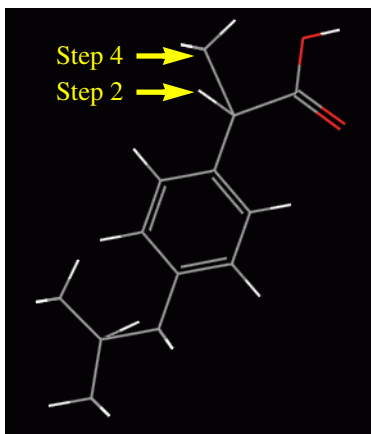


Figure 3.19. The hydrogen and carbon atoms to move.

Next, we'll move one of the methyl group out of the plane.

3. Choose +Z from the Move button menu on the Build panel toolbar.



4. Click twice on the methyl carbon atom in the $-\text{CH}(\text{CH}_3)\text{COOH}$ group.
5. Now rotate the structure (middle mouse button) to see what happened.

Notice that the methyl hydrogens moved with the carbon atoms. Groups that are attached to an atom move with the atom.

6. Click the Reset Workspace button after you have looked at the rotated structure.



Finally, we'll move one of the CH hydrogen atoms in plane.

7. Choose XY from the Move button menu on the Build panel toolbar.



8. Click on the CH hydrogen atom in the $-\text{CH}_2\text{CH}(\text{CH}_3)_2$ group (see [Figure 3.20](#)).

A purple cube appears around the atom, to show that it has been selected for movement in the xy-plane.



Figure 3.20. Moving the CH hydrogen in the xy plane.

9. Click somewhere between the two methyl groups.

The hydrogen atom moves to the new location.

Moving atoms with these tools is useful when you are drawing a structure, so that you get approximately the right geometry, but it doesn't give structures that are close to the "real" structure. In the next exercise, you'll see how to improve the structure.

3.17 Cleaning Up Structures

Maestro's fragment libraries allow you to build molecules out of fragments with optimized geometries, and the structures you create usually have reasonable geometries. When you have been using Draw mode, the structures usually need cleaning up. The Build panel allows you to "clean up" the structure using the UFF (universal force field) minimizer.

- Click the Clean up geometry button on the Build panel toolbar.



As the structure is refined, a dialog box appears that says "Cleanup in progress", and you can see the structure changing in the Workspace. When cleanup is finished, notice how the functional groups have moved (see [Figure 3.21](#)).

The UFF minimizer is suitable for refining manually built structures, but it offers less accuracy and flexibility than the minimizers and force fields present in software such as MacroModel, and it sometimes doesn't give good geometries for groups that should be planar. Always check your structures after using the UFF minimizer.

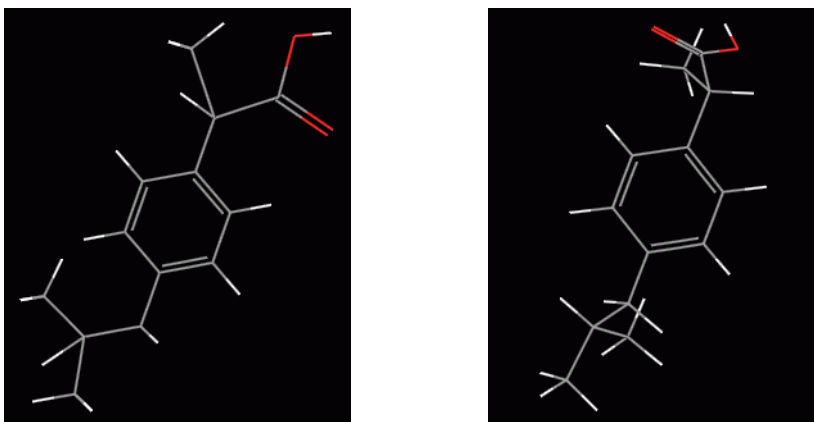


Figure 3.21. The ibuprofen structure before (left) and after (right) cleanup.

3.18 Changing the Chirality

Ibuprofen is a chiral molecule. The configuration at the chiral atom of the actual molecule (the CH carbon atom in the $-\text{CH}(\text{CH}_3)\text{COOH}$ group) is S. We need to check the chirality of this atom, and change it if it isn't correct. First, we'll label the atoms with their chirality labels.

1. Choose Chirality from the Label atoms button menu on the main toolbar.



A label appears on the chiral atom. For the structure that you built, it should be R. If you can't see the label clearly, rotate the molecule.

Now we'll change the chirality from R to S.

2. Choose Chirality from the Adjust button menu on the main toolbar (*not* the Build panel toolbar).



To change the chirality, you have to pick three atoms: the chiral atom, then the two that you want to stay where they are. The other two atoms are swapped, with their attached groups.

3. Click on the chiral carbon atom, then click on the ring carbon and the carboxylate carbon (see [Figure 3.22](#)).

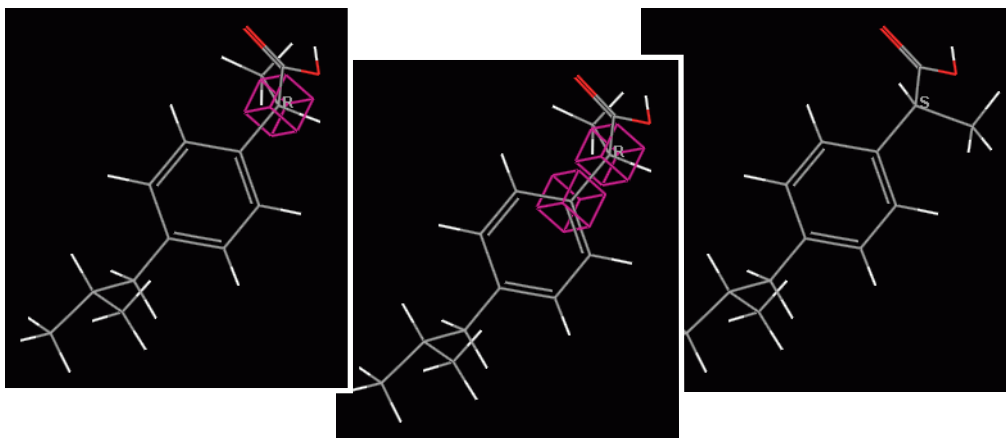


Figure 3.22. Changing the configuration of ibuprofen from R to S.

On the first and second clicks, a purple cube appears around the picked atoms. On the third click, the groups are moved, and the chirality label changes from R to S.

Congratulations! You have completed the transformation of aspirin to tylenol to ibuprofen, and have learnt how to use most of the tools in the Build panel. In the next few chapters, you'll see how to use some of the other capabilities of Maestro.

3.19 Exporting Structures to a File

Now that you've created these three molecules, you can share them by exporting them to a disk file. SD format is used by a number of other programs, so we'll export them in this format.

1. In the Project Table, click the aspirin entry, then shift-click the ibuprofen entry.

The three structures should be highlighted in yellow, which means they are selected.

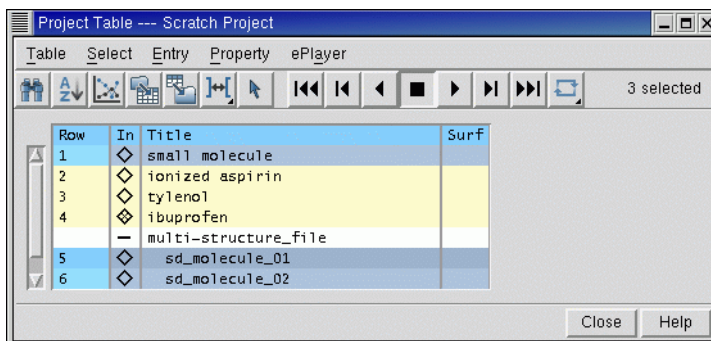


Figure 3.23. The Project Table panel with the three entries selected.

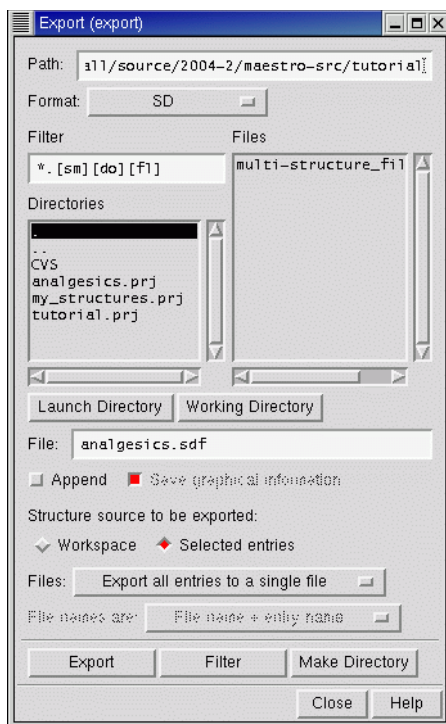


Figure 3.24. The Export panel.

2. Choose Export from the Table menu in the Project Table panel, to open the Export panel.
3. Choose SD from the Format option menu.
4. If you want to save the file in place other than your tutorial directory, navigate to the place you want to save the structure.
5. Enter a name in the File text box, such as `analgesics.sdf`.
6. Click Export.

The file should now be saved to disk.

7. Click Close to close the Export panel.

To learn more about to exporting your structures, see [Section 3.2](#) of the *Maestro User Manual*.

Changing the Appearance of Structures

Now that you know how to import and build structures using Maestro, you're ready to change the way structures appear in the Workspace. Using Maestro's visualization tools, you can color atoms, change the three-dimensional representation of atoms and bonds, and view molecular surfaces. This section of the tutorial will demonstrate how to make such changes, allowing you to draw attention to important structural features and create presentation-quality graphics.

4.1 Getting Ready for the Exercises in this Chapter

To perform the exercises in this chapter, you'll need to import the structure file `visualization_exercises.mae`. This file includes a protein-ligand complex originally taken from the PDB (Protein Data Bank) and "prepared" using Schrödinger's protein preparation tools.¹ This protein-ligand complex contains a protein structure, a ligand, and several active site water molecules. If you haven't already copied this file from the installation subdirectory `maestro-vversion/tutorial`, copy it now from this directory to your working directory.

To import the file into Maestro:

1. Click the Import structures toolbar button to open the Import panel.



2. Ensure that the file format is set to Maestro.
3. Select the file `visualization_exercises.mae` from the Files list.
4. Click Import to import the multi-structure Maestro file.
5. Close the Import panel by clicking the Close button in the bottom right-hand corner.


You're now ready to begin the exercises in this chapter.

1. PDB structures are experimental protein structures. The data in these structures can be incomplete: missing hydrogen atoms, formal charges, and so on. The protein preparation process ensures that all atoms are present and have the correct charges and bonding patterns, and can also do other tasks to "clean up" the structures.

4.2 Displaying and Undisplaying Atoms

One of the most common changes you might want to make to the molecular representation is to display or undisplay structural features. In order to draw attention to certain parts of a molecular structure, it is often helpful to display only a portion of the entire structure. Rather than doing so by deleting other parts of the structure, this can be done using Maestro's "undisplay" and "display" tools. These tools allow you to make visible parts of the structure invisible, and to make invisible parts of the structure visible again.

The Maestro toolbar includes several buttons for displaying and undisplaying atoms, shown below. This section of the tutorial will demonstrate how to use these tools. You will use them again in a later part of this chapter.

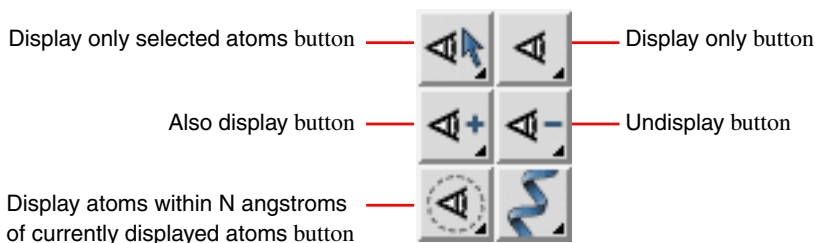

4.2.1 Using the “Display Only” Buttons

There are two buttons in the toolbar that allow you to display only a subset of the atoms in the Workspace. These buttons are the Display only selected atoms button, which allows you to select individual groups of atoms for display, and the Display only button, which allows you to choose predefined structural feature types for display.

This exercise demonstrates how to use these toolbar buttons.

When you imported the file in the previous exercise, it should have been automatically included in the Workspace. If it wasn't, open the Project Table panel and include the entry Protein-ligand complex. See [Section 2.6 on page 12](#) for more information on including and excluding entries.

1. Choose Atoms from the Display only selected atoms menu button (see [Figure 4.1](#)).



(If you don't remember how to do this, see the explanation on [page 7](#).) When you move the pointer into the Workspace, you should see a yellow square with the letter "A" next to it, indicating that you can select atoms in the Workspace.



Figure 4.1. Display only selected atoms *and* Display only **button menus**.

2. Hold down the left mouse button, drag the pointer over any set of atoms in the Workspace (creating a dotted yellow box), and let go of the left mouse button.

The atoms that were not within the yellow box are undisplayed—only the atoms that you selected remain visible.

3. Choose All from the Display only button menu (see [Figure 4.1](#)).



All atoms in the Workspace are displayed again.

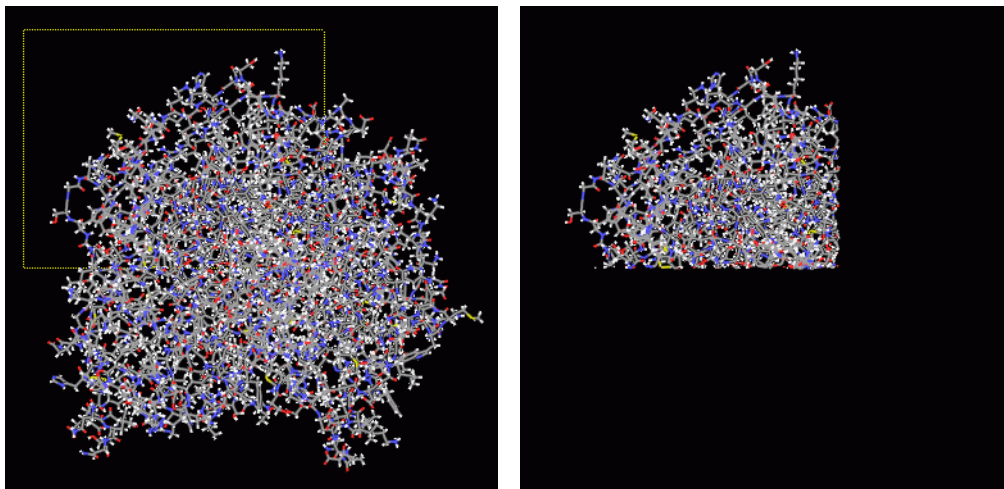


Figure 4.2. Atoms selected for display (left) and the result (right).

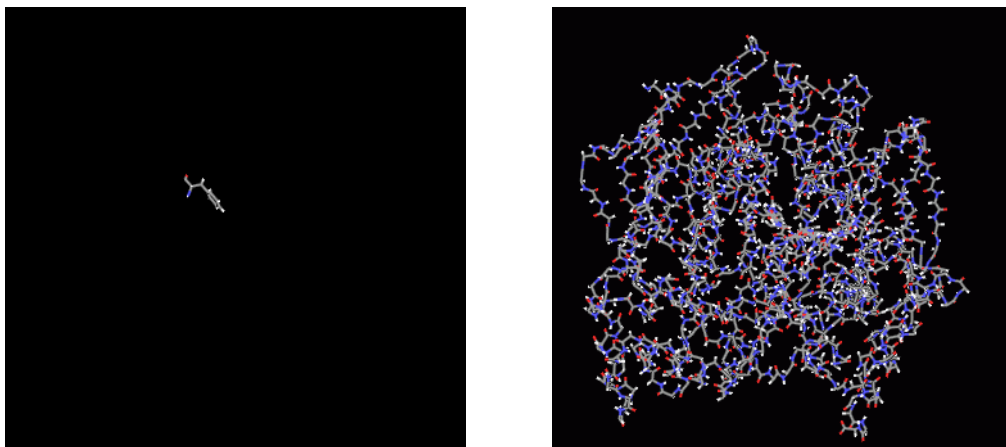


Figure 4.3. Residue displayed (left) and protein backbone displayed (right).

4. Choose Residues from the Display only selected atoms button menu.

In a protein, a “residue” typically refers to a single amino acid subunit, a ligand, a water molecule, or a cofactor. When you move the pointer into the Workspace, you should see a yellow square with the letter “R” next to it, indicating that you can select residues in the Workspace.

5. Click on any one atom in the Workspace.

Maestro displays only the residue that contains the atom you clicked on. All other atoms are undisplayed.

6. Choose Protein Backbone from the Display only button menu.

Now only the protein backbone is visible, and all other atoms are invisible.

7. Choose All from the Display only button menu, to redisplay all atoms.

The toolbar buttons behave similarly when other modes are selected from the button menus. For example, if you choose Molecules from the Display only selected atoms button menu, you can click on any molecule in the Workspace, and all other molecules will be undisplayed.

4.2.2 Using the Other Display Buttons

In the previous section we saw how to display only a certain set of atoms in a single operation. However, Maestro also allows you to display subsets of invisible atoms without undisplaying any currently visible atoms, and to selectively undisplay a subset of currently visible atoms. These actions can be performed using the Also display and Undisplay toolbar buttons.

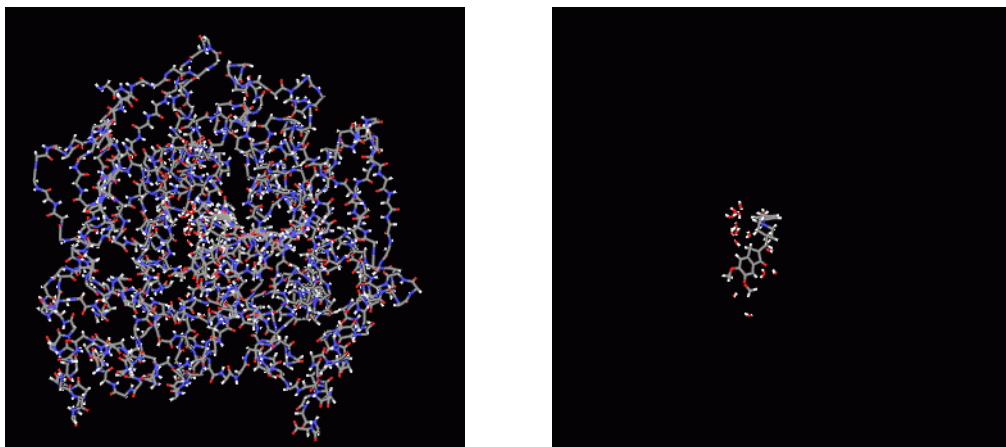


Figure 4.4. Protein side chains undisplayed (left) and entire protein undisplayed (right).

Maestro also allows you to display atoms within a given distance of any currently displayed atoms. This can be done using the Display residues within N angstroms of currently displayed atoms toolbar button (which we'll shorten to Display residues within N angstroms).

This exercise demonstrates how the Also display, Undisplay, and Display residues within N angstroms toolbar buttons can be used to selectively display the structural features that you are interested in.

1. Choose All from the Also display button menu.



This is to make sure that all atoms are displayed. Note that this operation will have no apparent effect if you displayed all atoms at the end of the previous step.

2. Choose Protein Side Chains from the Undisplay button menu.



Observe that the protein backbone, ligand, and water molecules are now visible, but the side chains have disappeared.

3. Choose Protein Backbone from the Undisplay button menu.

The protein backbone has disappeared, and only the water molecules and the ligand remain visible in the Workspace.

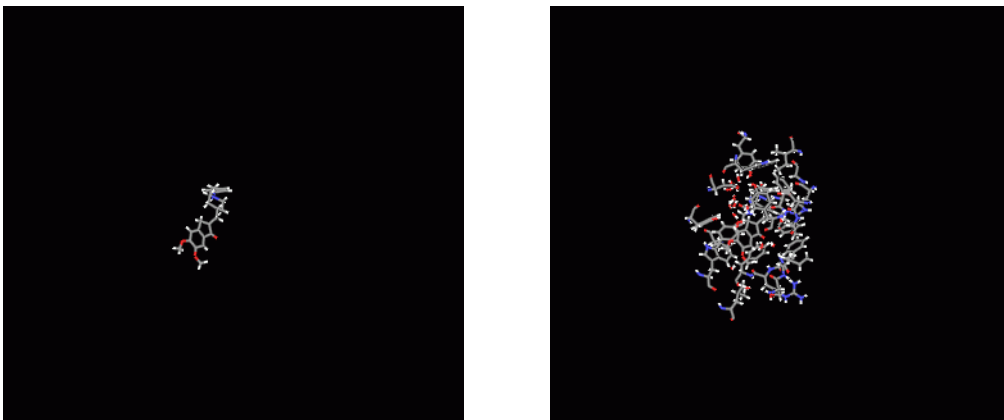


Figure 4.5. Ligand only (left) and ligand with residues within 4Å (right).

4. Choose Waters from the Undisplay button menu, to undisplay the water molecules.

Only the ligand remains visible in the Workspace.

Now that the ligand is the only molecule displayed in the Workspace, we can use the Display residues within N angstroms toolbar button to display the active site features.

5. Choose +4 Angstroms from the Display residues within N angstroms button menu



Residues that have any atoms within 4 Å of the currently displayed ligand are now visible. These residues are all part of the active site, and include some water molecules.

6. To display the remaining water molecules without undisplaying any of the currently visible atoms, choose Waters from the Also display button menu.

Observe that the ligand, active site residues, and waters are all visible. If you want to enlarge the view of the structure, click the Fit to screen button.



7. To display the protein backbone without undisplaying any of the currently visible atoms, choose Protein Backbone from the Also display button menu.

Observe that the protein backbone, ligand, active site residues, and water molecules are all visible. To view the entire structure, click the Fit to screen button again.

Now that you know how to operate the various tools that control atom visibility, you can apply the same basic steps when working with structures of your own. Keep in mind that there may be more than one way to display or undisplay any given set of atoms—feel free to use whatever method works best for you. We'll be using many of these controls again in a later exercise.

You may have noticed that most of the menus for the toolbar buttons described in the above exercises all had an option called **Select**. When you choose this option, Maestro opens the Atom Selection dialog box. This tool allows you to select combinations of atoms that are a bit more complicated than those in the tutorial examples. To learn more about the Atom Selection dialog box, see [Section 5.3](#) of the *Maestro User Manual*.

4.3 Changing the Atom Coloring Scheme


When working in a three-dimensional modeling environment, you can use color to show information about a structure. In most cases, Maestro uses an element-based coloring scheme by default. In this exercise, we'll use atom coloring schemes to visualize information about our protein structure. Take a moment to look over the protein structure in the Workspace. Right now, the color scheme represents various elements in the Workspace, but it can also tell you about calculated and experimental structural properties.

1. Choose Atom partial charge from the Color Scheme button menu.



The atoms are colored by their partial charge. The correspondence of the colors to the partial charges is shown in the table below. The partial charge is an estimate of the “real” charge on an atom, as opposed to the “formal” charge, which is always an integer and is derived from valency rules. These partial charges came from a MacroModel calculation.

Color	Charge Range
blue	$\geq +0.25$
aquamarine	+0.15 to +0.25
green	+0.05 to +0.15
white	-0.05 to +0.05
yellow	-0.15 to -0.05
orange	-0.25 to -0.15
red	≤ -0.25

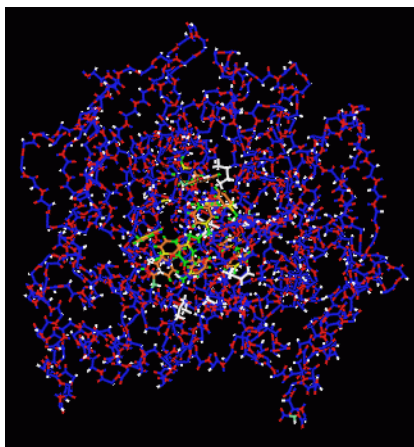


Another property that can be visualized using color is the “*B* factor” that was experimentally determined when the protein structure was resolved. The *B* factor is a measure of how mobile a particular region of a protein is. Side chains and loops, for example, tend to have higher *B* factors while the relatively rigid helix and sheet regions of the protein backbone tend to have lower *B* factors.

2. To view the *B* factors, choose Atom PDB B Factor (Temperature Factor) from the Color Scheme button menu.

The color scheme used by Maestro is shown in the table below.


Color	B Factor Range
gray	0.0
blue	0.001 – 20.0
green	20.001 – 40.0
yellow	40.001 – 60.0
orange	60.001 – 80.0
red	> 80.0

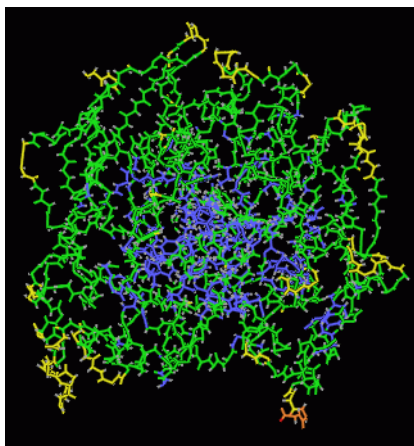

Color schemes can be useful for identifying parts of the displayed structure. In a big molecule like this protein-ligand complex, it's hard to tell what's the protein, what's the ligand, and what's something else. Each of these pieces is usually a separate molecule, so you can color them by molecule number to tell which is which.

3. Choose Molecule number from the Color Scheme button menu.
4. To return to the element-based coloring scheme, choose Element from the Color Scheme button menu.

Using the Atom Coloring panel, it's possible to apply custom color schemes to structures. For more information on using the Atom Coloring panel, see [Section 6.1](#) of the *Maestro User Manual*.

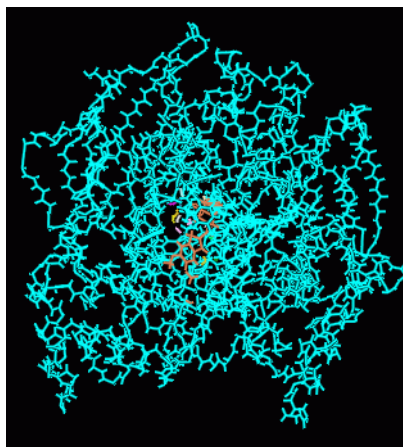

4.4 Changing the Appearance of Atoms and Bonds

Prior to this section, the tutorial operations were all performed using a simple “wire frame” representation for molecular structures. In this section, you’ll learn how to represent atoms and bonds as spheres and tubes. Not only can you make it easier to view the structure by changing the way that atoms and bonds appear, but you can also use various representation modes to differentiate substructures in the Maestro Workspace.

The Maestro toolbar includes several buttons used to control the representation of atoms and bonds in the Workspace. Structures can be represented as wires, space-filling CPK atoms, balls and sticks, or tubes.



To begin this exercise, you will need to ensure that you can see the ligand and the water molecules clearly. Follow the instructions below to set up the Workspace for this exercise.

1. Zoom in on the ligand until it occupies a large part of the Workspace. (Drag horizontally with the middle and right mouse buttons.)
2. Right-click on an atom in the middle of the ligand, to center the ligand in the Workspace.

Now we’ll start changing the representation. First, we’ll display the structure in Ball & Stick representation instead of wire frame, and then convert the protein back to wire frame.

3. Double-click the Draw atoms in Ball & Stick button menu ([page 54](#), upper left).





Double-clicking applies the action of the button to all atoms in the Workspace. You can use a double-click for all four representation buttons.

4. Choose Molecule from the Draw bonds in wire button menu, and click on a protein atom.

Observe that both the protein backbone and active site residues are displayed as wires, but the ligand and the water molecule remain unchanged ([page 54](#), upper right). The ligand and water molecules are now much easier to distinguish from the protein.

Next, we’ll change the representation of the ligand and the residues in the active site.

5. To represent the ligand using space-filling CPK atoms, choose Molecules from the Draw atoms in CPK button menu, and click on any atom in the ligand ([page 54](#), lower left).
6. To represent active site residues as tubes, choose Residues from the Draw bonds in tube button menu, and click on an atom in each active site residue ([page 54](#), lower right).

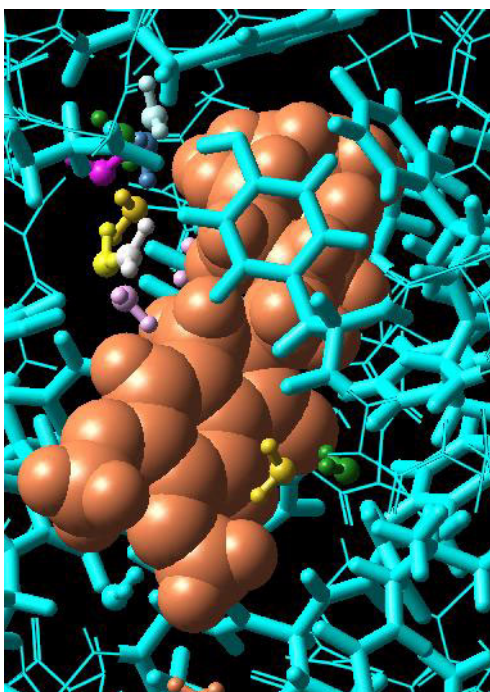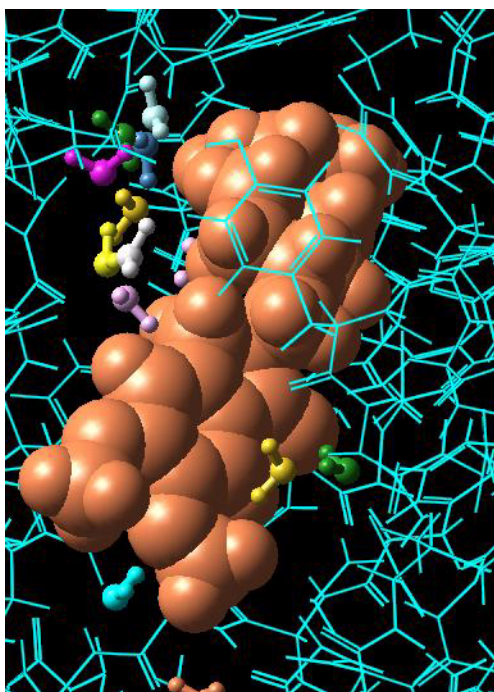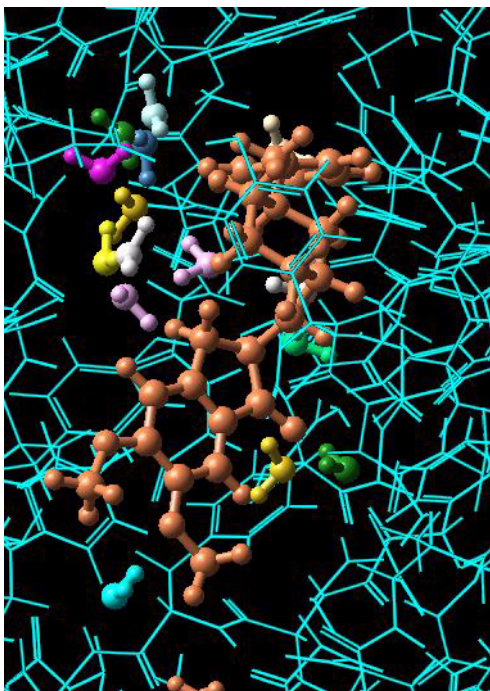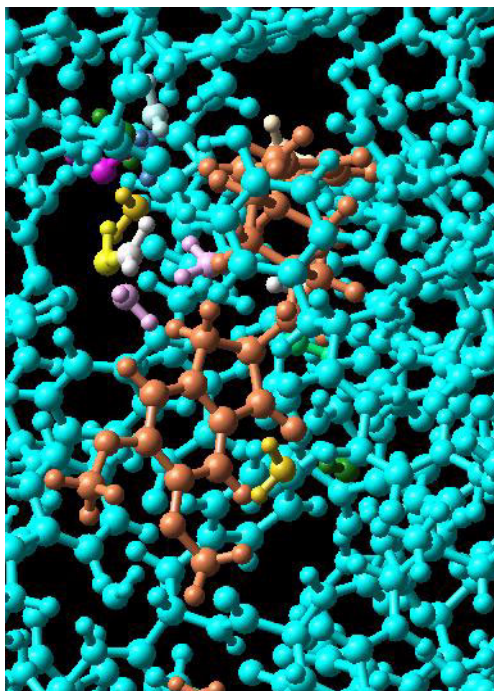

Beyond the options in the Maestro toolbar, it is also possible to customize settings such as the width of tubes, the size of spheres, and so on. For more information on changing these settings, refer to [Section 6.2](#) of the *Maestro User Manual*.

4.5 Displaying the Protein Backbone as a Ribbon

In addition to representing atoms and bonds using 3D features, Maestro can also represent protein backbones using three-dimensional “ribbons”. Ribbons can be turned on or off using the Show, hide, or color ribbons button on the Maestro toolbar:



In this exercise, you’ll represent the protein backbone using a ribbon.

1. Click the Fit to screen toolbar button, so you can see all of the protein.



2. Choose Show Ribbons from the Show, hide, or color ribbons button menu.

The protein is now represented by ribbons that trace the protein backbone instead of atoms, but the ligand and the water molecules remain in their previous representation.

3. Choose Show Ribbons and Atoms from the Show, hide, or color ribbons button menu.

The protein atoms reappear, with the ribbons displayed as well. Observe that the ribbon follows the protein backbone.

4. Choose Hide Ribbons from the Show, hide, or color ribbons button menu.

The ribbons disappear, leaving the atoms on display.

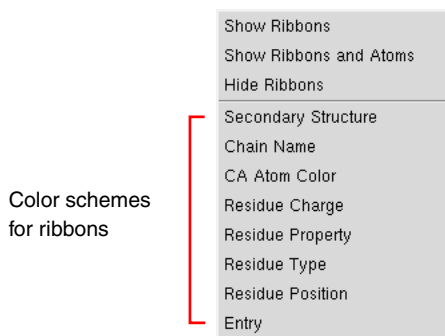


Figure 4.6. Show, hide, or color ribbons *button menu*.

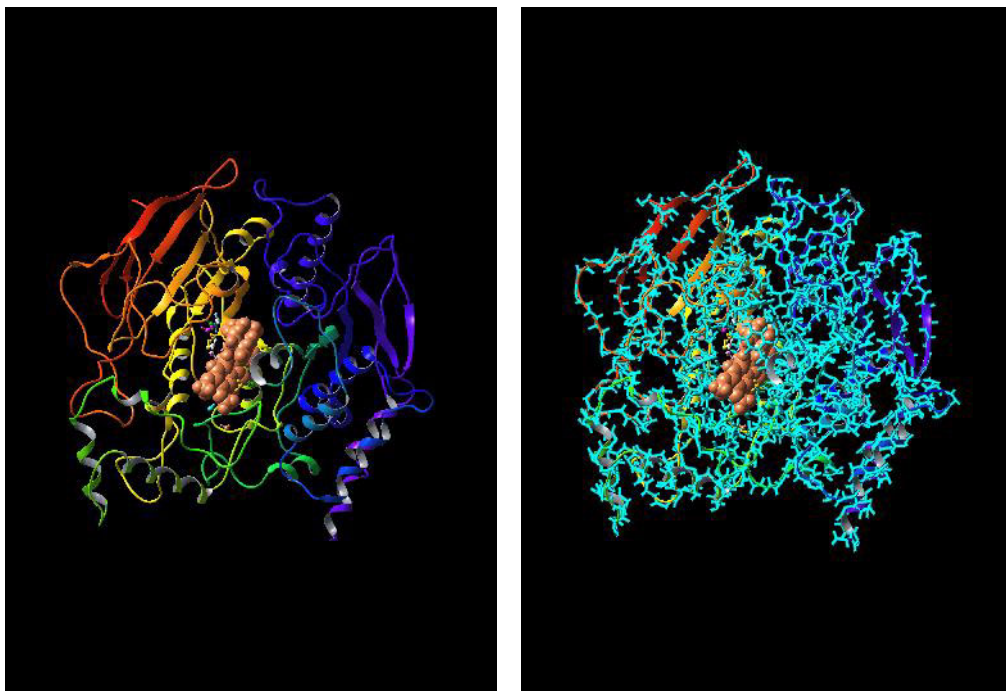


Figure 4.7. Ribbon representation of protein with atoms (right) and without (left).

The default coloring scheme for ribbons is to color according to residue position. As you follow the ribbon along, you'll see the colors progress through the spectrum as the residue number increases. The menu on the ribbons button includes alternate coloring schemes that you can select.

Once you have experience with Maestro's ribbon rendering options, you may wish to further customize the look and feel of protein ribbons. For more information on doing so, see [Section 6.2.4](#) of the *Maestro User Manual*.

4.6 Creating and Displaying a Molecular Surface

In addition to custom rendering for higher-level structural features like protein backbones, Maestro can also represent the “surfaces” that are created by groups of molecules. Surfaces can represent the volume of a molecule as approximated by van der Waals radii. Surfaces can also display properties such as the electron density created by programs like Jaguar. In this example, we’ll create and visualize a “molecular surface,” which is a smooth, contiguous surface that represents the solvent-accessible area presented by a structure. We’ll use the protein from the previous exercise, so leave the structure in the Workspace.

1. Choose Molecular Surface from the Surfaces submenu of the Display menu on the Maestro menu bar.

The Molecular Surface panel appears.

2. Delete the text in the Name text box at the top of the panel and type Tutorial Surface.

This text box contains the text Molecular Surface by default.

Below the Name text box is an area titled Atoms for surface display. This is used to define the structure for which Maestro creates a surface. For this exercise we’ll select the protein structure only:

3. Choose Molecules from the Pick option menu.

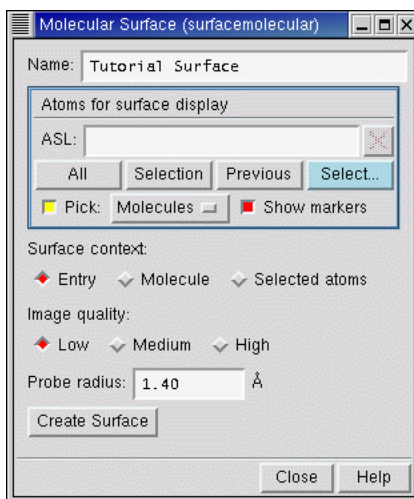


Figure 4.8. The Molecular Surface panel.

- Click on any atom in the protein to select the entire protein.

The wire bonds in the protein are temporarily colored a purple color, to show that they have been selected for the surface. Observe that the text in the ASL text box is automatically filled in once you click on the protein. ASL stands for “Atom Specification Language”, and is a way of describing which atoms have been selected. If you want to be more selective in the atoms you choose, you can click the **Select** button, which opens the Atom Selection dialog box. For more information on using this tool, see [Section 5.3](#) of the *Maestro User Manual*.

The rest of the settings can be left at their default values. For large molecules such as proteins, the default image quality setting of Low is usually appropriate. For small molecules, higher quality settings may be useful.

- Click **Create Surface** to generate the molecular surface.

Once the surface is created, the **Surface Table** panel automatically appears. Notice that you can see inside the surface where the ligand atoms are. This is because the surface doesn’t cover these atoms: we chose to put a surface on only the protein.

- Click **Close** to close the Molecular Surface panel.

Next, we’ll experiment with the color and transparency of this surface. These properties can be changed in the **Display Options** dialog box.

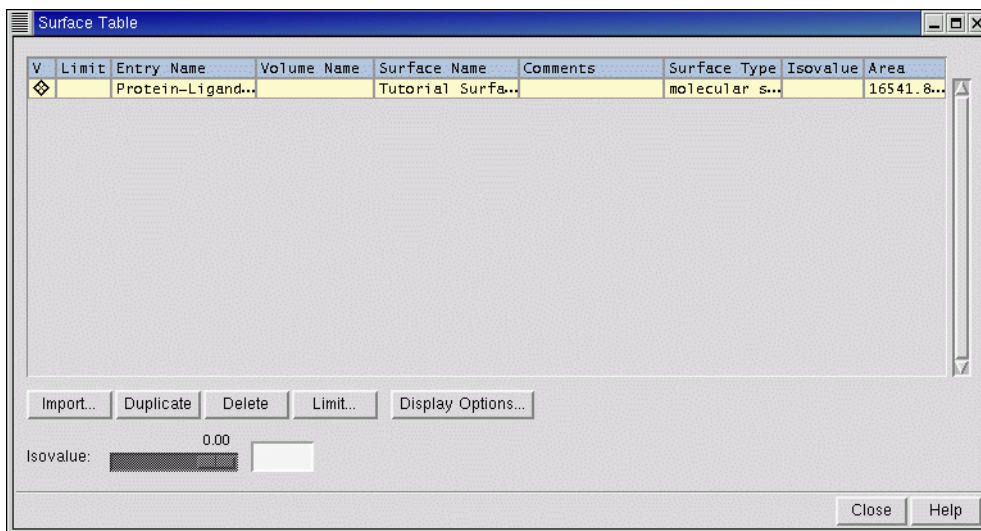


Figure 4.9. The Surface Table *panel*.

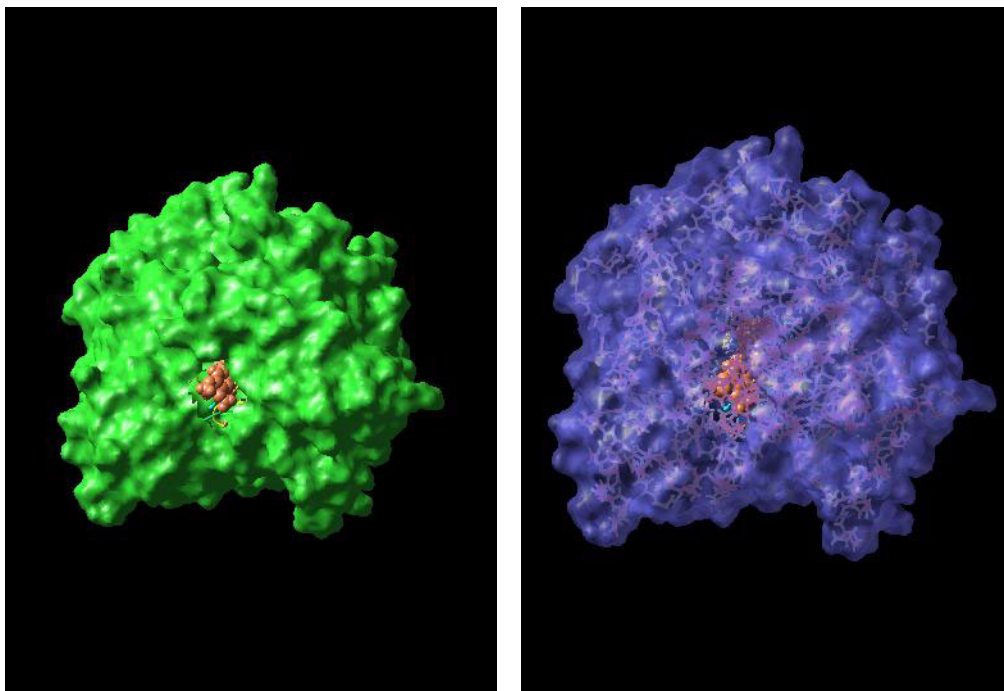


Figure 4.10. Molecular surface, green at 100% transparency (left) and light blue at 50% transparency (right).

7. Click the Display Options button near the bottom of the Surface Table panel.

The Display Options dialog box appears.

8. Use the Transparency slider to set the transparency to 50 (fifty percent transparent).
9. Select Light Blue from the Color option menu.
10. Click OK to accept and view the changes (see [Figure 4.10](#)).

The dialog box closes and the changes are applied. If you want to experiment more with the display options, click Apply instead. The dialog box remains open, and you can make whatever other changes you would like to make.

11. Once you've viewed the surface, control-click the diamond in the leftmost column of the Surface Table to undisplay the surface.
12. Close the Surface Table panel by clicking the Close button in the bottom right-hand corner.

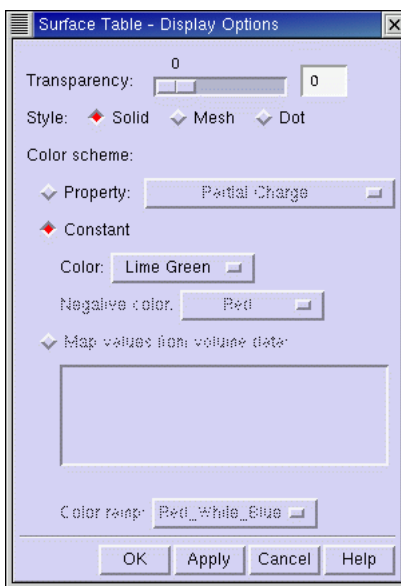


Figure 4.11. The Surface Table - Display Options *dialog box*.

Maestro is also capable of generating several other surface types. The interfaces used to create these surfaces are generally similar to the Molecular Surfaces panel, and all settings for the generated surfaces are controlled via the Surface Table panel. For more information on creating and viewing surfaces using Maestro, see [Chapter 11](#) of the *Maestro User Manual*.

Measuring, Analyzing, and Superimposing Structures

Once you're satisfied with the way that your structures appear in the Workspace, you can go on to examine the structures and their properties. For example, you can examine the geometric differences between two conformations of a molecule by superimposing them, or you can have Maestro display the hydrogen-bonding interactions between a ligand and a receptor. In this chapter, we'll start with basic operations like measuring distances and angles before performing more complicated structural analyses.

For more information on the full extent of Maestro's analysis capabilities, see [Chapter 9](#) of the *Maestro User Manual*.

5.1 Getting Ready for the Exercises in this Chapter

To perform the exercises in this chapter, you'll need to import the structure file `analysis_exercises.mae`. This is a multi-structure file that includes several small molecules, a set of conformations created by MacroModel, and a ligand that has been docked to a protein structure using Glide. If you haven't already copied this file from `$SCHRODINGER/maestro-vversion/tutorial`, copy it now from this directory to your working directory.

To import the file into Maestro:

1. Click the Import structures toolbar button to open the Import panel.
2. Ensure that the file format is set to Maestro.
3. Select the file `analysis_exercises.mae` from the Files list.
4. Click Import to import the multi-structure Maestro file.
5. Close the Import panel by clicking the Close button in the bottom right-hand corner.
6. If the Project Table panel isn't open, open it by clicking the Open/Close Project Table toolbar button.



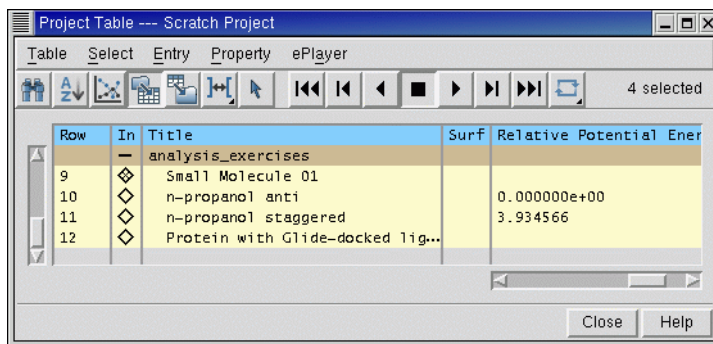


Figure 5.1. The Project Table panel with the new entries.

In the Project Table panel, you'll notice that several new rows have been added. These rows represent new entries in the Project Table. Some of these entries contain only one molecule, while other entries contain multiple molecules.

The Project Table has an easy-to-use set of controls that allow you to record structure-related properties, present the results of calculations, and include and exclude entries in the Workspace. You've already done some of these operations in earlier chapters of this tutorial. The exercises in this chapter will walk you through any necessary Project Table manipulations, but for more information on the Project Table, see [Chapter 8](#) of the *Maestro User Manual*.

We'll be using the Project Table quite a bit in these exercises, so we're going to leave it open. If it gets in the way, you can move it, or you can close it and re-open it when you need to use it again. To close the Project Table panel, click the Open/Close Project Table button in the toolbar (it should be indented), or click the Close button in the bottom right-hand corner of the Project Table panel.

You're now ready to begin the exercises in this chapter.

5.2 Measuring Distances

One simple analysis task is to determine the distance between two atoms. These can be any two atoms, bonded or non-bonded. In this example, we'll use the Measure distances, angles, or dihedrals button on Maestro's toolbar to determine bonded and non-bonded interatomic distances on a structure that has been optimized using Jaguar.

To begin, we'll make sure that the appropriate Project Table entry is included in the Workspace.

1. In the Project Table, find the entry named Small Molecule 01.

2. Click the diamond in the In column for the entry to ensure that this entry is included in the Workspace and all others are excluded.
3. Choose Distance from the Measure distances, angles, or dihedrals button menu in the main toolbar.



(If you don't remember how to do this or are starting with this chapter, see the explanation of how to use button menus on [page 7](#).)

4. Click on any atom in the Workspace.

A purple cube appears over the atom, to indicate that you have selected it as the starting point for your measurement.

5. Click on an atom that is bonded to the atom you've just selected.

A dashed marker appears between the atoms, and the interatomic distance (measured in angstroms) is displayed over the marker. If you rotate the molecule in the Workspace (by holding down the middle mouse button and moving the mouse), you'll see that the measurement marker moves with the structure.

Next, we'll measure a distance between two non-bonded atoms. Because the measurements tool is still active (the Measure distances, angles, or dihedrals button is indented to show this), with Distance selected, you don't need to repeat the process of selecting Distance from the toolbar button menu.

6. Click on the same atom that you clicked on in [Step 4](#).

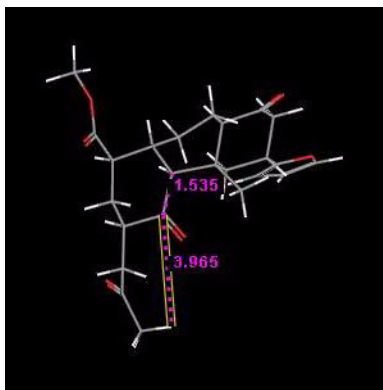


Figure 5.2. Distance measurements.

Note that the purple box appears again, indicating that you've selected this atom as the starting point for your measurement.

7. Now click on any atom in the Workspace that's not bonded to the atom that you clicked on in the previous step.

You'll see that a new distance marker appears between the two atoms that you just clicked on. You may also have noticed that the yellow lines disappeared from the original measurement when you made the new measurement. This is because the yellow marker on both sides of the dashed line simply indicates that it was the most recently measured distance.

Finally, after measuring the distances, we can remove the markers from the Workspace.

8. Choose **Delete Measurements** from the **Measure distances, angles, or dihedrals** button menu.

The structure remains in the Workspace unchanged, and the distance markers disappear.

5.3 Measuring Angles

Another easily performed analysis task is angle measurement. In this example, we'll use the structure from the previous exercise to see that angles can be measured in much the same way that distances can be measured.

1. Choose **Angle** from the **Measure distances, angles, or dihedrals** button menu.



2. Click the oxygen of the carbonyl near the middle of the molecule.

You'll notice that a purple cube appears to indicate that you've selected an atom to begin an angle measurement.

3. Click the carbonyl carbon atom.

You'll see a second purple cube appear, indicating that you've selected the clicked atom to be part of your angle measurement. This second atom clicked will always be the vertex of the angle that is measured.

4. Finally, click on a carbon atom that is bonded to the atom you picked in the last step.

After you click on the third atom, the purple cubes vanish, and a green angle marker appears.

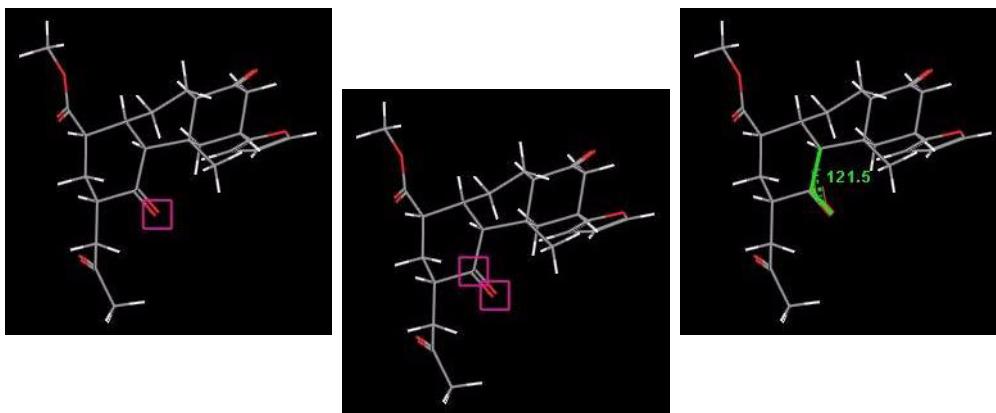


Figure 5.3. Measuring an angle.

If you like, you can define several other angle measurements on the structure. Just like distance measurements, multiple angle measurements can appear in the Workspace at any given time, and you can measure angles between non-bonded atoms.

5. Finally, delete the angle markers from the Workspace by choosing Delete Measurements from the Measure distances, angles, or dihedrals button menu.

The angle markers disappear, and the structure remains unchanged in the Workspace.

5.4 Measuring Dihedrals

Once you know how to measure distances and angles, it's a simple matter to measure a dihedral angle, the torsion defined by any four atoms. For this example, we'll measure the dihedral angles in two different conformations of n-propanol. We'll begin this exercise by placing one of the n-propanol structures in the Workspace.

1. Locate the entry n-propanol anti in the Project Table.
2. Include the n-propanol anti entry in the Workspace and exclude all other entries, by clicking the diamond to the left of the title.
3. Choose Dihedral from the Measure distances, angles, or dihedrals button menu.
4. Locate the terminal methyl group in the molecule and click on the carbon atom in that group.

A purple cube appears, indicating that you've selected this atom to be used as part of the dihedral measurement.

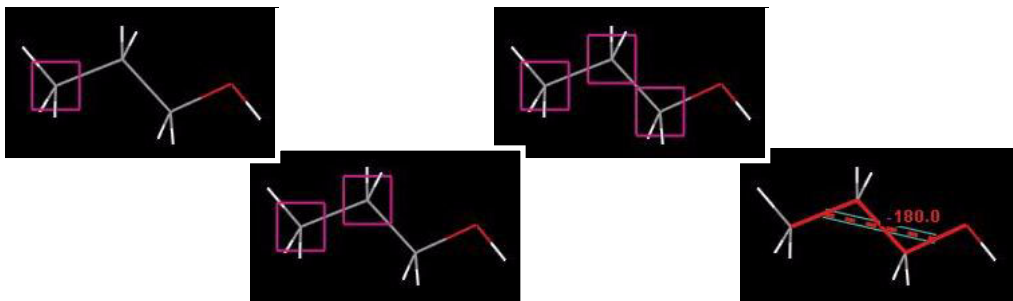


Figure 5.4. Measuring a dihedral angle.

5. Next, click on the carbon atom that is bonded to the one you clicked on in the previous step.

After doing so, another purple cube appears.

6. Click on the one remaining carbon atom in the structure.

A third purple cube appears.

7. Finally, click on the red oxygen atom.

A marker is superimposed over the structure, and the measurement (180°) is displayed over the marker.

After measuring the dihedral angle of the “anti” conformation, we’ll include the “staggered” conformation in the Workspace and measure the corresponding dihedral angle.

8. Locate the entry n-propanol staggered in the Project Table.
9. Control-click (click while holding down the CTRL key) the diamond next to the entry name to include this entry in the Workspace.

You’ll notice that control-clicking makes the “staggered” entry appear in the Workspace next to the “anti” entry rather than replacing it. This is because control-clicking has no effect on the inclusion of other entries: it only affects the entry you are working with.

Because the measurement tool is still active (indicated by the indented appearance of the Measure distances, angles, or dihedrals button in the toolbar), with Dihedral selected, you don’t need to click the button again to measure another dihedral angle.

To measure the dihedral, you click on four atoms, just as for the “anti” conformation.

10. Click the carbon atom in the terminal methyl group of the entry that you just added to in the Workspace, then on the carbon atom that’s bonded to it, then on the remaining carbon atom, and finally on the oxygen atom.

A purple cube appears for each of the first three atoms. When you click on the fourth atom, the purple cubes vanish, and a dihedral marker appears. The dihedral measurement (-63.9°) is visible near the marker.

11. After you've examined both dihedral measurements, delete the dihedral markers from the Workspace by choosing Delete Measurements from the Measure distances, angles, or dihedrals button menu.

As you saw, it's possible to view measurements on multiple Project Table entries at the same time. You can display as many different measurements at the same time as you like, and include angle and distance measurements along with dihedral measurements.

5.5 Superimposing Multiple Structures

Since we've already placed multiple conformations of a molecule in the workspace, we're ready to superimpose the two structures. Maestro's Superposition panel allows you to easily align two or more structures and calculate the root-mean-square deviation (RMSD) between atom pairs. This can be useful when examining the results of MacroModel conformation searches, comparing computational results to experimentally determined structures, or any other time you might want to gauge geometric differences between multiple structures. In this exercise, we'll open the Superposition panel and then go on to align the conformations used in the previous exercise.

1. If the two Project Table entries used in the previous exercise are not present in the Workspace, locate the entries in the project table, click on the diamond for n-propanol staggered, and then control-click on the diamond for n-propanol anti.
2. Open the Superposition panel by choosing Superposition from the Tools menu of the Maestro menu bar.

When the Superposition panel appears, you'll see that it has a number of options for defining molecules and substructures that you can align. Near the top of the panel are controls that allow you to pick individual atom pairs that you can align. However, we'll use the Superimpose by ASL section in this example.

3. Select the Create property for included entries option.
4. Click All to perform the RMSD calculation for all atoms.

You'll notice that one of the conformations is aligned on top of the other, and the RMSD value is again printed in the table at the bottom of the panel.

If you scroll the Project Table to the right, you'll see that a property has been added that contains the RMSD value for the staggered conformation.

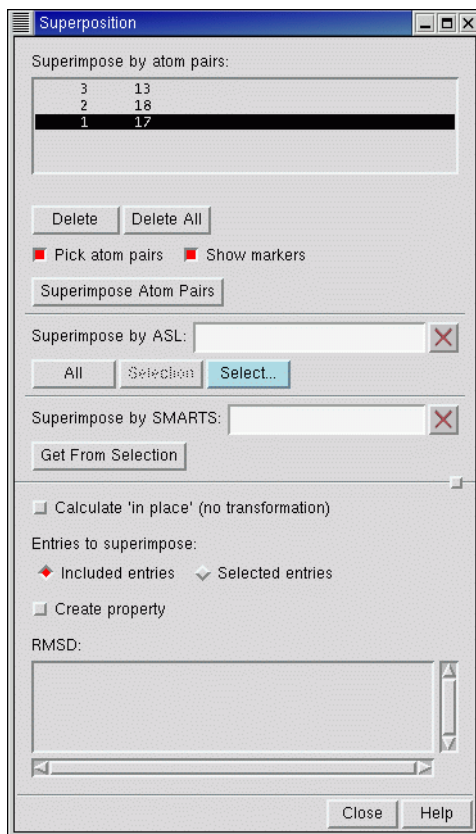


Figure 5.5. The Superposition panel.

Maestro can also calculate a root-mean-square deviation (RMSD) for structures without aligning them. This capability is useful when you want to compare two structures in the context of another structure, such as two ligands docked to a receptor. To calculate the RMSD without aligning the structures, select the Calculate 'in place' (no transformation) option.

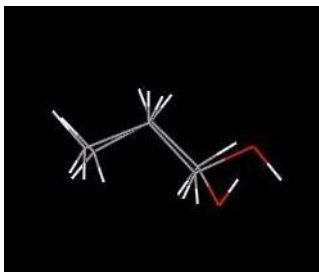


Figure 5.6. The superimposed propanol conformers.

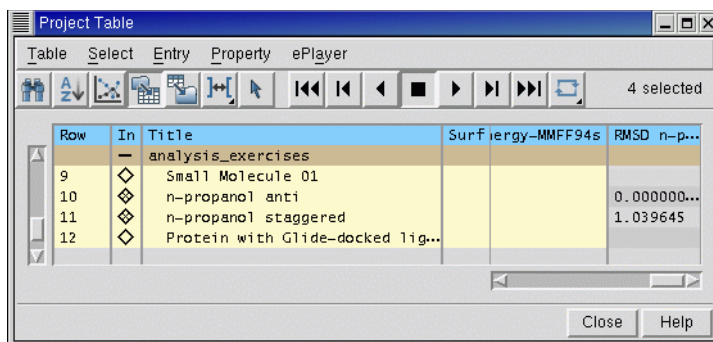


Figure 5.7. The Project Table panel with the RMSD property added.

In addition to the abovementioned tools for superimposing using atom pairs, you may have noticed that the Superposition panel contained a blue button labeled Select. This button opens the Atom Selection dialog box, an interface that allows you to select groups of atoms that can be used as input for most Maestro operations. For more information on the Atom Selection dialog box, see [Section 5.3](#) of the *Maestro User Manual*. For a more thorough description of the Superposition panel, see [Section 9.3](#) of the *Maestro User Manual*.

5.6 Displaying Hydrogen Bonds

In this exercise, we'll use the results of a Glide docking job to demonstrate one of Maestro's useful analysis tools, the ability to display hydrogen bonding interactions. By displaying hydrogen bonding interactions, it's possible to visualize the interactions that make ligand binding energetically favorable. We'll begin this exercise by including the protein-ligand complex in the Maestro Workspace, and then use the toolbar to display and undisplay hydrogen bonding interactions.

1. Clear the Workspace. (Click the Clear Workspace button on the toolbar.)
2. In the Project Table, click on the diamond next to the entry named Protein with Glide-docked ligand.

The protein complex appears in the Workspace. The protein is colored blue and the ligand is colored orange. You might want to zoom in on the ligand and center the ligand in the Workspace. To do this, drag horizontally with the middle and right mouse buttons, then right-click on a ligand atom.

3. To display hydrogen-bonding interactions between the ligand and the receptor, choose Inter H-bonds from the Display H-bonds button menu.



4. Click on any atom in the ligand (colored orange), to display intermolecular hydrogen bonds between that molecule and all other molecules in the Workspace.

You'll see that hydrogen-bonding interactions are marked with dashed yellow lines.

5. To display intramolecular hydrogen bonding interactions, choose Intra H-bonds from the Display H-bonds button menu.

6. Click on any atom in the protein structure (colored blue) to visualize intramolecular hydrogen bonds.

You'll see many dashed yellow markers appear—among these are the hydrogen-bonding interactions that form helices and sheets, resulting in the overall structure of the protein.

7. Undisplay the hydrogen bonds by choosing Delete H-bonds from the Display H-bonds button menu.

8. Finally, close the project table by clicking on the Open/Close Project Table button in the toolbar (it should be indented), or by clicking the Close button in the bottom right-hand corner of the Project Table panel.

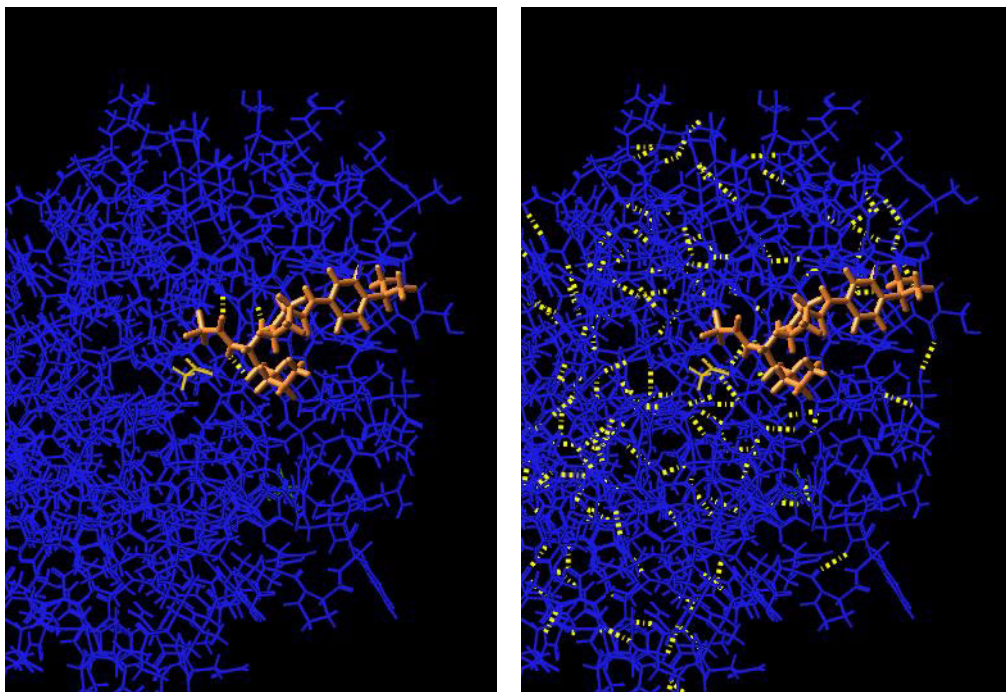


Figure 5.8. Protein-ligand complex showing hydrogen bonds between the ligand and the protein (left) and internal protein hydrogen bonds (right).

Congratulations, you've just completed a quick survey of Maestro's analysis capabilities. In addition to the analysis tasks covered in this chapter, Maestro can also display intermolecular steric contacts, add H-bond counts to the project table, and align protein structures. For more information on Maestro's analysis tasks, see [Chapter 9](#) of the *Maestro User Manual*.

In the next chapter, you'll see an overview of some of Maestro's more advanced visualization and analysis capabilities, and a brief description of other Schrödinger products.

Beyond the Basics

6.1 Advanced Maestro Features

Beyond the Maestro operations surveyed in the previous chapters, there are other advanced tasks that Maestro can perform. This section offers an overview of a few such tasks, explaining when they might come in handy and where to go for more information.

- **Plotting Data:** Using Maestro's plotting tools, it's possible to examine relationships between structure-related data that's saved in the Project Table. For more information on Maestro's plotting tools, see [Chapter 10](#) of the *Maestro User Manual*. For more advanced statistical analysis, the program Strike is integrated into Maestro—see the *Strike User Manual* for details.
- **Viewing vibration animations:** Both Jaguar and MacroModel allow chemists to visualize vibrational modes using Maestro. For more information on calculations, see [Section 4.11](#) of the *Jaguar User Manual*. For information on viewing MacroModel vibration animations, see [Section 20.4](#) of the *MacroModel User Manual*.
- **Viewing electronic structure:** Beyond the molecular structure representations described in this tutorial, Maestro also allows chemists to view the electronic structure of a molecule, as determined by Jaguar, using surfaces. For more information on viewing electronic structure information, see [Section 4.12](#) of the *Jaguar User Manual* and [Section 11.3](#) of the *Maestro User Manual*.
- **Visualizing active site properties:** In addition to the shape-based volumetric surfaces that Maestro can calculate, Maestro can also map regions of space that correspond to active site properties such as hydrophobicity. For more information on mapping active sites with Maestro, see [Section 11.2](#) of the *Maestro User Manual*.
- **Command Scripting:** You can automate repetitive tasks using Maestro's command scripting tools. Because all actions in the Maestro interface are automatically logged by the Command Script Editor, there's no need to learn the actual commands themselves—just perform the desired actions using the Maestro interface, and save the command log for future use. For more information, see [Section 13.3](#) of the *Maestro User Manual*.
- **Python Scripting:** Maestro's support for Python scripting gives users even more power than command scripting. Using Python, it's possible to connect discrete steps in the research pipeline and to create new interfaces that are suited to your exact needs. For

more information on Python scripting with Maestro, see [Section 13.2](#) of the *Maestro User Manual*.

6.2 Programs Accessible from Maestro

As described in the introduction to the first chapter, Maestro is more than a molecular visualization environment—it's also the interface to all of Schrödinger's chemical simulation software. The interfaces to these products, available from the Applications menu on the main Maestro toolbar, include:

- **CombiGlide:** CombiGlide is a structure-based virtual screening program for the design of optimal, focused combinatorial libraries. CombiGlide significantly accelerates lead discovery, and streamlines lead optimization efforts.
- **Epik:** Combining the proven reliability of Hammett and Taft methods with powerful tautomerization tools, Epik is a capable tool for adjusting ligand protonation states.
- **Glide:** Designed to offer the full of range of performance from extremely high throughput to extremely high accuracy, Glide is Schrödinger's ligand docking program. Through the Maestro interface, it's possible to screen millions of compounds against a target receptor, and retrieve and evaluate likely hits.
- **Jaguar:** Emphasizing speed and accuracy, Jaguar is Schrödinger's widely used software for ab initio quantum mechanical (QM) calculations. The Maestro interface to Jaguar allows users to easily set up calculations based on the desired task and level of theory.
- **Liaison:** The Maestro interface for Liaison facilitates the setup and applications of linear interaction (LIA) calculations to predict the binding energies of lead candidates.
- **MacroModel:** A general purpose molecular modeling (MM) program, MacroModel is equipped with a full range of proven computational methodologies. Using the Maestro interface, it's possible to set up and run anything from simple molecular minimizations to active site conformational searches.
- **Phase:** Phase is Schrödinger's full-featured tool for ligand-based drug design, and includes a step-by-step wizard-like interface that intuitively guides researchers through the process of creating a pharmacophore model, creating a 3D database, and screening databases for matches to the pharmacophore model.
- **Prime:** An accurate program for protein structure prediction, Prime boasts an intuitive, easily mastered step-by-step Maestro interface that leads researchers from sequence data to finished protein models.
- **QikProp:** The Maestro QikProp panel allows researchers to perform extremely rapid ADME property predictions, for hundreds of thousands of compounds per hour.

- **QSite:** Using the Maestro interface for QSite, researchers can easily set up and interpret the results of QM/MM calculations, where reactive centers are modeled using Jaguar's accurate ab initio methodologies.
- **SiteMap:** SiteMap characterizes the regions near a protein binding site so that researchers can identify sites that are likely to be good targets for potential new drugs.
- **Strike:** Strike is closely integrated within Maestro, and affords researchers a chemically aware package of statistics tools, facilitating QSAR and statistical modeling.

Getting Help

Schrödinger software is distributed with documentation in PDF format. If the documentation is not installed in `$SCHRODINGER/docs` on a computer that you have access to, you should install it or ask your system administrator to install it.

For help installing and setting up licenses for Schrödinger software and installing documentation, see the *Installation Guide*. For information on running jobs, see the *Job Control Guide*.

Maestro has automatic, context-sensitive help (Auto-Help and Balloon Help, or tooltips), and an online help system. To get help, follow the steps below.

- Check the Auto-Help text box, which is located at the foot of the main window. If help is available for the task you are performing, it is automatically displayed there. Auto-Help contains a single line of information. For more detailed information, use the online help.
- If you want information about a GUI element, such as a button or option, there may be Balloon Help for the item. Pause the cursor over the element. If the Balloon Help does not appear, check that Show Balloon Help is selected in the Help menu of the main window. If there is Balloon Help for the element, it appears within a few seconds.
- For information about a panel or the folder that is displayed in a panel, click the Help button in the panel. The Help panel is opened and a relevant help topic is displayed.
- For other information in the online help, open the Help panel and locate the topic by searching or by category. You can open the Help panel by choosing Help from the Help menu on the main menu bar or by pressing CTRL+H.

To view a list of all available Maestro-related help topics, choose Maestro from the Categories menu of the Categories tab. Double-click a topic title to view the topic.

If you do not find the information you need in the Maestro help system, check the following sources:

- *Maestro User Manual*, for detailed information on using Maestro
- *Maestro Command Reference Manual*, for information on Maestro commands
- *Scripting with Python*, for information on using Python scripts with Maestro
- Frequently Asked Questions pages, at https://www.schrodinger.com/Maestro_FAQ.html

The manuals are also available in PDF format from the Schrödinger [Support Center](#). Information on additions and corrections to the manuals is available from this web page.

If you have questions that are not answered from any of the above sources, contact Schrödinger using the information below.

E-mail: help@schrodinger.com

USPS: Schrödinger, 101 SW Main Street, Suite 1300, Portland, OR 97204

Phone: (503) 299-1150

Fax: (503) 299-4532

WWW: <http://www.schrodinger.com>

FTP: <ftp://ftp.schrodinger.com>

Generally, e-mail correspondence is best because you can send machine output, if necessary. When sending e-mail messages, please include the following information, most of which can be obtained by entering `$SCHRODINGER/machid` at a command prompt:

- All relevant user input and machine output
- Maestro purchaser (company, research institution, or individual)
- Primary Maestro user
- Computer platform type
- Operating system with version number
- Maestro version number
- mmshare version number

Copyright Notices

NCSA HDF5 Software Library and Utilities

Copyright Notice and Statement for NCSA HDF5 (Hierarchical Data Format 5) Software Library and Utilities

NCSA HDF5 (Hierarchical Data Format 5) Software Library and Utilities

Copyright 1998, 1999, 2000, 2001, 2002, 2003, 2004 by the Board of Trustees of the University of Illinois

All rights reserved.

Contributors: National Center for Supercomputing Applications (NCSA) at the University of Illinois at Urbana-Champaign (UIUC), Lawrence Livermore National Laboratory (LLNL), Sandia National Laboratories (SNL), Los Alamos National Laboratory (LANL), Jean-loup Gailly and Mark Adler (gzip library).

Redistribution and use in source and binary forms, with or without modification, are permitted for any purpose (including commercial purposes) provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions, and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions, and the following disclaimer in the documentation and/or materials provided with the distribution.
3. In addition, redistributions of modified forms of the source or binary code must carry prominent notices stating that the original code was changed and the date of the change.
4. All publications or advertising materials mentioning features or use of this software are asked, but not required, to acknowledge that it was developed by the National Center for Supercomputing Applications at the University of Illinois at Urbana-Champaign and to credit the contributors.
5. Neither the name of the University nor the names of the Contributors may be used to endorse or promote products derived from this software without specific prior written permission from the University or the Contributors, as appropriate for the name(s) to be used.

6. THIS SOFTWARE IS PROVIDED BY THE UNIVERSITY AND THE CONTRIBUTORS “AS IS” WITH NO WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED. In no event shall the University or the Contributors be liable for any damages suffered by the users arising out of the use of this software, even if advised of the possibility of such damage.

Portions of HDF5 were developed with support from the University of California, Lawrence Livermore National Laboratory (UC LLNL). The following statement applies to those portions of the product and must be retained in any redistribution of source code, binaries, documentation, and/or accompanying materials:

This work was partially produced at the University of California, Lawrence Livermore National Laboratory (UC LLNL) under contract no. W-7405-ENG-48 (Contract 48) between the U.S. Department of Energy (DOE) and The Regents of the University of California (University) for the operation of UC LLNL.

DISCLAIMER: This work was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor the University of California nor any of their employees, makes any warranty, express or implied, or assumes any liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately-owned rights. Reference herein to any specific commercial products, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or the University of California. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or the University of California, and shall not be used for advertising or product endorsement purposes.

C and C++ Libraries for Parsing PDB Records

The C and C++ libraries for parsing PDB records are Copyright (C) 1989 The Regents of the University of California. All rights reserved.

Redistribution and use in source and binary forms are permitted provided that the above copyright notice and this paragraph are duplicated in all such forms and that any documentation, advertising materials, and other materials related to such distribution and use acknowledge that the software was developed by the University of California, San Francisco. The name of the University may not be used to endorse or promote products derived from this software without specific prior written permission. THIS SOFTWARE IS PROVIDED “AS IS” AND WITHOUT ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

120 West 45th Street
32nd Floor
New York, NY 10036

101 SW Main Street
Suite 1300
Portland, OR 97204

3655 Nobel Drive
Suite 430
San Diego, CA 92122

Dynamostraße 13
68165 Mannheim
Germany

QuatroHouse, Frimley Road
Camberley GU16 7ER
United Kingdom

SCHRÖDINGER.