

Maestro 7.5

User Manual

Copyright © 2006 Schrödinger, LLC. All rights reserved. CombiGlide, Epik, Glide, Impact, Jaguar, Liaison, LigPrep, Maestro, Phase, Prime, QikProp, QikFit, QikSim, QSite, SiteMap, and Strike are trademarks of Schrödinger, LLC.

Schrödinger and MacroModel are registered trademarks of Schrödinger, LLC.

The Visualization Toolkit (VTK) is a copyrighted work (1993-2002) of Ken Martin, Will Schroeder, Bill Lorensen. All rights reserved.

Python is a copyrighted work of the Python Software Foundation. All rights reserved.

The C and C++ libraries for parsing PDB records are a copyrighted work (1989) of the Regents of the University of California. All rights reserved.

The NCSA HDF5 (Hierarchical Data Format 5) Software Library and Utilities is a copyrighted work (1998-2004) of the Board of Trustees of the University of Illinois. All rights reserved.

See the [Copyright Notices](#) for full copyright details.

To the maximum extent permitted by applicable law, this publication is provided “as is” without warranty of any kind. This publication may contain trademarks of other companies.

Please note that any third party programs (“Third Party Programs”) or third party Web sites (“Linked Sites”) referred to in this document may be subject to third party license agreements and fees. Schrödinger, LLC and its affiliates have no responsibility or liability, directly or indirectly, for the Third Party Programs or for the Linked Sites or for any damage or loss alleged to be caused by or in connection with use of or reliance thereon. Any warranties that we make regarding our own products and services do not apply to the Third Party Programs or Linked Sites, or to the interaction between, or interoperability of, our products and services and the Third Party Programs. Referrals and links to Third Party Programs and Linked Sites do not constitute an endorsement of such Third Party Programs or Linked Sites.

Revision A, April 2006

Contents

Document Conventions	xiii
Chapter 1: Maestro Overview	1
1.1 General Interface Design	1
1.1.1 Maestro Windows	1
1.1.2 Mouse Functions.....	2
1.1.3 Keyboard Shortcuts	2
1.2 Maestro Workflow	3
1.3 Maestro Projects.....	3
1.4 Job Launching and Incorporation.....	4
1.5 Citing Maestro in Publications.....	6
Chapter 2: The Maestro Main Window.....	7
2.1 Starting Maestro	7
2.1.1 Starting Maestro on a Windows Host.....	7
2.1.2 Starting Maestro on a UNIX Host	8
2.2 The Main Window	9
2.3 The Menu Bar	11
2.4 The Toolbar	11
2.5 The Status Bar	15
2.6 The Sequence Viewer	15
2.7 Shortcut Menus	17
2.8 Mouse Functions	19
2.9 Shortcut Keys	20
2.10 Undoing Workspace Operations	21
2.10.1 Undoing Operations on the Structures.....	21
2.10.2 Undoing Operations on the View	21
2.11 Ending a Maestro Session.....	22

Chapter 3: Importing and Exporting Structures and Data	23
3.1 Importing Structures	23
3.1.1 Selecting Files.....	24
3.1.2 Selecting Import Settings.....	26
3.1.3 Error Reporting for PDB Files	27
3.1.4 Entry Naming Conventions for Imported Structures	28
3.1.5 Reading Jaguar Input Files	29
3.2 Exporting Structures	30
3.2.1 Selecting Export Settings.....	30
3.2.2 Exporting Multiple Files.....	32
3.3 Exporting Data to a Spreadsheet	33
3.4 Importing Data From a Spreadsheet	34
Chapter 4: Building Structures	37
4.1 The Build Panel	37
4.2 Building a Structure From Fragments	39
4.2.1 Building Structures Using Place Mode.....	41
4.2.2 Building Structures Using Grow Mode	41
4.3 Building a Structure From Atoms	43
4.4 Changing Elements	44
4.5 Changing the Bond Order and Formal Charge	45
4.6 Fusing or Connecting Structures	46
4.7 Adjusting the Geometry	49
4.7.1 Adjusting Parameters.....	49
4.7.2 Cleaning up the Geometry	51
4.8 Applying a Hydrogen Treatment	52
4.9 Changing Atom Properties	53
4.10 Changing Residue Properties	57
4.11 Defining Dummy Atoms	59

4.12 Deleting Atoms	59
Chapter 5: Selecting Atoms	61
5.1 Toolbar Buttons	61
5.2 Picking Tools	61
5.3 The Atom Selection Dialog Box	63
5.3.1 Selecting Atoms by Property	65
5.3.2 Combining and Modifying Atom Selections	65
5.3.3 Selecting Atoms by Proximity	66
5.3.4 Editing and Storing Expressions	67
5.3.5 Examples of Atom Selection	67
Chapter 6: Displaying Structures	69
6.1 Changing Atom Color	69
6.1.1 Applying a Single Color	69
6.1.2 Using Predefined Schemes to Color Atoms	71
6.2 Changing Molecular Representations	72
6.2.1 Changing the Representation of Atoms	72
6.2.2 Changing the Representation of Bonds	73
6.2.3 Changing Representation Attributes	74
6.2.4 Rendering Proteins as Ribbons	76
6.3 Labeling Atoms	78
6.4 Displaying and Undisplaying Atoms	80
Chapter 7: Manipulating Structures	83
7.1 Global Transformations	83
7.2 Local Transformations	84
7.2.1 Selecting Atoms for Transformation	86
7.2.2 Selecting a Rotation Center	86
7.3 Tiling Multiple Entries	86

Chapter 8: Projects.....	87
8.1 Project Operations	87
8.1.1 Creating and Opening Projects.....	87
8.1.2 Saving and Closing Projects	88
8.1.3 Deleting Projects.....	89
8.1.4 The Project Selector Panels.....	89
8.1.5 The Project Summary Panel	90
8.2 Adding Entries to a Project	91
8.2.1 Importing Structures From a File	91
8.2.2 Creating Entries From Workspace Structures.....	91
8.2.3 Incorporating Entries From Job Output.....	92
8.2.4 Merging Entries From Another Project.....	92
8.3 The Project Table Panel.....	93
8.3.1 The Project Table Toolbar	94
8.3.2 The Project Table Menus	95
8.3.3 Mouse Functions in the Project Table	96
8.3.4 Project Table Shortcut Menus	97
8.3.5 Project Table Shortcut Keys	98
8.3.6 Configuring the Project Table.....	99
8.3.7 Finding Text in the Project Table	100
8.4 Selecting Project Entries	100
8.4.1 The Select Menu.....	101
8.4.2 Selecting Entries By Property	101
8.4.3 Selecting Entries Using the Plot Panel	103
8.4.4 Selection Examples	103
8.5 Operating on Selected Entries	104
8.5.1 Including, Excluding, and Fixing Entries	104
8.5.2 Limiting the Number of Entries Displayed.....	105
8.5.3 Renaming Entries	105
8.5.4 Duplicating Entries	105
8.5.5 Merging Multiple Entries	106
8.5.6 Splitting Entries by Molecule.....	106

8.5.7	Deleting Entries.....	106
8.5.8	Moving Entries	106
8.5.9	Sorting the Entries of a Project	107
8.6	The ePlayer	108
8.6.1	ePlayer Modes	108
8.6.2	ePlayer Options.....	109
8.7	Entry Properties.....	110
8.7.1	Creating New Properties	111
8.7.2	Changing Property Values	112
8.7.3	Setting the Number of Decimals in a Property.....	114
8.7.4	Renaming and Deleting Properties	114
8.7.5	Displaying Selected Properties	114
8.7.6	Moving and Resizing Property Columns.....	115
8.7.7	Exporting and Importing Data	116
8.8	Undoing Project Operations	116
Chapter 9: Tools		117
9.1	Displaying Markers	117
9.2	Making Measurements in the Workspace	117
9.2.1	Measuring Distances, Angles, and Dihedrals	117
9.2.2	Displaying Hydrogen Bonds and Contacts.....	119
9.2.3	Displaying H-H NMR Coupling Constants.....	121
9.2.4	Storing the Results of Measurements	121
9.2.5	Defining Dummy Atoms for Measurements	121
9.3	Superimposing Structures	122
9.4	Sets	124
9.4.1	Creating and Deleting Sets	125
9.4.2	Selecting Atoms for a Set.....	125
9.4.3	Reading and Writing Sets	125
9.5	Assigning and Aligning Protein Structures	126
9.6	Assigning Bond Orders.....	126

Chapter 10: Plotting	127
10.1 Creating Plots	127
10.2 Viewing and Manipulating Plots	128
10.2.1 Selecting Plots	129
10.2.2 Panning and Zooming	130
10.2.3 Labeling and Unlabeling Data Points	130
10.2.4 Selecting and Deselecting Project Entries	131
10.2.5 Including and Excluding Project Entries	131
10.3 Editing Plots and Plot Settings	131
10.3.1 Modifying Plot Attributes	132
10.3.2 Operating on Plot Series	133
10.3.3 Changing X and Y Axis Display	134
10.4 Updating Plots to Reflect Project Table Changes	135
10.5 Saving Plot Images	136
Chapter 11: Surfaces	137
11.1 Generating Surfaces	137
11.1.1 Surface Generation Controls	137
11.1.2 Molecular Structure Surfaces	139
11.1.3 Extended Radius Surfaces	139
11.1.4 van der Waals Surfaces	140
11.2 Hydrophobic/philic Maps	141
11.2.1 Background	141
11.2.2 Mapping Algorithm	143
11.2.3 Specifying a Structure to Be Mapped	143
11.2.4 Defining the Mapping Box	143
11.2.5 Running the Job	144
11.2.6 Viewing the Site Map Surface	144
11.3 Importing Surface or Volume Files	145
11.4 The Surface Table Panel	145
11.4.1 Importing, Modifying, and Deleting Surfaces	147

11.4.2 Changing the Appearance of the Surface	148
Chapter 12: Customizing Maestro	151
12.1 Changing and Saving Panel Layout.....	151
12.2 Setting Preferences	152
12.2.1 Changing Font Size for Interface Component Text.....	152
12.2.2 Setting Project Preferences	153
12.2.3 Specifying a Default Working Directory	155
12.2.4 Customizing Atom Label Appearance.....	155
12.2.5 Setting Mouse Preferences.....	156
12.2.6 Setting File Suffix Preferences.....	157
12.2.7 Setting Builder Preferences	157
12.2.8 Setting Feedback Preferences	160
12.3 Changing the Workspace Appearance.....	160
12.3.1 Showing or Hiding Controls and Displays	160
12.3.2 Setting the Background Color	160
12.4 Setting Up 3D Viewing.....	161
12.4.1 Specifying a Stereo Viewing Method	161
12.4.2 Enabling Perspective.....	162
12.4.3 Enabling Fog	163
12.4.4 Depth Cues	165
Chapter 13: Scripting in Maestro	167
13.1 The Scripts Menu	167
13.2 Customization Using Python Scripts	169
13.3 Customization Using Command Scripts.....	170
13.3.1 Command History	171
13.3.2 Building a Command Script	171
13.3.3 Saving a Command Script	172
13.3.4 Opening Existing Command Scripts	172
13.3.5 Running and Stopping Scripts	172
13.3.6 Customizations You Can Perform With a Script.....	172

13.4 Creating Command Aliases	174
13.5 Creating Macros	175
13.5.1 Named Macros.....	175
13.5.2 Function Key Macros	176
Chapter 14: Printing and Saving Workspace Images	177
14.1 Printing an Image (File or Printer)	177
14.2 Creating TIFF and JPEG Image Files	179
Chapter 15: Help	183
15.1 Maestro Online Help	183
15.1.1 Locating a Topic by Searching	183
15.1.2 Locating a Topic by Category.....	184
15.1.3 Searching Within a Topic	184
15.2 Context-Sensitive Help	185
15.3 Manuals	186
15.4 The Schrödinger Web Site	186
15.5 Technical Support	186
Appendix A: The Maestro File Format	189
A.1 Basic File Description	189
A.2 Data Blocks	189
A.3 Compressed Format	191
A.4 Data Item Names	191
A.5 Example Maestro File	191
Appendix B: Atom Types	195
Appendix C: Entry Suffixes	201

Appendix D: Utilities.....	203
D.1 Structure Conversion	203
D.1.1 Conversions to and From MacroModel Format	203
D.1.2 Conversions to and From Mol2 Format	204
D.1.3 Conversions to and From PDB Format: pdbconvert	205
D.1.4 Conversions to and From SD Format: sdconvert.....	207
D.2 Structure Extraction	208
D.3 Display of Properties	210
D.4 Structure Preparation	211
D.5 Unique Names	212
Appendix E: Customizing Colors.....	215
E.1 Color Definitions	215
E.2 Color Schemes	216
E.3 Color Ramps	217
Appendix F: Setting Environment Variables on Windows	219
Copyright Notices	221
Glossary.....	223
Index.....	225

Document Conventions

In addition to the use of italics for names of documents, the font conventions that are used in this document are summarized in the table below.

Font	Example	Use
Sans serif	Project Table	Names of GUI features, such as panels, menus, menu items, buttons, and labels
Monospace	<code>\$\$SCHRODINGER/maestro</code>	File names, directory names, commands, environment variables, and screen output
Italic	<i>filename</i>	Text that the user must replace with a value
Sans serif uppercase	CTRL+H	Keyboard keys

In descriptions of command syntax, the following UNIX conventions are used: braces { } enclose a choice of required items, square brackets [] enclose optional items, and the bar symbol | separates items in a list from which one item must be chosen. Lines of command syntax that wrap should be interpreted as a single command.

File names and paths are generally represented with the UNIX conventions. If you are running on Windows, replace the forward slash character / with the backslash character \.

In this document, to *type* text means to type the required text in the specified location, and to *enter* text means to type the required text, then press the ENTER key.

References to literature sources are given in square brackets, like this: [10].

Maestro Overview

Maestro™ is the graphical user interface (GUI) for all of Schrödinger's computational programs: CombiGlide™, Epik™, Glide™, Impact™, Jaguar™, Liaison™, LigPrep™, MacroModel®, Phase™, Prime™, QikProp™, QSite™, SiteMap™, and Strike™. It contains tools for building, displaying, and manipulating chemical structures; for organizing, loading and storing these structures and associated data; and for setting up, submitting, monitoring, and visualizing the results of calculations on these structures. Maestro's Job Control facility manages jobs submitted from Maestro and from the command line to both local and remote hosts. Maestro can be run locally and submit jobs to any host that you have access to.

The Maestro interface uses the OpenGL graphics tools, and can take advantage of hardware graphics capabilities, including stereo viewing capabilities. Maestro runs on SGI and Linux platforms. For information on Maestro requirements, software and hardware, see the *Installation Guide*. For information on configuring stereo viewing, see [Section 12.4 on page 161](#).

This manual contains an introduction to the Maestro graphical user interface (GUI) and a description of how to use Maestro's settings, panels, and features to build, import, and manipulate molecular structures. No prior knowledge of Maestro is assumed. If you cannot find the information you are looking for in this document, see the Maestro online help. For help with preparing and starting computations, see the user manual for the related product.

This chapter provides an overview of Maestro—the general interface design, Maestro projects, Maestro workflow, and running jobs from Maestro.

1.1 General Interface Design

The Maestro interface design follows most common interface conventions, but it has some characteristics that are unique. The general operation of the Maestro interface is described in this section. Specific features are described in the relevant sections or chapters.

1.1.1 Maestro Windows

Maestro is based on a main window, from which you can open a range of secondary windows for the performance of various tasks. The main window is described in detail in [Chapter 2](#).

Most of the windows that are opened from Maestro are amodal, and are called *panels*. More than one panel can be open at any given time, and a panel need not be closed for an action to be carried out.

Maestro also has dialog boxes, most of which are modal: you must click OK to carry out the action and you cannot open any of the panels while the dialog box is open. There are also some dialog boxes that allow you to carry out limited operations in the Workspace or in a panel, and allow you to open the help viewer.

1.1.2 Mouse Functions

Maestro supports common mouse functions:

- The left button is used for selecting: choosing menu items, clicking buttons and selecting objects. This button is also used for resizing and moving panels. In the description of mouse actions, “click” always means left-click.
- The right button is used for opening a shortcut (context-sensitive) menu, where these menus are available.
- Shift-click is used to select a contiguous range of items in a list, and control-click is used to select and deselect a single item in a list without affecting the selection of other items.
- Dragging operations are supported. For instance, in the Workspace, dragging selects multiple objects. In the Project Table, dragging selected entries allows you to reposition the entries; dragging a column heading moves the column; dragging with the middle mouse button on the boundary of a row or column resizes the row or column.

There are also specialized mouse functions in the Workspace, which are described in [Section 2.8 on page 19](#).

If you have the handedness on your mouse set to “left”, the mouse functions are the mirror image of those described: the right mouse button is used for picking, and the left button is used for shortcut menus.

If you have a two-button mouse, ensure that it is configured for three-button mouse simulation. Then the middle mouse button is simulated by pressing or holding down both buttons.

1.1.3 Keyboard Shortcuts

Maestro supports keyboard shortcuts in some of its panels, including the main window. These shortcuts are documented with the relevant panel. The shortcut keys do not work if you have Num Lock, Caps Lock, or Scroll Lock enabled.

1.2 Maestro Workflow

The main Maestro workflow¹ involves performing actions on the displayed structures, to change the appearance of the structures or to change the content of the structures. In many cases, you will want to change only a part of the structure, which involves selecting the atoms that you want to change.

Maestro has two ways of performing an action on a structure or a part of a structure:

- Choose an action first and then select the atoms to apply the action to.

Actions are usually chosen either from a toolbar button menu or by opening a panel and selecting options.

- Select atoms first and then choose an action to perform on them.

Actions can be performed by opening a panel, selecting options and clicking the Selection button, or by choosing an item from the selection shortcut menu. See [Section 2.7 on page 17](#) for more details on this shortcut menu.

Maestro's atom selection tools provide you with a range of options for selecting single atoms, entire structural units, or complex combinations of atom sets. These tools are described in detail in [Chapter 5](#).

1.3 Maestro Projects

When you use Maestro, you are always working within a *project*. A project is a collection of chemical structures and their associated data. These structures and their data are organized into *entries*, each of which can consist of multiple molecules and their properties. The project is represented in the Project Table, which displays an ordered list of entries and any associated data. You can open the Project Table panel by choosing Show Table from the Project menu, or by clicking the Open/Close Project Table button on the toolbar.



If you do not specify a project when you start Maestro, Maestro creates a *scratch* project. You can work in a scratch project, but you must save it in order to use it in future sessions.

Entries are represented by rows in the Project Table. Each row contains the row number, the title, the entry's Workspace inclusion state (the In column), a button to open the Surfaces panel if there are surfaces associated with the entry, the entry name, and any properties associated

1. This is to be distinguished from the workflow for a given product, which usually involves setting up and running jobs, and which is supported by Maestro panels.

with the entry. If there are no surfaces associated with any entry, the Surf column of the Project Table is empty.

You control which entries are displayed in the Workspace from the Project Table. You can use entries as input for most of the computational programs. You can select entries as input for the ePlayer, which displays the selected structures in sequence. You can duplicate, combine, rename, split, and sort entries; create properties for entries, import structures as entries, and export entries in various formats. Although the Project Table does not have spreadsheet capabilities, you can export properties to and import properties from a spreadsheet.

When you build molecules in the Workspace, they constitute a *scratch* entry until you choose to save the structures as project entries. The scratch entry is not saved with the project unless you explicitly add it to the project. However, you can use a scratch entry as input for some jobs, and the results can be incorporated.

For products that have wizard-based interfaces (CombiGlide, Phase, and Prime), data is stored inside the project, grouped as entities called “runs”. Output from each step is incorporated into the next step and stored in the project, but not via the Project Table. The data in the runs can include structural data, but do not make use of the Project Table directly. Instead, structures can be copied to and from the Project Table at various points in the run.

More information on projects is given in [Chapter 8](#).

1.4 Job Launching and Incorporation

Maestro’s job launching and incorporation capabilities are designed to make it easy to manage multiple structures and their passage through multiple computational programs, without having to deal with the details of the file system.

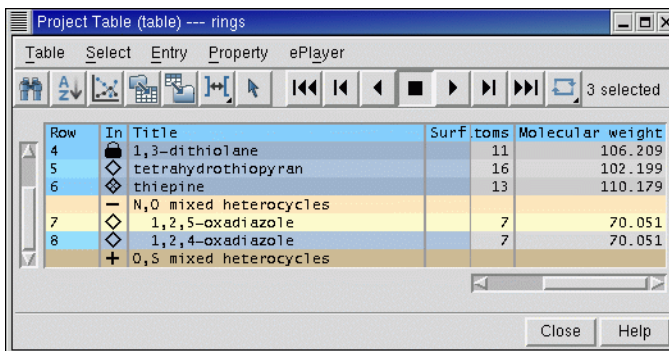


Figure 1.1. The Project Table panel.

Maestro has dedicated panels for each product for preparing and submitting jobs. To use these panels, choose the appropriate product and task from the Applications menu and its submenus. Set the appropriate options in the panel, then click **Start** to open the Start dialog box and set options for running the job. See [Section 4.2](#) of the *Job Control Guide* for a detailed description of the options in the Start dialog box. When you have finished, click **Start** to launch the job and open the Monitor panel.

Input files for the job are written to Maestro's current working directory unless otherwise specified (see [Section 12.2.3 on page 155](#)). This directory can be changed by choosing **Change Directory** from the Maestro menu and navigating to the new directory, or entering the command `cd path` in the command input area of the main window.

Once a job is launched, the Monitor panel is automatically displayed and job progress is posted to the monitor window. You may cease or resume monitoring at any time, and you can monitor any currently running job. If a job finishes while it is being monitored, its output is incorporated into the current project, according to the settings used to launch the job. If a job that is not currently being monitored finishes successfully, you must select it in the Monitor panel and click **Monitor** to incorporate the results. Multiple structures are incorporated as entry groups.

For the programs that do not incorporate results, the output files are placed in the file I/O directory that you specify using the Preferences panel.

For information on job control and incorporation, see the *Job Control Guide*.

Epik, Impact, Jaguar, LigPrep, MacroModel, QikProp, and QSite

These computational programs, which have been fully integrated with the Project Facility, can use Maestro structure files, project entries and structures displayed in the Workspace as input, and can incorporate job results into the project. You can specify whether the output structures will be appended to the project, replace the input entries, or be ignored when the job finishes. As the job progresses, any interim structures are shown in the Workspace.

Some of these computational programs have additional capabilities:

- Epik and LigPrep jobs can read structures from external files. This choice of structural input has no connection with a Maestro project, so the output is required to be written to an external file.
- Jaguar maintains its own input and output files and directory structure.

Glide, Liaison, and SiteMap

These computational programs have differing degrees of integration with the Project Facility. Glide can use selected entries from the Project Table as input for docking, but does not incorporate the results. The pose files that are generated are better viewed in the Glide Pose Viewer

than in the Project Table. Liaison takes input from structures in the Workspace or the Project Table, and the results are automatically appended to the project. Protein preparation uses the Workspace structure as input and the resulting structure is incorporated. SiteMap uses the Workspace as the source of input, and incorporates the results into the project.

Phase

The three stages of Phase workflow behave differently in relation to the Maestro project. Development of a pharmacophore model uses a wizard-based interface. It reads structures from the project and stores information in the project, but does not incorporate results into the Project Table. Database preparation has no connection to a Maestro project. The database search automatically adds the hits to the Project Table: there are no other incorporation options.

Prime

Prime's Protein Structure Prediction uses a wizard-based interface. The input is a sequence rather than a structure, though the sequence can be taken from a structure in the Workspace. At the final step, the structures can be copied into the Project Table. The standalone Refinement facility is similar to MacroModel and QikProp: it takes the structure in the Workspace as input and incorporates the results into the Project Table according to the selected option.

CombiGlide

Like Prime and Phase, CombiGlide adds data to a project in the form of "runs", but structures are not incorporated into the Project Table unless explicitly requested.

Strike

Strike is intimately linked with the Project Facility, uses descriptor data from single or multiple Project Table entries as input, and automatically incorporates job results into the project without overwriting existing data. Workspace and selected entry molecules are used to determine which molecules are used by Strike for a given calculation. Job output can be monitored in the Monitor panel. Jobs can be run on any machine listed within the `schrodinger.hosts` file. Models generated by Strike can be imported from one project to another or exported for use at the command line.

1.5 Citing Maestro in Publications

The use of this program should be acknowledged in publications as:

Maestro, version 7.5, Schrödinger, LLC, New York, NY, 2006.

The Maestro Main Window

This chapter describes how to start and end a Maestro session, and gives a detailed description of the main window and its functions.

2.1 Starting Maestro

The mechanism for starting Maestro is different for UNIX and Windows hosts. This mechanism is described for each of these host types in the following subsections.

When Maestro starts, it reads resource files that define your preferences from your preferences directory. For information on the location of this directory and on setting preferences, see [Chapter 12](#). Any resource files that it does not find in this location it reads from the resource subdirectory in the Schrödinger software installation, `maestro-vversion/data/res`, where *version* is the 5-digit Maestro version number.

The default project that Maestro starts with is a scratch project (see [Section 1.3 on page 3](#) for information on scratch projects). This project is stored in the directory specified by the `TEMP_PROJECT` environment variable. The default location is

UNIX: `$HOME/.schrodinger/tmp`
Windows: `%USERPROFILE%\schrodinger\tmp`

If space is limited on your home file system, consider setting `TEMP_PROJECT` to a directory in which you have plenty of space to avoid filling the home file system.

On UNIX, the directory from which Maestro was started is called the *launch directory*. On Windows, this directory is the directory specified by the `USERPROFILE` environment variable.¹ The launch directory becomes Maestro's *current working directory*. For information on changing the current working directory, see [Section 1.4 on page 4](#).

2.1.1 Starting Maestro on a Windows Host

To start Maestro, double-click the Maestro icon on the desktop, or use the Start menu. Maestro is in the Schrödinger submenu. This submenu can also be used to open the release notes and the documentation index in your browser.

1. See [Appendix F](#) for information on setting environment variables on Windows.

2.1.2 Starting Maestro on a UNIX Host

Before starting Maestro, you must set the `SCHRODINGER` environment variable to point to the installation directory, by entering the following command at a shell prompt:

```
cshtcsh:      setenv SCHRODINGER installation-directory
bash/ksh:    export SCHRODINGER=installation-directory
```

You might also need to set the `DISPLAY` environment variable, if it is not set automatically when you log in. To determine if you need to set this variable, enter the command:

```
echo $DISPLAY
```

If the response is a blank line, set the variable by entering the following command:

```
cshtcsh:      setenv DISPLAY display-machine-name:0.0
bash/ksh:    export DISPLAY=display-machine-name:0.0
```

After you set the `SCHRODINGER` and `DISPLAY` environment variables, you can start Maestro using the command:

```
$SCHRODINGER/maestro [options] [filename]
```

If you add the `$SCHRODINGER` directory to your path, you only need to enter the command `maestro`.

The options for the `maestro` command are listed in [Table 2.1](#), including two options for specific problems. The optional *filename* specifies a file in Maestro or MacroModel format.

Table 2.1. Options for the `maestro` command.

Option	Description
<code>-b PDB-file</code>	Read a PDB file. Cannot be used with <code>-m</code> .
<code>-c command-script</code>	Run the specified Maestro command script when Maestro starts.
<code>-h</code>	Print a usage message, but do not start Maestro.
<code>-m filename</code>	Load structures from a Maestro (<code>.mae</code>) or MacroModel (<code>.dat</code>) structure file. Cannot be used with <code>-b</code> . The <code>-m</code> can be omitted if this option is at the end of the option list.
<code>-p project</code>	Load the specified Maestro project.
<code>-s number</code>	Specify the first structure to load when a Maestro file is loaded.
<code>-t number</code>	Specify the total number of structures to load when a Maestro file is loaded.
<code>-v</code>	Display the Maestro version number, but do not start Maestro.

Table 2.1. Options for the `maestro` command. (Continued)

Option	Description
<code>-I</code>	Display extra information about X visual settings when Maestro starts.
<code>-SGL</code>	Use the Schrödinger-supplied OpenGL library. The default is to use the system OpenGL library.
<code>-nosoftaa</code>	SGI only. On some SGIs anti-aliasing (smoothing) for line widths greater than 1 is only possible in a software implementation. On certain machines there may be a performance penalty when using the software implementation. Use this option to disable software anti-aliasing.
<code>-nosplash</code>	Disable display of the splash screen when starting Maestro.
<code>-nosurfdisplaylist</code>	Don't use display lists for surface rendering (required to fix display problems on some hardware platforms).

2.2 The Main Window

The Maestro main window is shown in [Figure 2.1](#). The following main window components are always visible:

- **Title bar**—displays the project name and the current working directory
- **Auto-Help**—automatically displays context-sensitive help
- **Main menu bar**—provides access to panels
- **Workspace**—displays molecular structures and other 3D graphical objects

The following main window components can be displayed or hidden by choosing the component from the Display menu. Your choice of which main window components are displayed is persistent between Maestro sessions.

- **Clipping planes window**—displays a small, top view of the Workspace that shows the clipping planes and viewing volume indicators. This feature is hidden by default. To display it, choose Clipping Planes from the Display menu
- **Toolbar**—contains buttons for many common tasks and also provides tools for displaying and manipulating structures, as well as organizing the Workspace
- **Status bar**—displays information about a particular atom, or about structures in the Workspace
- **Sequence viewer**—shows the sequences for proteins displayed in the Workspace. This feature is hidden by default. To display it, choose Sequence Viewer from the Display menu

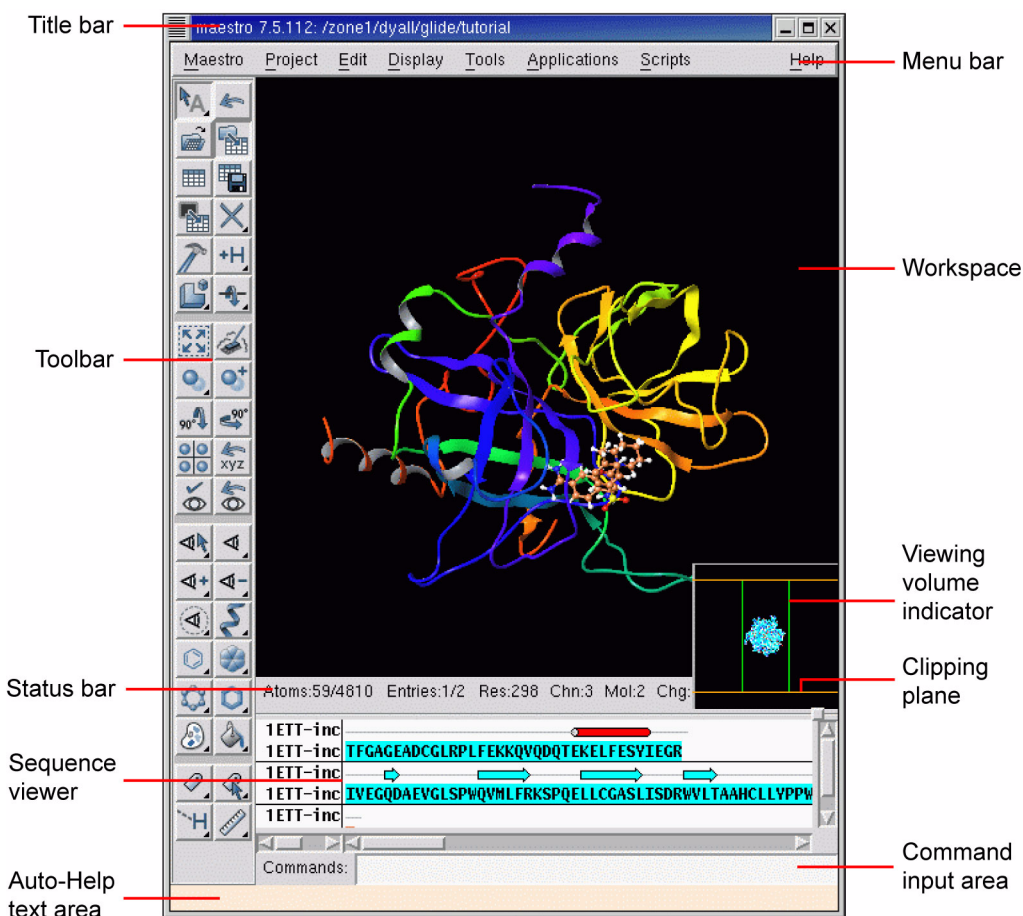


Figure 2.1. The Maestro main window.

- **Command input area**—provides a place to enter Maestro commands. This feature is hidden by default. To display it, choose Command Input Area from the Display menu

When a distinction between components in the main window and those in other panels is needed, the term *main* is applied to the main window components (e.g., main toolbar).

You can expand the Workspace to occupy the full screen, by pressing CTRL+=. All other components and panels are hidden. To return to the previous display, press CTRL+= again.

The following sections describe some of these main window components in more detail.

2.3 The Menu Bar

The menus on the main menu bar provide access to panels, allow you to enter commands, save or load a panel layout, set preferences, set up Maestro command aliases, and quit Maestro.

- **Maestro**—save or print images in the Workspace, enter system commands, save or load a panel layout, set preferences, set up Maestro command aliases, and quit Maestro.
- **Project**—open and close projects, import and export structures, make a snapshot, and annotate a project. These actions can also be performed from the Project Table panel. For more information on projects, see [Chapter 8](#).
- **Edit**—undo actions, build and modify structures, define command scripts and macros, and find atoms in the Workspace.
- **Display**—control the display of the contents of the Workspace, arrange panels, and display or hide main window components.
- **Tools**—group atoms; measure, align, and superimpose structures; assign structural features; and view and visualize data.
- **Applications**—set up, submit, and monitor jobs for Schrödinger’s computational programs and solutions. Some products have a submenu from which you can choose the task to be performed.
- **Scripts**—manage, install, and run scripts. For details, see [Section 13.1 on page 167](#).
- **Help**—open the Help panel, open the documentation index, display information about Maestro and other Schrödinger products, and display or hide Balloon Help (tooltips).

The functions of the items on the menus are described in detail in later chapters of this manual.

2.4 The Toolbar

The main toolbar contains three kinds of buttons for performing common tasks:



Action—Perform a simple task, like clearing the Workspace.



Display—Open or close a panel or open a dialog box, such as the Project Table panel.



Menu—Display a *button menu*. These buttons have a triangle in the lower right corner.



Figure 2.2. The Workspace selection *button menu* and the Adjust distances, angles or dihedrals *button menu*.

There are four types of items on button menus, and all four types can be on the same menu (see Figure 2.2):

- **Action**—Perform an action immediately.
- **Display**—Open a panel or dialog box.
- **Object types for selection**—Choose Atoms, Bonds, Residues, Chains, Molecules, or Entries, then click on an atom in the Workspace to perform the action on all the atoms in that structural unit.

The object type is marked on the menu with a red diamond and the button is indented to indicate the action to be performed.

- **Other setting**—Set a state, choose an attribute, or choose a parameter type and click on an atom in the Workspace to act upon that parameter.

The toolbar buttons are described below. Many of the features are described in later chapters.

Workspace selection

- Choose an object type for selecting
- Open the Atom Selection dialog box



Undo/Redo

Undo or redo the last action. Performs the same function as the Undo item on the Edit menu, and changes to an arrow pointing in the opposite direction when an Undo has been performed, indicating that its next action is Redo.

Open a project

Open the Open Project dialog box.



Import structures

Open the Import panel.

Open/Close Project Table

Open the Project Table panel or close it if it is open.



Save as

Open the Save Project As dialog box, to save the project with a new name.

Create entry from Workspace

Open a dialog box in which you can create an entry in the current project using the contents of the Workspace.

**Delete**

- Choose an object type for deleting
- Delete hydrogens and waters
- Open the Atom Selection dialog box
- Delete other items associated with the structures in the Workspace
- Double-click to delete all atoms

Open/Close Build panel

Open the Build panel or close it if it is open.

**Add hydrogens**

- Choose an object type for applying a hydrogen treatment
- Open the Atom Selection dialog box
- Double-click to add hydrogens to all structures in the Workspace

Local transformation

- Choose an object type for transforming
- Open the Advanced Transformations panel

**Adjust distances, angles or dihedrals**

- Choose a parameter for adjusting
- Delete adjustments

Fit to screen

Scale the displayed structure to fit into the Workspace and reset the center of rotation.

**Clear Workspace**

Clear all atoms from the Workspace

Set fog display state

Choose a fog state. Automatic means fog is on when there are more than 40 atoms in the Workspace, off otherwise.

**Enhance depth cues**

Optimize fogging and other depth cues based on what is in the Workspace.

Rotate around X axis by 90 degrees

Rotate the Workspace contents around the X axis by 90 degrees.

**Rotate around Y axis by 90 degrees**

Rotate the Workspace contents around the Y axis by 90 degrees.

Tile entries

Arrange entries in a rectangular grid in the Workspace.

**Reset Workspace**

Reset the rotation, translation, and zoom of the Workspace to the default state.

Save view

Save the current view of the Workspace: orientation, location, and zoom.

**Restore view**

Restore the last saved view of the Workspace: orientation, location, and zoom.

Display only selected atoms

- Choose an object type for displaying
- Double-click to display all atoms

**Display only**

- Choose a predefined atom category
- Open the Atom Selection dialog box

Also display

- Choose a predefined atom category
- Open the Atom Selection dialog box

**Undisplay**

- Choose a predefined atom category
- Open the Atom Selection dialog box

Display residues within N angstroms of currently displayed atoms

- Choose a radius
- Open a dialog box to set a value



Show, hide, or color ribbons

- Choose show or hide ribbons
- Choose a color scheme for coloring ribbons

Draw bonds in wire

- Choose an object type for drawing bonds in wire representation
- Open the Atom Selection dialog box
- Double-click to apply to all atoms



Draw atoms in CPK

- Choose an object type for drawing bonds in CPK representation
- Open the Atom Selection dialog box
- Double-click to apply to all atoms

Draw atoms in Ball & Stick

- Choose an object type for drawing bonds in Ball & Stick representation
- Open the Atom Selection dialog box
- Double-click to apply to all atoms



Draw bonds in tube

- Choose an object type for drawing bonds in tube representation
- Open the Atom Selection dialog box
- Double-click to apply to all atoms

Color all atoms by scheme

Choose a predefined color scheme



Color residue by constant color

- Choose a color for applying to residues
- Double-click to color all atoms

Label atoms

- Choose a predefined label type
- Delete labels



Label picked atoms

- Choose an object type for labeling atoms
- Open the Atom Selection dialog box
- Open the Composition tab of the Atom Labels panel
- Delete labels

Display H-bonds

- Choose bond type:

intra—displays H-bonds within the molecule you select

inter—displays H-bonds between the molecule you select and all other atoms in the Workspace.

- Delete H-bonds



Measure distances, angles or dihedrals

- Choose a parameter for displaying measurements
- Delete measurements

2.5 The Status Bar

The status bar displays information about the contents of the Workspace or about a particular atom, depending on where the pointer pauses. If the pointer is not over an atom, information about the contents of the Workspace is displayed, in the format

Atoms:*ma/na* Entries:*me/ne* Res:*nr* Chn:*nc* Mol:*nm* Chg:*n*

where the fields in the display have the following meanings:

Atoms	Number of atoms in the Workspace: <i>ma</i> is the number displayed, <i>na</i> is the total
Entries	Number of entries in the Workspace: <i>me</i> is the number displayed, <i>ne</i> is the total
Res	Number of residues <i>nr</i> in the Workspace
Chn	Number of chains <i>nc</i> in the Workspace
Mol	Number of molecules <i>nm</i> in the Workspace
Chg	Total charge <i>n</i> of the Workspace contents

If the pointer is over an atom, and the preference to show feedback on atoms is set, the status bar displays the chain, residue number, element, PDB atom name, formal charge, and a user-selected property, in the format

chain:residue number : (element) PDB :charge: "userchoice"

where the fields in the display have the following meanings:

<i>chain</i>	Chain name
<i>residue</i>	Residue name
<i>number</i>	Residue number and insertion code
<i>element</i>	Element
<i>PDB</i>	PDB atom name
<i>charge</i>	Formal charge of the atom
<i>userchoice</i>	Selected Project Table property

The preference to show feedback and the user-selected property can be set in the Feedback tab of the Preferences panel—see [Section 12.2.8 on page 160](#).

2.6 The Sequence Viewer

The sequence viewer displays the sequence of the proteins that are included in the Workspace. Each chain is displayed on a separate line. The left column of the sequence viewer lists the entry name for the sequence, with the chain name appended. The right column displays the actual sequence as a list of single-letter codes, and also displays the secondary structure assignment, if one is available.

The background colors of the sequences reflect the current coloring scheme for the alpha carbon in the Workspace. Each residue in the sequence also has a tooltip that gives the residue name, number, and insertion code.

If the protein has a secondary structure assignment, the assignment is displayed above the sequence, using the following symbols:

- Helices are displayed as tubes, colored red
- Strands are displayed as arrows, colored cyan
- Loops are displayed as a gray line.

You can select residues in the sequence and perform actions on these residues.

- To select a range of residues, either click the first and shift-click the last, or drag over the residues.
- To select a single residue, click the residue.
- To add residues to the selection or remove residues from the selection, control-click the residues.

The selected residues are highlighted as white letters on a black background, and are marked in yellow in the Workspace.

The sequence viewer and the Workspace are synchronized: selection of atoms in the Workspace selects the corresponding atoms in the sequence viewer, and selection of residues in the sequence viewer selects them in the Workspace. Likewise, changing the color scheme in the Workspace changes the color scheme in the sequence viewer, and vice versa. Structural changes, such as mutating a residue, are also reflected in the sequence viewer.

The sequence viewer has a shortcut menu (right-click and hold), which you can use to perform various actions on the sequences or the display. If you do not have any residues selected in the sequence viewer, the sequence viewer shortcut menu is displayed. The items on this menu are described in [Table 2.2](#). If you have residues selected, the selection shortcut menu is displayed (see [Table 2.3](#)). In this context there is one extra item on the Color Scheme submenu: Proximity. If you choose this item, the residues within a given distance of the selected residues are colored orange, and the selected residues are colored red. All others are colored gray. The distance is specified in the Sequence Viewer - Preferences dialog box, which you open from the sequence viewer shortcut menu.

Table 2.2. Sequence viewer shortcut menu items.

Menu Item	Action
Hide	Hide the sequence viewer. To redisplay it, use the Display menu.
Color Scheme	Color the sequence with the scheme selected from the submenu. Choose from None, Color by Chain, Residue Type, Residue Property, Secondary Structure, Template ID.
Font Size	Set the font size for the sequence viewer. Choose Small, Medium, Large, or Huge from the submenu.
View SSA	Show or hide the secondary structure assignment for the sequence.
Exclude Entry	Exclude the entry containing this chain from the Workspace. The chain is removed from the sequence viewer. If the entry contains multiple chains, all the chains in the entry are removed from the sequence viewer.
Legend	Open the Sequence Viewer Legend dialog box, which explains the sequence color schemes.
Preferences	Open the Sequence Viewer - Preferences dialog box, in which you can set the proximity cutoff.

2.7 Shortcut Menus

Several of the Workspace components have shortcut menus, which you open by right-clicking and holding over the component. The shortcut menus for the status bar and the command input area only have one item, **Hide**, which hides the component. In addition to **Hide**, the toolbar shortcut menu has **Move to Left** and **Move to Right**, which move the toolbar to the left or right side of the Workspace.

The Workspace has three shortcut menus: one for the Workspace itself (the *Workspace shortcut menu*), one for the atoms in the current selection (the *selection shortcut menu*), and one for bonds (the *bond shortcut menu*). The selection shortcut menu is also available from the sequence viewer. The items on these menus are described in [Table 2.3](#), [Table 2.4](#) and [Table 2.5](#). When you right-click and hold in the Workspace:

- On an unselected atom, the atom becomes the selection and the selection shortcut menu is displayed.
- On a selected atom, the selection shortcut menu is displayed.
- On a bond, the two atoms in the bond are selected and the bond shortcut menu is displayed.
- On empty space, if there is an atom selection, the selection shortcut menu is displayed, otherwise the Workspace shortcut menu is displayed.

Table 2.3. Selection shortcut menu items.

Menu Item	Action
Invert Selection	Select the unselected atoms and deselect the selected atoms.
Representation	Set the molecular representation to wire frame, ball and stick, CPK, or tube.
Color Scheme	Apply the color scheme from the submenu to the selected atoms.
Labels	Apply the label from the submenu to the selected atoms.
Display	Display the entire molecule or residue containing the selected atoms. This action is useful if part of the molecule or residue was previously undisplayed.
Display Only	Display only the selected atoms.
Undisplay	Undisplay the selected atoms.
Delete	Delete the selected atoms.
Set Element	Change the selected atoms to atoms of the element chosen from the submenu.
Mutate Residue	Change the residues that contain the selected atoms to the type chosen from the submenu.
Exclude Entry	Exclude the entries that contain the selected atoms from the Workspace.
Select Entry in Table	Select the entries that contain the selected atoms in the Project Table.
Fit to Screen	Fit the displayed structure to the Workspace and reset the center of rotation.

Table 2.4. Workspace shortcut menu items.

Menu Item	Action
Background Color	Set the Workspace background color to black, white, gray, or blue.
Print	Open the Print panel to print an image of the Workspace or create a postscript file of the Workspace.
Save Image	Open the Save Image panel to save an image of the Workspace in TIFF or JPEG format.
Save View	Save the current view of the Workspace: orientation, location, and zoom.
Restore View	Restore the last saved view of the Workspace: orientation, location, and zoom.
Fit to Screen	Fit the displayed structure to the Workspace and reset the center of rotation.
Create Entry from Workspace	Create a Project Table entry from the Workspace.
Clear Workspace	Clear the Workspace.
Delete	Delete labels, measurements, or adjustments.

Table 2.5. Bond shortcut menu items.

Menu Item	Action
Order	Set the bond order to zero, one, two, or three.
Measure	Measure the bond length and display it. The measurement is added to the list in the Measurements panel.
Adjust	Adjust the bond length or the dihedral angle around this bond.
Delete	Delete the bond.

2.8 Mouse Functions

The Workspace has special uses for the middle and right mouse buttons. These are used on their own and in combination with the SHIFT and CTRL keys to perform common operations such as rotation, translation, centering, and zooming. See the table below for a summary and [Chapter 7](#) for more detailed information.

Table 2.6. Mapping of Workspace operations to mouse actions.

Mouse Button	Keyboard	Motion	Action
Left		click, drag	Select
Left	SHIFT	click, drag	Toggle the selection
Middle		drag	Rotate about X and Y axes
Middle	SHIFT	drag vertically	Rotate about X axis
Middle	SHIFT	drag horizontally	Rotate about Y axis
Middle	CTRL	drag horizontally	Rotate about Z axis
Middle	SHIFT + CTRL	drag horizontally	Zoom
Middle		drag highlighted parameter	Adjust bond, angle or dihedral (see Section 4.7.1 on page 49)
Right		click	Spot-center on selection Toggle selection
Right		click and hold	Display shortcut menu
Right		drag	Translate in the X-Y plane
Right	SHIFT	drag vertically	Translate along the X axis

Table 2.6. Mapping of Workspace operations to mouse actions. (Continued)

Mouse Button	Keyboard	Motion	Action
Right	SHIFT	drag horizontally	Translate along the Y axis
Right	CTRL	drag horizontally	Translate along the Z axis
Middle & Right		drag horizontally	Zoom

2.9 Shortcut Keys

Some frequently used operations in the main window have been assigned shortcut key combinations. The shortcuts, their functions, and their menu equivalents are listed in [Table 2.7](#). Shortcut keys in the Project Table panel and the Plot panel are listed in the chapters that describe those panels.

Note: If Num Lock, Caps Lock, or Scroll Lock is enabled, the shortcut keys do not work.

Table 2.7. Shortcut keys in the Maestro main window.

Keys	Action	Equivalent Menu Choice
CTRL+B	Open Build panel	Edit > Build
CTRL+C	Create entry	Project > Create Entry From Workspace
CTRL+E	Open Command Script Editor panel	Edit > Command Script Editor
CTRL+H	Open Help panel	Help > Help
CTRL+I	Open Import panel	Project > Import Structures
CTRL+M	Open Measurements panel	Analysis > Measurement
CTRL+N	New project	Project > New Project
CTRL+O	Open project	Project > Open Project
CTRL+P	Print	Maestro > Print
CTRL+Q	Quit	Maestro > Quit
CTRL+S	Open Sets panel	Analysis > Sets
CTRL+T	Open Project Table panel	Project > Show Table
CTRL+W	Close project	Project > Close Project
CTRL+Z	Undo/Redo last command	Edit > Undo/Redo
CTRL+=	Enter and exit full screen mode (Workspace occupies full screen)	None

2.10 Undoing Workspace Operations

There are two classes of operations in the Workspace: operations that act on the structure, and operations that act on the view of the molecule. The former include local transformations, other coordinate changes, and adding or removing atoms from the display. The latter include global transformations: translation, rotation, and magnification (zoom). These operations do not change the coordinates recorded in the project entry, but store a local copy that is translated, rotated, or magnified.

2.10.1 Undoing Operations on the Structures

Single operations in the Workspace can be undone by choosing *Undo action* from the Edit menu. This item is dynamic: its name includes a description of the last undoable action, and changes to *Redo* after it has been used. Selecting *Redo action* reimplements the action that was undone. This menu item is mapped to the Undo/Redo button on the toolbar.

To undo the last operation, do one of the following:

- Click the Undo/Redo button on the toolbar.



- Choose *Undo action* from the Edit menu.
- Press CTRL+Z.

2.10.2 Undoing Operations on the View

Global molecular transformations, such as rotations and translations, cannot be undone with the Undo option on the Edit menu or in Undo Blocks. To undo such a transformation, you must save the orientation of the structure before you do the global transformation. You can also revert to the default orientation and zoom. The default orientation is defined by the coordinates of the structure. The default zoom value is the value you get by clicking the Fit to screen button on the toolbar.



To save an orientation:

- Click the Save view button in the toolbar.



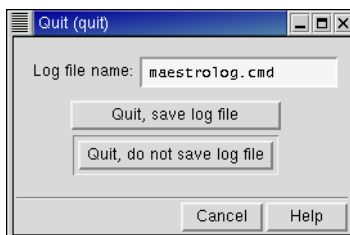


Figure 2.3. The Quit panel.

To restore the saved orientation:

- Click the Restore view button in the toolbar.



To revert to the default orientation and zoom:

- Click the Reset workspace button in the toolbar.



2.11 Ending a Maestro Session

To end a Maestro session, choose Quit from the Maestro menu. If you want to save a log file with a record of all operations performed in the current session, click Quit, save log file on the Quit panel. This information can be useful to Schrödinger support staff when responding to any problem you report. If you do not want to save a log file, click Quit, do not save log file. If you are working in a scratch project and have not saved the project, you are prompted to do so. If you decide that you do not want to end the session, click Cancel.

Importing and Exporting Structures and Data

The primary place for reading and writing data from Maestro is a project. When you open a project, the structures and data from that project are read into Maestro. As you modify structures and data they are automatically saved in the project (unless you have selected an option not to do so). Projects are described in detail in [Chapter 8](#). In addition to reading and writing to projects, Maestro can also read and write structural information in a variety of formats, including its own native Maestro format, MacroModel, PDB, Sybyl Mol2, and MDL SD format. Through the Jaguar panel, structural information can be obtained from the input and output files of a wide variety of quantum chemical programs.

Obtaining structures from a file in Maestro is called *importing* rather than *reading* because the structures are imported into the current project. Any structures that do not exist as entries in the current project are external to the project, and therefore must be imported. Similarly, the term *exporting* is used instead of *writing* because a file is created for use outside the current project.

Data that is associated with structures can also be imported and exported independently to spreadsheet files. This capability enables properties that are created by Schrödinger software to be read by a spreadsheet program, manipulated, and re-imported into the project.

3.1 Importing Structures

Structures are imported using the Import panel. You can import data from files in Maestro, MacroModel, PDB, Sybyl Mol2, MDL SD, and a range of other formats. If the files contain multiple structures, such as the output of a conformational search, you can select the structures to import. When importing structures, Maestro places each structure into a separate, new entry in the current project. For the Maestro file type, molecular representation information can be imported along with a structure. You can undo the import operation, but the changes made to the Project Table remain.

To open the Import panel, do one of the following:

- Choose Import Structures from the Project menu on the main menu bar.
- Choose Structures from the Import submenu of the Table menu in the Project Table panel.
- Click the Import structures button on the main toolbar.



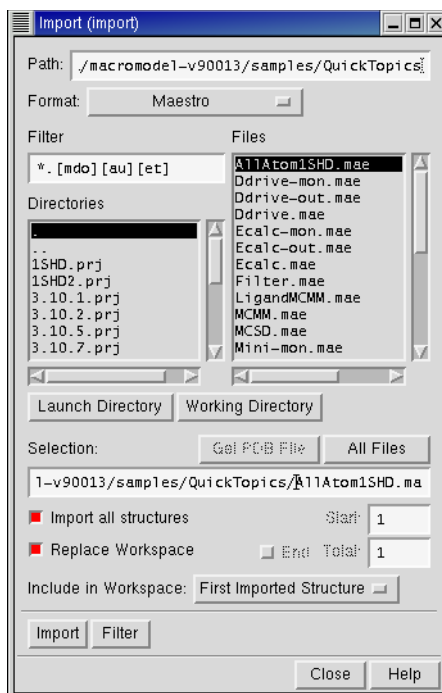


Figure 3.1. The Import panel.

3.1.1 Selecting Files

Selecting a file type: Before you choose files, you must specify the file type by choosing a file format from the Format menu. This choice determines the file extensions to be used, and thus the filter displayed in the Filter text box and the files displayed in the Files list.

The menu options are described below:

Maestro

Import from Maestro structure files, which have a .mae extension, or MacroModel input or output structure files, which have .dat or .out extensions. Properties are imported as well as structures.

PDB

Import PDB files, which have either a .pdb or a .ent extension. Where multiple atomic coordinates exist for a single PDB entry, Maestro chooses the atoms with the highest occupancy ratio. Properties are imported as well as structures.

Mol2

Import Sybyl Mol2 files, which have a `.mol2` extension.

SD

Import MDL SD files, which have either a `.sdf` extension for multiple structure files or a `.mol` extension for files with a single structure. Properties are imported as well as structures.

Babel

Import structures from a range of other programs. The format conversion is performed by the program Babel. The list of formats provided here is a subset of the formats that Babel can convert. If you want to import structures in other formats using Babel, you can use Jaguar to read the files. See [Section 11.1.5](#) of the *Jaguar User Manual* for more information, including the full list of file types read by Babel.

On UNIX hosts you can use Babel from the command line to convert structures to a format that Maestro can import directly. To do so, you must first set the environment variable `BABEL_DIR`:

```
csh/tcsh:    setenv BABEL_DIR $SCHRODINGER/mmshare-vversion/data/babel
bash/ksh:    export BABEL_DIR=$SCHRODINGER/mmshare-vversion/data/babel
```

To obtain usage information, enter the `babel` command without arguments:

```
$SCHRODINGER/mmshare-vversion/bin/platform/babel
```

Here, *version* is the 5-digit mmshare version number, and *platform* indicates the platform on which you are running, e.g. `linux-x86`.

Any

Import a file of any type. The file type is determined from the extension of the selected file or the file name you enter in the Selection text box. If the extension is not unique, Maestro does not import the file. The formats for which Maestro recognizes the extension are Maestro, Mol and Mol2, PDB, SD, XYZ, Mopac Cartesian and Internal, and Spartan. This option is useful if you want to locate any of these kinds of structure files.

Selecting a directory: To specify the directory from which to import, you can enter it in the Path text box, navigate to it using the Directories list, or click Launch Directory or Current Directory to select the directory from which you launched Maestro or the current working directory. The path and the directory and file lists are updated for the chosen directory.

Selecting a file: To select a file, you can type its name (and path) in the Selection text box, or select it from the Files list.

If you selected one of the four specific file types from the first section of the menu, or some of the Babel file types, the file that you select must have the correct extension. If you type the name of the file in the Selection text box, you must include the correct extension. In the Files list, only files with the correct extension are displayed.

If the file does not appear in the list, check that you have chosen the correct file type in the Format section, or select All in the Format section to view all files. You can also view files with a different extension by altering or deleting the filter and pressing RETURN or clicking Filter.

Selecting multiple files: You can import structures from more than one file at a time. To select multiple input files, use the SHIFT and CTRL keys in combination with mouse clicks to select as many files as you wish in the Files list, or click All Files to select all files in the Files list.

Importing PDB files: If you have the PDB database and Prime installed, you can import a structure directly from the Import panel:

1. In the Import panel, choose PDB from the Format option.
2. Click Get PDB File and enter the PDB ID and chain name in the dialog box

The structure is copied to your local directory so you can import it into Maestro.

Note: PDB files do not support multiple structures, only multiple sets of coordinates for the same molecule, e.g. different NMR refinements. Certain PDB structures do not follow the PDB standard and have varying number of atoms in different MODELs. Do not import these PDB structures, as they will cause Maestro to fail.

3.1.2 Selecting Import Settings

Once you have selected the files to import structures from, you must choose which structures to import and how to display them.

Selecting structures to import: Maestro can import multiple structures from multiple files. If you do not want to import all structures from each file, deselect Import all structures. You can then enter the index of the first structure in the Start text box, and the total number of structures to import in the Total text box. You can import to the end of the file by selecting End instead of specifying the total number. These settings are applied to all files: they cannot be set for each file in a multiple file selection. To import different numbers of structures from a set of files, select each file in turn and set the number of structures to import.

Selecting structures to display: When you import structures from a file, you can set options to control which structures to display, and whether the structures already in the workspace are cleared or not. If Replace workspace is selected (the default setting), Maestro automatically clears existing structures from the Workspace. To import new structures without removing the Workspace contents, deselect Replace workspace.

To specify which structures are to be included in the workspace, choose one of the options from the Include in workspace option menu: No Imported Structures, First Imported Structure, or All Imported Structures. The default is First Imported Structure. If you choose No Imported Structures and Replace Workspace is selected, the Workspace is cleared.

Importing graphical information: If you import structures from a Maestro file, you can also import graphical information for the structures, such as the molecular representation. Maestro's default behavior is to read in and use this information when importing structures from a file. If you do not want imported structures to be displayed using the saved representation information, clear the Use graphical info option.

3.1.3 Error Reporting for PDB Files

When Maestro imports structures from PDB files, the color scheme it uses reflects the success or failure of the identification of atoms and the placement of bonds. Maestro attempts to place bonds according to a set of standard residue templates, the CONECT records of the PDB file, and geometry.

A few amino acids require that choices be made as to where to place the double bonds and formal charges. Schrödinger's conventions for these choices are given in [Table 3.1](#).

Table 3.1. Placement of double bonds and formal charges for imported amino acids.

Amino Acid Pair	Double Bond	Formal Charge
ARN/ARG	between CZ and NH1	for ARG, +1 formal charge on NH1
ASP/ASH	between CG and OD1	for ASP, -1 formal charge on OD2
GLU/GLH	between CD and OE1	for GLU, -1 formal charge on OE2

If the placement of bonds or identification of atoms fails, the regions of the structure in which problems were identified are colored according to the type of error. The color scheme is given in [Table 3.2](#). This color scheme is available as a standard color scheme in Maestro, so you can make use of it later on.

Maestro imports PDB files by running the utility `pdbconvert`, which is described in [Appendix D](#). The conversion creates a temporary Maestro-formatted file, and imports the structures from this file. In addition to the color scheme, which is set up in the temporary file, Maestro displays warning messages to indicate the status of the PDB conversion, which are listed in [Table 3.3](#). The color schemes are stored in the project, and can be applied later as a standard color scheme, or used to select residues in the Atom Selection dialog box.

Table 3.2. Coloring scheme for structures imported from PDB files.

Color	Description
orange	Non-standard residues connected by geometry and/or CONECT records. Unless duplicate CONECT records are specified, only single bonds are assigned, so you will need to add multiple bonds.
red	Standard residue, but the input PDB file had missing atoms.
blue	A standard residue that has some atom names unknown to the standard connection template. The unrecognized atom names were connected by geometry.
green	A residue with an alternate location indicator. The first one listed in the PDB file is the one used to generate the converted structure.
gray	Standard residues connected by standard templates. Confidence in bond orders assigned to these residues is high.

Table 3.3. PDB conversion warning messages.

Status	Message and Explanation
OK	No message. Conversion appears successful.
WARNING	“Problems found while converting file.” A temporary PDB output file was generated successfully, but the assignments were not all successful.
ERROR	“Could not convert file to <i>formatname</i> format.” The conversion failed, usually because the temporary output file could not be created. There may be a permissions problem, or the disk may be full.

3.1.4 Entry Naming Conventions for Imported Structures

The names of entries in a project are required to be unique. Maestro enforces the uniqueness by adding suffixes to entry names to resolve any duplications. When structures are imported from a file, they are assigned entry names according to the following procedure.

- If the file already contains entry names, these names are read and used as the base of the entry name. If the file does not contain entry names, the stem of the file name (the portion before the extension) is used as the base of the entry names.
- The suffix *-impn* is added to the base entry name, where *n* is a number beginning at 1, and is incremented for each import operation. Consequently, entries imported from different files automatically have different entry names.
- If the resulting entry name is not unique, the suffix *.m* is appended to the entry name, where *m* is a number beginning at 1, and is incremented for each import operation.

For example, if you imported three conformers of decane from a file in which each conformer had the name decane, the entry names would be decane-imp1.1, decane-imp1.2, and decane-imp1.3.

The same naming mechanism is used when calculation results are incorporated into a project, except that the suffix `-incj` is added to the entry name. See [Appendix C](#) for information on entry names.

3.1.5 Reading Jaguar Input Files

Jaguar maintains its own input files, which include calculation options in addition to the molecular geometry, and may also include molecular orbitals. Therefore, Jaguar has a separate file reading mechanism. When you read input files through Jaguar, the structures are displayed in the Workspace and Maestro-formatted files are created.

To read input files through Jaguar, choose Jaguar from the Applications menu, then click Read in the Jaguar panel. The Read File file selector is displayed. Navigate to the directory containing the file, and select the file. Choose the file type from the File Format option menu, and change the Read as setting if necessary. Click OK to read the file.

For more information, see the *Jaguar User Manual*.

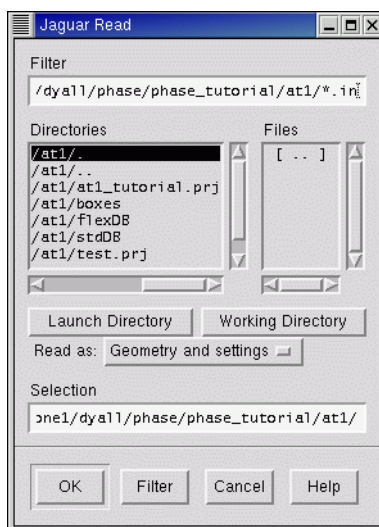


Figure 3.2. The Jaguar Read File file selector.

3.2 Exporting Structures

Using the Export panel, you can write out structure files in Maestro, MacroModel, PDB, Sybyl Mol2, and MDL SD formats, and in various other formats available through the program Babel. You can use the Workspace contents or selected project entries as input, and you can append structures to an existing file and retain additional information when the file format allows it. When exporting project entries, Maestro provides options for exporting each selected entry to its own file, or placing all selected entries in the same file.

Project data is automatically saved, unless you are running in a non-automatic synchronization mode. The structures from a project are saved in a Maestro-formatted file, along with other data. If you want to save just the structures and their associated data, or if you want to save structures in other file formats, you can export the structures.

To open the Export panel, choose Export Structures from the Project menu on the main menu bar, or choose Structures from the Export submenu of the Table menu in the Project Table panel.

When you have made all the necessary selections, click Export. A dialog box confirms that the files were written. If a file by the specified name already exists, the Exporting-Overwrite? dialog box is displayed and you must choose an action.

3.2.1 Selecting Export Settings

The Export panel is laid out in a similar fashion to the Import panel, with the file and file type selection tools at the top, and various options under them.

Selecting a file format: To create a file, you must select the type of file you wish to write from the types given in the Format option menu, which include Maestro, MMOD (MacroModel), PDB, Mol2, SD, and a range of file types converted by Babel. The file extension in most cases is determined by the file type. Default file extensions for the files in the first section of the menu can be set in the File Suffix tab of the Preferences panel. Properties from the Project Table are exported as well as structures to Maestro and SD files, and any relevant properties are exported to PDB files.

Selecting a directory: To select a directory, you can enter it in the Path text box, navigate to it in the Directories list, or click Launch Directory or Current Directory to select the directory from which you launched Maestro or the current working directory. The path and the directory and file lists are updated for the chosen directory.

Selecting a file: If the file exists, you can select it in the Files list or enter its name in the File text box. If you want to export to a new file, enter the name in the File text box. You can specify an absolute or relative path in the File text box in addition to the file name.

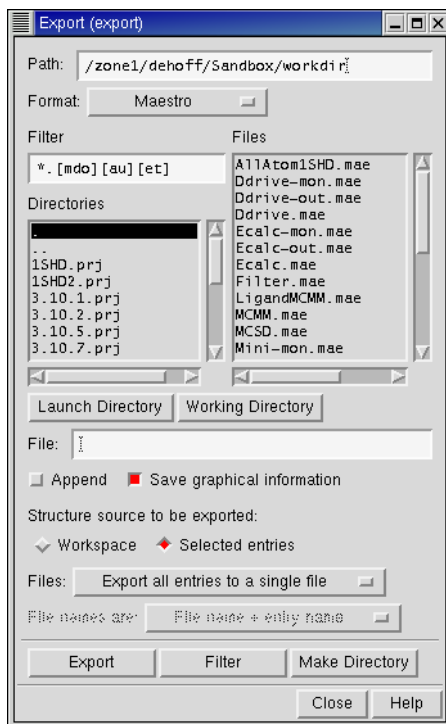


Figure 3.3. The Export panel.

Appending or overwriting files: To append a structure to an existing file, select Append. If the file exists and you do not select Append, it is overwritten.

Saving graphical information: If you are writing a Maestro file, you can select Save Graphical Information to store the molecular representation information for the structure you are saving. This option is ignored for other file types.

Specifying the source of the structures: The Export panel can be used to save structures directly from the Workspace or to write files containing structures associated with one or more project entries. To write out the Workspace to a file, select Workspace in the Structure Source To Be Exported section. This is equivalent to writing the included entries. To write out the entries that are selected in the Project Table to the exported file, select Selected Entries.

Exporting to single or multiple files: You can export entries to a single, multi-structure file, or to multiple files that each contain a single entry. From the Files option menu, choose either Export All Entries To A Single File or Export Each Entry As An Individual File. If you export structures using one of the Babel-supported formats, you can only export a single structure to each file. In this case, you should not choose Export All Entries To A Single File.

3.2.2 Exporting Multiple Files

If you want to save each entry to a separate file, you must specify how to name the files by choosing one of three naming schemes from the File Names Are option menu. The naming schemes are described below.

File Name + Entry Name

If you choose File Name + Entry Name, you must select a file or type in a file stem. The entry name and the file suffix are appended to the file stem. If you choose a file name, the suffix is the suffix of that file name. For example, if you are exporting the entries `combination-3`, `combination-4`, and `combination-5` in Maestro format, and you select the file `structure.dat`, the new files are named

```
structure_combination-3.dat
structure_combination-4.dat
structure_combination-5.dat
```

If you typed in the file stem `HIV`, the following files would be written:

```
HIV_combination-3.mae
HIV_combination-4.mae
HIV_combination-5.mae
```

File Name + Automatic Number

If you select File Name + Automatic Number, you must specify a filename. The files are given the same name and extension as the selected file, but numbers are added in sequence to the base file name. For the same example as above, if you select the file `structure.dat`, the new files are named

```
structure_1.dat
structure_2.dat
structure_3.dat
```

Just Entry Names

If you choose Just Entry Names, the file name is constructed from the entry name. For the above example, the new files are named

```
combination-3.mae
combination-4.mae
combination-5.mae
```

and a dialog box confirms that the files were written to the desired directory.

3.3 Exporting Data to a Spreadsheet

The Project Table does not have spreadsheet capabilities, but you can export data to a spreadsheet file, manipulate it in a spreadsheet program, then re-import it into Maestro. Exporting data to a spreadsheet file is done from the Export Spreadsheet panel. The data that you can export are the properties in the Project Table. To export data, you must first open Project Table panel, which you can do by choosing Show Table from the Project menu in the main window or clicking the Show/Close project table button in the toolbar.



When the Project Table panel is displayed, choose Spreadsheet from the Export submenu of the Table menu. The Export Spreadsheet panel is displayed. The layout of this panel is similar to that of the Export and Import panels, with file selection tools at the top and action buttons at the bottom.

Selecting a directory: To select a directory, you can enter it in the Path text box, navigate to it in the Directories list, or click Launch Directory or Current Directory to select the directory from which you launched Maestro or the current working directory. The path and the directory and file lists are updated for the chosen directory.



Figure 3.4. The Export Spreadsheet panel.

Selecting a file: If the file exists, you can select it in the Files list or enter its name in the File text box. If you want to export to a new file, enter the name in the File text box. You can specify an absolute or relative path in the File text box in addition to the file name. Existing files will be overwritten.

Selecting entries to export data from: The Entries option menu allows you to choose different classes of entries to export data from: All, Shown, Selected, and Included. For more information on showing, selecting, and including entries, see [Section 8.4 on page 100](#) and [Section 8.5 on page 104](#).

Selecting properties to export: You can choose to export all properties or only the shown properties, using the Properties option menu. The row number and the inclusion state are not exported, but all other properties are. The entry name is always in the first column.

Selecting a delimiter: You can choose from the two standard delimiters, a comma and a tab character, or choose other, then enter your own delimiter in the text box.

Once you have made your selection, click Export. The file is written and its name is displayed in the Files list.

3.4 Importing Data From a Spreadsheet

Once you have exported data from Maestro and manipulated it in a spreadsheet, you can re-import the data into Maestro. You can also import data that did not originate from Maestro. Importing data from a spreadsheet file is done from the Import Spreadsheet panel, and the data is added to entries in the Project Table. The data that you import must therefore match entries in the project. For data that you exported from Maestro, matching should be no problem, but for data that did not originate from Maestro, you must ensure that there is a proper match to the entries in the Project Table.

To import data, you must first open the Project Table panel, which you can do by choosing Show Table from the Project menu in the main window, or clicking the Show/Close project table button in the toolbar.



When the Project Table panel is displayed, choose Spreadsheet from the Import submenu of the Table menu. The Import Spreadsheet panel is displayed. The layout of this panel is similar to that of the Export and Import panels, with file selection tools at the top and action buttons at the bottom.

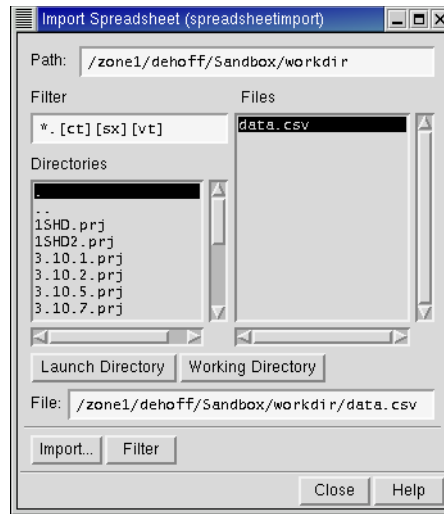


Figure 3.5. The Import Spreadsheet panel.

To select a file to import data from, you can navigate to it with the file selection tools, or enter its name in the File text box.

When you have selected the file, click Import. The Choose Import Key dialog box is displayed. In this dialog box, you determine how the rows of the spreadsheet are mapped to the entries in the Project Table, and select the delimiter. The File key column lists the properties in the spreadsheet file. The values of the key that you choose in this column are matched to the values of the property that you choose in the Match to property in Project Table column. For example, if both the spreadsheet and the Project Table have a Title property, you can match the rows by title. If the file key value for a given row in the spreadsheet matches more than one entry, the properties for that row are copied to all entries; if several file key values match a single entry, the row for the last file key value that matches is used to import data.

To complete the import, choose the delimiter that is used in the spreadsheet file, then click OK. The values in the spreadsheet are imported into the Project Table.

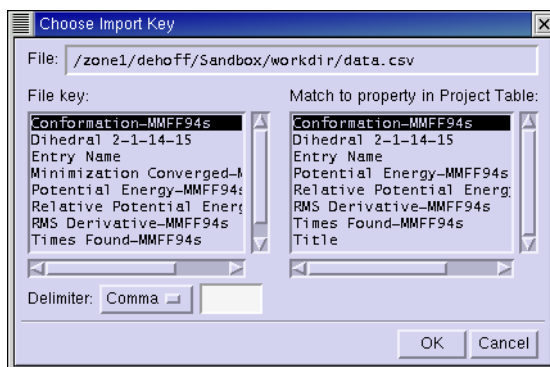


Figure 3.6. The Choose Import Key *dialog box*.

Building Structures

Maestro provides a variety of tools for building and modifying structures. In the Build panel, you can create molecular structures from predefined fragments or from individual atoms, and modify structures and atom and residue properties. Other panels enable you to add hydrogen atoms, adjust the geometry, connect or fuse molecules, and define centroids of groups of atoms. You can also set preferences for the behavior when certain building tasks are performed in the Preferences panel.

Some of the operations for building and modifying structures involve selection of atoms in the Workspace. For information on selecting atoms, see [Chapter 5](#).

4.1 The Build Panel

The Build panel consists of a toolbar and three tabs. The primary building tools are contained in the Fragments tab and the toolbar. The Fragments tab contains several fragment libraries from which you can select a fragment and add it to the Workspace, either by positioning it manually or by using predefined grow rules. The other two tabs, Atom Properties and Residue Properties, contain tools for changing the properties of atoms or residues.

To open the Build panel, choose Build from the Edit menu, or click the Open/Close Build panel button on the main toolbar.



The Fragments tab is displayed by default.

The toolbar of the Build panel provides quick access to tools for drawing and modifying structures and labeling atoms. See [Section 2.4 on page 11](#) for a description of the types of toolbar buttons. The toolbar buttons are described below, and their use is described later.



Draw structures

Choose an element for drawing structures freehand in the Workspace (default C). Each click in the Workspace places an atom and connects it to the previous atom.



Delete

Choose an object for deleting. Same as the [Delete](#) button on the main toolbar, see [page 13](#).



Set element

Choose an element for changing atoms in the Workspace (default C). Click an atom to change it to the selected element.



Increment bond order

Select a bond to increase its bond order by one, to a maximum of 3.



Decrement bond order

Select a bond to decrease its bond order by one, to a minimum of 0.



Increment formal charge

Select an atom to increase its formal charge by one.



Decrement formal charge

Select an atom to decrease its formal charge by one.



Move

Choose a direction for moving atoms, then click the atom to be moved. Moves in the XY plane are made by clicking the new location. Moves in the Z direction are made in 0.5 Å increments.



Label

Apply heteroatom labels as you build a structure. The label consists of the element name and formal charge, and is applied to atoms other than C and H.



Display Connect/Fuse panel

Open the Connect & Fuse panel so you can connect structures (create bonds between structures) or fuse structures (replace atoms of one structure with those of another).



Display Adjust panel

Open the Adjust panel so you can change bond lengths, bond angles, dihedral angles, or atom chiralities.



Add hydrogens

Choose an atom type for applying the current hydrogen treatment. Same as the [Add hydrogens](#) button on the main toolbar, see [page 13](#).



Symmetrize Workspace

Open the Symmetrize Workspace panel for symmetrizing the geometry of the structure in the Workspace.



Clean up geometry

Cleans up the geometry of the structure in the Workspace.

4.2 Building a Structure From Fragments

Whenever possible, you should use fragments to build structures. Building from fragments ensures the correct geometry within each fragment, and thereby reduces the work needed to minimize the structure. The following fragment libraries are available with Maestro:

- Organic
- Rings
- N-heterocycles
- O-heterocycles
- S-heterocycles
- Mixed heterocycles
- Carbohydrate furanose
- Carbohydrate pyranose
- L carbohydrate furanose
- L carbohydrate pyranose
- Metal centers
- Metal ligands
- Amino acids
- Nonstandard amino acids
- Protein capping groups
- Deoxyribonucleic acids
- Ribonucleic acids
- Double-stranded A DNA
- Double-stranded B DNA
- Double-stranded A RNA
- Double-stranded B RNA

These libraries are described briefly below.

- **Organic**—This library contains the common organic functional groups, a selection of simple ring structures, some common ions such as hydroxide, sulfate, nitrate, ammonium, and a hydrogen molecule.
- **Rings**—This library contains some more complex ring structures than those in the organic fragments library, including some fused ring structures, a basic steroid structure, porphyrin and phthalocyanine, buckminsterfullerene, and dodecaborane.
- **Heterocycles**—The four heterocycle fragment libraries provide a range of ring structures that contain nitrogen, oxygen, and sulfur atoms.
- **Furanose and pyranose structures**—Sugars can assume cyclized forms that include either five or six members. The five-membered ring form is called the furanose form, and the six-membered version is called the pyranose form. The natural state of most sugars is the pyranose form, but some derivatives assume a furanose form. Maestro provides both options in both D and L forms.
- **Metal centers and ligands**—These libraries contain common ligands used in inorganic chemistry and metal centers with various coordination numbers, but with the specific metal undefined. To change the metal type, retype the center once it is placed in the Workspace, using the Elements menu option in the Atom Properties tab.

- **Amino acids**—The Amino acids library contains the 22 naturally-occurring amino acids. The Nonstandard amino acids library includes the D forms of the 22 naturally-occurring amino acids and various other residues.
- **Nucleic acids**—The Deoxyribonucleic acids library contains Adenosine, Cytosine, Guanine, and Thymine, while the Ribonucleic acids library contains Adenosine, Cytosine, Guanine, and Uracil. The double-stranded DNA libraries contain deoxyribonucleic and ribonucleic acid base pairs in the A and B helix formations.
- **Protein capping groups**—This library contains some common molecules used to cap peptides.

To select a fragment library, choose the library from the Fragments option menu. The area below this menu contains a grid of buttons that you use to select fragments. Each button displays either the name of the fragment, an abbreviation for the fragment name, or the fragment structure. If a structure is displayed, the name of the fragment is given in the Balloon help for the button.

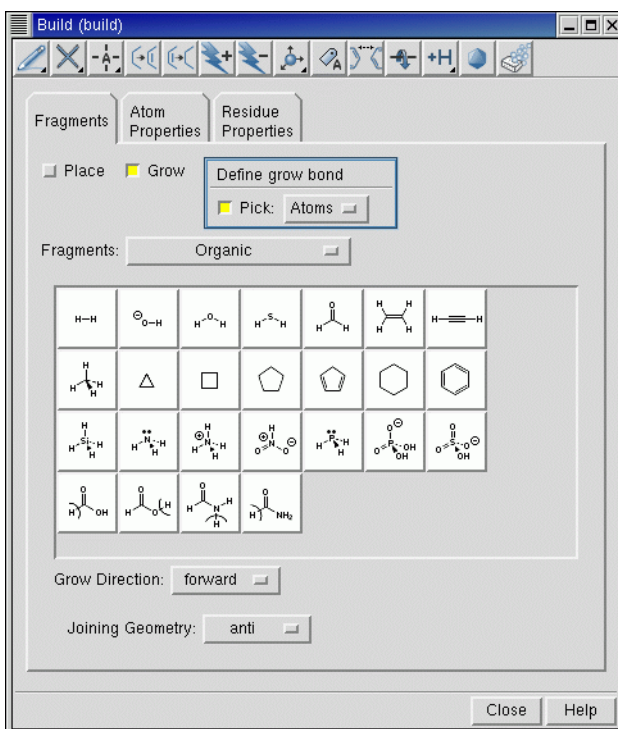


Figure 4.1. The Fragments *tab* of the Build *panel*, showing the Organic *library*.

If the particular fragment you want to use is not in the supplied libraries, you can build it from other fragments, or select a fragment with a geometry similar to the one you want to use, then replace atoms with those of different elements from the Atom Properties tab.

The fragment libraries and grow rules are read and interpreted at run time. If you are interested in creating your own fragments or fragment libraries, please write to help@schrodinger.com to obtain the appropriate instructions.

4.2.1 Building Structures Using Place Mode

You can use Place mode to create multiple molecular structures in the Workspace, to assemble structures from multiple fragments, and to change a fragment into another fragment within a structure. To enter Place mode, select Place in the Fragments tab.

To create multiple structures in the Workspace in Place mode, choose a library from the Fragments option menu, select a fragment, and click in the Workspace. Repeat this process, clicking in a different place in the Workspace for each new fragment you want to place. The fragments are placed in a predefined orientation at the location you select, and are completed with hydrogen atoms to form molecules.

If you want to join these structures to form a single molecule, you can use the tools in the Connect/Fuse panel. See [Section 4.6 on page 46](#) for details.

You can build a single structure from multiple fragments in Place mode. Once you have placed the first fragment, you can join a second fragment to it by selecting the fragment from the fragments list and clicking in the Workspace on the atom that you want the new fragment to replace. If you click on a terminal atom, the fragment replaces the terminal atom. If you click on a non-terminal atom, that atom and its attached hydrogen atoms are replaced by the new fragment. For example, clicking on a secondary amine nitrogen with a methyl fragment gives a disubstituted methyl center, and clicking on a disubstituted methyl carbon with an amine fragment produces a secondary amine center.

You can also mutate one residue into another in a protein. When you select Amino acids from the Fragments menu, a Mutate button appears beside the Define grow bond picking tools. To mutate a residue, select Mutate, select the new residue from the list and pick an atom in the residue you want to mutate. The new residue replaces the old residue.

4.2.2 Building Structures Using Grow Mode

In Grow mode you can take advantage of the predefined “grow” rules that dictate where and how fragments are connected to build structures. To enter Grow mode, select Grow in the Fragments tab.

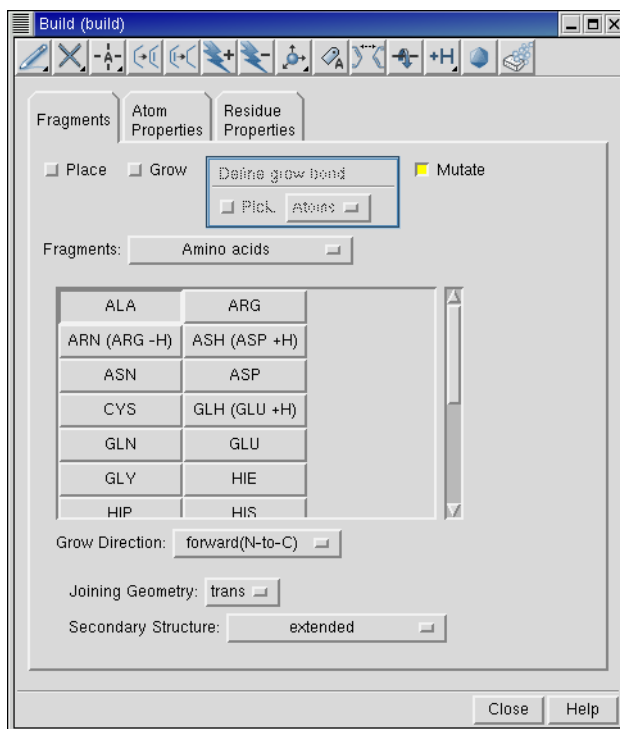


Figure 4.2. The Fragments *tab of the Build panel, showing the Mutate option.*

To create a structure in Grow mode, click the fragments in the fragment library that you want to assemble in the Workspace. At each click, the selected fragment is added to the structure, and a green arrow marks the location of the *grow bond*—the bond that will be replaced when the next fragment is added.

To grow a structure from a location other than that indicated by the current grow bond, you can choose a new grow bond using one of the following methods:

- Choose Atom from the Pick option menu of the Define Grow Bond section, and pick two atoms in the order that defines the grow direction.
- Choose Bond from the Pick option menu of the Define Grow Bond section, and pick a bond. The grow direction arrow points to the atom closest to where you click.
- Choose backward from the Grow Direction option menu. The grow bond moves to the opposite end of the molecule.

In addition to the fragments, the Fragments tab of the Build panel contains controls for setting the joining geometry, and, where applicable, secondary structure and helix geometry. The

Joining Geometry option menu provides context-sensitive options for attaching new fragments to the existing structure. For example, under the Deoxyribonucleic acid fragment type, the joining geometry options are Standard Helix and User Defined.

If you choose to define your own joining geometry, you must specify values for the torsions that define applicable secondary and tertiary structures. When you select the User Defined option and click on the name of the fragment, a dialog box prompting you to enter an angle value is displayed. Enter a value and click OK.

4.3 Building a Structure From Atoms

The toolbar of the Build panel contains tools for drawing structures in the Workspace atom by atom, and the Atom Properties tab and the toolbar contain tools for retyping atoms within a structure. Maestro does not automatically correct geometries of structures created in Draw mode. You can adjust the geometry after building the structure from the Adjust panel. You should also minimize the energy of the completed structure.

When using Draw mode to create a molecular structure, only a few elements—H, C, N, O, P, and S—are available. If you want to include different elements, use these elements to draw the structure, and then change the desired atoms—see [Section 4.4](#). Maestro allows you to violate the valency of an atom in Draw mode so that you can change it to an atom of some other element.

The basic operations of Draw mode are described below.

- To begin drawing with the currently selected element, click the Draw structures button.



The Draw structures button appears indented, indicating that you are in Draw mode. The default element is carbon.

- To place an atom, click in the Workspace.

The atom appears as an asterisk enclosed in a purple cube. The cube indicates that the atom is active, which means that it will be automatically connected to the atom that is placed next. Each subsequent click in the Workspace places another atom and bonds it to the previous atom with a single bond. To finish drawing, click the active atom. The cube vanishes.

- To start drawing from another location in the structure, deactivate the active atom by clicking on it, then click on the atom that you want to start drawing from. This atom becomes the active atom.

- To connect the active atom to an existing atom, click the existing atom.
- To create a multiple bond between the active atom and an adjacent atom, click the adjacent atom. Another bond is added between the two atoms.
- To change the element that you are drawing with, select the element from the Draw button menu. The next new atom you place will be an atom of the new element. If you connect existing atoms of a different element, the element of the existing atom is not changed: thus, it does not matter which element is selected when you draw a bond between two existing atoms.

The atoms that you draw are placed in the x-y plane. If you want to create a 3D structure, you can rotate the structure in the Workspace with the middle mouse button, then continue drawing atoms. Rotation does not change the active atom or exit Draw mode. You can also move atoms that you have placed using the Move button menu.



To move an atom in the x-y plane, choose X-Y from the button menu, click the atom you want to move, then click the new location. If the atom is in a ring, only the atom is moved. If the atom is not in a ring, the atom and its attached atoms are moved. To move an atom in the z direction, choose +Z or -Z from the menu, and click the atom to be moved. The atom is moved 0.5 Å in the chosen direction.

4.4 Changing Elements

Whether you build a molecule in Place mode, Grow mode, or Draw mode, the range of elements available to you is limited. To include other elements in your structure, you must first build with existing elements, then change the element to another element in the appropriate locations in the structure. You do not have to complete the structure before you make the change, you only have to place an atom at the desired location.

Changing an atom from one element to another can be done from the Build panel toolbar or from the Atom Properties tab. From the Set Element button menu in the toolbar



you can choose an element from H, C, N, O, P, or S, then click on an atom in the Workspace to change it to the selected element.

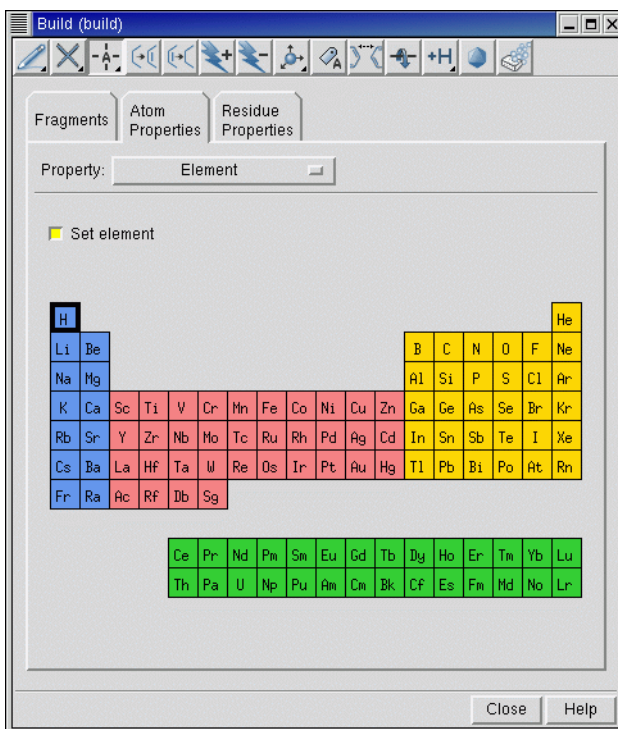


Figure 4.3. The Atom Properties *tab* showing tools for the Element *property*.

If you want to change an atom to some other element, choose Element from the Property menu in the Atom Properties tab. A periodic table is displayed, from which you can select an element, then click on an atom in the Workspace to change it to that element. You do not need to select Set Element: this is done automatically when you click on an element. However, you can select Set Element if it is unselected and you want to use the currently selected element.

4.5 Changing the Bond Order and Formal Charge

After creating a structure using the tools in the Build panel, you may wish to change some of the bond orders or formal charges. You can make these changes with the tools on the Build panel toolbar (as shown below) or by right-clicking a bond and choosing the appropriate item from the shortcut menu¹.

1. By default, right-clicking displays a context-sensitive shortcut menu. You can turn off this option by choosing Preferences from the Maestro menu and selecting the appropriate option in the Mouse tab.

- To increase the bond order, click the Increment bond order toolbar button



then click the bonds whose order you want to increase. Each click increments the bond order by one, up to a maximum of three. You can also increase the bond order in Draw mode by drawing another bond between the two atoms.

- To decrease the bond order, click the Decrement bond order toolbar button



then click the bonds whose order you want to decrease. Each click decrements the bond order by one, to a minimum of 0. Zero-order bonds are displayed with a green line.

When you change the bond order, the valence states of the atoms can be violated. Maestro gives no warning of this violation, so it is up to you to ensure that your structure is valid.

- To increase the formal charge, click the Increment formal charge toolbar button



then click the atoms whose formal charge you want to increase. You can only increase the formal charge one atom at a time. However, the formal charge is also set when you change the MacroModel atom type.

- To decrease the formal charge, click the Decrement formal charge toolbar button



then click the atoms whose formal charge you want to decrease. Like increasing the formal charge, you can only decrease the formal charge one atom at a time, and the formal charge is set when you change the MacroModel atom type.

4.6 Fusing or Connecting Structures

To build some structures, it is sometimes convenient to connect or to fuse two existing structures. Connecting two structures creates a new bond between two selected atoms, with deletion of appropriate hydrogens. An example is the creation of bipyridyl from two pyridine rings, by connecting the 2-carbon on one ring with the 5-carbon on the other, and deleting the hydrogen atoms. Fusing two structures deletes atoms. An example is the fusion of two benzene rings, in

which two carbon and two hydrogen atoms are deleted from one ring, and bonds formed to the other ring with deletion of two hydrogen atoms.

The tools for fusing and connecting structures are available from the Connect & Fuse panel. In the connection and fusion operations, an attempt is made to create structures with reasonable geometries, and a check is made that valency is not violated. This is in contrast to Draw mode, which does not evaluate bonding or geometry.

For both connect and fuse operations, you must first select atom pairs to be joined. The atoms you select must not be hydrogen atoms. You can select as many pairs as you like. The first atoms in all selected pairs must all come from the same molecule, and the second atoms in the pairs must come from the same molecule. The molecule for the set of first atoms set need not be different from the molecule for the set of second atoms. If you select a pair of atoms in the same molecule, a new bond is drawn between the atoms.

To select atoms for a connect or fuse operation, ensure that Pick to define atom pairs is selected and click on the Workspace atoms to define the pairs. Atom numbers of the selected pairs appear in the Connect or fuse atoms list. Also, if Show Markers is selected, each pair is marked with a dashed green line. For the pair that is selected in the list, the green line is flanked by solid yellow lines. The lines disappear once the connect or fuse operation is completed.

When fusing structures, the order in which you select the atom pair is important. The atom selected first replaced by the second selected atom. For example, when fusing an amino nitrogen and an ethyl carbon, selecting the nitrogen first results in its deletion, followed by attachment of the remaining atoms of the amine to the ethyl carbon. Selecting the ethyl carbon first results in deletion of the carbon and retention of the nitrogen.

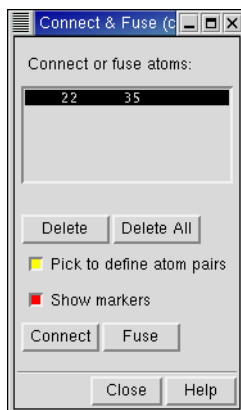


Figure 4.4. The Connect & Fuse panel.

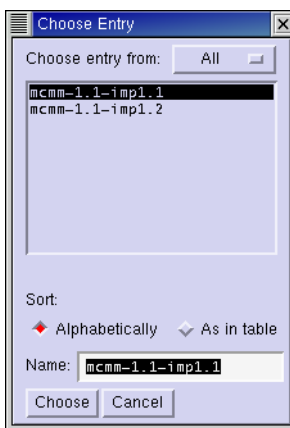


Figure 4.5. The Choose Entry *dialog box*.

To delete an atom pair from the Connect or Fuse Atoms list and cancel the connection or fusion of those atoms, select the pair in the list and then click Delete. To clear the atom pair list, click Delete All.

When you have selected all the pairs of atoms that you want to connect or fuse, click the appropriate button—Connect or Fuse.

If you are performing a connection or fusion between structures belonging to different project entries (the source entries), the result must belong to a single entry. When you click Connect or Fuse, the Choose Entry dialog box is displayed, and you must specify the entry with which the new structure should be associated. Subsequent actions depend on which entry you choose:

- If you choose one of the two source entries, the other entry is deleted.
- If you choose an existing entry other than the two source entries, the new structure is added to the existing entry and the source entries are deleted.

The Choose Entry dialog box provides some tools for selecting the destination entry: a menu to limit the entries listed, and some sort options.

- If you provide a new entry name in the Name text box, the source entries are deleted, and a new entry is created that contains the fused or connected structure.

You can undo a connect or fuse operation, and the changes to the entries are also undone.

4.7 Adjusting the Geometry

Before you use a structure for a calculation, you might want to adjust the geometry, especially if you used Draw mode, or if the Grow mode orientation was incorrect.

4.7.1 Adjusting Parameters

In addition to the adjustments you can make with the Move button, described on [page 38](#), you can also change bond lengths, bond angles, dihedral angles, and atom chirality. You can make these adjustments either in the Workspace directly or by using the tools in the Adjust panel.

In the Workspace

The basic process for adjusting parameters directly in the Workspace is the same for bond lengths, bond angles, dihedral angles, and chirality. This description uses bond length as an example:

1. Click the Adjust distances, angles or dihedrals button in the main toolbar:



and choose Distance.

2. In the Workspace, click on two atoms to define the bond.

The bond is highlighted with two yellow solid lines, a purple dotted line, and an arrow indicating the direction of adjustment. (Angles are highlighted in green and orange, and dihedrals are highlighted in red and turquoise.)

3. With the middle mouse button, click on the defined bond and drag to the right to lengthen the bond or to the left to shorten it.

You can also choose bond lengths or dihedrals for adjustment by right-clicking on a bond in the Workspace and choosing Distance or Dihedral from the shortcut menu.

From the Adjust Panel

You can open the Adjust panel from the Edit menu in the main window, or with the Display Adjust panel button on the Build panel toolbar.



The tools for adjusting bond lengths, bond angles, and dihedral angles are similar, so the description given here is generic.

To define the geometric parameter that you want to adjust, choose Atoms or Bonds from the Pick option menu in the Define atoms for *operation* section, and pick the atoms or bonds. For a bond length, you must pick two atoms or one bond; for a bond angle, you must pick three atoms or two bonds; for a dihedral angle you must pick four atoms or one bond—the bond about which the rotation is performed. When you pick this bond, the groups attached to the atom nearest the point you click are the ones that are moved.

In addition to the standard picking tools just described, you can also specify dihedrals in proteins using standard names (psi, phi, omega, etc.). Click Select in the Dihedral tab and select dihedrals by the residues in each entry (or chain), or by standard names.

The atom numbers of the selected atoms are listed in the table in the center of the tab. You can pick as many sets of atoms as you like before doing any adjustments. The last atom that you pick is the “target atom”—the atom that is moved in the adjustment process.

Once you have chosen the atoms, you must decide which of the attached atoms to move along with the target atom. You can choose from All attached atoms, All attached terminal atoms, or Only single atom. If the target atom is part of a ring, the first option is unavailable.

To make the adjustment, you can enter a new value for the geometric parameter in the text box, or you can use the mouse. If you enter a new value, the value is changed and the Workspace updated when you press ENTER. To adjust the geometry using the mouse, select Allow adjusting then drag with the middle mouse button to adjust the geometry.

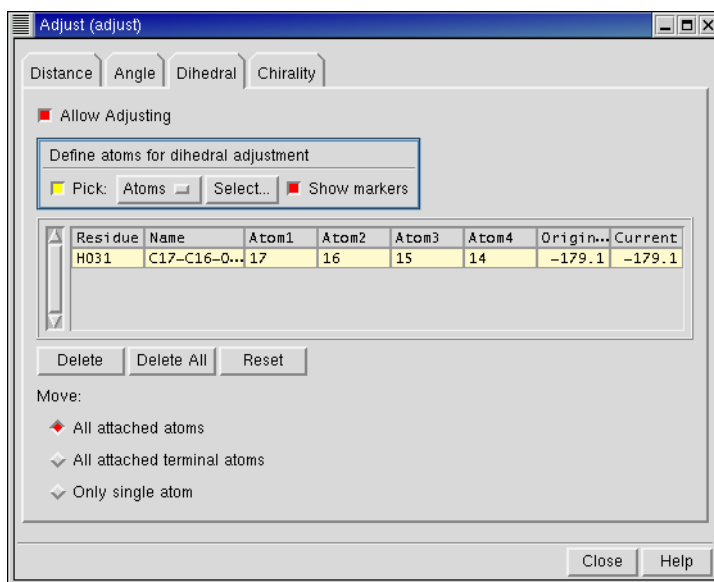


Figure 4.6. The Adjust panel, showing the Dihedral tab.

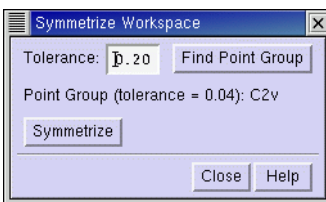


Figure 4.7. The Symmetrize Workspace panel.

To adjust the stereochemistry at a chiral center, use the picking tools in the Chirality tab. You can invert the chirality of selected chiral centers or of all chiral centers. This tab also has buttons for displaying the R and S labels for the configuration at chiral centers, and for clearing the labels. If you click Label Chiral Centers, all other labels are cleared.

4.7.2 Cleaning up the Geometry

Symmetrizing structures

The Geometry Symmetrizer allows you to symmetrize the geometry of the structure in the Workspace, subject to a given tolerance.

To use the Geometry Symmetrizer:

1. Choose Symmetrize Workspace from the Edit menu, or click the Symmetrize Workspace toolbar button in the Build panel.



2. Click Find Point Group and check to see whether the desired point group is found for the structure in the Workspace. If not, increase the tolerance and click Find Point Group again.
3. When the desired point group is found, click Symmetrize.

Cleaning up the Geometry

To clean up the structure in the Workspace, the Geometry Cleanup utility minimizes the energy of the structure using a Universal Force Field (UFF). It realigns atoms that are too crowded and corrects unphysical geometry that can occur when building molecules from the Build panel. In general, it eliminates distortions that can cause jobs to fail, such as MacroModel or Jaguar jobs.

Note: Geometry Cleanup only works for structures with less than 300 atoms, and can be slow for structures with more than 200 atoms.

To run this utility, choose Clean Up Geometry from the Edit menu, or click the Clean up geometry toolbar button in the Build panel.



4.8 Applying a Hydrogen Treatment

When building molecules or preparing proteins, you might want to apply a hydrogen treatment that is compatible with the force field you intend to use in calculations on those structures. You can also use a hydrogen treatment to add hydrogen atoms to a structure that is used for any program that requires explicit hydrogen atoms, such as QikProp or Jaguar. Hydrogen treatments can treat hydrogen atoms explicitly or implicitly, as part of another atom. To select and apply a hydrogen treatment, you can use the controls in the Hydrogen Treatment panel, which you open from the Edit menu.

To select a hydrogen treatment, choose an option from the Treatment option menu. A description of the treatment is displayed in the Description text box, including the force fields for which the treatment is valid. This information is also given in [Table 4.1](#).

Table 4.1. Appropriate hydrogen treatments for Maestro-supported force fields.

Hydrogen Treatment Method	MM3	MM2	AMBER	AMBER94	MMFF	MMFFS	OPLS-AA, OPLS_200x	OPLS
All-atom with No-Lp	X	X		X	X	X	X	
All-atom with Osp3-Lp		X						
Csp3 United-atom with S-Lp ^a			X					
Csp3 United-atom with No-Lp ^a								X
All atom with S-Lp			X					
Csp2/sp3 United atom with No-Lp ^a								X

- a. In a united atom representation, hydrogen atoms are incorporated into the dimensions of the heavy atom to which they are attached. That is, they are implicit.

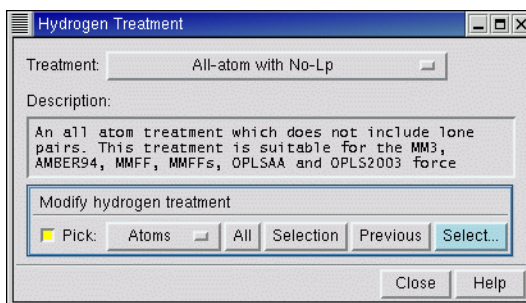


Figure 4.8. The Hydrogen Treatment panel.

You can apply the selected treatment to the entire structure or to selected parts of the structure. To apply the treatment to the entire structure click **All** in the Modify Hydrogen Treatment section. To apply the treatment to selected atoms, use the picking tools in the Modify Hydrogen Treatment section or click **Select** to open the Atom Selection dialog box and define the selected atoms. If the structure already complies with the treatment you selected, no changes will be made when you apply the treatment.

You can also apply the current hydrogen treatment to selected atoms, molecules, residues, chains, or entries from the main toolbar or the Build panel toolbar. Select one of these entities from the Add hydrogens button menu,



then click on atoms in the Workspace to apply the treatment. Double-click this button to apply the treatment to all atoms in the Workspace.

Each time you select a force field in the Potential section of a MacroModel or Impact energy panel, the hydrogen treatment is automatically changed to comply with that force field, and the hydrogen treatment for any atoms you subsequently draw in the Workspace will have the new hydrogen treatment. If you attempt to start a MacroModel, Impact, or Jaguar job using an input structure with an incompatible hydrogen treatment, the job fails.

4.9 Changing Atom Properties

In addition to changing the element and the formal charge, you can change a variety of other atom and residue properties. The tools for making changes are in the Atom Properties tab. Each choice of property from the Property menu displays a set of tools that you can use to change the property. Apart from the Element property, the properties are listed below with a description of the tools. Before you change atom properties, it may be helpful to display atom labels. For atom labeling information, see [Section 6.3 on page 78](#).

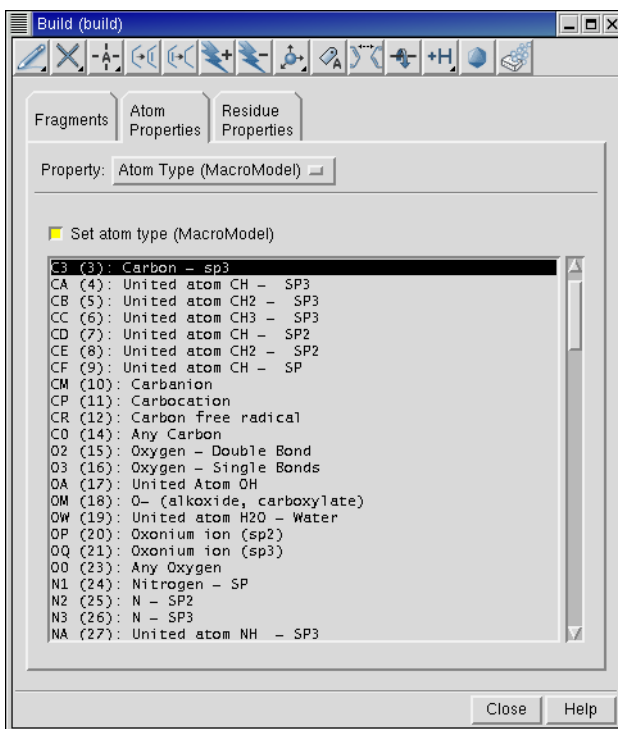


Figure 4.9. The Atom Type (MacroModel) tools in the Atom Properties tab.

Atom type (MacroModel)

A list of atom type symbols with a brief description is displayed (see [Figure 4.9](#)). The atom type symbols represent atoms with specific properties, such as sp^3 hybridization for a C atom with type symbol C3, or a free radical for a C atom with type symbol CR. The use of atom types is provided primarily for MacroModel users, though Maestro always assigns an atom type to each atom. For information on atom types, see [Appendix B](#). To retype an atom, choose the atom type and click on the atoms in the Workspace that you want to change.

Partial Charge

To apply new partial charges to atoms of interest, enter new charges in the Partial Charge and Solvation Charge text boxes, and use the picking tools to apply the partial charges to the atoms of interest. Both of these charges are used in MacroModel calculations. The solvation charge is used in MacroModel GB/SA solvation calculations. Generally, these two values should be the same. Default values for these charges are supplied with MacroModel force fields. You can also use partial charges from Jaguar electrostatic potential calculations, which are read by default when the results Jaguar calculations are incorporated or imported.

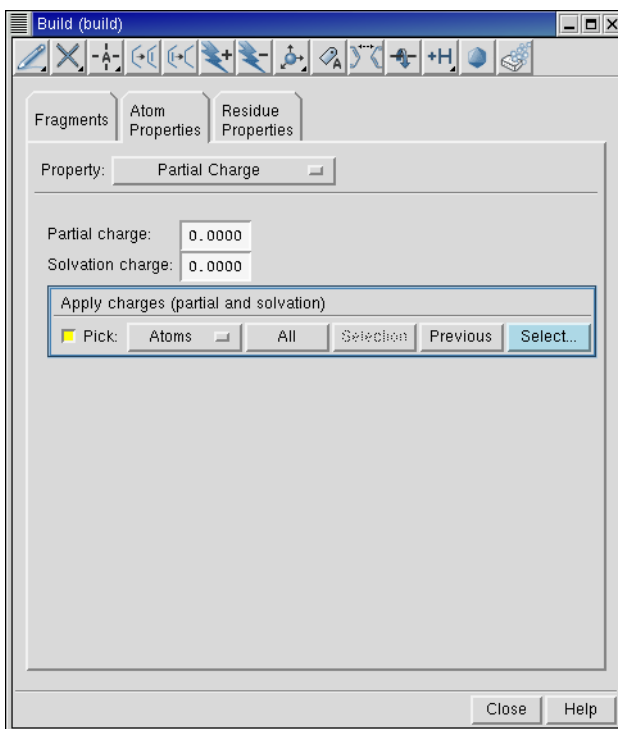


Figure 4.10. The Partial Charge tools in the Atom Properties tab.

PDB Atom Name

A PDB atom name is a combination of a chemical symbol and a connection indicator. The first two characters contain the element symbol, right-justified; the third character is a remoteness indicator (A=alpha, B=beta, G=gamma, and so on), and the fourth, if it is needed, is a numeric branch designator. For instance, the name “_CA_” represents an alpha carbon, “_CB_” a beta carbon, and “_HB_” a hydrogen on a beta carbon. In these descriptions, the underscore characters represent spaces, which are required for compliance with PDB semantics: “CA__” is a calcium atom, and is different from “_CA_”.

Occasionally Maestro is not able to assign PDB names to all the atoms. To manually apply missing PDB atom names, enter the desired name in the PDB atom name text box and use the Apply PDB atom name picking tools to select the atoms to which you want to apply the name.

For some purposes, such as for ligands in Prime induced fit calculations, it is useful to have unique atom names within residues. To make the atom names unique, use the Set unique PDB atom names within residues picking tools. The unique names are composed of the element symbol, right justified in the first two characters, followed by a numeric index for that element.

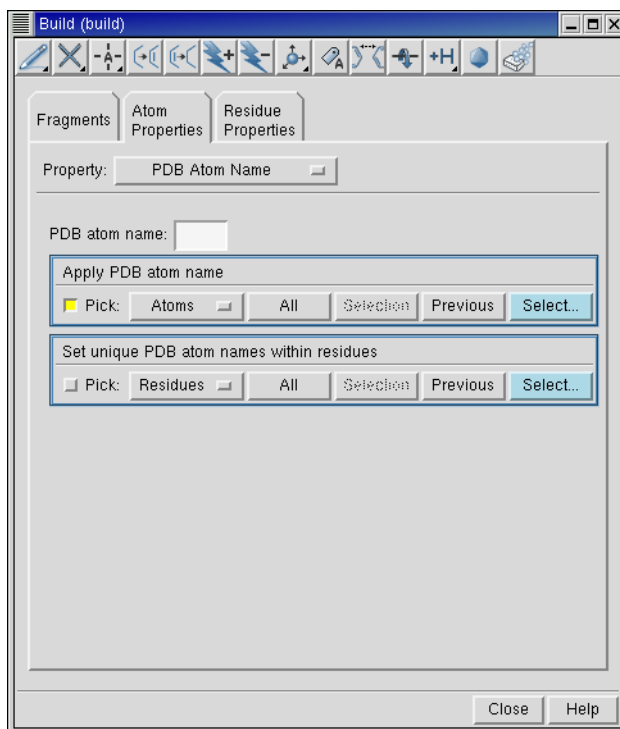


Figure 4.11. The PDB Atom Name *tools in the Atom Properties tab.*

Grow Name

This name is the name used internally for a fragment when growing a structure. The intended use of this property is labeling new fragments for a fragment library.

Atom Name

You can assign one of three types of names to atoms in the Workspace: user-specified, canonical, or unique. User-specified names can be any length. To assign a user-specified name, type the desired name in the Atom name text box of the Set atom name section, and use the picking tools to apply the name. You must display the atom name labels before you can see the new name.

A canonical atom name is an automatically assigned combination of atom property pairs: Element & Atom Number or Residue Number & PDB Atom Name. To assign canonical names, select the desired pair of properties from the Composition option menu, specify a separator symbol in the Separator text box, and click Apply To All Atoms. You cannot selectively apply canonical names.

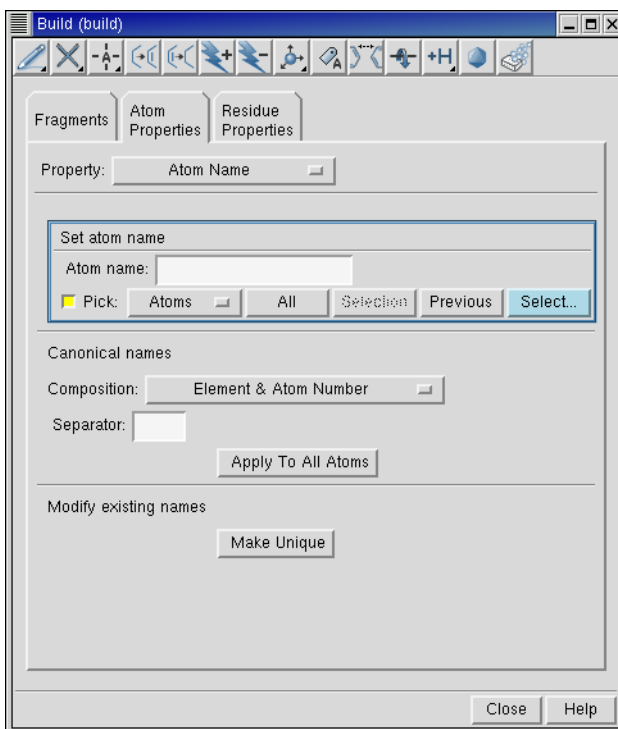


Figure 4.12. The Atom Name tools in the Atom Properties tab.

Unique names are generated using the element and atom number of each atom. If this naming scheme produces atoms with identical names, Maestro adds numeric suffixes, followed by tilde (~) symbols to distinguish the names. To apply unique names to all atoms, click Make Unique in the Modify existing names section. As for canonical names, you cannot selectively apply unique names.

4.10 Changing Residue Properties

Residue properties can be changed in the Residue Properties tab. This tab is laid out in a similar fashion to the Atom Properties tab, with a Property menu and a set of tools displayed for each choice. The residue properties and their tools are described below. In order to view the results, you must first display the appropriate atom labels.

Residue Number

To renumber residues, in the Renumber Residues section enter the starting value in the Renumber From text box, then use the picking tools to select the residues to renumber. Renum-

bering a single residue is the same as assigning its number. To apply a specific number to one or more residues, in the Apply residue number section enter the number you want to apply in the Residue number text box, then use the picking tools to select the residues to apply the number to.

Residue Name

To name residues, enter the residue PDB name in the Residue PDB name text box, and use the picking tools to select the residues to name. PDB residue names are three characters long and are right-justified. For example, the residue name G3 must be entered with a leading space, as “ G3”, whereas the residue name C19 should have no leading space, “C19”.

Chain Name

You can assign a single-letter name to chains. The definition of a chain is not restricted to the items labeled Chain in the pick option menus: you can define a chain in any way that you like, and assign it a name. To assign the name, enter a character in the Chain Name text box and use the picking tools to define the chain.

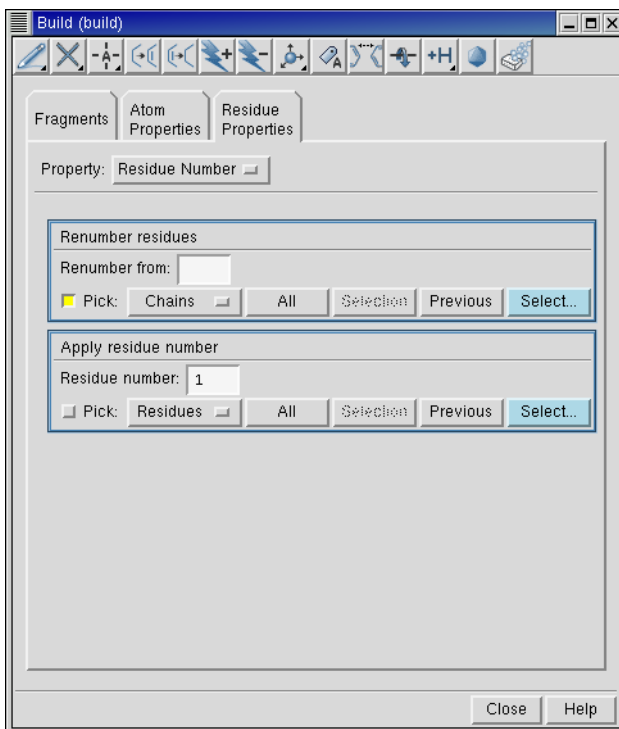


Figure 4.13. The Residue Number tools in the Residue Properties tab.

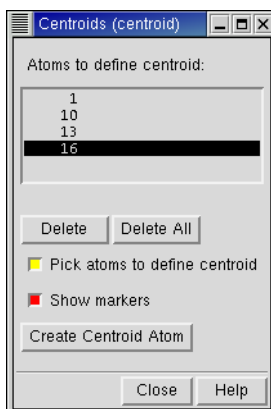


Figure 4.14. The Centroids *panel*.

4.11 Defining Dummy Atoms

Dummy atoms are used by several of the computational programs for defining geometries and geometric parameters. You can retype an atom to a dummy atom type using the MacroModel Types atom property described earlier. You can also place a dummy atom at the centroid (mean position) of a group of atoms, such as a ring. You can then change the dummy atom to a real atom, and thus place real atoms at the centroid of some system of interest.

The centroid capability is available in the Centroids panel, which you open from the Edit menu. To define a centroid, select Pick atoms to define centroid, pick the atoms in the Workspace. The atoms are listed in the Atoms to define centroid list. If you have Show markers selected, the atoms are marked with a salmon-colored star, and the atom that is selected in the list is marked with a lavender-colored star. You can delete atoms from the list by selecting the list item and clicking Delete. Once you have picked all the atoms, click Create Centroid Atom. The To define another centroid, click Delete All before picking atoms for the new centroid.

4.12 Deleting Atoms

If you need to delete atoms from a structure, you can use the Delete button menu on either the main toolbar or the Build panel toolbar.



This button menu offers choices for the type of object to delete, and allows you to open the Atom Selection dialog box to select atoms for deletion based on their properties.

When you delete atoms, the terminal atoms that are attached to this atom can be deleted, or the terminal atoms can be left in place. Terminal atoms are atoms that are bonded only to that atom, such as hydrogen or a halogen. Deleting terminal atoms would not, for example, delete a methyl group attached to the deleted atom. This action is controlled by an option that you can set in the Builder tab of the Preferences panel.

If you want to delete all atoms in the Workspace, use the Clear Workspace button, which is located on the main toolbar.

Selecting Atoms

Maestro has a powerful set of tools for selecting atoms in a structure: toolbar buttons, picking tools in panels, and the Atom Selection dialog box. These tools allow you to select atoms in two ways:

- Select atoms first and apply an action to them
- Choose an action first and then select atoms for that action

5.1 Toolbar Buttons

The small triangle in the lower right corner of a toolbar button indicates that the button contains a menu. Many of these buttons allow you to choose an object type for selecting. Choose Atoms, Bonds, Residues, Chains, Molecules, or Entries, then click on an atom in the Workspace to perform the action on all the atoms in that structural unit.

For example, to select atoms in residues with the Workspace selection button:

1. Choose Residues from the Workspace selection toolbar button:



The button changes to:



2. Click on an atom in a residue in the Workspace to select all the atoms in that residue.

You can also drag to select multiple residues.

5.2 Picking Tools

The picking tools are embedded in each panel in which you need to select atoms to apply an operation. The picking tools in a panel can include one or more of the following:

- Pick option menu—Allows you to choose an object type for selecting. [The Pick Option Menu](#) is described on [page 63](#).
- All button—Performs the action on all atoms in the Workspace.

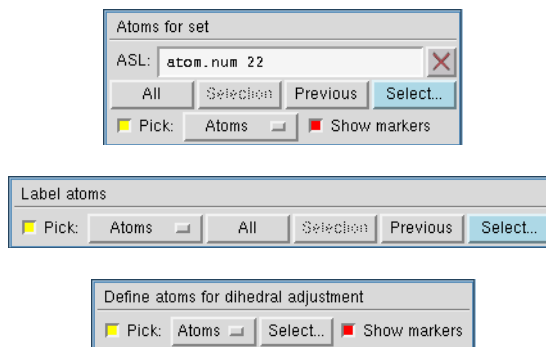


Figure 5.1. Examples of picking tools. Top: from the Sets panel; Center: from the Atom Labels panel; Bottom: from the Adjust panel.

- Selection button—Performs the action on any atoms already selected in the Workspace.
- Previous button—Performs the action on the most recent atom selection defined in the Atom Selection dialog box (see [Section 5.3 on page 64](#)).
- Select button—Opens the Atom Selection dialog box (see [Section 5.3 on page 64](#)).
- ASL text box—Allows you to type in an ASL expression for selecting atoms.

ASL stands for Atom Specification Language, and is described in detail in the [Maestro Command Reference Manual](#).

- Clear button—Clears the current selection.



- Show markers option—Marks the selected atoms in the Workspace

The range of picking tools available in a panel varies according to the types of structural units that are needed for the operation and whether picking causes immediate action or requires an action to be specified later. The ASL text box and Clear button are included only when picking merely defines a set of atoms for later use. For example, the Label Atoms section of the Atom Labels panel has a Pick option menu, an All button and a Select button, and the action is immediate. The Sets panel has a full set of picking tools, and stores the atoms in a defined set that can be used later.

The Pick Option Menu

With the Pick option menu, you can choose an object type (Atoms, Bonds, Residues, Chains, Molecules, or Entries) and click on an atom or a bond in the Workspace to select all the atoms

in that structural unit. For instance, choose **Atoms** and click on an atom in the Workspace to select just that atom. Choose **Residues** and click on an atom to select all the atoms in that residue. Choose **Bonds** and click on a bond to select both atoms that define that bond.

When you select an object type, the cursor in the Workspace appears as a box with a letter: A, B, R, C, M, and E represent Atoms, Bonds, Residues, Chains, Molecules, and Entries, respectively. Successive clicks add atoms to the selection.

The number of object types available in the menu depends on the operation. For instance, when you define a dihedral angle, you can pick either four atoms or three bonds. The Pick option menu for this operation therefore only contains Atoms and Bonds.

You do not need to select the Pick option before choosing an item from the Pick option menu; it is automatically selected when you choose a menu item.

5.3 The Atom Selection Dialog Box

If you want to make atom selections based on more complex criteria, such as selecting all the carbon atoms in a protein backbone, you can use the Atom Selection dialog box. The Atom Selection dialog box allows you to select atoms by defining complicated ASL expressions without the need to understand the ASL syntax.

To open the Atom Selection dialog box, choose **Select** from a button menu or click **Select** in the picking tools of a panel. The Atom Selection dialog box is modal: you must dismiss it before you can perform any other operations, with the exception of rotation and translation of structures in the Workspace, and selecting atoms.

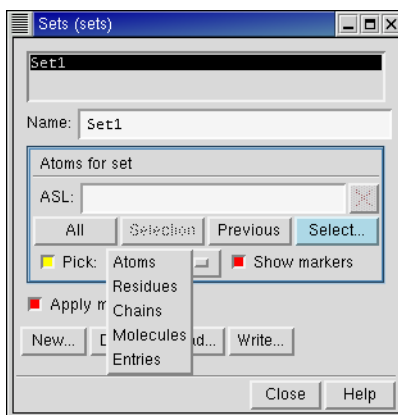


Figure 5.2. The Sets panel, showing the Pick option menu.

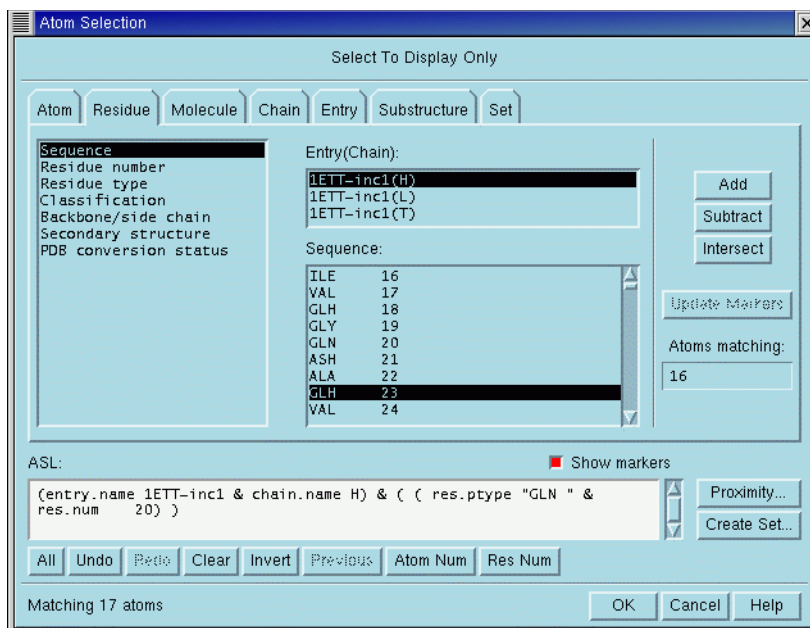


Figure 5.3. The Atom Selection *dialog box*.

When you close the Atom Selection dialog box, the ASL expression is automatically stored and can be recalled with the Previous button on any set of picking tools (see [Section 5.2 on page 61](#)).

To select atoms

1. Create an atom group by selecting a property and value (see [Section 5.3.1 on page 66](#)).
2. Click Add or Subtract to define the current selection (see [Section 5.3.2 on page 66](#)).
3. Create more atom groups as needed and click Add, Subtract, or Intersect to modify the current selection (see [Section 5.3.2 on page 66](#)).
4. When you are satisfied with the current selection, click OK to apply the operation to the selection.

Some operations, such as deleting atoms, take effect immediately. Others merely define a set of atoms to be used in a subsequent task, such as selecting atoms for creation of a surface.

The features of the Atom Selection dialog box are described in the following sections.

5.3.1 Selecting Atoms by Property

You can create an atom group by selecting a property and a value for that property. These properties are grouped into a number of categories, each in their own tab: Atom, Residue, Molecule, Chain, Entry, Substructure and Set. For each category there is a list of properties shown on the left of the tab, from which you can select one item. When you select an item, tools for that property are displayed in the center of the tab. The three types of tools are:

- Text box. You can specify the allowed values for the property by using lists of numbers, closed ranges (e.g. 12-42) and open ranges (e.g. ≥ 30).
- List. You can choose one or more entries from the list by clicking, in combination with the SHIFT and CTRL keys. Only items that are currently in the Workspace are displayed in the list.
- Set of options, which can be nested. You can choose only one of the possible combinations by selecting the options.

With text boxes and lists, you can also choose items by selecting atoms in the Workspace. For example, if you choose the Element property in the Atom category, then select a carbon atom in the Workspace, the “C” item in the list of elements is selected.

The current atom group is represented in two ways: the atoms are marked in purple in the Workspace (if Show markers is selected) and the value in the Atoms Matching field indicates how many atoms are matched by the current property specification.

5.3.2 Combining and Modifying Atom Selections

Once you have created an atom group, you can combine it with the current selection. The Add, Subtract, and Intersect buttons combine the ASL expression for the atom group with the existing ASL expression in the ASL text box (the current selection), as follows:

- The Add button includes the atom group in the current selection by performing a boolean OR operation on the ASL expressions.
- The Subtract button excludes the atom group from the current selection by performing a boolean AND NOT operation on the ASL expressions.
- The Intersect button restricts the selection to the atoms that are in both the atom group and the current selection by performing a boolean AND operation on the ASL expressions.

Note: Intersect is not available for the first atom group that is defined.

In addition to combining ASL expressions to modify the current selection, you can use the following buttons:

- All—select all atoms in the Workspace
- Undo—undo the last operation
- Red—redo the last operation
- Clear—clear the ASL expression
- Invert—negate an entire expression
- Previous—recall the atom selection defined the last time the Atom Selection dialog box was used
- Atom Num—convert the ASL expression to list each atom explicitly by atom number
- Res Num—convert the ASL expression to list each residue explicitly by residue number

If Show Markers is selected, the atoms corresponding to the current selection (the expression in the ASL text box) are highlighted with light blue markers in the Workspace.

5.3.3 Selecting Atoms by Proximity

Once you have made an initial selection, you can select atoms based on their distance from the atoms defined by the current ASL expression (the source atoms) in the Proximity dialog box, see [Figure 5.4 on page 68](#). To open this dialog box, click Proximity.

You can specify the distance in angstroms or in terms of the number of bonds. If the boundary for inclusion intersects a molecule or a residue, you can choose to include the entire molecule or residue using the Fill buttons. You can either include or exclude the source atoms in the selection. For example, you might want to include residues within a certain distance of a metal atom, but not include the atom itself. The Update Markers button updates the markers for the proximity selection so that you can view the selection in the workspace before accepting it.

5.3.4 Editing and Storing Expressions

You can undo multiple operations with repeated clicks on the Undo button. The Redo button reverses undo operations, and can also be used multiple times. (These are the buttons in the dialog box, not the toolbar button.)

You can also type ASL expressions directly into the ASL text box. The hand-edited expressions are considered part of the last operation performed with the mouse and can be undone and redone as part of that operation, but cannot be undone or redone separately.

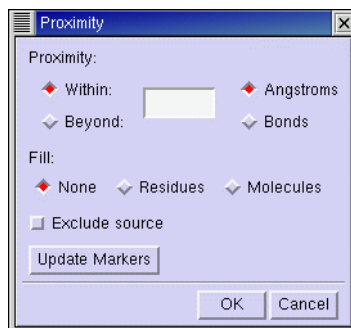


Figure 5.4. The Proximity *dialog box*.

The Create Set button opens a dialog box that enables you to name and store the existing ASL expression for use as a set definition. The sets created in this manner are the same as sets created with the Sets panel, and are listed in the Set tab. You can therefore create sets from the Atom Selection dialog box, and use them for atom selection later. For more information on creating and saving sets, see [Section 9.4 on page 124](#).

5.3.5 Examples of Atom Selection

Below are some examples of the use of the Atom Selection dialog box. It is assumed that you have the Atom Selection dialog box open. To apply your selection and dismiss the dialog box, click OK. These examples assume that the protein and the ligand are in separate entries in the Workspace.

To select all carbon atoms in the backbone of a protein:

1. In the Atom tab, select Element from the list of properties.
2. Select “C” from the list of elements.
3. Click Add to add the atom group to the ASL expression.
4. Select Backbone/Side Chain from the list of properties.
5. Select the Backbone option.
6. Click Intersect to restrict the current ASL expression to be applied only to protein backbone atoms.
7. If you want to continue to the next example, click Clear to clear the current ASL expression.

To select all side chains within 5.0 Å of a ligand in the active site of a ligand/receptor complex:

1. In the Entry tab, select the entry name for the ligand from the list of entries.
2. Click Add to make this entry name the current ASL expression.
3. Click Proximity to open the Proximity dialog box
4. Enter 5.0 in the text box, click Residues in the Fill section, select Exclude source, then click OK.
5. In the Residue tab, select the Backbone/Side Chain property.
6. Select the Sidechain option.
7. Click Intersect to restrict the current ASL expression to include only atoms in the protein side chain.
8. If you want to continue to the next example, click Clear to clear the current ASL expression.

To select a range of residues but exclude proline residues:

1. In the Residue tab, select Sequence from the property list.
2. Select the appropriate chain from the list of chains in the top scrolling list.
3. Select a range of residues in that chain from the lower scrolling list.
4. Click Add to make this residue range the current ASL expression.
5. Select the Residue Type property.
6. Select “PRO” from the list of residue types.
7. Click Subtract to exclude proline residues from the current ASL expression.

Displaying Structures

Maestro provides a wide selection of options for displaying structures. For example, you can change the color of some or all of the atoms, change the molecular representation in which the atoms are depicted, label atoms by their properties, or selectively undisplay certain atoms so that you can better examine small subsections of a structure.

6.1 Changing Atom Color

Changing the color of a structure can help you to visualize its properties and its three dimensional structure. The Atom Coloring panel controls allow you to apply a new color to all or some of the Workspace atoms, or to use predefined or user-defined schemes to color specific atom groups. You can also use the shortcuts available in the toolbar to change the atom coloring.

To open the Atom Coloring panel, choose Atom Coloring from the Display menu.

6.1.1 Applying a Single Color

To apply a single color scheme to some or all atoms in the Workspace, you can use the tools in the Atom Color tab. This tab contains a color palette in the By Palette subtab (see [Figure 6.1 on page 70](#)), and a list of color names in the By Name subtab [Figure 6.2 on page 70](#)). You can use either of these to select a color. Once you have selected a color, use the picking tools in the subtab to apply the color to selected atoms, or click the All button to apply the color to all atoms in the Workspace.

You can also apply a color from the Color residues by constant color button menu in the toolbar.



This button displays a menu from which you can choose from a selection of colors. The color choice is persistent. You can then pick residues to apply the color, or double-click the button to apply the color to all atoms. If the color you want to use is not on the menu, you can use the Atom Color tab to select the color. For information on adding custom colors, see [Appendix E](#).

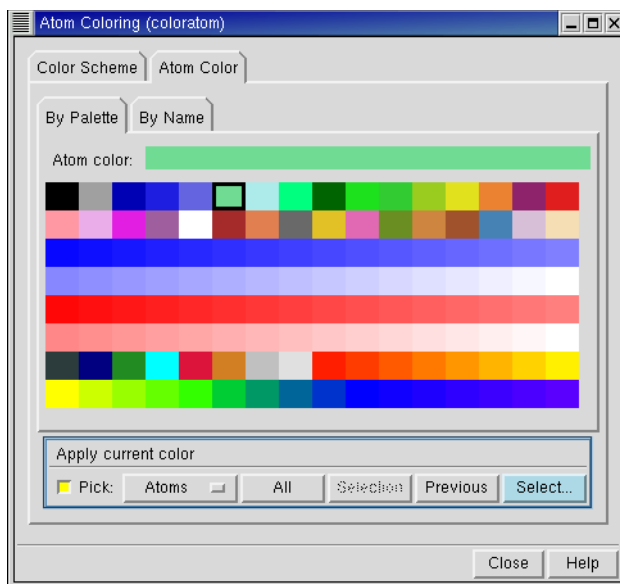


Figure 6.1. The By Palette *tab of the* Atom Coloring *panel.*

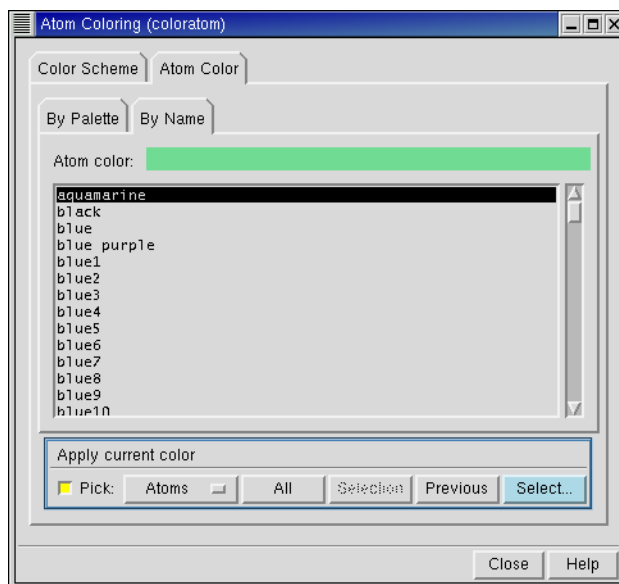


Figure 6.2. The By Name *tab of the* Atom Coloring *panel.*

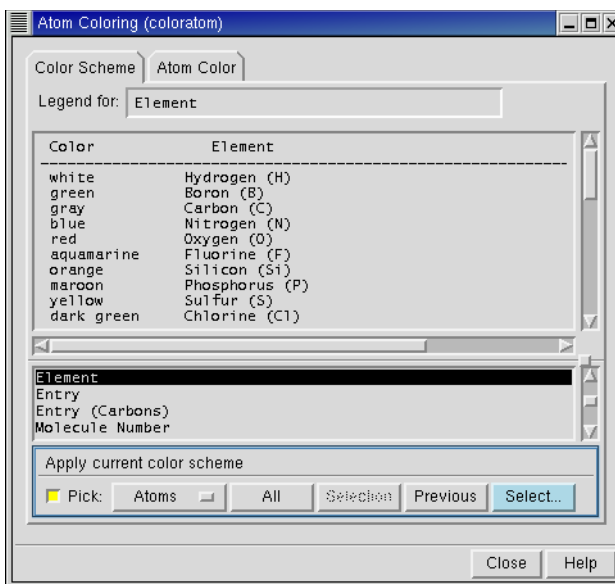


Figure 6.3. The Color Scheme *tab of the* Atom Coloring *panel.*

6.1.2 Using Predefined Schemes to Color Atoms

To apply coloring to some or all Workspace atoms using one of Maestro's predefined color schemes, you can use the tools in the Color Scheme tab (see [Figure 6.3](#)). Click on the desired scheme name in the list in the lower portion of the panel. The legend at the top of the panel describes the colors that will be applied and the criteria used to apply them. Once you have selected a scheme, use the picking tools to select the atoms you wish to color, or click the All button to apply the selected color scheme to all atoms in the Workspace.

You can also apply a color scheme to all atoms in the Workspace from the Color all atoms by scheme button menu in the toolbar.



When you click and hold on this icon a menu is displayed from which you can choose the scheme. The color scheme is applied when you make your menu choice. The schemes are described in the Color Scheme tab.

For information on creating custom color schemes, see [Appendix E](#).

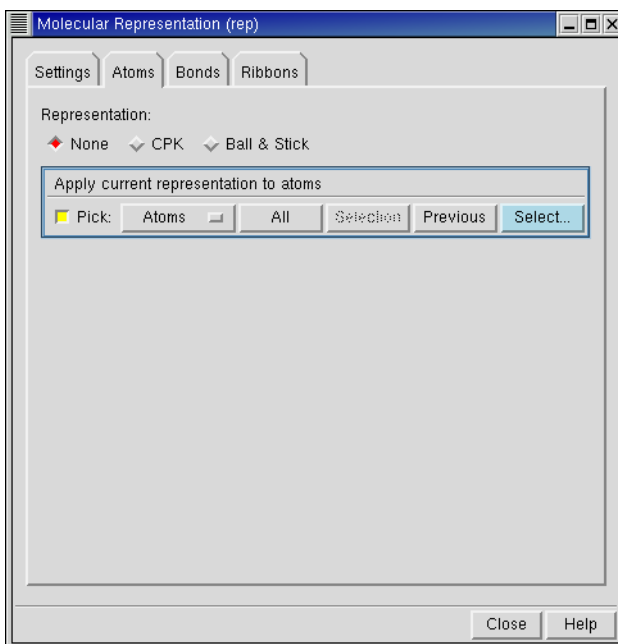


Figure 6.4. The Atoms *tab of the* Molecular Representation *panel*.

6.2 Changing Molecular Representations

The Molecular Representation panel contains four tabs, Settings, Atoms, Bonds, and Ribbons. The Ribbons tab contains controls for representing the structure of proteins. The first three tabs contain controls for setting molecular representation for the atoms and bonds in the Workspace and for setting the molecular representation in which subsequently drawn atoms and bonds are displayed. Using the controls in these tabs, you can choose to change the molecular representation for all atoms and bonds in the Workspace, for only the atoms currently in the Workspace, or for only the bonds currently in the Workspace. You can also use the shortcuts available in the toolbar to change the molecular representation.

To open the Molecular Representation panel, choose Molecular Representation from the Display menu.

6.2.1 Changing the Representation of Atoms

In the Atoms tab, you can choose between None, CPK and Ball & Stick for the representation of atoms. The None option means that only the bonds are depicted. The bond representation is set in the Bonds tab. To change the representation of the atoms, use the picking tools or the Select button in the Apply current representation to atoms section to select atoms to apply the repre-

sentation, or click All to apply the representation to all atoms. You can mix representations within a single molecule by picking.

The Ball & Stick option sets a representation for both atoms and bonds. When you pick an individual atom in this representation, half of each bond attached to the picked atom is displayed in stick form.

You can apply a molecular representation using the Draw atoms in CPK representation and Draw atoms in Ball & Stick representation buttons on the toolbar. Click either button to pick the objects to apply the representation to. Double-click either button to apply the representation to all atoms. You can choose the object to pick from the button menus.



6.2.2 Changing the Representation of Bonds

In the Bonds tab you can choose between None, Wire, and Tube for the representation of bonds. The None option means that only the atoms are displayed. To change the bond representation for selected atoms, use the picking tools or the Select button in the Apply current representation to atoms section. Click All to apply the representation to all bonds. You can mix representations within a single molecule by picking.

If you choose the Wire representation, several other options are also available in the Wire Settings section:

- **Style:** You can select Split to display distinctly colored half-bonds, or Blend to display bonds as color gradients. Selecting Blend can improve graphics performance.
- **Represent bond orders with:** You can choose Multiple Lines to display multiple bonds as multiple lines, or Thickness to display multiple bonds as thicker lines, as in earlier versions of Maestro.
- **Show Bond Orders:** You can choose when to display multiple bonds. Automatically means that bond orders are displayed at high zoom but are not displayed at low zoom. The threshold for the display of bond orders can be changed with the Maestro command `repall borderscale=value`. The default for *value* is 30 Å, and represents the width of the portion of the structure visible in the Workspace.
- **Smoothing:** You can smooth the lines that represent bonds. Smoothing Wire bonds especially enhances the viewing of larger structures. However, smoothing can prevent multiple bonds from being displayed correctly in the Thickness representation. If you do not see thick bonds, deselect this option.

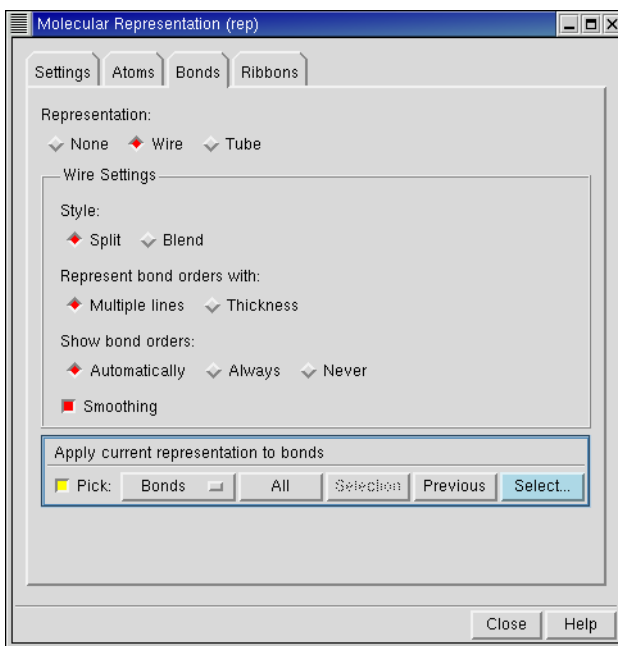


Figure 6.5. The Bonds *tab of the* Molecular Representation *panel.*

To apply the new representation, use the picking tools or the Select button to select bonds, or click All to select all bonds.

You can also apply a molecular representation using the Draw bonds in wire and Draw bonds in tube buttons in the toolbar.



Click either of these buttons to pick the object to apply the representation to. Double-click either button to apply the representation to all atoms. You can choose the object to pick from the button menus.

6.2.3 Changing Representation Attributes

The controls in the Settings tab allow you to customize the appearance of the Wire, CPK, Tube, and Ball & Stick representations. It also provides a quick means of applying a representation to all atoms to view the changes made.

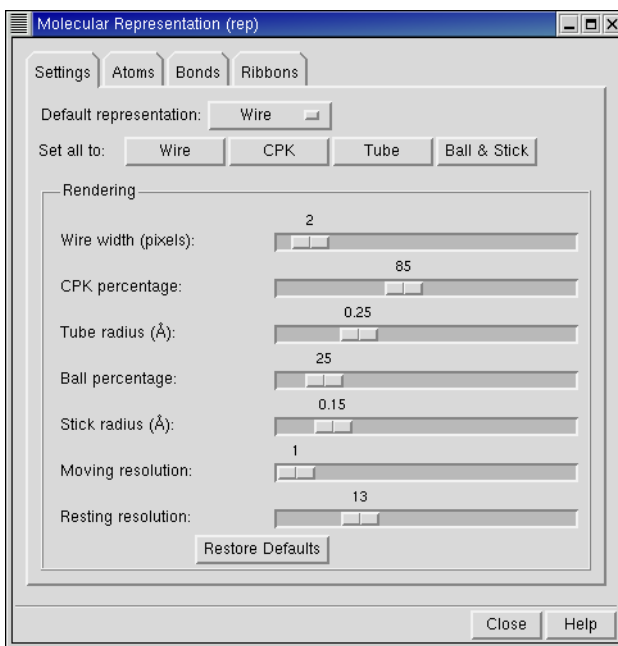


Figure 6.6. The Settings *tab* of the Molecular Representation *panel*.

To change the default representation, so that any subsequently placed atoms use the new representation, select a representation from the **Default Representation** to menu at the top of the tab. Changing the default representation does not affect the representation of any atoms or bonds currently in the Workspace.

You can adjust the parameters of the representations using the sliders in the **Rendering** section and apply the representation using the buttons. If the representation is already selected, changes made using the sliders take effect immediately.

The seven slider controls have the following functions:

- **Wire Width (pixels)**—adjust the width of the lines used to represent bonds in the Wire representation.
- **CPK Percentage**—change the volume of the spheres in CPK mode as a percentage of the van der Waals radius.
- **Tube Radius (Å)**—change the radius of the tubes in Tube mode.
- **Ball Percentage**—change the volume of the spheres in Ball & Stick mode as a percentage of the van der Waals radius.

- Stick Radius (Å)—change the radius of the sticks in Ball & Stick mode.
- Moving Resolution—change the resolution of the ball rendering during translation and rotation. A low resolution can speed up translations and rotations.
- Resting Resolution—change the resolution of the ball rendering when the structure is not being moved.

6.2.4 Rendering Proteins as Ribbons

The controls in the Ribbons tab enable you to display protein backbones as a ribbon diagram. Displaying ribbons can simplify the analysis of the protein's secondary and tertiary structures. You can choose from seven ribbon styles and nine color schemes, and apply the ribbon style by picking or using the Atom Selection dialog box. You can mix representations within a single protein by picking.

You can open the Molecular Representation panel with the Ribbons tab displayed by choosing Ribbons from the Display menu on the main menu bar.

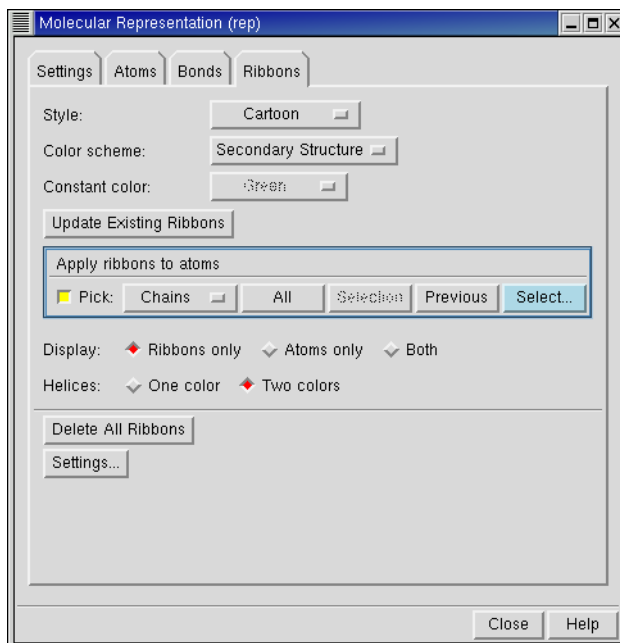


Figure 6.7. The Ribbons *tab of the* Molecular Representation *panel*.

Ribbons are defined by a style and a color scheme. The choices of style are:

- Cartoon—display loops as thin tubes, sheets as arrows and helices as ribbons.
- Ribbon—thin, curved flat surface that follows the backbone.
- Tube—curved tube that follows the backbone.
- Thin Tube—thin curved tube that follows the backbone.
- Curved Line—curved line that follows the backbone.
- CA Trace—straight line that connects alpha carbon atoms.
- CA Trace Tube—straight tube that connects alpha carbon atoms.

With the exception of the CA Atom Color scheme, the color schemes are described in the Color Scheme tab of the Atom Coloring panel. In the CA Atom Color scheme, the ribbon is colored with the current color of the alpha carbon atoms. Thus, if you change the color of certain residues and reapply this color scheme, the section of the ribbon for which you changed the residue color changes color to match. You can also select a two-color representation of helices, using the Helices options in the lower part of the panel. The chosen color scheme is applied to the outside of the helix, and the inside of the helix is gray.

To apply the chosen style and color scheme, you can use the picking tools in the Apply Ribbons To Atoms section, or click Select to use the Atom Selection dialog box, or click All. You can apply the current scheme to all existing ribbons by clicking Update Existing Ribbons. If you want to change the scheme for a part of the protein or for one existing ribbon rather than all ribbons, use the tools in the Apply Ribbons to Atoms section.

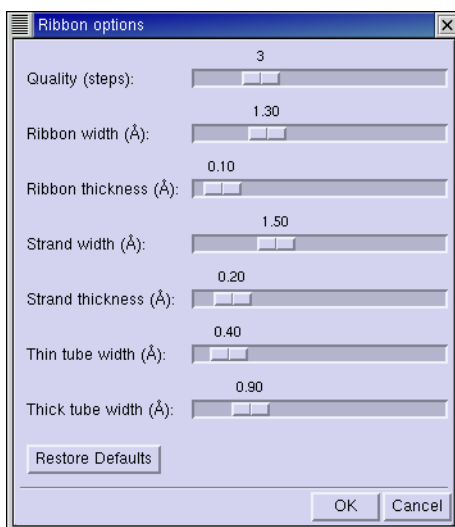


Figure 6.8. The Ribbon options *dialog box*.

The display options in the lower part of the tab allow you to display the ribbons only, the atoms only, or both. The atom display operates only on the already displayed atoms, and only on the atoms for which a ribbon is defined.

You can apply the color schemes and the display options using the Show, hide, or color ribbons toolbar button menu.



If you want to change the quality of the ribbons, or the width or thickness of the ribbons, strands, or tubes, click Settings. The Ribbon options dialog box is displayed. This dialog box has slider controls for quality, ribbon width, ribbon thickness, strand width, strand thickness, thin tube width, and thick tube width, and a button to restore the default values.

6.3 Labeling Atoms

Labels are useful tools for building, identifying, and analyzing molecular structures since they can be used to display any combination of predefined information or user-specified text. For instance, you might want to display atom numbers so that you can retype a specific atom, or to pick atom pairs for superposition, inclusion in a set, or measurement.

Atom labels are persistent. The labels are stored as part of the entry, and are restored when an entry is excluded from then included in the Workspace. However, labels that are dependent on the numbering in the Workspace, such as Molecule Number, may change because they are relative to the Workspace, not the entry.

To open the Atom Labels panel, choose Atoms Labels from the Display menu.

Options for labeling atoms with a variety of atom, molecule, residue, chain, and entry properties are provided in the Composition tab of the Atom Labels panel. To compose a label, select as many properties as you like. The properties are separated by spaces in the label. If you select User Text, you can enter text in the text box to include your own information in the labels.

The color of the labels is the color of the labeled atom, by default. To use a fixed color for the label, deselect Color label by atom in the Color tab, and select a color in the By Palette or By Name subtab. These subtabs are identical to those in the Atom Coloring panel.

Instead of choosing your own labels, you can choose from two fixed labeling schemes that label non-carbon atoms with the element name. In the Hetatom Labels tab, you can select Labels on hetero and hydrogen atoms to label all non-carbon atoms with the element label, or Labels on hetero atoms to label all atoms other than carbon and hydrogen with the element label. When you select either of these schemes, the labels are applied immediately to all rele-

vant atoms. The Labels on hetero atoms scheme is also available from the Build panel toolbar, and is applied as atoms are added to the Workspace.

To select the atoms for labeling, use the tools in the Label atoms section. You can use the picking tools, click All to label all atoms, or click Select and use the Atom Selection dialog box to choose the atoms to label. When you modify the label contents or appearance, you can update all labels with the new definitions by clicking Update Existing Labels.

If you want to remove some or all labels, select Remove in the Mode section. The Label atoms section is relabeled Clear Labels. Use the tools in this section to remove labels. When you select Remove, the label composition is no longer relevant.

You can also label atoms from the toolbar. To apply a single property from a subset of the full range of properties to all atoms, use the Label atoms button menu:

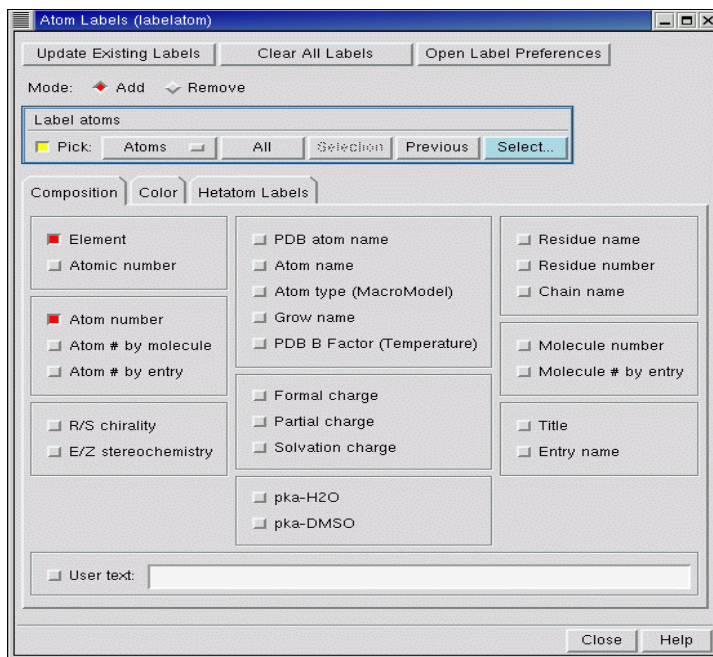


Figure 6.9. The Atom Labels *panel*.

To apply the current labeling scheme by picking, use the Label picked atoms button menu:



Double-click this button to apply the current labeling scheme to all atoms. From the button menu you can choose Composition to open the Atom Labels panel with the Composition tab displayed, you can choose Select to open the Atom Selection dialog box to select the atoms to label, and you can choose Delete Labels to delete all labels.

If you want to customize the appearance, position, and behavior of the label, you can do so in the Label tab of the Preferences panel. See [Section 12.2.4 on page 155](#) for information.

6.4 Displaying and Undisplaying Atoms

The Display/Undisplay Atoms panel allows you to select the atoms that are displayed or are not displayed in the Workspace. For example, you can hide from view all atoms, user-defined atom subsets, or predefined subsets such as protein backbones, sidechains, water molecules, or hydrogen atoms. If you want to display and undisplay project entries, use the In setting in the Project Table panel. See [Section 8.5.1 on page 104](#) for more information about operating on project entries.

To open the Display/Undisplay Atoms panel, choose Display/Undisplay Atoms from the Display menu (see [Figure 6.10](#)).

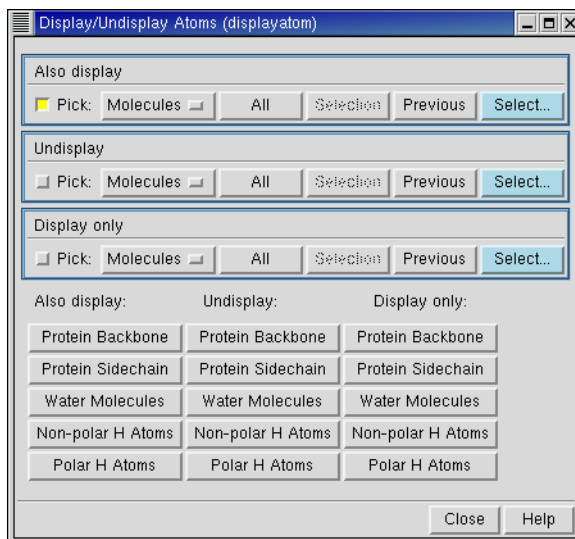


Figure 6.10. The Display/Undisplay Atoms panel.

The Display/Undisplay Atoms panel is divided into two sections. Each section contains control for the following three operations.

- Also Display—add the selected atoms to the atoms that are already displayed.
- Undisplay—do not display the selected atoms.
- Display Only—display only the selected atoms.

In the upper section, you can choose the atoms for the operation using picking tools, including a Pick option menu, an All button and the Atom Selection dialog box, which is opened using the Select button. You can also pick residues to display using the Display only this residue button in the toolbar:



In the lower section, you can apply the operation to the following structural units using the buttons under the operation heading:

- Protein Backbone
- Protein Sidechains
- Water Molecules
- Non-polar H Atoms
- Polar H Atoms

The tools in the lower section are also available from the toolbar as the Display Only, Also Display, and Undisplay button menus:



From these menu you can choose any of the structural units listed above, and you can choose All, All Hydrogens, or Select to open the Atom Selection dialog box and select the atoms to display or undisplay.

Manipulating Structures

Maestro provides a range of tools for manipulating structures in the Workspace. In addition to coordinate transformations such as rotation and translation, there are tools for centering a structure on a particular atom, zooming in on a structure, and arranging structures in a rectangular grid, which is called *tiling*.

Transformations can be carried out globally or locally. Global transformations move all of the atoms in the Workspace together as a unit. They do not change atom coordinates, only the orientation and location relative to the viewer. Local transformations move a user-defined subset of atoms relative to the rest of the atoms in the Workspace, and therefore must change the atom coordinates. Tiling also changes the atom coordinates.

The coordinate system in Maestro is oriented so that the positive z-axis points out of the screen, the x-axis is the horizontal axis, and the y-axis is the vertical axis.

When a project is closed, the current rotation, translation, and zoom settings are saved with the project and restored when the project is opened.

7.1 Global Transformations

Global transformations are defined as rotation and translation of the entire structure or set of structures in the Workspace. The axes for rotation and translation are assigned to combinations of mouse buttons with the SHIFT and CTRL keys. These assignments are given in [Table 7.1](#).

Global transformations do not actually change the coordinates of the structures in the Workspace. Rather, they change the view of the structures. You can save and restore a particular view of the structures from the toolbar—see [Section 2.10.2 on page 21](#) for details. However, if you want to change the coordinates of the structures to reflect the new view, you can do so by choosing **Update Coordinates** from the **Build** menu.

To perform a global transformation, position the cursor in the Workspace, and hold down the relevant key and mouse buttons while you move the mouse. You should take care that you do not click instead of drag, because clicking might have a particular action associated with it in the context in which you are working, and change the structure.

To select a new center that defines the axes for rotation, you can *spot-center* the molecule on an atom by right-clicking on the atom.

Table 7.1. Mapping of Workspace operations to mouse actions.

Operation	Action
Rotate about the x - and y -axes	Drag with middle mouse button
Rotate about the x -axis only	Drag vertically with SHIFT and middle mouse button
Rotate about the y -axis only	Drag horizontally with SHIFT and middle mouse button
Rotate about the z -axis	Drag horizontally with CTRL and middle mouse button
Spot-center on an atom	Right-click
Translate in the x - y plane	Drag with right mouse button
Translate along the y -axis	Drag vertically with SHIFT and right mouse button
Translate along the x -axis	Drag horizontally with SHIFT and right mouse button
Translate about the z -axis	Drag horizontally with CTRL and right mouse button
Zoom	Drag horizontally with middle and right mouse buttons or with SHIFT+CTRL and middle mouse button

The other operation that does not change the coordinates of the structure is zooming, which is also performed with the mouse (see Table 7.1). You can also fit the displayed size of the structure to the Workspace by clicking the Fit to screen button on the toolbar.



To restore the view to the default view, in which the view axes are aligned with the coordinate axes of the structure, click the Reset workspace button in the toolbar



7.2 Local Transformations

A local transformation in Maestro is the movement of a user-defined atom subset¹ with respect to the other atoms in the Workspace, such as the movement of one molecule with respect to another. If you want to perform a translation, you need only select the atoms for transformation. If you want to perform a rotation, a rotation center is required. Maestro provides a default, but if this is not satisfactory, you must choose a rotation center.

1. The subset can in fact include all atoms in the Workspace, so that the rotation or translation actually changes the coordinates of the entire contents of the workspace.

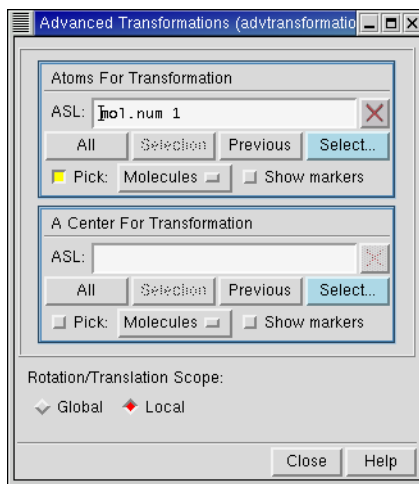


Figure 7.1. The Advanced Transformations panel.

The controls for local transformations are available from the Advanced Transformations panel, which you open by choosing Global/Local from the Local transformation button menu in the toolbar.



The picking tools for defining atom sets are also available directly from this button menu. The selections you make in the Advanced Transformations panel persist, even when you alter the structures in the Workspace. To clear them, click Clear.

Local transformations are performed with the same combination of keys and mouse movements as global transformations. When you begin a local transformation, Maestro automatically begins an undo block. As a result, selecting Undo *action* from the Edit menu reverses all operations that occurred between the time that the mouse button was pressed, and the time that it was released.

You can switch between global and local transformations by selecting Global or Local in the Rotation/Translation Scope section of the Advanced Transformations panel, by using the shortcut keys CTRL+G and CTRL+L in the Workspace, or by clicking on the Global/Local button on the toolbar. When the button is indented, local transformation is in effect.

While you are in local transformation mode, you can switch to global transformation mode temporarily by holding down the spacebar while you perform a transformation.

7.2.1 Selecting Atoms for Transformation

To select the atoms that you want to move, use the tools in the Atoms For Transformation section. You can pick atoms in the Workspace, click **Select** to select atoms using the Atom Selection dialog box, or enter an ASL expression in the ASL text box. You can also use the Local Transformation button menu to select an object and pick it in the Workspace.

When you use the picking tools for the Atoms for Transformation in the Advanced Transformations panel, the last ASL expression is replaced each time you make a new pick. If you want to pick more than a single structural unit, click **Select** and use the Atom Selection dialog box.

If **Show Markers** is selected, the atoms are marked with yellow circles connected by yellow lines. Once you have started to select the atoms in the Atoms For Transformation section, the **Local** option in the Rotation/Translation Scope section is selected.

In order to see the atoms you want to pick for a local transformation, it may be necessary to rotate the structure. If this is the case, select **Global** in the Rotation/Translation Scope section, rotate the structure, then reselect **Local** to define atoms for local transformation.

If you want to translate atoms, you need not define a center for transformation. Once you have selected the atoms for transformation, you can perform the desired translation.

7.2.2 Selecting a Rotation Center

To perform local rotations, a center for rotation must be defined. By default, Maestro selects the centroid of the atoms selected for transformation as the rotation center. If **Show Markers** is selected, the rotation center is displayed in the Workspace as a purple ball. If the assigned centroid is not satisfactory, select a new center using the tools in the A Center For Transformation section. Displaying atom number labels may be helpful for determining which atom you would like to specify as the rotation center.

7.3 Tiling Multiple Entries

To arrange multiple included entries into a grid pattern in the Workspace, click the **Tile Entries** button in the toolbar.



This operation changes the coordinates of the structures in the entries. The entries are arranged according to the order they appear in the Project Table. To reverse this operation, choose **Undo action** from the Edit menu or click the **Undo/Redo** button in the toolbar.

Projects

When you use Maestro, you are always working within a *project*. Projects are made up of *entries*, which can include multiple structures as well as properties of those structures. Project entries and their properties are displayed in the Project Table. You can perform a variety of operations on entries through the Project Table. By default, the Project Table panel is opened automatically when you open a project. You can change this behavior in the Preferences panel. To open the Project Table panel, choose Show Table from the Project menu, or click the Open/Close Project Table button in the toolbar.



The project data is stored in a directory whose basename has the extension `.prj`. This directory is created in Maestro's current working directory. If you have not explicitly opened an existing project or created a new project, a scratch project is created in the directory specified by the environment variable `TEMP_PROJECT`. The default location is:

UNIX: `$HOME/.schrodinger/tmp`

Windows: `%USERPROFILE%\schrodinger\tmp`

This chapter discusses in some detail how to use the Maestro Project Facility.

8.1 Project Operations

The Project menu of the main window enables you to perform operations on the project level. To perform operations on entries in the current project, use the menus in the Project Table panel. This section describes the operations that you can perform on projects.

8.1.1 Creating and Opening Projects

When you start Maestro without specifying a project to open, a scratch project is created. You can work with a scratch project, and perform nearly all the tasks that you can perform with a named project. If you want to start a new project and save its data, you can create a new named project by choosing New from the Project menu. The main difference between a new project and a scratch project is that the new project has a name and is automatically saved. When you choose New, a project selector panel is displayed. This panel enables you to choose the location of the new project, create a directory, and type in the project name.

To open an existing project, choose **Open** from the **Project** menu, or click the **Open a project** toolbar button.



A project selector panel is displayed that enables you to select a project to open. The project selector panels are described in more detail in [Section 8.1.4 on page 89](#). You can only have one project open at any time.

8.1.2 Saving and Closing Projects

The data in a project is normally saved automatically when you make changes to the project. Changes to entries that are made in the **Project Table** are saved automatically. The saving of changes to entries that are displayed in the **Workspace** depends on the preferences you set in the **Preferences** panel. By default these changes are saved automatically.

Because the data in a named project is automatically saved you do not need to explicitly save the project. When you have finished working in a project you need only close the project, by choosing **Close** from the **Project** menu. If you want to work on a different project, you do not need to close the old project before opening the new project. When you close a project without opening a new project, a scratch project is created.

If you are working in a scratch project and decide you want to save it as a named project, or if you want to save a named project with a different name, choose **Save As** from the **Project** menu or click the **Save As** toolbar button.



A project selector panel is displayed that enables you to enter the name and choose the location of the project, and if necessary, create a directory. The panel is described in [Section 8.1.4 on page 89](#). If you choose **Close** when you are working in a scratch project, you are prompted to save the scratch project.

When a project is saved, all view transformations (translation, rotation, and zoom) are saved with it and restored when the project is opened. In addition, all user-defined sets of atoms are saved with the project and included when the project is copied or merged with another project.

8.1.3 Deleting Projects

If you do not want to keep the current project you can remove it by choosing Delete from the Project menu. You will be prompted to confirm the deletion of the project. When you delete a project, any jobs associated with that project are killed and purged from the job database.

You can also delete the project from a UNIX shell by entering the command

```
rm -rf project_name
```

When you delete a project using `rm`, job files that were created using the project are not deleted unless they are located within the project directory.

8.1.4 The Project Selector Panels

Many of the project operations—creating, opening, saving—involve selection of a project location on disk. Each of these operations has a project selector panel. The project selector panels operate like file selectors, but project selector filters are configured so that Maestro project directories (those with `.prj` extensions) are the only visible items. Most of the features of these panels are common to all operations. The panel features are described below.

Project Availability Indicators: Only one instance of a given project can be open at a time. To indicate whether a project is in use, an asterisk (*) is placed in front of each available project and an exclamation point (!) in front of each locked project. A project is locked if you or someone else already has the project open.

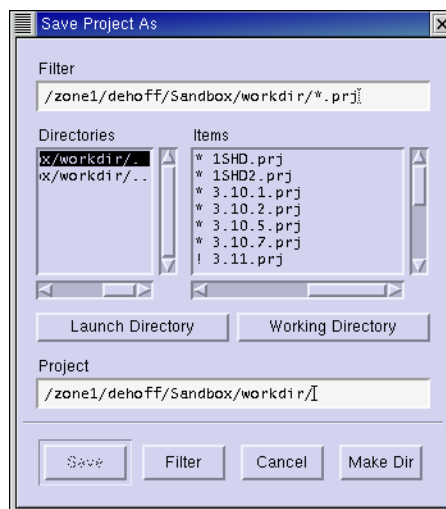


Figure 8.1. The Save Project as *project selector*.

Directory selection buttons: In addition to the usual directory and file lists, there are two buttons for selecting directories: Launch Directory and Working Directory. The first selects the directory from which Maestro was launched, and the second selects Maestro's current working directory.

Options: When you merge projects you can choose to include any user-created files by selecting Include User Created Files.

Project Text Box: The Project text box is automatically populated with the name of the current directory. When you navigate to locate a directory, the text box is updated to reflect the current selection. If you are creating a new project, saving a scratch project or copying a project, enter the project name at the end of the text. You do not need to add the `.prj` extension, which Maestro supplies if it is absent.

Making New Directories: Some of the project selectors have been equipped with the ability to make new directories. The Make Dir button in these project selectors opens the Make Directory directory selector. To create a directory, navigate to the desired location and type the new directory name at the end of the path information in the Project text box. You can only create one new directory level at a time. To create multiple directory levels you must use this procedure for each level you want to create.

8.1.5 The Project Summary Panel

Projects can be large, diverse collections of data. To help you keep track of the project, Maestro provides a Project Summary panel. This panel displays information about the current project, including annotation of the project, information on jobs run from the project, and information on runs saved from wizard-based panels. To open this panel, choose Summary from the Project menu. The panel has three sections, an Annotation section, a Project section, and a Jobs section.

Note: The Project Summary panel is not available under Windows.

Annotation section: This section displays the project annotation. Because projects are frequently modified, it may be helpful to annotate your project. The project creation date is automatically added to the project annotation.

You can edit the text in the annotation section, including the automatically generated text. The standard mouse actions (highlight to copy, middle-click to paste) can be used to copy and paste text. The main key strokes for editing text are:

- CTRL+W to cut the selection
- ALT+W to copy the selection
- CTRL+Y to paste the selection

These and other conventions are the same as those used by Emacs. For more information, visit the web site <http://www.tcl.tk/man/tcl8.4/TkCmd/text.htm#M152>¹.

Project section: This section displays statistics on the project, in two subsections, Entries and Runs. Clicking in the Entries subsection opens the Project Table panel. Clicking in the area for a product in the Runs subsection opens the wizard-based panel for that product.

Jobs section: This section displays the jobs in the job database for the current project. Clicking Refresh updates the display for all jobs, and is equivalent to clicking Monitor in the Monitor panel. Clicking Delete deletes the selected jobs from the database. You can select jobs with the usual click, shift-click and control-click key strokes.

8.2 Adding Entries to a Project

Projects are made up of entries. Entries can contain multiple structures and their properties. Each entry has both a name, by which it is identified internally and which can be modified by Maestro, and a title, which can be set or changed by the user but is otherwise not modified by Maestro. The title is displayed in the Project Table by default, but the entry name is not.

You can add entries to a project by importing structures from a file, by creating entries from structures you have built in the Workspace, or by merging entries from a different project. In addition, you can create new entries by splitting or duplicating an existing entry, or by incorporating the results of a job into the Project Table.

8.2.1 Importing Structures From a File

When Maestro imports structures from a file, the structures are copied, translated into entries, and incorporated into the current project. Structures can be imported from Maestro, Macro-Model, PDB, Mol2, SD, and other files using the Import panel. Multiple molecular structures can be imported simultaneously from multiple files, and you can choose which structures, if any, are included in the Workspace upon import. Entry names and titles are automatically assigned if they do not exist, and are read if they do exist. See [Chapter 3](#) for more information on importing structures.

8.2.2 Creating Entries From Workspace Structures

When you build structures in the Workspace from scratch, these structures constitute a scratch entry, and are not saved in the project. To save the scratch entry in the current project, choose Create Entry From Workspace from the Project menu, and type an entry name in the dialog box

1. Please see the [notice](#) regarding third party programs and third party Web sites on the copyright page at the front of this manual.

that is displayed. If you choose an existing entry name, Maestro appends `-wsen` to the entry name to ensure that the name is unique, where n is a number starting from 1. All structures in the Workspace are placed in the entry that you create.

You can also create an entry using the Create entry from Workspace button in the toolbar.



8.2.3 Incorporating Entries From Job Output

When you use selected or included entries from a project as input for most kinds of jobs, you have the option to incorporate the job results as entries in the Project Table. You can choose to append the new entries, to replace the input entries with the new entries, or you can choose to ignore the output and not incorporate them into the table at all. Incorporated entries have `-incn` appended to the entry name, where n is a number starting from 1.

If the job results contain multiple structures, an entry group is created for them on incorporation if one does not exist. See [Section 8.3.6 on page 99](#) for more information on entry groups.

For more information on incorporating job output into the Project Table, see [Section 4.8](#) of the *Job Control Guide*.

8.2.4 Merging Entries From Another Project

You can merge the entries from a specified project into the current project. There is no facility for merging a subset of entries. When the merge is performed, the entries are incorporated at the end of the Project Table list in the current project.

To copy the entries from a source project into the currently open project, choose **Merge** from the Project menu. In the Merge From Project project selector, navigate to the desired source project. If you want to include user-created files associated with the merge source project, select **Include user created files**.

If none of the incoming entries has a name identical to an entry in the current project, the entry names remain intact. If there are duplicate entry names in the projects to be merged, the Entry Name Resolution panel is displayed, and you must choose a resolution scheme in the Retain Entry Names In section—Current, Incoming, or User Specified. If you select Current, the names in the current project are kept, and the conflicting incoming entries are renamed. Likewise, if you select Incoming, the incoming entries retain their names, and the conflicting entries in the current project are renamed. If you want to choose which entries to rename, select User Specified, then click the Specify button. The User Specified Conflict Resolution panel is displayed, and you can choose to rename either the current or the incoming entry.

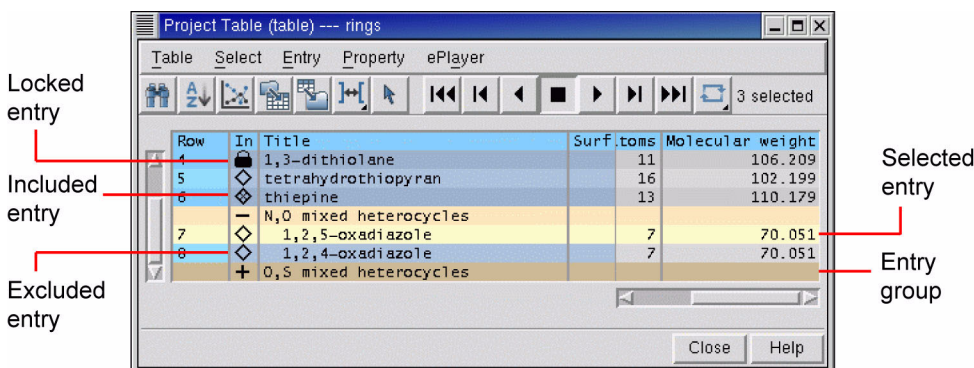


Figure 8.2. The Project Table *panel*.

The resolution scheme that Maestro offers may not produce especially meaningful names. For example, if in both the current and the merge source project, you have entries named *entry*, and you already have other entries in the current project named *entry-1*, *entry-2*, and *entry-3*, Maestro renames the incoming entry as *entry-4*. Since renaming entries in this fashion may destroy any location cues that you were using, you may want to go into each project separately and resolve any naming conflicts before attempting to merge projects. However, the titles of the entries are preserved, and these can be used to distinguish entries when the entry name is changed by Maestro.

8.3 The Project Table Panel

The Project Table panel provides access to a range of tools for managing the data in the project and operating on that data. The panel contains a menu bar, a toolbar, and the Project Table itself. This section provides a brief description of the basic table operations, the tools, and the use of the mouse and shortcut keys for performing various operations.

Each project entry is represented by a row in the Project Table. Each row contains the row number, an icon in the *In* column indicating whether the entry is displayed in the Workspace, the title, a button to open the Surfaces panel if there are surfaces associated with the entry, the entry name, and any properties associated with the entry. The first four columns of the project table are always visible: the row number, the *In* column, the entry name or the title (you can exchange these two from the *Table* menu), and the *Surf* column. The remainder of the columns appear in the scrolling section of the table.

Project entries have two states that represent how they are being used: the *selection* state and the *inclusion* state. Selected entries are highlighted in the Project Table and are available for operations such as copying, deleting, or use as input to a computational program. Selection of

entries is discussed in the next section. Included entries are displayed in the Workspace, and are available for manipulation, modification, and other actions that can be performed from the Workspace. The inclusion state is indicated in the In column. Excluded entries—entries that are not displayed in the Workspace—are marked with an empty diamond. Included entries are marked with a cross in the diamond. Fixed entries—entries whose display is not affected by the inclusion or exclusion of other entries—are marked with a padlock icon. Fixed entries are not frozen: you can perform operations on fixed entries in the Workspace, such as rotation and translation, or adding or deleting atoms. To undisplay a fixed entry, you must exclude it from the Workspace. The selection and inclusion states are stored in the project, so they are restored when you open the project.

8.3.1 The Project Table Toolbar

The toolbar contains two groups of buttons and a status display. The first set of buttons opens various panels that allow you to perform functions on the entries in the project table. The second set of buttons controls the ePlayer, which “plays through” the selected structures: each structure is displayed in the Workspace in sequence, at a given time interval. A description of the toolbar buttons is provided below. The use of these buttons is also discussed in later sections of this chapter.

**Find**

Open the Find panel for locating alphanumeric text in any column of the Project Table, except for the row number.

**Sort**

Open the Sort Project Table panel for sorting entries by up to three properties.

**Plot**

Open the Plot XY panel for plotting entry properties.

**Import structures**

Open the Import panel for importing structures into the project.

**Export structures**

Open the Export panel for exporting structures to a file.

**Columns menu**

Choose an option for adjusting the column widths.

**Entry selection**

Open the Entry Selection dialog box for selecting entries based on criteria for entry properties.



Go to start
Display the first selected structure.



Previous
Display the previous structure in the list of selected structures.



Play backward
Display the selected structures in sequence, moving toward the first.



Stop
Stop the ePlayer.



Play forward
Display the selected structures in sequence, moving toward the last.



Next
Display the next structure in the list of selected structures.



Go to end
Display the last selected structure.



ePlayer loop menu
Choose an option for repeating the display of the structures. **Single direction** displays structures in a single direction, then repeats. **Oscillate** reverses direction each time the beginning or end of the list is reached.

8.3.2 The Project Table Menus

The menus on the Project Table menu bar are described briefly below. The use of the menu items is discussed later in this chapter.

- **Table**—find text, sort entries, plot properties, import and export structures, and configure the Project Table.
- **Select**—select all entries, none, invert your selection, make a random selection, or select classes of entries using the Entry Selection dialog box and the Filter panel.
- **Entry**—include or exclude entries from the Workspace, perform various operations on the selected entries, and group or ungroup entries.
- **Property**—display and manipulate entry properties in the Project Table.
- **ePlayer**—view entries in succession, stop, reverse and set the ePlayer options.

8.3.3 Mouse Functions in the Project Table

The standard use of shift-click and control-click to select objects is supported in the Project Table. Middle-click is supported as an alternative to control-click. This behavior applies to selection of entries and inclusion of entries in the Workspace. Dragging to resize columns and to move rows is also supported. You can drag a set of non-contiguous entries to reposition them in the Project Table—see [Section 8.5.8 on page 106](#) for more information. You can edit cells by clicking in the cell and changing the text, then pressing ENTER or clicking elsewhere in the Project Table.

A summary of project-based mouse functions is provided in [Table 8.1](#).

Table 8.1. Mouse operations in the Project Table.

Task	Mouse Operation
Change a Boolean property value.	Click repeatedly in a cell to cycle through the possible values (On, Off, Clear).
Display the Entry menu for an entry.	Right-click anywhere in the entry. If the entry is not selected, it becomes the selected entry. The chosen action is applied to all selected entries.
Display a version of the Property menu for a property.	Right-click in the column header.
Edit the text or the value in a table cell.	Click in the cell and edit the text or value. If the entry containing the cell is not selected, it is selected. The cursor is placed at the end of the data by default.
Include an entry, exclude all others.	Click the In column of the entry.
Move selected entries.	Drag the entries.
Paste text into a table cell.	Middle-click.
Resize columns.	Drag the boundary with the middle mouse button.
Select an entry, deselect all others.	For an unselected entry, click anywhere in the row except the In column; for a selected entry, click the row number.
Select or include multiple entries.	Click the first entry then shift-click the last entry.
Toggle the entry selection or inclusion state.	Control-click or middle-click the entry or the In column.

8.3.4 Project Table Shortcut Menus

In addition to the menus on the menu bar, the Project Table panel has shortcut menus, which are displayed when you right-click and hold in the table. There are three shortcut menus: one for properties, one for entries, and one for entry groups.

- **Property shortcut menu:** This menu duplicates some items on the Property menu, and has some additional items. It is displayed when you right-click and hold over a column header, and the actions are applied to the property in the column. The items on this menu are described in [Table 8.2](#).
- **Entry shortcut menu:** This menu duplicates the Entry menu, and in addition has an Export Structures item, which opens the Export panel for export of the selected entries. It is displayed when you right-click and hold over any entry in the table. If the entry is part of the selection, the menu is displayed immediately. If the entry is not part of the selection, the entry is selected and all others deselected before the menu is displayed.
- **Entry group shortcut menu:** This menu has no counterpart on the menu bar. It is displayed when you right-click and hold over an entry group row. The items on this menu are described in [Table 8.3](#). There are four groups of items: for entry selection, for entry inclusion, for group management, and for actions on group members.

Table 8.2. Property shortcut menu items.

Menu item	Action
Hide	Hide (undisplay) the property.
Sort	Sort all rows in the Project Table by the values of this property.
Move	Move display to first or last column. The submenu contains First or Last.
Copy	Copy values from one property to another. Opens the Copy Property panel, in which you can specify the target and destination property. Properties can be copied for the selected entries or for all entries.
Rename	Rename properties. Opens the Rename Property panel, in which you can rename any property except those whose names are reserved.
Set Decimal Places	Set the decimal precision of real entries. Opens the Set Decimal Places dialog box, in which you can select properties and set the precision.
Clear Values	Clear the values for the selected property.
Delete	Delete the selected property.
Align	Align the data in the column. The submenu offers the choice of Left, Center, or Right.

Table 8.3. Entry group shortcut menu items.

Menu item	Action
Select All	Select all entries in the group.
Deselect All	Deselect all entries in the group.
Include All	Include in the Workspace the entries in the group, along with the entries that are already included.
Include Only All	Include in the Workspace only the entries in the group, and exclude all others. Does not affect fixed entries.
Include First	Include in the Workspace the first entry in the group, along with the entries that are already included.
Exclude All	Exclude the entries in the group from the Workspace.
Collapse or Expand	Show or hide all the group members. If the group members are visible, this item is Collapse, and the group members are hidden by this action (the group is collapsed). If the group members are not visible, this item is Expand, and the group members are shown by this action (the group is expanded).
Expand Only	Expand only this group, and collapse all others. Only the members of the selected group are shown.
Duplicate	Duplicate the group. The duplicate group is placed below the source group. Duplication makes copies of all the entries in the source group and places the copies in the new group.
Ungroup	Move the entries out of the group, and delete the group. The entries are placed at the end of the ungrouped entries.
Delete	Delete the group, and all the entries in it. A confirmation dialog is displayed when you choose this item.
Edit Title	Edit the title of the group. Only displayed when the pointer is in the Title column.
ePlay	Play through the entries in the group with the ePlayer.
Pose Viewer	Play through the entries in the group with the ePlayer, adding hydrogen bonds and counts.

8.3.5 Project Table Shortcut Keys

Some frequently used operations have been assigned shortcut key combinations. The shortcuts, their functions, and their menu equivalents are listed in [Table 8.4](#).

Note: If you have Num Lock, Caps Lock, or Scroll Lock enabled, the shortcut keys do not work.

Table 8.4. Shortcut keys in the project facility.

Keys	Action	Menu Equivalent
CTRL+A	Select all entries	Select > All
CTRL+F	Fix entry in Workspace	Entry > Fix
CTRL+I	Show import panel	Table > Import Structures
CTRL+N	Include only selected entries	Entry > Include Only
CTRL+U	Deselect all entries	Select > None
CTRL+X	Exclude selected entries	Entry > Exclude
CTRL+Z	Undo/Redo last command	Edit > Undo/Redo in main window

8.3.6 Configuring the Project Table

The Project Table panel provides tools for controlling what is displayed in the table. Some of these tools are also discussed later in this chapter.

Grouping and ungrouping entries. You can group the selected entries by choosing Group from the Entry menu. A new entry group is created at the end of the table, and the entries are moved to the new location. You can display or hide the members of the group by clicking the In column for the group.

To ungroup an entry group, select the group and choose Ungroup from the Entry menu. The group entries are placed at the end of the ungrouped entries, immediately above the first group.

You can add the selected entries to an existing group by choosing Move to Group from the Entry menu, and selecting the group in the dialog box that is displayed. You can move entries out of a group by selecting them and choosing Move to Row from the Entry menu.

Displaying a subset of properties. You can display a subset of properties by choosing Show from the Property menu, then using the Show Property dialog box to select the properties to show or hide. You can choose individual properties or a family of properties. Property families are the properties generated by an application, such as QikProp.

Resizing columns. You can resize columns by dragging the border of the column using the middle mouse button. You can also set the width of the columns to the width of the widest data item, the width the widest of the data and the heading, or to a fixed width, by choosing from the Columns button menu on the toolbar or the Columns submenu of the Table menu.

Aligning column data. You can align the data in a column on the left boundary, the right boundary, or the center of the column, by right-clicking and holding in the column header and choosing the alignment mode from the Align submenu of the Property shortcut menu.

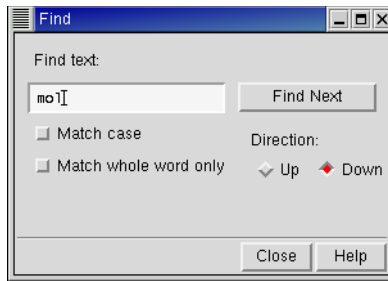


Figure 8.3. The Find panel.

Displaying the Entry Name. You can display the Entry Name in the scrolling area of the table by choosing Show Entry Name from the Property menu. This option can also be set as a preference in the Project tab of the Preferences panel.

Setting the precision of real values. You can set the precision of real-valued properties in the Set Precision dialog box, which you open from the Property menu.

8.3.7 Finding Text in the Project Table

If you want to locate text in the Project Table, whether it is a string or a numerical value, choose Find from the Table menu. The Find panel is displayed. Enter the text you want to find, select Down to search down from the first entry or Up to search up from the last entry, and click Find Next. The Project Table scrolls to the found entry, but the entry is not selected.

To increment the search, click Find Next again. The Project Table scrolls to the next instance of the search string. If the string is not found, a beep sounds when the end (or beginning) of the table is reached, and the Project Table scrolls to the original found entry.

Note: If you click on the Project Table panel, then click Find Next, the search starts over, searching from the beginning of the table if Down is selected, or from the end of the table if Up is selected.

8.4 Selecting Project Entries

Selection means choosing entry rows in the Project Table, which become the target for an action. The rows are highlighted when they are selected. Selection differs from “inclusion,” which is the display of the structures belonging to entries in the Workspace. Inclusion is described on [page 93](#). Selection is no different in principle from the selection of items in a list in any graphical interface. In addition to the usual mouse operations, summarized in [Table 8.1](#), Maestro provides additional capabilities for selecting entries based on their properties. Selection can include one or more entry groups.

Once entries are selected, they can be operated on using the mouse or the tools in the Entry menu. Selected entries can be renamed, duplicated, combined, split, deleted, exported, and included in the Workspace, or a movie of the structures associated with entries can be made with the ePlayer. Some operations on selected entries are mapped to shortcut keys, which are described in [Table 8.4 on page 99](#).

Selection states are persistent across Maestro sessions: entries that are selected when a project is closed will still be selected the next time that project is opened. When new entries are imported or incorporated into the Project Table, they are automatically selected.

8.4.1 The Select Menu

The Select menu provides tools for selecting entries. The first three items, All, None, and Invert, select all entries, no entries (deselect all entries), and invert the selection (select the unselected entries, and deselect the selected entries). Only, Add, and Deselect open the Entry Selection dialog box, so that you can select entries by a combination of properties. Random opens the Random Selection dialog box, in which you can randomly select a percentage of the selected entries or of all entries. The Filter submenu allows you to choose a named *filter* that is used to select entries, and to edit named filters. Filters are defined in terms of Entry Selection Language (ESL) expressions, and do not constitute a fixed set of entries. Whenever the filter is used, the ESL expression is applied to the current contents of the Project Table. For more information on ESL, see the [Maestro Command Reference Manual](#).

8.4.2 Selecting Entries By Property

The Entry Selection dialog box ([Figure 8.4](#)) enables you to select entries in the Project Table based on a combination of criteria for property values. Once you have defined a set of criteria for selecting entries, you can store this set as a named filter, and use the filter directly from the Select menu. You can open the Entry Selection dialog box from the Select menu or using the Select Only toolbar button.



The Entry Selection dialog box is similar in design to the Atom Selection dialog box. It allows you to select entries using GUI tools rather than entering an ESL expression.

At the top of the entry selection dialog box is a title bar, which defines the action to be taken when you click OK. In the center are two tabs, Properties and Filters, that you use to define property-based criteria for entry selection or choose named filters to combine with the entry selection that is being created (the *current* entry selection). On the right are buttons for combining the criterion or named filter with the current entry selection. These tools become

available when you have selected a property and a condition for that property, or a named filter. In the lower portion of the panel is an ESL text box that displays the current entry selection and a set of buttons for editing the current entry selection.

To define a property criterion, choose a property from the list in the Properties tab, then select a condition from the right side of the tab. The available conditions depend on the property type. For some conditions, you must enter a value or a range of values. The available conditions are:

- The property is defined.
- For Boolean properties, the value is true or false.
- For numerical (real or integer) properties, the value lies in a given range, for example, ≥ 5.0 or $1 - 7$. You can use any of the six operators $=$, \neq , $<$, $>$, \leq , and \geq to define the range of property values.
- For textual properties, the text matches a given string. “Match” means that the string given is contained in the property value. The string can include the wildcard characters $*$, which matches any number of characters, $?$, which matches any single character, and $\#$, which matches any single digit.

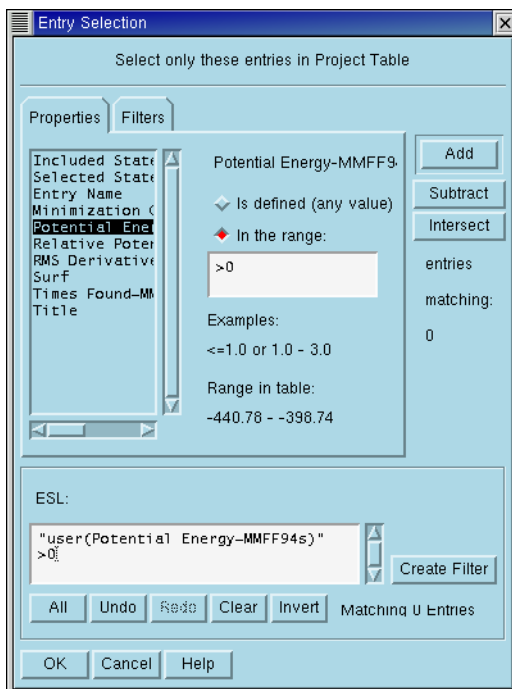


Figure 8.4. The Entry Selection dialog box.

In addition, the inclusion state and the selection state of the entry can be used as criteria.

To combine this criterion with the current entry selection, click one of the buttons to the right of the panel.

- **Add** adds entries that meet the new criterion to the current entry selection, using a logical OR operation on the ESL expressions. This is the default action, which is performed if you press ENTER.
- **Subtract** excludes entries that meet the new criterion from the current entry selection, using a logical AND NOT operation on the ESL expressions. If this is the first operation, the default current entry selection is `all`.
- **Intersect** limits the selection of entries to those in the current entry selection that also meet the new criterion, using a logical AND operation on the ESL expressions.

In addition to these three buttons, the **Invert** button at the bottom of the panel inverts the entire current entry selection, using a logical NOT operation. All entries that were selected are unselected, and all entries that were not selected are selected.

You can undo multiple operations with repeated clicks on the **Undo** button. The **Redo** button reverses undo operations, and can also be used multiple times. You can also type ESL expressions directly into the ESL text box. The hand-edited expressions are considered part of the last operation performed with the mouse and can be undone and redone as part of that operation, but cannot be undone or redone separately.

Once you have defined the criteria for selecting entries, the current entry selection can be stored as a named filter. To create a filter, click **Create Filter**. The **New Filter** dialog box opens, in which you can provide a name for the filter. The newly created filter appears in the **Filters** tab of the **Entry Selection** dialog box and on the **Filter** submenu of the **Select** menu.

8.4.3 Selecting Entries Using the Plot Panel

The **Plot** panel provides a means of graphically selecting entries and including them in the **Workspace**. By plotting a property value against the entry name, for example, you can view the properties and select the entries on the plot that you are interested in. See [Section 10.2.4 on page 131](#) for more information.

8.4.4 Selection Examples

There are generally multiple ways to perform each selection operation: from an option on the **Select** menu, using the **Entry Selection** dialog box, using the mouse, and using ESL. A few examples of the use of the **Entry Selection** dialog box are provided here.

To select a range of entries with the same base name:

1. Choose Only from the Select menu.
2. Select Entry name from the Properties list.
3. Select Matches in the conditions section and enter

*entry_name**

in the text box, where *entry_name* is the base name of the entries you want to select, and * is a wildcard character that represents any string.

4. Click Add, then click OK.

To add entries from the Workspace to the current selection:

1. Choose Add from the Select menu.
2. Select Included State from the Properties list.
3. Select Included (or fixed) in the conditions section.
4. Click Add, then click OK.

8.5 Operating on Selected Entries

The Entry menu provides tools to perform various operations on the selected table entries. These operations are discussed in the following sections.

8.5.1 Including, Excluding, and Fixing Entries

In addition to using mouse operations, you can include and exclude entries in the Workspace by selecting the entries and choosing Include, Include Only, or Exclude from the Entry menu. The last two of these menu items are mapped to the shortcut keys CTRL+N and CTRL+X. You can also include and exclude entries by property value from a plot in the Plot XY panel. See [Section 10.2.5 on page 131](#) for more information.

If you want to keep some entries in the Workspace while including and excluding other entries, such as a receptor and a set of ligands, you can select the entries you want to keep in the Workspace and choose Fix from the Entry menu, or press CTRL+F. These entries remain in the Workspace until you explicitly exclude them.

8.5.2 Limiting the Number of Entries Displayed

In large projects you may want to limit the number of entries displayed in the Project Table. To do this, select the entries, then choose **Show Only** from the Entry menu. Likewise you can hide the selected entries by choosing **Hide** from the Entry menu. To redisplay all entries, choose **Show All** from the Entry menu. The display of entries does not affect either the subsequent selection or inclusion of entries—though of course you can only use the mouse to select or include the entries that are displayed in the Project Table.

8.5.3 Renaming Entries

When importing entries from files, incorporating entries from jobs, merging entries with identical names, duplicating entries, or disassociating entries, Maestro automatically assigns new entry names. You do not have to accept the automatically assigned names; you may rename the entries using one of the methods described below.

To change the name of a single entry:

1. Double-click the entry name if the entry is not selected, or single-click the name if the entry is selected.
2. Type in the new name.
3. Press ENTER.

To rename multiple entries with the same base name:

1. Select the entries to be renamed.
2. Click the entry name for one of the entries.
3. Type in the new name.
4. Press CTRL+ENTER.

Maestro automatically adds incremental numbers to the name you assign, to make the entry names unique.

8.5.4 Duplicating Entries

To duplicate entries, select the entries and choose **Duplicate** from the Entry menu. Maestro places each duplicated entry after its source entry. The new entries retain the same base name as the originals, but have `-dup n` appended. The number n indexes the duplicates. The title is preserved when the copy is made.

8.5.5 Merging Multiple Entries

Maestro allows you to merge two or more selected entries, while leaving the original entries unmodified. To merge entries, select the entries and choose **Merge** from the **Entry** menu. The first new entry in a session is automatically named `merge`. Subsequent combination entries are named `merge.mrg1`, `merge.mrg2`, etc. When you merge entries, Maestro places the new entry after the last entry that was included in the merge operation. The new entry inherits the title from the first selected entry.

8.5.6 Splitting Entries by Molecule

Entries that contain multiple molecules can be split into entries that contain individual molecules. For example, you can split a protein-ligand complex into separate entries for the protein and the ligand. To split an entry, select it and choose **Split by Molecule** from the **Entry** menu. When an entry is split, a copy of the original entry is kept, and a new entry is created for each molecule in the original entry. The exception is that water molecules are collected in a single entry, with the suffix `-waters` appended to the original entry name. The new entries retain the same base name as the original entry, with `-disn` appended. The number *n* indexes the new entries.

8.5.7 Deleting Entries

You can delete entries from the Project Table by selecting the entries, then choosing **Delete** from the **Entry** menu or pressing **DELETE**. The entries disappear immediately: you are not prompted to confirm your choice. However, this function is undoable by selecting **Undo entry-delete** from the **Edit** menu of the Maestro main window.

8.5.8 Moving Entries

Entries in the Project Table can be moved to another location in the list. Multiple entries can be moved simultaneously, and entries need not be contiguous. To move entries, select them and then drag them to the desired new location. The entries are placed after the first unselected entry that precedes this location. For example, if you select entries 2, 4, and 6, and release the mouse button on entry 3, these three entries are placed after entry 1, because entry 1 is the first unselected entry that precedes entry 3. To move entries to the top of the table, drag them above the top of the table; to move entries to the end of the table, drag them below the end of the table.

You can also move the rows that you have selected by choosing **Move Rows** from the **Entry** menu. In the **Move selected entries to row** dialog box, you can specify the starting row number, or choose to move the rows to the beginning or the end of the table.

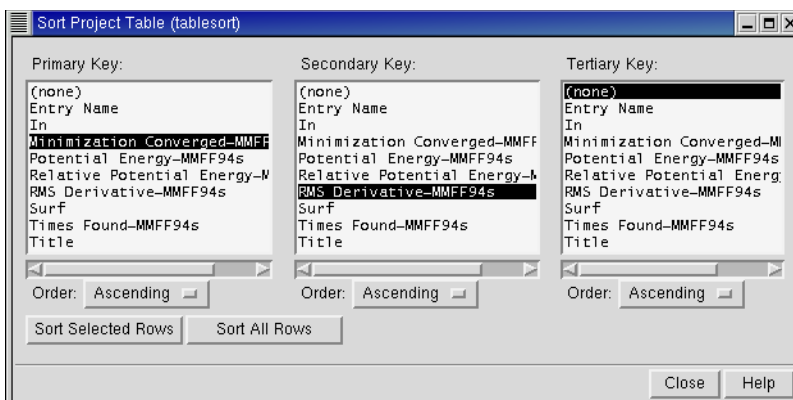


Figure 8.5. The Sort Project Table panel.

8.5.9 Sorting the Entries of a Project

You can sort all entries or the selected entries by property. For example, if you sort project entries in descending order based on their inclusion state, entries included in the Workspace are placed at the top of the list, and excluded entries are placed at the end. Sorting these entries in ascending order places included entries at the end of the list.

To sort entries by a single property, right-click and hold in the property header and choose Sort from the menu that is displayed.

To sort entries by more than one property, you can use the Sort Project Table panel. In this panel you can choose from three sort keys to sort the entries, and for each key you can choose to sort in ascending or descending order.

To sort entries:

1. Select the desired entries.
2. Choose Sort Selected Entries from the Table menu or click the Sort toolbar button.



The Sort Project Table panel is displayed.

3. Select a property from the Primary Key list.
4. If you want to sort entries whose primary key is the same by another property, select a property from the Secondary Key list, and likewise for the Tertiary Key list.
5. Select Ascending or Descending for each key for which you made a selection.

6. Click Sort Selected Rows or Sort All Rows.

If three keys is not sufficient to sort the table in the way you want, you can use the `tablesorthfields` command in the Command Input Area to choose as many keys as you want. See the *Maestro Command Reference Manual* for the command syntax.

8.6 The ePlayer

The ePlayer allows you to view in succession the structures associated with selected entries. By carefully selecting entries to be included in the Workspace, and then selecting an entry to be used as a reference, you can use the ePlayer to visualize structure relationships such as ligand-receptor docking. To use the ePlayer, you must have a project open with at least one entry selected. The ePlayer can be operated from the toolbar or from the ePlayer menu.

8.6.1 ePlayer Modes

The ePlayer has two modes, Continuous mode and Step mode. In Step mode, you manually advance the ePlayer to show each new entry structure. You can use the Next and Previous buttons on the toolbar, or the Next and Previous items on the ePlayer menu.



To go to the beginning or the end of the entry list, use the Go to Start or Go to End toolbar buttons or ePlayer menu items.



In Continuous mode the ePlayer automatically displays all selected entry structures in sequence. You can play through the structures once, or you can loop over the series by repeating the sequence of structures in the same direction, or by reversing the direction. These options can be selected from the Loop submenu of the ePlayer menu or the Loop button menu on the toolbar.



If you choose None, the ePlayer displays the selected structures once in the direction you choose and stops.

If you choose **Single Direction**, the ePlayer displays all selected structures in sequence, then loops back to the first structure and redisplay the entire series until you click another ePlayer directional button.

If you select **Oscillate**, the ePlayer displays all of the selected structures in the forward direction, then in the reverse direction, and repeats this sequence until you click another ePlayer directional button.

8.6.2 ePlayer Options

The ePlayer Options panel provides access to controls for the ePlayer and extra actions that can be performed during play. To open the ePlayer Options panel, choose **Options** from the ePlayer menu. The panel has three sections: general options, options for superimposing structures, and options for running a script.

General

You can adjust the speed at which the ePlayer moves through the specified structures by entering a value in seconds in the **Speed** text box. This is the time that each structure is displayed in the Workspace. The default is zero seconds per entry, which means that the speed of play is limited only by the speed of rendering of the structures.

In **Continuous** play mode, the structures are only temporarily displayed in the Workspace. By default, any changes in the structure are also temporary. However, you can save changes during play by selecting **Save workspace at each step**.

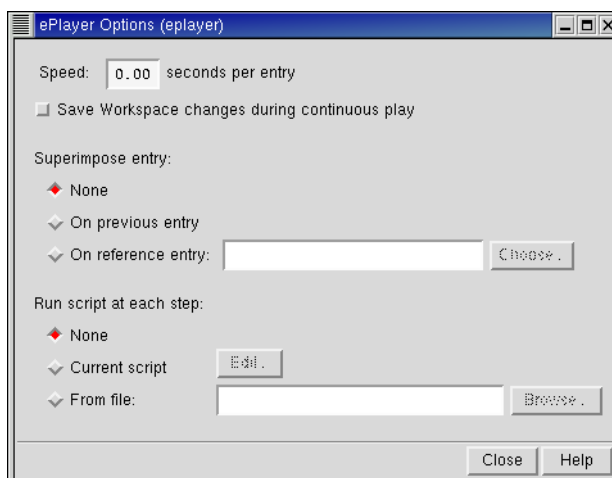


Figure 8.6. The ePlayer Options panel.

Superimpose Entry

One instance in which you may find the ePlayer very helpful is in superimposing entries onto each other or onto a reference structure. For instance, when viewing ligands, you might want to view them superimposed on a receptor. You can do this by fixing the receptor in the Workspace but selecting only the ligands.

If you choose **On Previous Entry**, each entry is superimposed on the previous entry in the series. If you choose **On Reference Entry**, you must specify a reference entry, either by typing its name in the Reference text box, or clicking the **Choose** button and choosing an entry from the list. The reference entry can be any entry, but to view it during play, you must fix it in the Workspace.

If you choose to superimpose structures in **Step** mode, the project synchronization mode that you have selected in the Preferences panel is applied. If you have chosen **Automatically** for the synchronization mode, superimposing structures changes the coordinates of the structures involved. If you do not want to save the coordinates that result from using the ePlayer, choose **Only When Specifically Requested** for the synchronization mode.

In **Continuous** mode, coordinate changes that result from superposition are saved to the project only if you have selected **Save workspace** at each step in the **General** section.

Run script at each step

In addition to superimposing structures, you can execute a command script when each structure is displayed. If you choose **Current Script**, the script that is currently loaded in the Command Script Editor panel is executed for each structure. You can open this panel and edit the script by clicking **Edit**. If you choose **From File**, you can choose the file that contains the command script you want to execute. This script is not loaded into the Command Script Editor panel, but is stored temporarily.

8.7 Entry Properties

Maestro supports incorporation and display of properties and calculated data for project entries. You can display, delete, sort, and filter entries based on property values. You can select a subset of properties to view, and you can add, delete, copy, and rename properties. Some properties of structures can be calculated by Maestro, such as the molecular weight, number of residues, and statistics of structural features. You can export properties to and import properties from a spreadsheet file. See [Section 3.3 on page 33](#) and [Section 3.4 on page 34](#) for information on exporting and importing properties.

Most of the operations you can perform on properties are available from the **Property** menu. Others, such as plotting, importing, and exporting, are available from the **Table** menu.

8.7.1 Creating New Properties

Apart from incorporating results from a calculation and importing data from a spreadsheet, there are some ways in which you can create new properties. One is to add a property column and insert values into it manually. Another is to copy an existing property and modify the property values. A third is to calculate one of a selection of new properties.

To add a property column to the Project Table, choose Add from the Property menu. In the Add Property panel, you can specify the name, the type, and the initial value of the new property (see Figure 8.7). When you add the property, all selected entries are assigned this property value. If you don't specify a value, the property is created but no values are assigned. For real-valued properties, you can set the precision to which the value is displayed.

To create a new property by copying an existing property, choose Copy from the Properties menu to open the Copy Property panel. In this panel, select the property in the From list, select <Add New Property> in the To list, provide a property name in the New Property Name text box, then click All Entries or Selected Entries. A new property column is created, and the values are copied for all entries or for the selected entries. See the next section for more information on copying properties.

Maestro provides a selection of properties that can be calculated for the entry from the structure itself. These calculations are done from the Calculate Properties panel, which you open by selecting Calculate from the Property menu. The available properties, given in the Properties list at the top of the panel, are:

- Number of atoms
- Number of residues
- Number of molecules
- Molecular weight
- Molecular formula
- Molecular charge
- Secondary structure content
- Substructure

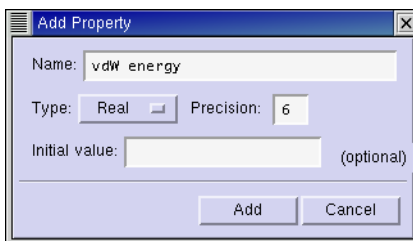


Figure 8.7. The Add Property panel.

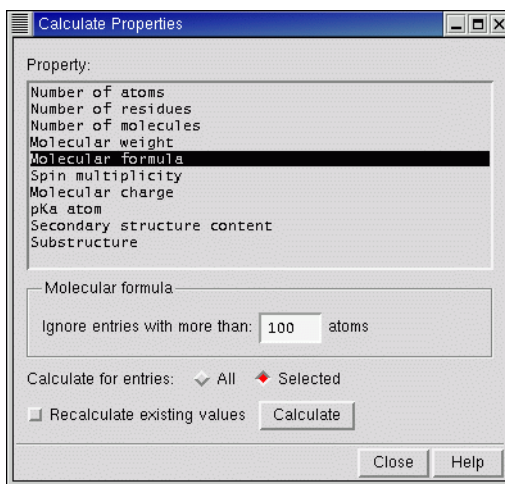


Figure 8.8. The Calculate Properties *panel*.

The Secondary structure content property is actually three properties: percentage of helices, percentage of beta strands, and percentage of random coils. The Substructure property is the number of occurrences in each entry of the substructure that you provide. Information on Maestro substructure notation is contained in the file `$SCHRODINGER/services-vversion/doc/mmsubs_syntax.txt`.

For some properties, there are options or settings that must be made in order to calculate the property. These are available in the area below the property list. To calculate the property, select it in the Property list, select any options and enter any data in the section below the list, select a Calculate For Entries option (All or Selected), then click Calculate.

If you have previously calculated the property for some entries and want to calculate it for other entries, you can choose to recalculate the existing values by selecting Recalculate existing values.

8.7.2 Changing Property Values

Maestro provides several means of changing the values of properties. You can edit individual cells and the selected cells of a single column, and you can copy values from another property.

To edit the contents of an individual Project Table cell, click the cell, edit the contents, then press ENTER or click elsewhere in the Project Table. If the cell is not in a selected entry, the entry it is in is selected and all others are deselected.

To set the value of a property to a particular value for the selected entries, click one of the cells, edit the contents, then press CTRL+ENTER. The new value is copied to all selected entries.

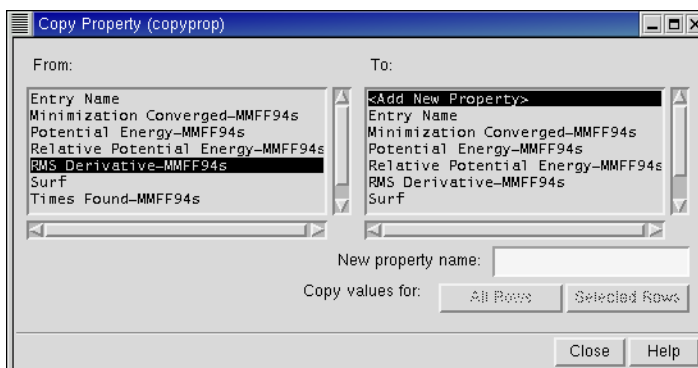


Figure 8.9. The Copy Property panel.

You can copy and paste property values using the mouse. To copy a value from the Project Table, click in the cell, then drag over the value you want to copy. To paste the value, click in the cell you want to copy the value to, then middle-click to paste the value into the cell. The value you copy does not need to come from Maestro—you can obtain it from any location.

You can copy values from one property to another in the Copy Property panel, which you open by selecting Copy from the Property menu. To copy values, specify the source property (From) and the destination property (To), then click All Rows or Selected Rows to copy the values.

The source and destination properties do not need to have the same data type. The result of the copy operation depends on the characteristics of the data. An entry must always have values for Entry Name, In, Surf, and Selected. These property values cannot be deleted. The Surf property cannot be overwritten. The Entry Name can only be replaced by a legal entry name. If the name is legal but not unique, it is made unique and a prompt to accept the modified name is displayed. If the source value is missing or the conversion fails when copying to the In or Selected properties, these properties are set to false. In addition, the Title property is required by Maestro, so if you delete it, Maestro creates a replacement title property.

The following conditions determine what happens when the destination is not one of the four special properties:

- If the destination property does not exist, it is created using the data type of the source property and the source data is copied into it.
- If the destination property does exist and is of the same type as the source, the destination property values are replaced.
- If the destination property exists and is of a different type, an attempt is made to convert the data type to that of the destination. If the conversion fails, the destination property is deleted. If it succeeds, the converted data replaces the original data.

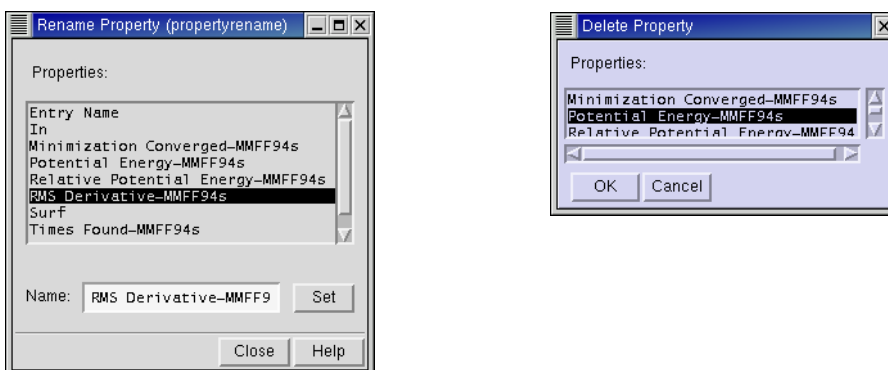


Figure 8.10. The Rename Property *panel* and the Delete Property *dialog box*.

- If fields in the source property column are empty, the corresponding fields in the destination property column are deleted.

Note that the row number is not a property of the entry.

8.7.3 Setting the Number of Decimals in a Property

You can change the number of decimal places displayed for a given property, by choosing Decimal Places from the Property menu. The Set Decimal Places dialog box is displayed, in which you can select one or more properties and set the number of decimal places to display.

8.7.4 Renaming and Deleting Properties

You can rename and delete properties other than the fixed properties, Entry Name, In, Surf, Title, and Selected. To rename a property, you can edit the table cell directly by clicking in the cell, changing the text, and pressing ENTER or clicking elsewhere in the Project Table. You can also rename properties by choosing Rename from the Property menu, and editing the name in the Rename Property panel.

To delete a property, choose Delete from the Property menu to open the Delete Property dialog box, select the properties in the list and click OK. You can choose multiple properties by shift-clicking and control-clicking the properties.

8.7.5 Displaying Selected Properties

If you want to limit the number of properties displayed, you can select the properties to be displayed in the Show Property dialog box. To open this dialog box, choose Show from the Property menu. You can select multiple properties in the Properties to show list by shift-clicking and control-clicking the list items.

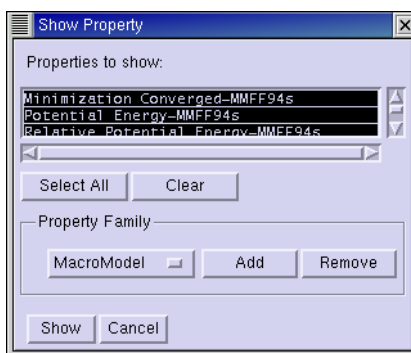


Figure 8.11. The Show Property dialog box.

Below this list is a set of tools for selecting properties. Select All selects all properties in the list and Clear clears the list. In the Property Family section, you can select properties that were generated by one of the computational programs, such as QikProp, or user-generated properties—properties that you added. Click Add to add properties of that family to the selection. Click Remove to remove properties of that family from the selection.

When you have chosen the properties, click Show to change the Project Table display. You can redisplay all properties by choosing Show All from the Property menu on the Project Table.

As well as reducing the amount of data displayed in the Project Table, limiting the displayed properties can be used to limit the properties that are exported to a spreadsheet file.

The Entry Name is treated separately from other properties. It is the unique identifier for an entry that is used and maintained by Maestro. You can show or hide the Entry Name by choosing Show Entry Name from the Property menu.

8.7.6 Moving and Resizing Property Columns

You can move and resize individual columns in the project table with the mouse. Dragging a column border with the middle mouse button resizes the column. Dragging the column heading moves the column. You can only move columns that are in the scrollable area. By default, the Title property is in the fixed area and the Entry Name property is in the scrollable area. You can swap these two columns by choosing Swap Title/Entry from the Columns submenu of the Table menu.

You can also resize all columns to fit the width of the data, to fit the width of the data and the heading, or to a fixed width by choosing the appropriate item from the Columns button menu or the Columns submenu of the Table menu.

8.7.7 Exporting and Importing Data

Maestro does not provide spreadsheet capabilities. However, you can export data to a spreadsheet file, operate on the properties, and re-import the modified data. These capabilities are available from the Import and Export submenus of the Table menu, and are discussed in detail in [Section 3.3 on page 33](#) and [Section 3.4 on page 34](#).

To specify a filename for the exported data, type a name in the File text box. To export the file to Maestro's current working directory, click Export. To export to a different directory, click Export To and navigate to another directory.

8.8 Undoing Project Operations

All changes made directly to the Project Table are saved immediately, regardless of your synchronization preference setting. The synchronization setting affects only the changes made in the Workspace. (For more information about synchronization preferences, see [Section 12.2.2 on page 153](#).) You can save a copy of your project using Maestro's Snapshot feature before performing operations that you might want to undo. To generate a snapshot of the current project, choose Make Snapshot from the Project menu.

If the synchronization mode you have selected is When approved by user, and you have modified entry structures in the Workspace but have not synchronized the changes, you may be prompted to do so. Otherwise, the Make Snapshot command is executed immediately and you receive no notification other than the temporary display of the cursor as a stopwatch.

Only one snapshot can be preserved at a time. Any previous snapshots are overwritten when you issue the Make Snapshot command.

To revert to the snapshot you made, choose Revert to Snapshot from the Project menu.

If you want to discard all work done to a particular project during a Maestro session, you can revert to the opening state of a project, which is saved by Maestro. Like reverting to a snapshot, reverting to the opening state of a project discards all work performed on that project since it was last opened in the current Maestro session. To scrap the project in its present state and begin again, choose Revert to Opening State from the Project menu. Since the reversion is not undoable, you will be prompted to confirm your choice.

Tools

Maestro is equipped with tools for investigating and changing the properties of structures in the Workspace. While some of these tools, such as plotting ([Chapter 10](#)) are discussed elsewhere, most of the tools are discussed here.

9.1 Displaying Markers

Maestro uses markers to mark atom selections in the Workspace. These markers provide visual confirmation that atoms have been successfully selected, or that an analysis task or operation has been performed.

For most operations, you can choose whether or not to display markers. To display the markers for a particular operation, select **Show Markers** in the panel used to perform the operation. When you are no longer using a property represented by a set of markers, you can hide the markers so that the Workspace does not become cluttered. The atom sets remain defined, so that you can redisplay the markers later.

Most panels display markers by default, and undisplay their associated markers when you close the panel. The exception is measurement markers, which are displayed even if the **Measurements** panel is hidden.

9.2 Making Measurements in the Workspace

The **Measurements** panel (see [Figure 9.1 on page 118](#)) provides tools to determine the geometric relationships between atoms in the Workspace, and other geometry-dependent properties. These tools are contained in six tabs—Distances, Angles, Dihedrals, H-Bonds, Contacts, and Coupling.

To open the **Measurements** panel, choose **Measurements** from the **Tools** menu.

9.2.1 Measuring Distances, Angles, and Dihedrals

The controls for Distances, Angles, and Dihedrals are very similar, as are the actual measurement procedures for these quantities. To measure a distance, you need to specify only two atoms or one bond, for an angle you must specify three atoms or two bonds, and for dihedrals you must specify four atoms or three bonds.

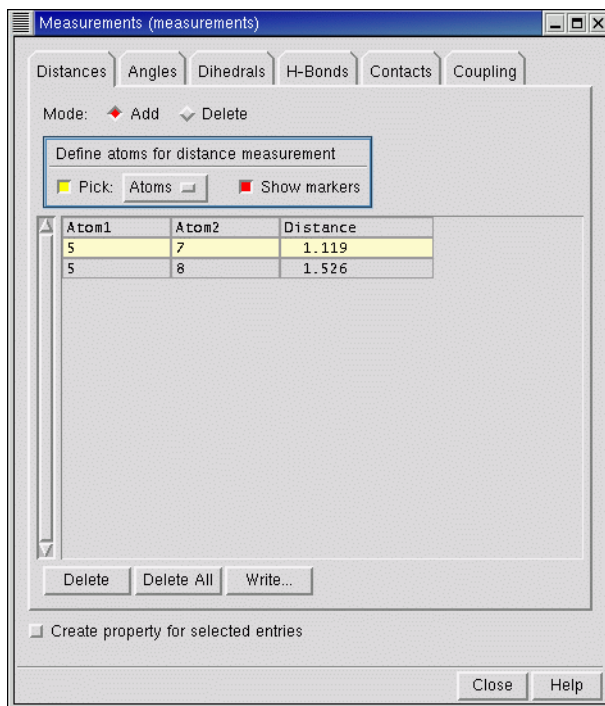


Figure 9.1. The Distances *tab of the* Measurements *panel*.

To select atoms for measurement, select the Add mode, then choose Atoms or Bonds from the Pick option menu in the Define atoms for measurement section, and select the atoms or bonds in the Workspace. When you have selected the requisite number of atoms, the measurement is made and entered into the measurement list. If Show markers is selected, a dashed line indicating the atom locations and the numerical value obtained from the measurement is displayed in the Workspace. Distances are marked in purple, angles in green, and dihedrals in red.

Once you have defined atoms for a measurement, those atoms remain defined until you delete the atom set or remove the structure from the Workspace. If you only temporarily want to undisplay the results of your measurements, deselect Show markers. However, if you no longer need the measurement value itself, you can remove it by selecting the corresponding list item and clicking Delete. To delete all of the atom sets currently defined for a particular type of measurement, click Delete All. You can also select the Delete option under Mode and click the defined atoms in the Workspace to delete a measurement.

You can make measurements with the Measure distances, angles or dihedrals toolbar button.



To make a measurement, choose Distance, Angle, or Dihedral from the button menu, then click the atoms in the Workspace. The measurement is made, and any options that are set in the Measurements panel are applied: for example, if Create property for selected entries is selected, the property is created.

9.2.2 Displaying Hydrogen Bonds and Contacts

The procedures for displaying H-Bonds and Contacts are similar. Unlike distance, angle, and dihedral measurements, you can display hydrogen bonds and contacts for an unlimited number of atoms, and you define the desired atoms by placing them in sets. Two sets can be defined at any given time, and you can choose to display the hydrogen bonds or contacts between the atoms of a set or between the two separate sets.

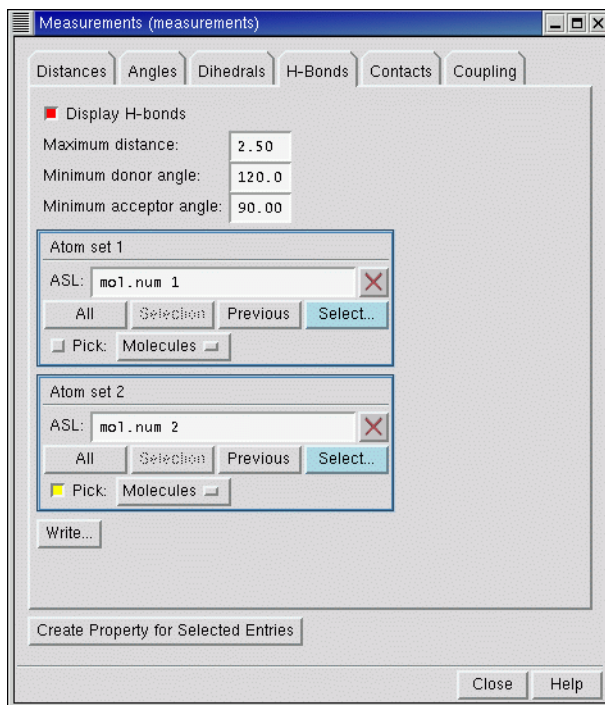


Figure 9.2. The H-Bonds *tab of the* Measurements *panel.*

To see hydrogen bonds or contacts marked in the Workspace, select Display H-bonds or Display Contacts. For contacts, you must also select the type of contact—Good, Bad, or Ugly—that you want to measure.

Default criteria for defining hydrogen bonds or contacts are given in the text boxes at the top of the tab. To change these criteria, enter a new value in the relevant text box. For hydrogen bonds you can set the Maximum distance, Minimum donor angle and Minimum acceptor angle. For contacts, you can alter the contact cutoff ratios in the Good contacts, Bad contacts, or Ugly contacts text boxes.

To define an atom set, use the picking tools, click Select to open the Atom Selection dialog box, or enter an ASL expression in the ASL text box. If you define only Atom set 1, the H-bonds or contacts are displayed within that set. If you define both Atom set 1 and Atom set 2, the H-bonds or contacts are displayed between the two sets.

When you have finished measuring, click Clear in each of the atom definition sections.

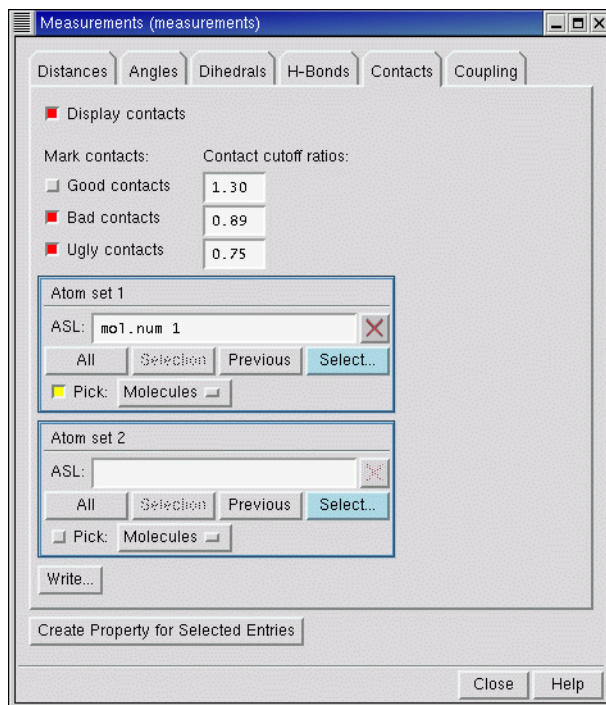


Figure 9.3. The Contacts tab of the Measurements panel.

You can also display hydrogen bonds from the toolbar. First, choose Inter H-bonds or Intra H-bonds from the Display H-bonds button menu:



Then select a molecule in the Workspace. If you choose Inter H-bonds, hydrogen bonds between the selected molecule and the rest of the atoms in the Workspace are displayed. If you choose Intra H-bonds, hydrogen bonds within the selected molecule are displayed.

9.2.3 Displaying H-H NMR Coupling Constants

NMR spin-spin coupling constants can be estimated for vicinal hydrogen atoms using a modified Karplus-type calculation. To display the estimated coupling constants, you must select two hydrogen atoms that are bonded to two carbon atoms, and these carbon atoms must be bonded to each other. The controls are the same as for distance measurements.

9.2.4 Storing the Results of Measurements

You can include measurements as properties in the Project Table, and you can export measurements to a tab-delimited or comma-delimited file.

To create properties for measurements, select **Create property for selected entries** before you define the measurements. The measured properties are added to the selected entries in the Project Table. Each distance, angle, or dihedral becomes a separate property, and is labeled with the name Distance, Angle, or Dihedral and the atom numbers that define the geometric quantity. For hydrogen bonds and contacts, the property is the number of H-bonds or contacts.

Each tab in the Measurements panel has a **Write** button, which allows you to write a file containing the measurements. The **Write** button opens a panel in which you can specify a delimiter for the file, the numbering of the atoms, and a file name. If you specify a name or a relative path in the File text box, the file is written relative to the Maestro working directory. If you click **Write To**, a dialog box is displayed in which you can navigate to the desired directory and select a file or enter a file name.

9.2.5 Defining Dummy Atoms for Measurements

Dummy atoms can be useful for making measurements. For example, you can measure the distance between the centroids of two ring systems by placing dummy atoms at the centroids. You can place dummy atoms in two ways:

1. Use the **Centroid** panel to define a centroid atom.

A dummy atom is placed at the mean position of the picked atoms.

2. Use Draw mode in the Build panel to place an atom in the Workspace, then retype it with the appropriate MacroModel atom type.

Several kinds of dummy atoms are defined in the Atom Type (MacroModel) list in the Atom Properties tab of the Build panel. If you attach the dummy atom to your structure, you can adjust its position using the Adjust panel.

9.3 Superimposing Structures

If you want to determine the similarity of structures such as a set of conformers or a family of ligands docked by Glide, you can display the structures in the Workspace and then use the tools in the Superposition panel to evaluate the differences in the atom positions. Each structure must be in a separate project entry.

The superposition scheme works by aligning the first selected atom in each structure, then attempting to align the second selected atom, and so on. The atoms are aligned to the first entry in the Project Table that is being used in the superposition (the *reference* entry). After superimposing atoms, the results of the operation are displayed in the RMSD text box near the bottom of the panel. The results include the molecule numbers, the RMS deviation of the atoms from those of the first structure, and the maximum difference between superimposed atom positions.

To open the Superposition panel, choose Superposition from the Tools menu.

There are two choices for the entries that are used for the superposition: the included entries (those in the Workspace), or the selected entries (those selected in the Project table). This choice is available under Entries to superimpose. If you choose Selected entries, you must have at least one entry displayed in the Workspace to be able to pick atoms for superposition.

To superimpose all the displayed structures, click All in the Superimpose by ASL section. The superposition is done using the atom numbers. This option is useful if each entry has the same number of atoms and is numbered the same way, which usually requires that the entries be conformers. If the structures you want to compare are not conformers, or if you only want to compare part of the structure, you can select the atoms to superimpose.

To compare two structures, you can select the atoms to superimpose in each structure in pairs. To define the atom pairs, select Pick atom pairs, then pick the atom pairs. You must pick one atom from each entry, and you must pick the first atom from the same entry each time. The atom pairs are displayed in the Superimpose by atom pairs text box as you pick them. If Show markers is selected, the atom pairs are connected by dashed lines. Once you have picked the first pair, Superimpose Atom Pairs, Delete, and Delete All are activated. If you need to delete a pair, select the pair in the list and click Delete. To delete all pairs, click Delete All. After all the atom pairs are selected, click Superimpose Atom Pairs to superimpose the two structures.

If you want to compare more than two structures or if picking individual pairs is not practical, you can superimpose atoms defined by an ASL expression or by a SMARTS pattern.

To superimpose using an ASL expression, click **Select** and construct the expression using the tools in the Atom Selection dialog box, or enter the expression in the Superimpose by ASL text box. For example, to use only carbon atoms for the superposition operation, type `atom.e C` and press ENTER. For more information on the Atom Specification Language, see the [Maestro Command Reference Manual](#).

To superimpose atoms defined by a SMARTS pattern, you can enter the pattern in the Superimpose by SMARTS text box, or you can create a SMARTS pattern from the atoms that are selected in the Workspace.

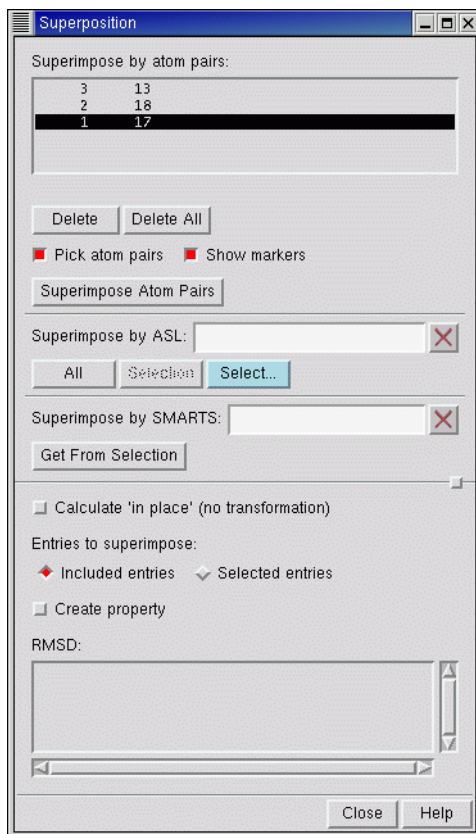


Figure 9.4. The Superposition panel.

To select the Workspace atoms and create a SMARTS pattern from them, first use the Workspace selection toolbar button to select the atoms.



Once you have selected the desired atoms, click **Get From Selection**. The SMARTS pattern corresponding to the selected atoms is entered in the **Superimpose by SMARTS** text box.

The default action when superimposing structures is to apply a transformation to the entries, then calculate the RMS difference in coordinates. If you want to preserve the coordinates of the superimposed structures, select **Calculate 'in place'** (no transformation). The RMS difference is then calculated without moving the structures. This is a convenient way of obtaining RMS differences for ligand poses in a receptor, for example.

If you do not use the calculate in place option, you can separate the superimposed entries by clicking the **Tile Entries** button in the toolbar.



The entries are aligned in a grid, but they are not returned to their original coordinates.

You can store the RMS deviation as a property for each entry that was superimposed. To create a property, select one of the options under **Create property** for before you perform the superposition. The property is labeled **RMSD *entry-name*** or **RMSD In Place *entry-name***, depending on how you performed the superposition. The reference entry has a value of zero for this property. The reference entry name is used for *entry-name* in the label. If you change the details of the superposition but do not change the reference entry, the property values are replaced. You can always rename the property if you want to keep it and perform other superpositions with the same reference entry. See [Section 8.7.4 on page 114](#) for information on renaming properties.

9.4 Sets

Sets are groups of atoms that can be used in a variety of Maestro operations. The **Sets Panel** is used to define these atom groups for subsequent operations, and can be used to save and read set definitions. Set definitions are automatically saved with the current project and are included when the project is copied or merged with another project.

To open the **Sets** panel, choose **Sets** from the **Tools** menu.

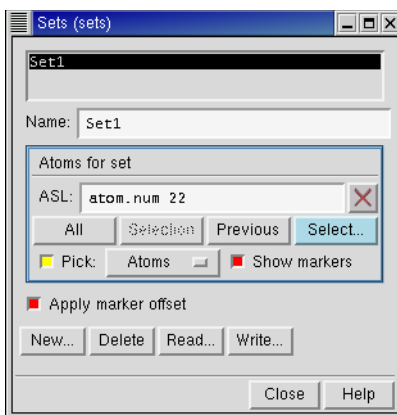


Figure 9.5. The Sets panel.

9.4.1 Creating and Deleting Sets

To create a new set, click the New button. In the dialog box that appears, enter a name for the set and click OK. This name is used to identify the set in later operations. The name is added to the list at the top of the Sets panel.

To delete a set, select the set name in the set list and click Delete.

9.4.2 Selecting Atoms for a Set

To define the atoms that make up a set, select the set from the set list, then use the tools in the Atoms For Set section to select the atoms. You can pick the atoms using the picking tools, enter an ASL expression in the ASL text box, or click Select to select atoms using the Atom Selection dialog box. The Atom Selection dialog box has a Create Set button, which you can use to create a set at any time, not just from the Sets panel. You can also select sets in the Atom Selection dialog box and perform logical operations on them to construct a new set.

If you select Show Markers, the atoms in the set are marked with green markers. If the markers obscure part of the structure that you want to see, you can move them slightly by selecting Apply marker offset. This option is selected by default.

For more information on ASL, see the *Maestro Command Reference Manual*, Chapter 3. For more information on the Atom Selection dialog box, see [Section 5.3 on page 63](#).

9.4.3 Reading and Writing Sets

You can save set definitions to disk for use with multiple projects. Click Write to open the Write Set Files panel. In this panel, navigate to the directory in which you wish to save the set defini-

tions, and enter a name for your sets definition file. The extension `.set` is automatically appended to the file name. Click **Write** to save the file and dismiss the panel.

To read an existing sets definition file, click **Read** in the **Sets** panel, then navigate to the desired directory and select the file in the **Read Set Files** panel, and click **Read**.

9.5 Assigning and Aligning Protein Structures

Maestro has tools for aligning protein structures and for assigning secondary structures in proteins.

If you want to bring two or more proteins into a common frame of reference by structurally aligning (superimposing) them, you can use the **Protein Structure Alignment** panel, which you open from the **Tools** menu. To use this panel, you must have a Prime license. For information on this panel, see [Chapter 12](#) of the *Prime User Manual* or the online help for the panel.

The **Assign Secondary Structure** command on the **Tools** menu provides a way of assigning secondary structure labels to a protein. The assignment occurs in two stages. First, all inter-backbone hydrogen bonds are identified. Second, a search is made for hydrogen-bonding patterns representing helices and strands. The assignment is used for the display of ribbons, the use of some color schemes, and for the **Workspace** sequence viewer.

9.6 Assigning Bond Orders

The **Assign Bond Orders** command on the **Tools** menu assigns bond orders to structures that do not have the correct bond orders. If you have a **Workspace** selection, bond orders are assigned to the selected atoms. Otherwise, the whole structure in the **Workspace** is assigned. Standard residues in a protein are not checked in the assignment, because they have templates for the bond orders, and the bond orders are usually assigned when the structure is imported. For molecules of less than 200 atoms, all bonds are assigned.

The assignment is based on a range of factors, including connectivity, bond length, bond angles and dihedral angles. Occasionally the bond order cannot be assigned, or is assigned incorrectly. You should always check your structure after using this command, especially for bonds that are unusual in any way. If you find any errors, you can correct them with the toolbar buttons in the **Build** panel—see [Section 4.5 on page 45](#).

If you have built a structure using **Draw** mode in the **Build** panel, you should add bond manually, because the structure is unlikely to meet the criteria for bond lengths, angles and dihedrals required for successful bond order assignment.

Assignment of bond orders is a necessary preliminary to hydrogen addition. If the bond orders are incorrect, hydrogen addition will place too many or too few hydrogen atoms.

Plotting

Maestro's Plot Facility allows you to generate XY plots of any two columns of data in the Project Table. Multiple plots are supported, and plot data is stored inside the project. When you have created a plot, you can select data points in the plot and display the corresponding entries in the Workspace.

10.1 Creating Plots

Plots are created from property data in the Project Table. To plot data, you must therefore have at least one entry in the Project Table. To create a new plot, first select the entries in the Project Table whose data you want to plot, then choose Plot from the Table menu or click the Plot button on the toolbar.



The Plot XY panel is displayed with a blank plot area. Now choose New Plot from the Plot menu. The New Plot dialog box is displayed.

To assign a name to the plot, enter a name in the Plot Name text field. The name appears in lists used to select plots, in the caption name underneath the plot, and in any TIFF or JPEG images that you make of the plot.

You can plot more than one set of data in a single plot. To distinguish the data sets, or *series*, enter a name in the Series Name text field. To create new series, however, you must use the New Series dialog box. See [Section 10.3.2 on page 133](#) for more information.

When you have assigned names, choose the property to be plotted on each axis from the X Property and Y Property lists.

The default plot type is a line plot. If you want to generate a plot of frequency of occurrence of a range of values against the midpoint of the value range, select Frequency Plot. An extra section is displayed. In this section you can choose to display the frequency on the x-axis or the y-axis. You can only choose one property for a frequency plot. The property that you select is divided evenly into ranges based on the maximum and minimum values of the property and the number of buckets you enter in the Buckets text box.

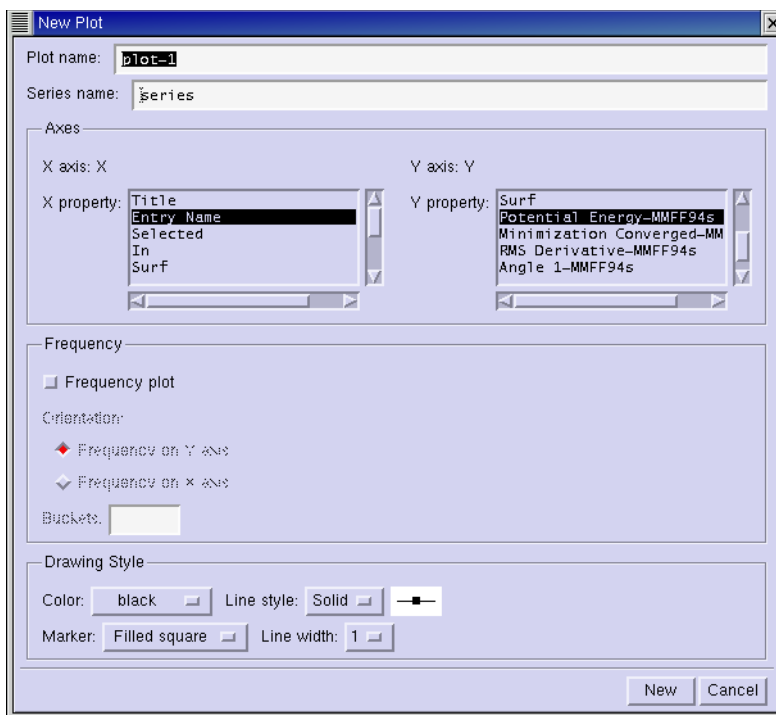


Figure 10.1. The New Plot *dialog box*.

You can modify the appearance of the data point markers and the plot line using the option menus in the Style of Drawing section: Color, Marker, Line Style, and Line Width.

When you have finished characterizing the plot, click New. The Plot XY panel is updated with the new plot.

10.2 Viewing and Manipulating Plots

The Plot XY panel is used to add, edit, view, and delete XY plots. To open the Plot XY panel, choose Show Plots from the Selection menu of the Project Table panel.

Plots are displayed in the main window of this panel. Plots can be tiled, or can be arranged either horizontally or vertically in the window. You can display as many plots as you like in the Plot XY panel, but with the addition of each plot, all plots are reduced in size.

The zoom bar is the area at the right of the Plot XY panel. It contains controls for zooming and panning the selected plots. You can display it or hide it by selecting Plot Zoom Bar from the Settings menu.

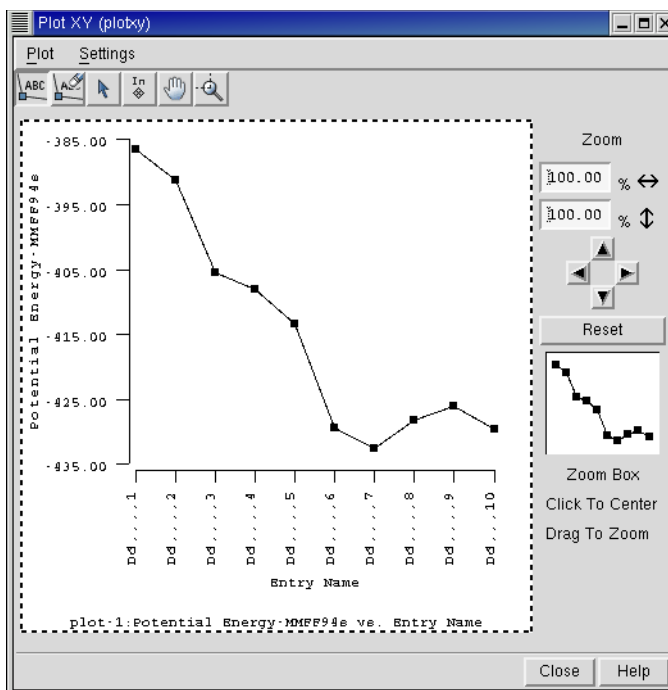


Figure 10.2. The Plot XY panel.

The Plot toolbar contains buttons for operations on the selected plot, such as labeling and unlabeled data points, selecting and deselecting entries associated with data points, including and excluding entries associated with data points, panning the plot, and zooming in on the plot.

10.2.1 Selecting Plots

You can select a plot by clicking in it. The selected plot is outlined with a dotted line. You can select more than one plot by clicking in the first plot and shift-clicking in additional plots. Once you have multiple plots selected, you can control-click in a plot to remove it from the selection. To select all plots, press CTRL+A; to deselect all plots, press CTRL+U.

You can also select plots from the Plot menu or from the list of plots in the Plot Settings panel. To open this panel, choose Plot Settings from the Settings menu or press CTRL+P.

As with entries in the Project Table, you can select plots and display plots independently. When you have selected a group of plots, you can display only those plots, add them to the display, or remove them from the display by using the commands on the Plot menu.

10.2.2 Panning and Zooming

The zoom bar controls operate on one or more selected plots. If you do not have any plots selected, then none of the zoom bar controls are active. If you have one plot selected, then the zoom controls display the plot's current zoom. If you have more than one plot selected, and all of the plots have the same zoom, that zoom is displayed. Otherwise, the zoom controls are blank. When you have multiple plots selected, the panning and zoom controls operate within each plot to pan the plotted data with respect to the axes, or zoom in on the plotted data.

The first items in the zoom bar are edit controls for directly setting the zoom for both horizontal and vertical directions. Below the edit controls are a series of arrow buttons for panning the displayed region of the selected plots. The **Reset** button restores the selected plots to their initial full-scale setting of 100% zoom and no panning. The final section of the zoom bar is the **Zoom Box**—a control for directly setting the zoom and pan relative to an image of the plot's data. This control is only active if you have a single plot selected. If you have no plots selected, or more than one plot selected, this control is blank. If the **Zoom Box** is active, you can use it to display only a portion of the selected plot. To change the current view of the plot, click inside the **Zoom Box** to center the view on that location. To zoom in, drag the mouse. The selected plot is zoomed to show only the rectangle defined by the mouse.

If you want to pan or zoom in on a single plot using the mouse, you can use the **Pan plot** and **Zoom box** buttons in the tool bar:



To pan a plot, click the **Pan plot** button, then drag the mouse in the plot in the desired direction. The data follows the mouse movement. To zoom in on a plot, click the **Zoom box** button, then drag the mouse over the region you want to zoom in on. When you release the mouse button, the plot zooms in to the defined region.

10.2.3 Labeling and Unlabeling Data Points

You can add labels to data points, or delete them, by clicking either the **Label Points** button or the **Unlabel Points** button, which are shown below.



To label points, click the **Label Points** button, then click a data point. The **Label Points** dialog box opens. In this dialog box you can label the point with any combination of: the corresponding entry name, the x-axis property value, and the y-axis property value. To remove data point labels, click **Unlabel Points**, then click on the labeled point.

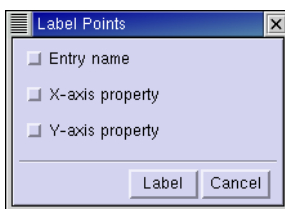


Figure 10.3. *The Label Points dialog box.*

10.2.4 Selecting and Deselecting Project Entries

If you click the Select entry button, shown below, clicking on a data point in the selected plot selects or unselects the corresponding entry in the Project Table. The usual rules for selection apply to selecting data points: clicking on its own includes one entry and excludes all others, shift-clicking includes an entry with the other included entries, and control-clicking includes or excludes the entry. You can also select multiple data points and their corresponding entries by dragging. All the data points in the bounding box have their entries selected or unselected in the Project Table.



10.2.5 Including and Excluding Project Entries

If you click the Include entry button, shown below, clicking on a data point includes the corresponding entry in the Workspace. The usual rules for selection apply to selecting data points: clicking on its own includes one entry and excludes all others, shift-clicking includes an entry with the other included entries, and control-clicking includes or excludes the entry. You can also select multiple data points and their corresponding entries by dragging. All the data points in the bounding box have their entries included in the Workspace.



10.3 Editing Plots and Plot Settings

The Plot Settings panel is used to select plots, edit plot properties, edit X and Y axes, and add, edit, and delete data series. To open the Plot Settings panel, choose Plot Settings from the Settings menu of the Plot XY panel or press CTRL+P. You can select multiple plots in the list of plots, but you can only change the settings when you have a single plot selected.

10.3.1 Modifying Plot Attributes

The tools in the Plot tab allow you to modify plot attributes other than the plot line.

Displaying and Undisplaying a Plot: To display a plot, select Plot displayed. To hide a plot, deselect Plot displayed.

Modifying a Plot Caption: When you create a new plot, a caption is displayed underneath the plot in the Plot XY panel by default. The default caption contains the assigned caption text as well as the plot name. To modify the caption content, type the new caption in the Caption text box. To change the contents or placement of the caption, or to undisplay the caption, select the appropriate options.

Displaying Selection and Inclusion Markers: By default, points corresponding to selected and included entries are not marked. To display the markers, select Display selected entry markers or Display included entry markers.

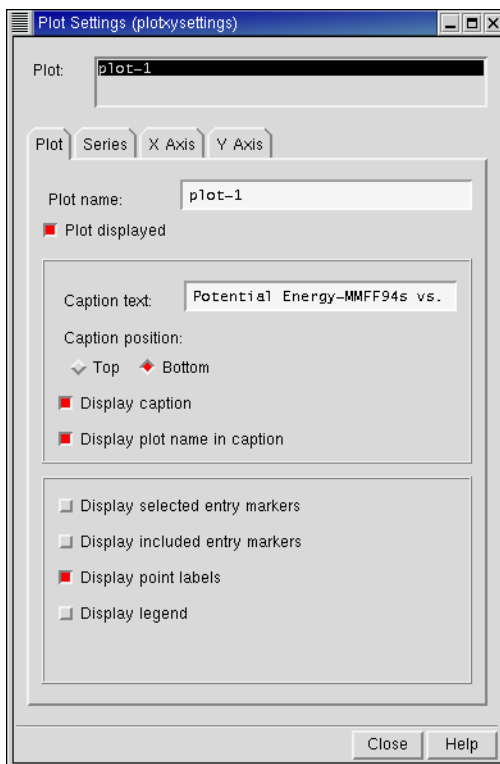


Figure 10.4. The Plot tab of the Plot Settings panel.

Displaying Data Point Labels: If Display point labels is selected, labels are displayed for data points that have labels. For frequency plots, currently, only the exact position of the point on the x-axis is available.

Displaying a Plot Legend: To display a plot legend that shows which lines and points correspond to each data series, select Display legend.

10.3.2 Operating on Plot Series

The tools in the Series tab allow you to manipulate attributes of the actual plot line, and to modify the Workspace selection based on the entries included in a particular plot series.

Selecting a Series: The series available for the selected plot are listed at the top of the tab. To perform operations on a series, select the series in the list. Information for the selected series is displayed in the lower portion of the Series tab.

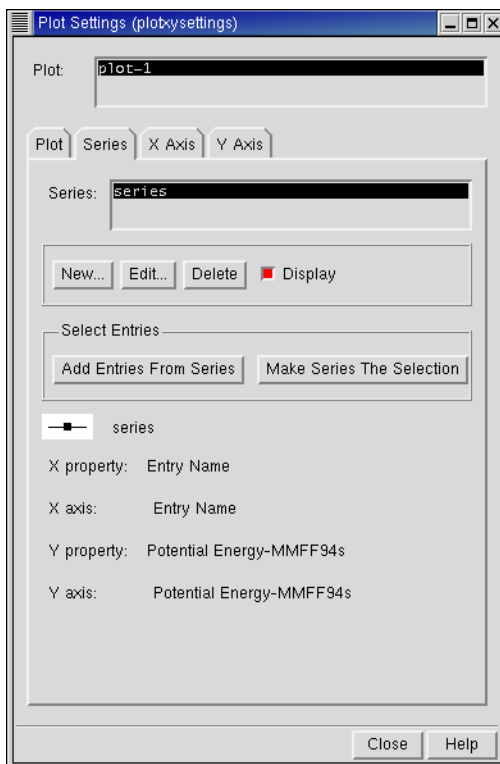


Figure 10.5. The Series *tab* of the Plot Settings *panel*.

Adding a Series: To add another data series to the selected plot, click **New**. The **New Series** dialog box is displayed. Select the desired values, and then click **New** to generate a plot line for the series. The controls are the same as in the **New Plot** dialog box—see [Section 10.1 on page 127](#) for more information.

Editing an Existing Series: To change the settings for a series, click **Edit**. The **Modify XY Data Series** dialog box is displayed. Modify the settings and click **Edit** to implement your changes. The controls are the same as in the **New Plot** dialog box—see [Section 10.1 on page 127](#) for more information. You can change all attributes of the series, including the name and the data to be displayed.

Deleting a Series: To delete the selected series click the **Delete** button.

Displaying or Undisplaying a Series: To display or undisplay the selected series, select or clear the **Display** button.

Adding Series Entries to the Project Table Selection: Click **Add Entries From Series** to add entries corresponding to the selected series to the entry selection in the **Project Table**.

Making a Series the Project Table Selection: To replace the entry selection in the **Project Table** with the entries corresponding to the selected series, click **Make Series The Selection**.

10.3.3 Changing X and Y Axis Display

The **X Axis** and **Y Axis** tabs contain tools that allow you change the names, scales, and display status of the x- and y-axes, respectively. You can display more than one x-axis or y-axis on each plot.

Selecting an Axis: The **X Axis** and **Y Axis** lists display all of the axes that exist in the selected plot. You can select an axis from the list and change its characteristics. Information on the selected axis is displayed at the bottom of the tab.

Creating an Additional Axis: To generate additional axes for the selected plot, click the **New** button and set the controls in the dialog box that is displayed. Click **New** to generate the new axis.

Editing Axis Definitions: To edit the title, minimum and maximum values, and number of scale markers for the X and Y axes, click the **Edit** button, set the controls in the dialog box that is displayed and click **Edit** to implement your changes.

Deleting an Axis: To delete an X or Y axis, select the axis from list and click **Delete**.

Displaying and Undisplaying an Axis: To display or undisplay an X and Y axis, select the axis from the list select or clear the **Display** button.

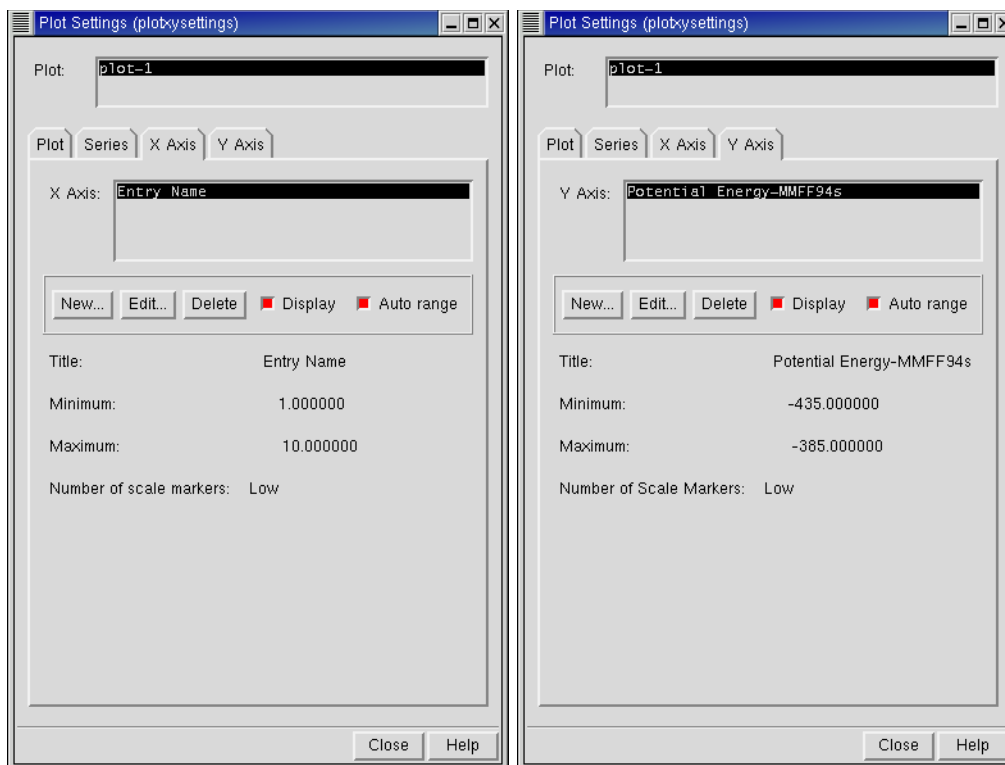


Figure 10.6. The X Axis and Y Axis tabs of the Plot Settings panel.

Setting the Axis Range: To set the axis range to the combined range of all the data series, select Auto range. To set the range manually, deselect Auto range, click Edit, and set the range.

10.4 Updating Plots to Reflect Project Table Changes

If the data values in the Project Table that are used for a plot have changed, you can update the plot with the new values by choosing Update Selected Plots from the Plot menu. For example, if you run a job with Replace selected, some cells in the table may acquire new values. If you want these values to replace old values in a plot, select the plot, then use the Update Selected Plots command.

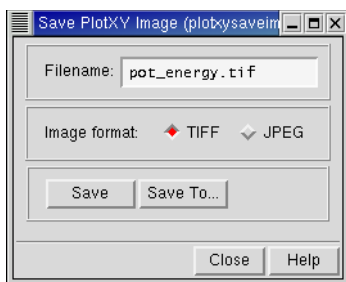


Figure 10.7. The Save PlotXY Image *panel*.

10.5 Saving Plot Images

The Save PlotXY Image panel is used to save an image of the plotting panel to a file. To open the Save PlotXY Image panel, choose Save PlotXY Image from the Plot menu.

Saved images have the same dimensions in pixels as the Plot XY panel. For images to save properly, other windows must not overlap the Plot XY panel. Two image formats can be selected from the panel: TIFF and JPEG.

TIFF

When Maestro saves an image as a TIFF file, it uses no compression. This means that the image is “lossless” (all of the information in the original image is retained), but the file can be quite large. The TIFF format is recommended for saving high quality Maestro images.

JPEG

JPEG is actually a compression method, not an image format. The JPEG file format specifies a means of compressing natural images that doesn't severely impair the perceived image quality. The JPEG method is almost always “lossy,” meaning that the reconstituted pixel data differs from that of the original.

In Maestro, the JPEG Quality parameter is fixed at 75, and the most common JPEG file format, JPEG-File Interchange Format (JFIF), is used.

Surfaces

The Maestro surface facility allows you to create and visualize molecular structure surfaces, extended radius surfaces, van der Waals surfaces, and maps of active sites in receptor molecules; to import surfaces; and to display electrostatic potential surfaces, density surfaces, and molecular orbitals generated by Jaguar.

Some surfaces are constructed directly as a set of surface elements. Other surfaces are derived from data for a function on a three-dimensional grid. This data is called a *volume*. The surface is created by selecting a value for the function and displaying the surface that corresponds to that value. The value is called an *isovalue*, and the surface is called an *isosurface*.

Generated surfaces are incorporated into the Project Table and are displayed in the Workspace by default when the associated entry is included. Surface data for any entry can be viewed using the Surface Table panel. From this panel, you can change surface display properties, such as transparency, display style, and color.

11.1 Generating Surfaces

From the Surfaces submenu of the Display menu in the main menu bar, you can generate three kinds of surfaces that have similar characteristics: molecular surfaces (Connolly surfaces), van der Waals surfaces, and extended radius surfaces.

Surfaces can only be generated for project entries. If you try to generate a surface for the scratch entry, you will be prompted to create a project entry from the scratch entry. Each surface is associated with only one entry, but an entry may have zero or more surfaces.

11.1.1 Surface Generation Controls

Most of the surface-generation settings are common to the Molecular Surface panel, the Extended Radius Surface panel, and the VDW Surface panel. All the settings are discussed below.

Name

You can name the surface by entering the name in the Name text box. If an entry has more than one surface, each surface name for that entry must be unique.

Atoms For Surface Display

The surface is displayed only for the atoms specified in the Atoms For Surface Display section of the panel. You can select atoms by picking in the Workspace, by entering an ASL expression in the ASL text box, by clicking All, or by clicking Select and using the Atom Selection dialog box. If you want to reset the atom selection, click Clear. If you have Show Markers selected, the atoms you select are marked in the Workspace.

Surface Context

The surface is actually generated for the atoms specified in the Surface context section, but only the portion of the surface defined by the Atoms for surface display is displayed. If you choose Entry or Molecule, the surface is generated for the entries or molecules of which the atoms selected for surface display are a part. If you choose Selected Atoms, the surface is generated only for the atoms selected for surface display, and appears as a complete surface.

You can also limit the surface dynamically once it is created, using the Limit Surface dialog box. See [Section 11.4.1 on page 147](#) for details.

Image Quality

In the Molecular Surface panel you can set the quality of the image to Low, Medium, or High. These controls set the grid spacing for surface generation. The high quality surface is smoother but takes longer to generate. The van der Waals and extended radius surfaces do not use grids, so there is no corresponding image quality section.

Probe Radius

This setting appears only on the Molecular Surface panel and the Extended Radius Surface panel. Since the van der Waals surface is, by definition, built directly on the van der Waals spheres, a probe radius setting is not needed for this panel.

The value in the Probe Radius text box controls the radius of the probe sphere that is “rolled over” the van der Waals surface of the molecule. The default probe radius of 1.4 Å corresponds to water.

Create Surface

To start the molecular surface generation, click Create Surface. During this time Maestro is inoperative and no user interaction is possible. For large surfaces and high quality, surface generation can be time consuming.

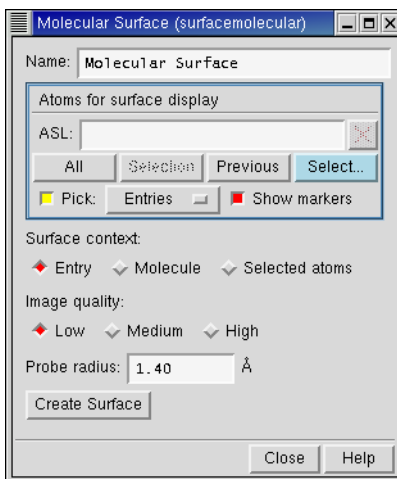


Figure 11.1. The Molecular Surface panel.

11.1.2 Molecular Structure Surfaces

In the Molecular Surface panel, you can generate a surface that follows the molecular surface, but is smoother than the CPK surface. The best known type of molecular surface is the Connolly surface. Imagine rolling a ball, representing a solvent molecule, over the assembly of spheres constituting the van der Waals surface. Where the outside of the ball contacts the van der Waals spheres, the surfaces coincide. Where two van der Waals spheres intersect, the Connolly surface follows the surface of the solvent sphere, producing fillets rather than cusps. Changing the solvent sphere or probe radius affects the size and shape of the surface. Reducing the grid spacing improves the accuracy of the surface, but increases the calculation time.

To open the Molecular Surface pane, choose Molecular Surface from the Surfaces submenu of the Display menu. Instructions for using the settings in the Molecular Surface panel are described in [Section 11.1.1](#) on page 137.

11.1.3 Extended Radius Surfaces

In the Extended Radius Surface panel, you can generate a surface constructed on van der Waals spheres in a similar way to a Connolly surface. In an extended radius surface, the point that corresponds to the center of the probe is used to generate the surface, rather than the point representing the closest contact of the probe with the surface. An extended radius surface therefore shows roughly the surface that the centroid of a solvent or a ligand cannot penetrate. Thus, looking at bound ligands with wire bond representation, the ligands would “hug” but not penetrate the extended-radius surface.

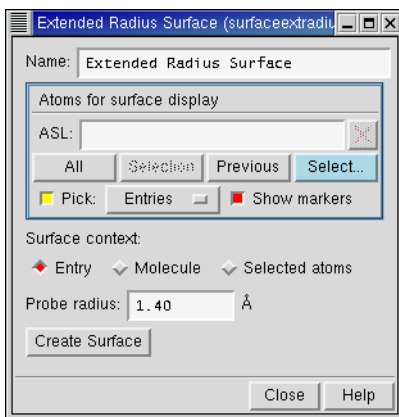


Figure 11.2. The Extended Radius Surface *panel*.

To open the Extended Radius Surface panel, choose Extended Radius from the Surfaces submenu of the Display menu. Instructions for using the settings in the Extended Radius Surface panel are described in [Section 11.1.1 on page 137](#).

11.1.4 van der Waals Surfaces

The VDW Surface panel is used to generate a van der Waals surface. This surface is identical to the CPK molecular representation with the percentage set to 100.

To open the VDW Surface panel, choose van der Waals from the Surfaces submenu of the Display menu. Instructions for using the settings in the VDW Surface panel are described in [Section 11.1.1 on page 137](#).

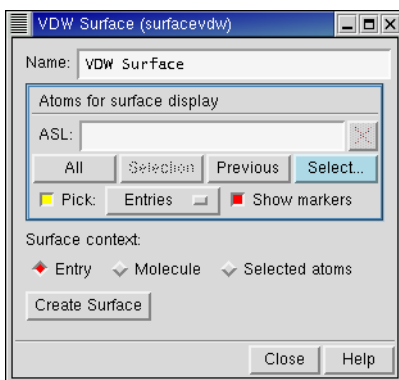


Figure 11.3. The VDW Surface *panel*.

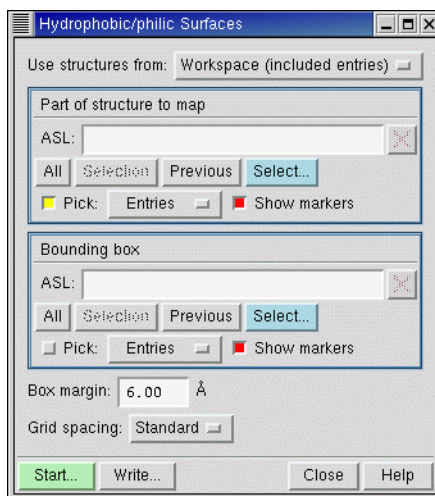


Figure 11.4. The Hydrophobic/philic Map *panel*.

11.2 Hydrophobic/philic Maps

Maestro provides a tool that is designed to help you to visualize the preferred locations of ligand atoms in a receptor site. Using the Hydrophobic/philic Map panel, you can generate a surface that displays hydrophobic and hydrophilic regions within the active site of a protein. To open the Hydrophobic/philic Map panel, choose Hydrophobic/philic Map from the Surfaces submenu of the Display menu.

11.2.1 Background

Given a receptor structure, the accessible space in the active site is partitioned into three types of regions:

1. **Hydrophobic**—regions that are favorable for occupancy by hydrophobic ligand groups
2. **Hydrophilic**—regions that are favorable for occupancy by hydrophilic ligand groups
3. **Neither hydrophobic nor hydrophilic**—regions that are of mixed character or that are far enough from the receptor surface to be similar to bulk water, where to a first approximation any group could reside with little effect on the binding affinity.

Hydrophobic and hydrophilic regions are marked by surface contours that enclose the region. The “neither” regions are implicit; these are simply regions that are accessible to the ligand but are not marked as being either hydrophobic or hydrophilic.

By revealing “targets of opportunity”—e.g., hydrophobic regions that have room to accommodate a larger hydrophobic group—active site maps can aid in the design of new ligands. Alternatively, by showing the degree to which poses produced by a program like Glide display or violate proper complementarity to the receptor site, site maps can aid in the evaluation of docking hits. The “neither” regions are also important because they are regions in which the physical properties of the ligand can be changed, for example to make the ligand more or less soluble, with minimal expected effect on the binding affinity.

The active-site mapping procedure operates in a manner analogous to Goodford’s GRID algorithm (P. Goodford, *J. Med. Chem.* **1985**, 28, 849). Hydrophilic and hydrophobic regions are defined in a way that considers both spatial proximity to the receptor and suitability for occupancy by solvent water. A putative van der Waals energy and the magnitude and direction of the electric field (calculated using a distance-dependent dielectric formulation) are computed for a probe centered at each grid point by considering interactions with all atoms of the receptor site within a defined cutoff distance. In contrast to techniques that color-code the receptor surface to represent hydrophilicity or hydrophobicity, these maps depend on more than the character of the nearest receptor atom. Moreover, the maps explicitly show the shape and extent of hydrophilic and hydrophobic regions, something a receptor-surface display cannot do. The maps behave rather like an “extraradius” surface in that the atoms of a stick figure representation of the ligand can approach but should not penetrate the map surface (except in hydrogen bonding regions, where internuclear distances are expected to be smaller than normal contact van der Waals distances).

Hydrophilic Map: The hydrophilicity is defined by adding an “electric-field reward” term to the vdW energy:

$$\text{Grid_philic} = \text{vdW_energy} + \text{oriented-dipole_energy} \quad (1)$$

where the oriented-dipole energy is necessarily negative. Hydrophilic regions are those within which the sum of the two terms is sufficiently negative, and are displayed by contouring the hydrophilic grid at a prescribed negative isovalue, typically -6 kcal mol^{-1} .

Hydrophobic Map: The quantity representing hydrophobicity is constructed by adding an oppositely signed (positive) “electric-field penalty” term to the vdW term:

$$\text{Grid_phobic} = \text{vdW_energy} - 0.15 * \text{oriented-dipole_energy} \quad (2)$$

Hydrophobic regions are those for which the favorable van der Waals term is not too strongly degraded by the positive electric-field penalty. Qualitatively, hydrophobic regions are those that lie suitably close to the surface of the receptor but for which the water-dipole-orienting electric field produced by the receptor is sufficiently small. In short, these are regions where

something would like to be, but water would not. Hydrophobic regions are displayed by contouring the hydrophobic grid at a suitably negative isovalue, typically $-0.5 \text{ kcal mol}^{-1}$.

11.2.2 Mapping Algorithm

The first step is to define a rectilinear box that contains an active site and to define grid points with a typical grid spacing of 1 \AA within the box.

Next, van der Waals energies and x , y , and z electrostatic field components are computed at each of the grid points. Receptor atom partial charges and van der Waals parameters are taken from the OPLS-AA force field. A probe is represented by a van der Waals sphere of radius 1.5 \AA , well depth $0.2 \text{ kcal mol}^{-1}$, and point dipole moment 2.3 Debye . The probe's point dipole is oriented along the electric field to give a minimum (negative) electrostatic energy and is offset from the vdW body (which is located at the grid point and represents the oxygen atom of a water molecule) toward the center of an optimally oriented O–H bond. A smoothing procedure is applied during the calculation of the receptor's electric field to avoid artificial singularities. The hydrophilic and hydrophobic grid values are then determined from Eqs. 1 and 2.

Finally, the hydrophilic and hydrophobic grids are read by Maestro, which contours the grids at the empirically selected default isovalues cited above (-6 and $-0.5 \text{ kcal mol}^{-1}$) and displays the corresponding hydrophilic and hydrophobic volumes as surfaces.

11.2.3 Specifying a Structure to Be Mapped

To create a map of a receptor structure, you must first include the structure in the Workspace. The structure must belong to an entry in the current project: you cannot map structures from more than one entry or from a scratch entry.

Once you have displayed the structure, specify the part of the structure to be mapped with the tools in the Part Of Structure To Map section. You can pick atoms in the Workspace, enter an ASL expression, click All, or click Select to open the Atom Selection dialog box and select atoms. If Show Markers is selected, a dark pink trace is superimposed on the selected atoms and bonds.

To clear the atom selection, click Clear. Clearing the selection does not delete existing maps.

11.2.4 Defining the Mapping Box

In addition to specifying the structure for which a map is to be generated, you must also set the location and size of the “box” of atoms to be mapped, using the tools in the Bounding Box section. You can pick atoms in the Workspace, enter an ASL expression, click All, or click Select to open the Atom Selection dialog box and select atoms. Mapping box atoms must be

associated with a project entry; they cannot be part of the scratch entry. If Show Markers is selected, an orange box is displayed in the Workspace to mark the mapping box region.

The box margin provides a buffer around the minimal defined mapping box enclosure. Adding a buffer region may be useful if you are investigating modification of the ligand to obtain a better fit in the active site, and would like to generate a surface large enough to show how modified or added atoms would interact.

To alter the box margin, enter a value into the Box Margin text box. If you have defined the mapping box region using `fillres` in an ASL expression, you may want to set the box margin value to zero.

11.2.5 Running the Job

In contrast to the generation of other types of surfaces, mapping of an active site is run as a job. After you have configured the Hydrophobic/philic Map panel settings, you can choose to either run the job or to simply write out the job files so that you can run the calculation later.

To write the job files without running the job, click the Write button. You can use these files to subsequently launch the calculation from the command line.

To set up and run the calculation, click Start. In the Start dialog box, enter a unique name for the job, enter your user name and host (if they are different from the machine on which you are starting the job), and click Start to launch the job. The Monitor panel is displayed and the job's progress is printed to the window in the upper portion of the panel.

If the job finishes while it is being monitored, the results are automatically incorporated and the site map is displayed. If the job is not being monitored when it finishes, you can open the Monitor panel from the Applications menu in the main window and monitor the completed job to incorporate the results.

Once the map for a particular entry has been generated, the surface information is incorporated into the Project Table and marked in the Surf column for the entry with a button bearing the letter S. Click this button to open the Surface Table panel and view the data associated with this surface.

11.2.6 Viewing the Site Map Surface

The site map is displayed as a pair of isosurfaces generated from the hydrophilic and hydrophobic grids. The default isovalues can be changed in the Surfaces panel if more expansive or more restrictive philic and phobic volumes are appropriate in a given application. If you need to see inside these surfaces, you can display them as Mesh surfaces, or if you display them as Solid surfaces (the default), you can increase the transparency.

11.3 Importing Surface or Volume Files

You can import Maestro `.vis` files or Jaguar `.plt` files, associate this data with a project entry, and display the data as a surface.

To import a file, choose Import from the Surfaces submenu of the Display menu, or click the Import button in the Surface Table panel. Browse for the desired `.plt` or `.vis` file and click Import. In the Choose Entry dialog box, choose the entry to which the surface or volume data belongs. Once you have associated the data, it becomes a part of the Maestro project. It can be displayed by clicking the S button in the Surf column of the Project Table or by clicking the diamond in the V column of the surface table in the Surface Table panel.

Jaguar `.plt` files always contain volume data—3D grids with data values at each point. They may contain orbital, electron density, spin density, or electrostatic potential information. These files are created using the appropriate keywords in the **gen** section in the Jaguar input file.

A Maestro `.vis` file contains volumes, surfaces, or both. A volume contains all of the data values for each of the points on the grid. Once the file is imported into Maestro, one or more surfaces can be created. In addition, if the `.vis` file still contains the original volume data used to generate an imported surface, you can change the isovalue in Maestro and generate new surfaces using different isovalues. If a Maestro `.vis` file contains only surfaces, then it contains only a description of the surface and not the original volume data from which it came. You cannot change the isovalue for such a surface in Maestro.

You can create `.plt` or `.vis` files using the Jaguar panel. When you run a plot job from the Jaguar panel, the surface and volume data is automatically incorporated into the project and displayed. See [Section 4.12](#) of the *Jaguar User Manual* for details.

11.4 The Surface Table Panel

The Surface Table panel is used to import, modify, and view surfaces and their properties. The panel consists of a table of surfaces for all entries that are included in the Workspace, and a set of controls below the table. The columns of the surface table are described in [Table 11.1](#). Scratch entries are not part of a project and thus cannot be used for surfaces.

To open the Surface Table panel, choose Surface Table from the Surfaces submenu of the Display menu.

- To display a surface in the Workspace, click the diamond in the V (Visible) field for the desired surface.
- To change the name of a surface, enter a new name in the Surface Name field.

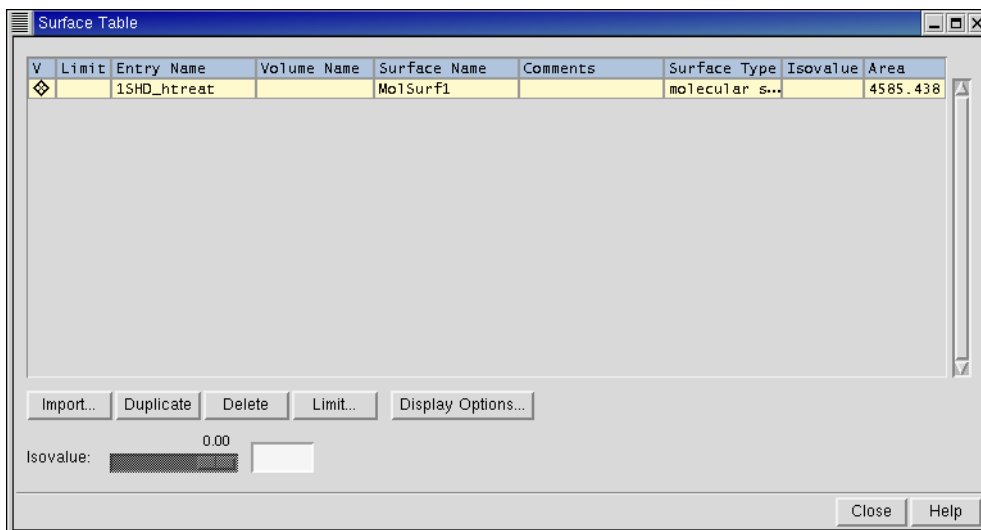


Figure 11.5. The Surface Table panel.

Table 11.1. Surface table columns.

Column	Purpose
V	Displays or undisplay the surface for that entry in the Workspace.
Limit	If limits are set for this entry, a checkbox appears. You can then check or uncheck the box to turn limits on or off. To add or edit limits, click the Limit button.
Entry Name	Lists the name of the entry that contains the surface.
Volume Name	Lists the name of the volume from which the isosurface was created.
Surface Name	Edit this field to create a unique surface name within each entry.
Comments	Edit this field to add your own comments. Jaguar fills this field in with molecular orbitals.
Surface Type	Lists the type of surface to be displayed: molecular surface, Van der Waals, extended radius, density, potential, orbital, spin, hydrophobic, or hydrophilic.
Isovalue	Lists the value that was used to create the surface. To change this value, use the Isovalue slider at the bottom of the panel or enter a value in the text box next to it. For regular surfaces, this field is blank.
Area	Lists the surface area of the limited surface (summation of all polygons).

At the bottom of the panel are buttons for importing, duplicating, deleting, and limiting surfaces, as well as controlling the display of surfaces. You can also adjust the isovalue of a volume with the slider or the text box. These controls and buttons are described in the following sections.

11.4.1 Importing, Modifying, and Deleting Surfaces

Importing Surfaces

To import a `.vis` or `.plt` file from the Surface Table panel, click Import, browse to select the desired surface file, and click Import. In the Choose Entry dialog box, select an entry with which to associate the imported data and click Choose. You can also import surfaces by choosing Import from the Surfaces submenu of the Display menu in the main window.

When you import a surface or volume file, the surfaces and volumes are already named. To change the name of the displayed surface or volume, edit the name in the Surface Name field in the surfaces table.

Changing isovalues

The isovalue is the current value of the property that was used to create the surface. Surfaces generated by Jaguar and hydrophilic/phobic surfaces have volume data for which an isovalue is specified. Molecular surfaces, van der Waals surfaces, and extended radius surfaces do not have volume data.

You can change the isovalue for a surface that has an associated volume: either enter a new number in the Isovvalue text box at the bottom of the panel or use the slider to adjust the isovalue. The surface is redisplayed with the new isovalue.

To create a new isovalue for a surface, select the surface in the table and click the Duplicate button at the bottom of the panel. Then you can change the isovalue for the new surface.

Limiting surfaces

You can limit the visible portion of the surface to the portion that is within a specified distance from a selection of atoms. Click Limit to open the Limit panel. Select the atoms using the picking tools in the Limit to ASL section, enter a value in the Distance text box, and click Apply. The limit is applied to the surface. You can change the atom selection, or click OK to accept the limit. To remove the limit, click the square in the Limit column of the surfaces table.

Deleting surfaces

To delete a surface, select it in the surfaces table and click Delete.

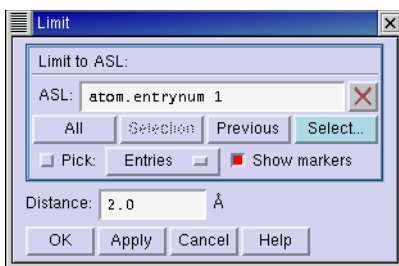


Figure 11.6. *The Limit panel.*

11.4.2 Changing the Appearance of the Surface

You can adjust various characteristics of the surface's appearance using the controls in the Display Options panel.

Style and Transparency

Surfaces can be displayed in three styles: Solid, Mesh and Dot. The transparency of Solid surfaces can be changed using the Transparency slider. Mesh and Dot surfaces are opaque unless the Transparency slider is set to 100%, when they become invisible.

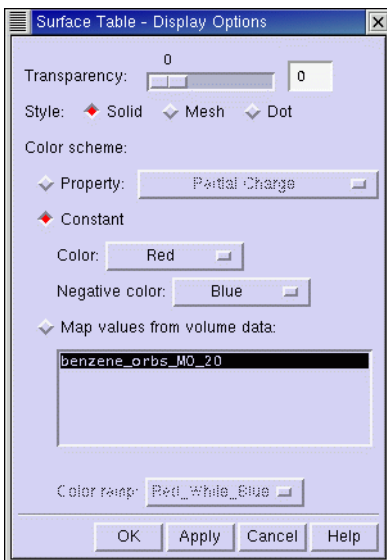


Figure 11.7. *The Surface Table - Display Options dialog box.*

Color Scheme

You can color a surface using the following options in the Color Scheme section of the panel:

- **Property**—each part of the surface is colored according to the value of the property on the nearest atom to that part of the surface. Choose a property from the Property option menu.
- **Constant**—the entire surface is colored with the same color, which you can choose from the Color option menu.

When a surface has both positive and negative portions (as with orbitals), the **Negative** color option menu is active and you can select the color in which the negative surface is displayed.

- **Map values from volume data**—the surface is colored based on the value of the volume on the surface, which you select in the list.

The color fades between the colors selected in the Color ramp option menu according to the value of the volume property on the surface. The range of colors is determined by the range of values in the volume. This scheme can be used, for example, to display the values of the electrostatic potential obtained from a Jaguar calculation on the van der Waals surface. You can also create your own color ramp (see [Section E.3 on page 217](#) for instructions).

Customizing Maestro

Many aspects of the Maestro interface can be customized from within Maestro. This chapter describes customizing the panel layout, the Maestro program preferences, the Workspace appearance and 3D viewing. Other customization tasks that can be performed using scripts are described in [Chapter 13](#), and various color-related customizations that can be performed by editing resource files are described in [Appendix E](#).

The changes described in this chapter can be saved explicitly, so that they are preserved in subsequent Maestro sessions. The changes are stored in the following directory, which is referred to hereafter as the *preferences* directory:

UNIX: `$HOME/.schrodinger/maestrovv`

Windows: `%USERPROFILE%\schrodinger\maestrovv`

Here, *vv* is the two-digit Maestro version number (75 for this release), and *USERPROFILE* is the user profile environment variable on Windows, which is normally set to the user Documents and Settings directory. For information on setting environment variables on Windows, see [Appendix F](#).

Other changes can be preserved only by creating a startup script that issues the Maestro commands necessary to generate the desired changes when Maestro is launched. For more information on startup scripts, see [Section 13.3 on page 170](#).

12.1 Changing and Saving Panel Layout

Although effort has been made to optimize the default window positions and dimensions, you may want to change them. Resizing and repositioning windows is done in the normal way by dragging a corner, a side, or the title bar of the panel. Here, “windows” means panels, dialog boxes, and the main window. You can also close all panels (except the main window), restore the panels you had open, and arrange the panels in a cascade, from upper right to lower left, by choosing the corresponding items from the Display menu.

You can save the size and location changes in a layout file, by choosing **Save Panel Layout** from the **Maestro** menu in the main window. Maestro immediately saves your panel settings to the file named `layout.res` in your preferences directory. If you want to return the panels to their default sizes and locations, delete the `layout.res` file. If you change the panel size and

location during a session and want to revert to the last saved layout, choose **Restore Saved Layout** from the **Maestro** menu.

When you upgrade to a new version of Maestro, the panel layouts you saved in a previous version are not accessible in the new version. The reason for this is that, in the new release, panels are often upgraded and their elements rearranged. Using an old panel layout may hide important features of the panel. Once you have installed the new version of Maestro, start the program, arrange the panels the way you want them, and save the panel layouts. Maestro will use the new panel layouts for all subsequent sessions.

12.2 Setting Preferences

The Maestro Preferences panel has seven tabs, each of which contains controls for modifying Maestro's default behavior according to your personal preference. Changes made to interface settings in the Preferences panel are persistent: they define Maestro's default behavior until you modify them. A file called `prefer.cmd`, which contains all of your Maestro preferences, is created and saved to your preferences directory.

To open the Preferences panel, choose **Preferences** from the **Maestro** menu.

12.2.1 Changing Font Size for Interface Component Text

The **Font Size** tab of the Preferences panel allows you to change the font size of the text that is displayed in the GUI. Changing the font size may change the amount of space each Maestro panel or window requires.

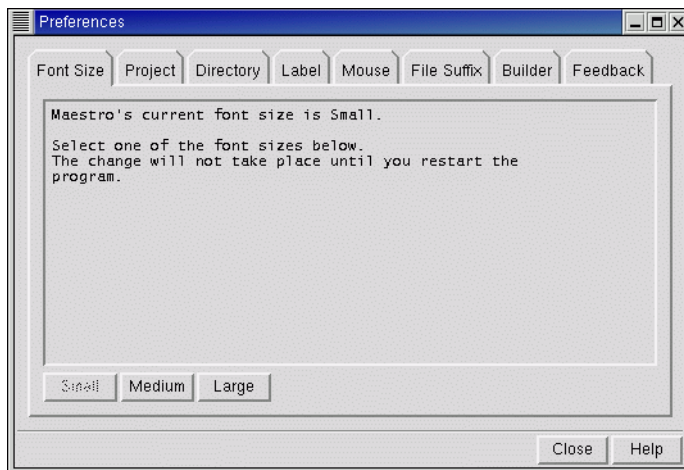


Figure 12.1. The **Font Size** tab of the Preferences panel.

To select a new default font size, click Small, Medium, or Large. The actual point size depends on the size of your monitor screen. The default is Small. After you have made your selection, restart Maestro to implement the change.

12.2.2 Setting Project Preferences

The Project tab of the Preferences panel contains a number of options for saving modifications of a Workspace structure to the corresponding project entry, saving projects, and opening projects.

Synchronize workspace changes

When you include an entry in the Workspace, you are working with a copy of the structures in that entry. By default, Maestro saves, or *synchronizes*, changes in the Workspace version of the structure to the entry with which it is associated when you make the changes. The other two synchronization preferences allow you to dictate what changes are saved to the project. If you select When approved, you are prompted to save Workspace changes when you attempt to exclude that entry from the Workspace. If you choose not to synchronize, your changes are lost. If you select Only when specifically requested, Maestro does not prompt you to save changes made to Workspace structures: you must click the exclamation point icon in the status bar to save any changes.

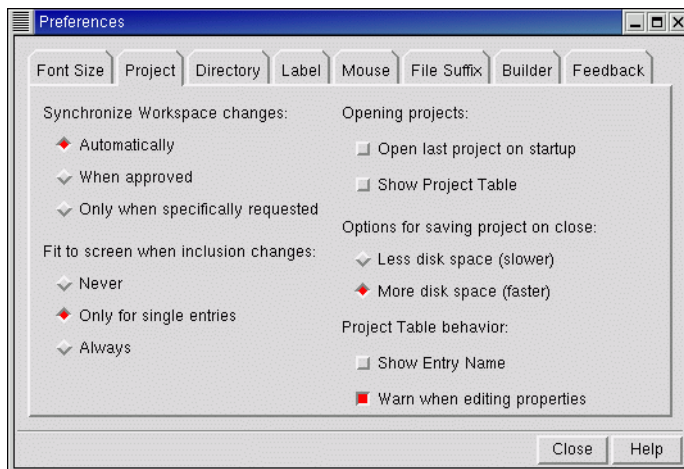


Figure 12.2. The Project *tab* of the Preferences *panel*.

Fit to screen when inclusion changes

In this section you can set preferences for when an automatic fit-to-screen is performed. The action only applies when you include or exclude an entry from the Workspace. You can choose between **Never**, which turns off automatic fitting, **Only for single entries**, which performs a fit only when there is a single entry in the Workspace, and **Always**, which performs a fit whenever an entry is included or excluded.

Opening projects

You can choose to open the last project you were working on when you start Maestro. This feature is particularly useful if you are working on one project over a period of time. You can always open a specific project when you start Maestro by using the `-p` option. See [Section 2.1 on page 7](#) for more information on Maestro options.

You can also choose to display the Project Table when you open a named project. This is the default behavior, but it does not apply to the scratch project.

Options for saving project on close

When you open a project, the initial state is stored in compressed format, and the working copy is stored in uncompressed format. When you close a project, the initial state is replaced with a compressed version of the current state of the project. You can choose to keep or to remove the uncompressed current state by selecting one of the following options:

- **Less Disk Space (slower):** The uncompressed version is removed. When you open the project again, the compressed initial state must be uncompressed before it is read. This option is useful if storage space is an issue.
- **More disk space (faster):** The uncompressed version is kept. When you open the project again, the uncompressed version is read.

Prior to the Maestro 5.1 release, projects were only stored in uncompressed format. When you open a project from a previous release, a copy of the project is compressed and stored as the initial state.

Project Table behavior

The entry name can be displayed as a column in the scrolling area of the Project Table. It is hidden by default. To display it select **Show Entry Name**.

When properties are edited, a warning message is displayed by default, and you must confirm the change to the property. If you want to turn off this warning you can deselect **Warn when editing properties**.

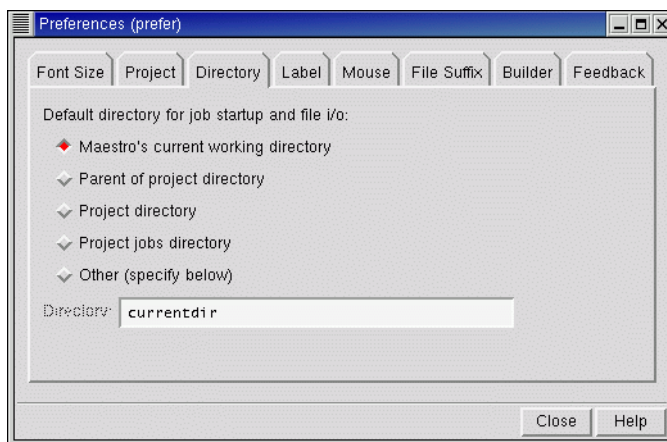


Figure 12.3. The Directory tab of the Preferences panel.

12.2.3 Specifying a Default Working Directory

The Maestro working directory is the directory to which file and project selectors open, and to which job files are written. The default initial working directory is the Maestro startup directory. By changing the settings in the Directory tab of the Preferences panel, you can specify a default working directory.

In addition to the option Maestro's current working directory, you can choose from Project directory, Parent of project directory, Project jobs directory, and Other. If you choose any of the project-related options, the working directory is automatically updated when you begin working in a different project. If you choose Other, you must enter a directory in the Directory text box. The Other option may be useful if you keep all Maestro-related files in one area of your home directory.

If you change the working directory using the `cd` command in the command input area, the selection made in the Directory tab is superseded. However, if the selection is project-related, it is reapplied if you open a new project. If you use the filter in a file or project selector, the filter is retained until you change it or change the working directory.

12.2.4 Customizing Atom Label Appearance

The Label tab of the Preferences panel contains controls for defining the separator between properties in atom labels, the font, and the location of the label. You can also choose whether to display labels for hidden atoms, and display the property name in the label as well as the property. Changes made in the label preferences are immediately applied to existing labels.

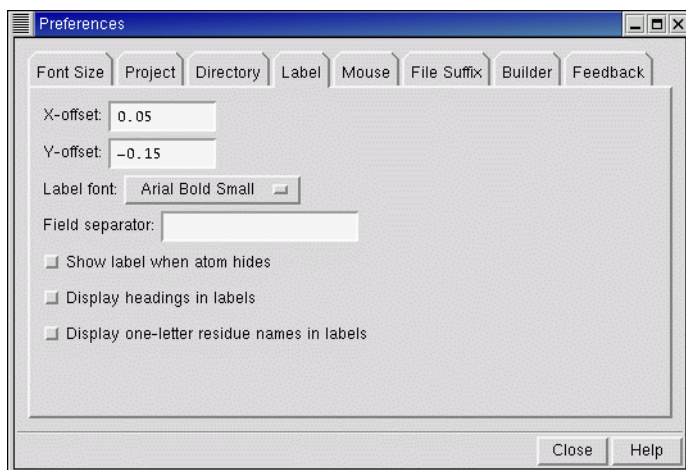


Figure 12.4. The Label tab of the Preferences panel.

- **X-Offset and Y-Offset:** To reposition labels, enter new values in these text boxes. The offsets define the location of the lower left corner of the label relative to the atom position.
- **Label Font:** To change the font, choose a font style from this option menu.
- **Field Separator:** By default, the atom label fields are separated by a space. To add a separator such as a colon, semicolon, or vertical bar, type the character in this text box.
- **Show label when atom hides:** You can display labels for atoms that are hidden by selecting this option. Labels for hidden (undisplayed) atoms are shown in the Workspace, but do not appear to be associated with an atom. Undisplaying atoms is not the same as excluding entries from the Workspace. Labels, like ribbons, can only be shown for included atoms; however, these atoms need not be displayed.
- **Display headings in labels:** If you select this option, the label includes the descriptions of the label fields. For example, if you selected Molecule Number, the label is `molnum=n` instead of simply `n`, where `n` is the number of the molecule.
- **Display one-letter residue names in labels:** If you select this option, the one-letter residue names are displayed instead of the usual three-letter or four-letter names. Only the common (naturally occurring) residues have one-letter names; the rest are labeled “X”.

12.2.5 Setting Mouse Preferences

If you want to change the sensitivity of rotation of Workspace structures to mouse movement, use the Sensitivity in rotation slider in the Mouse tab. Increasing the sensitivity means that moving the mouse results in greater movement of the molecule in the Workspace.

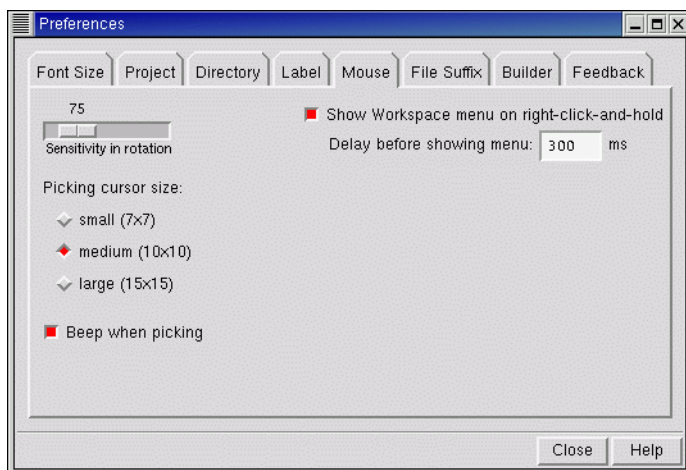


Figure 12.5. The Mouse tab of the Preferences panel.

You can select the picking cursor size by choosing small, medium or large.

You can choose whether or not audio feedback is given when you use the mouse to pick atoms in the Workspace. By default, Beep when picking is selected.

You can also choose whether to display a shortcut menu when you right-click in the Workspace, and set the delay before this menu appears.

12.2.6 Setting File Suffix Preferences

You can change the default file suffixes that are used for the file types that Maestro recognizes in the File Suffix tab. These suffixes are appended to file names when you export structures from the Export panel. They do not include suffixes for file formats converted by Babel.

12.2.7 Setting Builder Preferences

In the Builder tab you can set options to control the use of united atoms, automatic adjustment of bond lengths, removal of hydrogens to maintain valence during build operations, and deletion of terminal atoms. The options are persistent between Maestro sessions, and are selected by default.

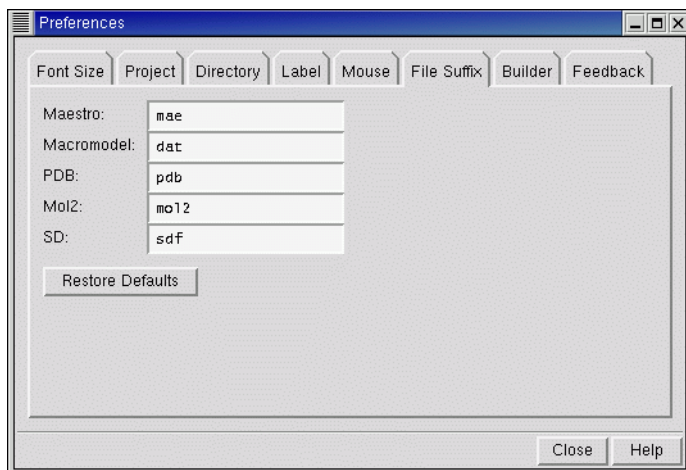


Figure 12.6. The File Suffix *tab of the* Preferences *panel.*

Allow united atom types while building

The use of united atom types can cause problems for users (especially Jaguar users) who do not want implicit hydrogen atoms in their structures. United atom types with implicit hydrogens are assigned by default when hydrogens are deleted from heavy atoms. If this option is deselected, united atom types are never assigned when atoms and bonds are explicitly deleted, atoms are retyped, or structures are drawn freehand. Instead, the atom types remain unchanged, or explicit radical types are assigned if they are available.

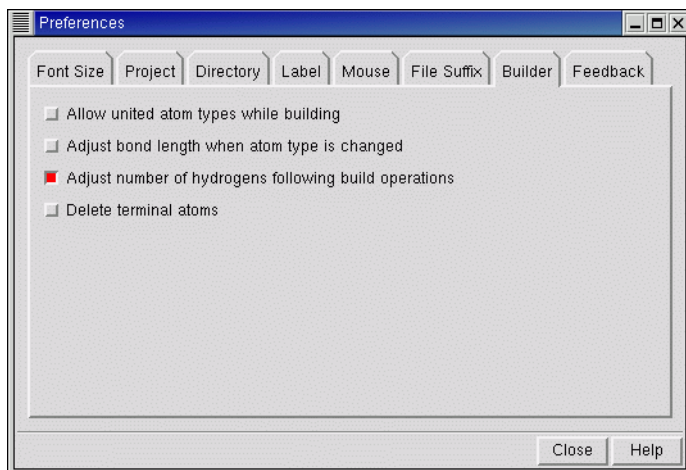


Figure 12.7. The Builder *tab of the* Preferences *panel.*

United atom types are used regardless of this option if a hydrogen treatment that requires united atom types is applied, or when reading in a structure file (such as a PDB file) that does not have explicit hydrogens present.

Adjust bond length when atom type is changed

When atom types are changed, bond lengths are automatically adjusted to give a “reasonable” bond length. This behavior is undesirable if the geometry is already correct (such as from a crystal structure) and it is only necessary to change the atom types. Deselecting this option prevents the adjustment of bond lengths.

Adjust number of hydrogens during build operations

The number of hydrogen atoms is adjusted by default to maintain a normal valence during formal charge, bond order, or retyping changes. This behavior is not always desirable, and can be prevented by deselecting this option.

Delete terminal atoms

When you delete an atom, the terminal atoms that are attached to this atom can be deleted, or they can be left in place. Terminal atoms are atoms that are bonded only to that atom, such as hydrogen or a halogen. Deleting terminal atoms would not, for example, delete a methyl group attached to the deleted atom.

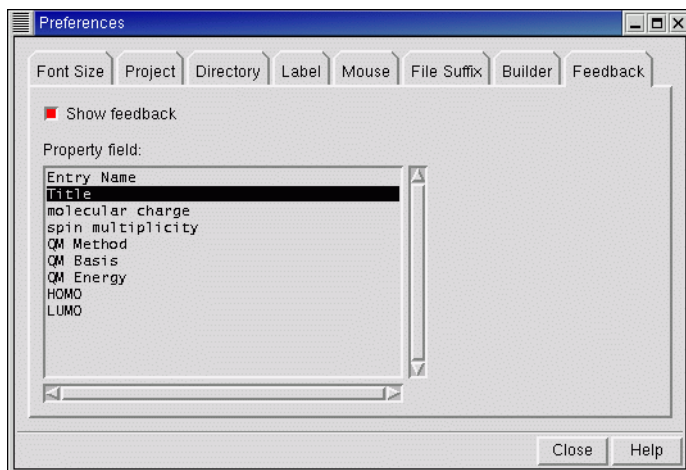


Figure 12.8. The Feedback tab of the Preferences panel.

12.2.8 Setting Feedback Preferences

In the Feedback tab, you can set the option to display information in the status bar when the pointer pauses over an atom in the Workspace. This tab provides a list of properties, from which you can choose one property for the property field of the atom information. See [Section 2.5 on page 15](#) for more information about what is displayed.

12.3 Changing the Workspace Appearance

You can adjust the appearance of the Workspace to fit the way you work best. You can hide or display many of the main window controls and change the background color of the main window.

12.3.1 Showing or Hiding Controls and Displays

The following five items in the Display menu hide or show main window controls and displays:

- Toolbar
- Clipping Planes Window
- Command Input Area
- Status Bar
- Sequence Viewer

You can also right-click on any of these controls or displays and choose Hide from the shortcut menu.

12.3.2 Setting the Background Color

The default background color of the main window is black, but the color can be changed to adjust the contrast between the background and the Workspace atoms. The background color can be retained for printing.

To change the background color, choose Display Options from the Display menu of the main window and select the Background tab. Choose a color in the By Palette subtab, or select a color by name in the By Name subtab. When a color is selected, a sample of the color is displayed in the Background color box.

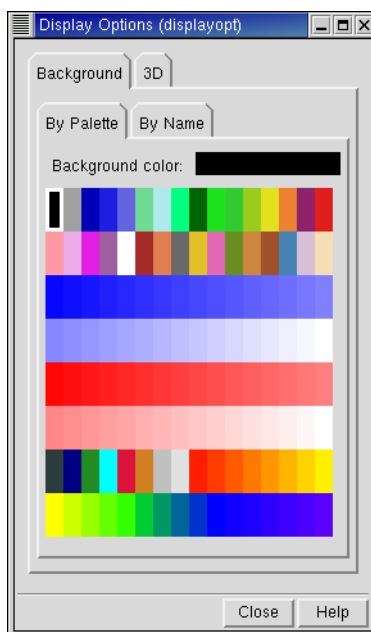


Figure 12.9. The Background *tab* of the Display Options *panel*.

12.4 Setting Up 3D Viewing

Maestro provides several tools to enhance viewing of complex structures. These tools include stereo viewing, perspective, fogging distant atoms, and depth cues and are described in the sections that follow.

12.4.1 Specifying a Stereo Viewing Method

Maestro's default stereo viewing method is wall-eyed. In this method of stereo viewing, the left eye looks at the left image, and the right eye looks at the right. The other option is cross-eyed viewing, in which the left eye looks at the right image and the right eye looks at the left.

Maestro also supports hardware stereo display. On both Linux and SGI platforms, stereo-in-a-window is supported, and full-screen stereo is supported on SGI platforms. For details, see the [Installation Guide](#).

To set the stereo viewing options, choose Display Options from the Display menu in the main window and select the Stereo subtab of the 3D tab. To select a stereo viewing method, choose Wall-eyed, Cross-eyed or Hardware from the Method option menu. To view the Workspace in stereo, select Display in stereo.

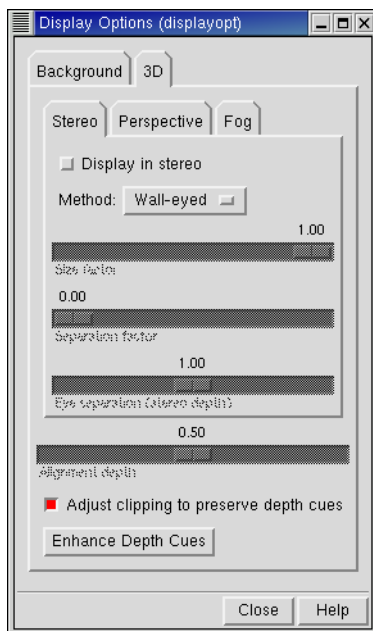


Figure 12.10. *The Stereo tab of the Display Options panel.*

If you are using the cross-eyed or wall-eyed method and the displayed stereo pairs are more than 2.5 inches from each other, adjust the Size factor, Separation factor, and Eye separation controls to bring the structures closer. This will make stereo viewing easier.

With any stereo method, you can adjust the eye separation and the alignment depth. Adjusting the eye separation changes the amount that the stereo image appears to extend into and out of the screen—the stereo depth—and adjusting the alignment depth makes the stereo image appear to move into or out of the screen.

12.4.2 Enabling Perspective

Maestro's Perspective feature enhances 3D graphical representation by scaling items in the Workspace in proportion to their distance from the viewing plane ("camera").

The perspective controls are in the Perspective subtab of the 3D tab. You can adjust the scaling using the Perspective scale slider. The green lines in the Top View window move to reflect the viewing volume available at a given value on the Perspective scale.

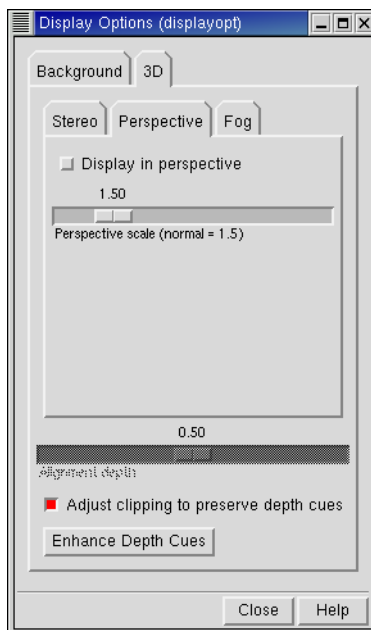


Figure 12.11. The Perspective *tab of the* Display Options *panel.*

12.4.3 Enabling Fog

Maestro's Fog feature causes atom colors in the displayed structure to fade toward that of the background as atom coordinates grow increasingly negative in the *z* direction. This feature is particularly useful for depth-cueing when viewing large structures, such as proteins, and when using the Wire molecular representation.

The fog controls are in the Fog subtab of the 3D tab. Fogging is turned off by default. To enable fogging, select On or Automatic from the Display fog option menu. Automatic fogging means that fogging is turned on when a molecule is bigger than a given number of atoms. Below this threshold, fogging is turned off. The threshold is 40 atoms, and can be set with the Maestro command `displayopt fogcutoff=number`. You can independently apply fogging to atom labels using the Display fog on atom labels option. You can also turn fogging on or off with the Set fog display state button menu in the main toolbar.



The default method for fogging employs linear fall-off for the brightness of the displayed structural elements. To exaggerate the fogging effect, you can select either the Exponential

method or Exponential squared method. When you are displaying a structure in perspective, you should use linear fall-off.

The points at which fogging begins and ends can be controlled with the Linear start and Linear end sliders. Adjusting Linear start moves the boundary at which fogging starts between the front (0.0) and the back (1.0) of the displayed structure. Adjusting Linear end moves the fog boundary that is behind the structure further back, to a maximum of twice the depth of the displayed portion of the structure. By adjusting these boundaries you can control the amount by which the structure appears to fade into the background.

For exponential and exponential squared fogging methods, the Density slider allows you to change the contrast of the colors used in fogging. Higher density makes closer portions of the structure appear brighter and more distant ones fade more quickly into the background.

Fog options are preserved between Maestro sessions.

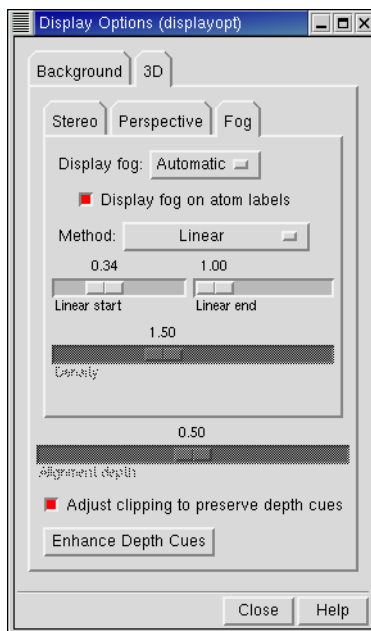


Figure 12.12. The Fog tab of the Display Options panel.

12.4.4 Depth Cues

When you are viewing in stereo or in perspective and you zoom in on a structure, the structure can appear to “flatten”. This happens because the entire depth of the structure remains in the field of view when you zoom in: it is not clipped at all, and the depth cues (eye separation or perspective scale) are not preserved. You can preserve depth cues by selecting the Adjust clipping to preserve depth cues option. When you do this, the front part of the structure is clipped as you zoom in, so that you see only the parts of the structure that are “in front of you”, not the parts that are “behind you”.

The depth cues are affected by the location of the clipping planes: the closer they are to the structure, the better the depth cues. To move the clipping planes in to fit the displayed structure, click Enhance Depth Cues in the 3D tab, or click the Enhance depth cues button on the main toolbar.



Scripting in Maestro

Maestro provides several facilities for creating scripts. You can generate scripts using Python or the Macros panel and install them on the Scripts menu. You can also create scripts and assign them to one of the 12 function keys on your keyboard.

13.1 The Scripts Menu

The Scripts menu in the main window is configurable and allows you to install, manage, and run Python scripts or named macros directly from Maestro.

To install a Python script:

1. Choose Manage from the Scripts menu.

The Manage Scripts panel opens (see [Figure 13.1](#)).

2. Click Install Scripts.

The Install Scripts panel opens (see [Figure 13.2](#)). You can install the scripts that came with the distribution or browse to a directory containing your own scripts.

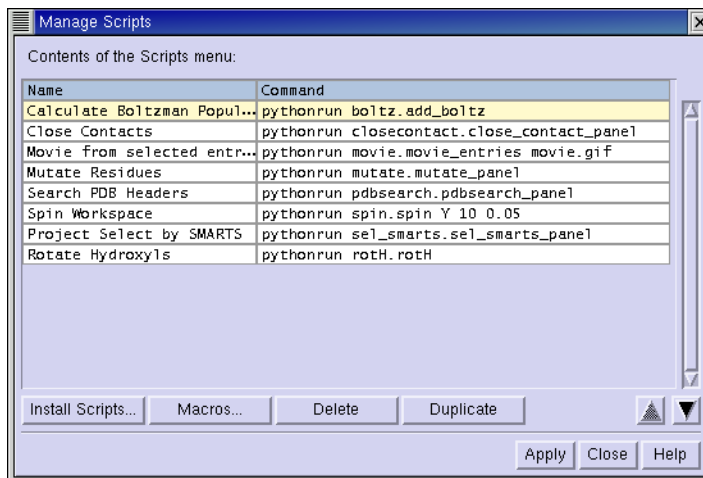


Figure 13.1. The Manage Scripts panel.

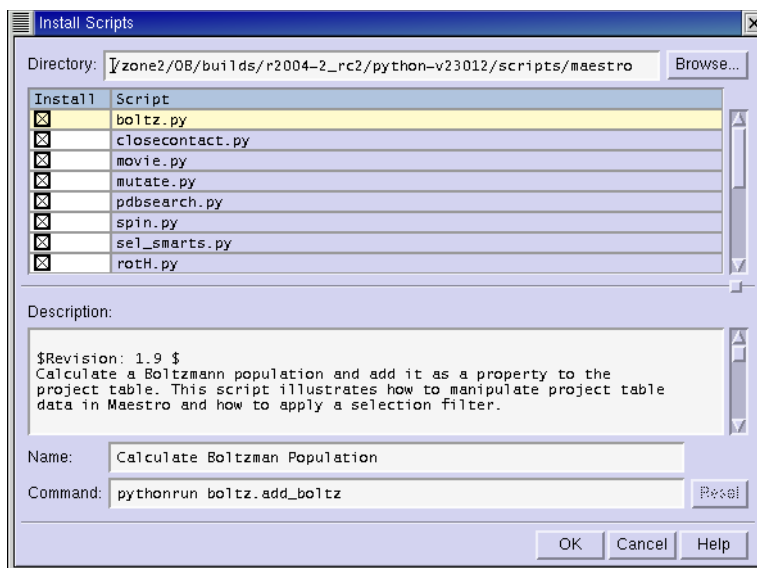


Figure 13.2. The Install Scripts panel.

3. Click in the Install column to select a script.

The script's name, the command that runs it, and a description (if available) appear in the text boxes at the bottom of the panel. The script name is what appears in the Scripts menu.

4. Click OK.

The script is added to the Scripts menu and the panel closes. If the script is already on the menu, you are asked if you want to overwrite it.

Note: The Install Scripts panel only displays files with a .py suffix, the standard suffix for Python scripts.

To install a macro:

1. Choose Manage from the Scripts menu.

The Manage Scripts panel opens (see [Figure 13.1](#)).

2. Click Macros.

The Macros panel opens.

3. Select a macro from the list.
4. Click Install.

To change the name or command of a script or macro, edit these fields directly in the Contents of the Scripts menu table.

You can delete or duplicate scripts and macros by using the Delete and Duplicate buttons. You can change the order in which the items appear in the Scripts menu by selecting a script and clicking the up and down arrow buttons.

To create submenus in the Scripts menu, edit the script names and insert a colon to separate the submenu name from the menu item. For example:

```
Color Atoms:Green
```

creates the submenu Color Atoms with the menu item Green.

To see the changes you have made to script names and locations, click Apply, then click on the Scripts menu. When you are satisfied with the changes, click Close to close the Manage Scripts panel.

13.2 Customization Using Python Scripts

Python is a full-featured scripting language that has been embedded in Maestro to extend scripting facilities. The Python capabilities within Maestro include access to Maestro functionality for dealing with chemical structures, projects, and Maestro files.

The two main Python commands used in Maestro are:

- `pythonrun`—executes a Python module. You can also use the alias `pyrun`. The syntax is:

```
pythonrun module.function
```

- `pythonimport`—rereads a Python file so that the next time you use the `pythonrun` command, it uses the updated version of the module. (You can also use the alias `pyimp`.)

The distribution includes a number of Python scripts designed for interaction with the structure in the Workspace. You can use these scripts, as well as create your own. With the Maestro Scripts menu, you can install, manage, and run your Python scripts directly from Maestro. For more information on the Scripts menu, see [Section 13.1 on page 167](#).

For more information on using Python with Maestro, see [Scripting with Python](#).

13.3 Customization Using Command Scripts

Some Maestro interface customization options can be set for the present session, but are not automatically retained for future sessions. To use these settings each time you launch a new Maestro session, you can create a startup script that will make the settings for you.

You can create startup scripts even if you are not familiar with the Maestro Command Language because most of the commands that you may want to include in a startup script can be extracted from the Command History log in the Command Script Editor panel.

To create a startup script, use the controls on the Command Script Editor panel, and then save the script in your preferences directory under the file name `maestro.cmd`. Each time you subsequently launch Maestro, this script is used to override the system-wide startup settings with your personal preferences.

To open the Command Script Editor panel, choose Command Script Editor from the Edit menu.

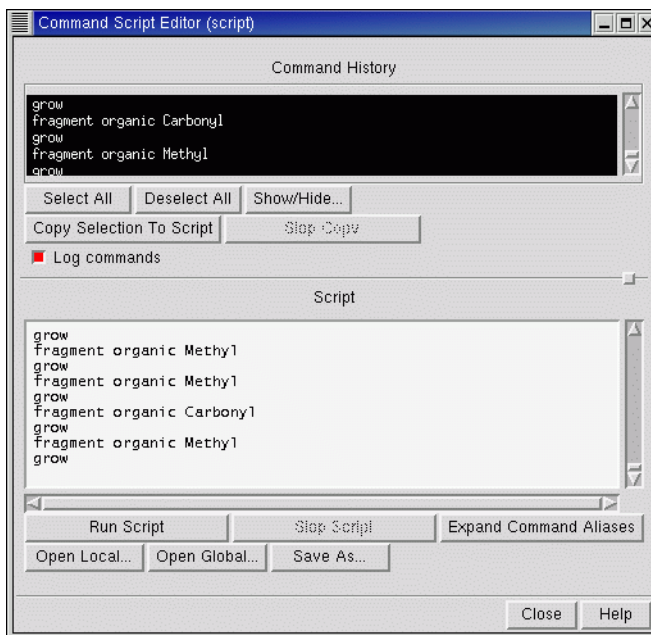


Figure 13.3. The Command Script Editor *panel*.

13.3.1 Command History

By default, most of the internal Maestro commands generated when you perform actions are logged in the Command History list in the Command Script Editor panel (see [Figure 13.3](#)). This written record of your commands is provided mainly to make the scripting process simpler. Since you can perform an operation, then cut and paste the resulting command into the Script list and save the script, it is not necessary to memorize the Maestro Command Language.

You can control which types of commands are displayed by using the Show/Hide Commands panel (see [Figure 13.4](#)). To open this panel, click Show/Hide in the Command Script Editor panel. If a command is not logged, open the relevant tab and select the button corresponding to the command, click Rebuild Command History, and then click Close. The command of interest should now appear in the list.

13.3.2 Building a Command Script

To build a command script from commands logged in the Command History list, you need to select and copy the commands of interest. Individual commands can be selected or deselected by clicking. A group of commands can be selected by dragging. To select all of the commands in the list, click Select All. To deselect all the commands in the list, click Deselect All.

Once you have selected the desired commands, copy them to the Script text area by clicking Copy Selection To Script. Copying many commands may take some time. To interrupt the copying process, click Stop Copy. This button is only active during a copy operation.

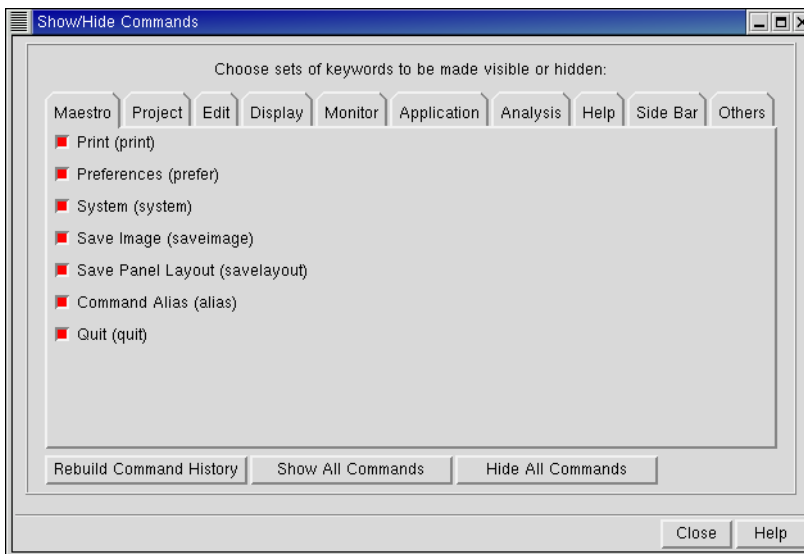


Figure 13.4. The Show/Hide Commands *panel*.

The commands are copied to the location of the cursor, which is placed at the end of the text after a copy operation. Therefore, repeated copy operations append the new command selection to the commands in the Script text area. If you want to create a new script, you must first delete the commands from the Script text area.

The commands in the Command History list remain selected after the copy operation, and must be explicitly deselected.

You can also directly enter commands in the Script text area.

13.3.3 Saving a Command Script

To save the command script displayed in the Script text area, click **Save As** and navigate to the directory into which you want to place the file. Then type the name of the file to be created at the end of the line of text in the Selection box, and click **Write**. You must include the `.cmd` extension when specifying the file name.

13.3.4 Opening Existing Command Scripts

The Command Script Editor panel provides two options for locating existing command scripts:

- Click **Open Local** to browse in your working directory.
- Click **Open Global** to browse in the `maestro-vversion/data/res` subdirectory of the Schrödinger software installation.

You can then navigate to the desired directory, select a command script, and click **Open** to load the script into the Script text area. If there are any commands currently in the Script text area, you are asked whether you want to replace them with the script you are about to open.

13.3.5 Running and Stopping Scripts

After you load a script into the Script text area, you can execute it by clicking **Run Script**. If **Log commands** is selected, the script commands are added to the Command History list, even if the command script doesn't perform any operations that are visible in the Workspace.

To halt the execution of a script, click **Stop Script**.

13.3.6 Customizations You Can Perform With a Script

As mentioned above, some customization settings in Maestro are not persistent across sessions. You can configure Maestro on startup according to your personal preferences by performing any of the desired operations described below, making a script containing the resulting commands, and then saving the script as `maestro.cmd` in the preferences directory.

Before you create scripts to perform the following customizations, open the Command Script Editor panel, ensure that Log commands is selected and that all commands are shown.

If you already have a `maestro.cmd` file in your preferences directory and you want to save the commands in this file, open it in the Command Script Editor panel before you perform the new customizations. If you no longer want to override the system-wide settings with your personal startup preferences, delete your `maestro.cmd` file, or move it out of this directory.

After performing each customization, select the relevant commands in the Command Script Editor panel and copy the resulting commands into the Script text area. When you have completed all customizations, save the script.

To set the Workspace background color:

1. Choose Display Options from the Display menu of the main window.
2. Click the Background tab.
3. Choose a color in the By Palette tab, or select a color name in the By Name tab.

To specify a stereo viewing method:

1. Choose Display Options from the Display menu of the main window.
2. In the 3D tab, click the Stereo tab.
3. Choose a stereo display technique from the Method option menu.
4. Adjust the Size Factor and Scale Factor if necessary.

To enable perspective:

1. Choose Display Options from the Display menu of the main window.
2. In the 3D tab, click the Perspective tab.
3. Adjust the Perspective Scale, if needed.

To hide the toolbar:

- Choose Toolbar from the Display menu of the main window.

Alternatively, you can type the following command into the Scripts text area, or add it directly to your `maestro.cmd` file.

```
hidepanel toolbar
```

To hide the command input area:

- Choose Command Input Area from the Display menu of the main window.

Alternatively, you can type the following command into the Scripts text area, or add it directly to your `maestro.cmd` file.

```
hidepanel cominput
```

To hide the clipping planes window:

- Choose Clipping Planes Window from the Display menu of the main window.

Alternatively, you can type the following command into the Scripts text area, or add it directly to your `maestro.cmd` file.

```
hidepanel topview
```

13.4 Creating Command Aliases

If you choose to use command scripts to perform operations in Maestro, you may find that there are many that you use frequently. The controls on the Command Alias panel allow you to assign shorter names to these commands. For example, you might want to define the alias label for the `labelatom` command. You can then use these aliases in commands typed in the command input area.

To open the Command Aliases panel, choose Command Aliases from the Maestro menu.

Defining a New Alias

To enter a new command alias, click New, type the abbreviation that you wish to use for the command in the Alias name dialog box, and click OK. In the Command Aliases panel, enter the command in the Alias definition box for the abbreviation you entered in the previous step.

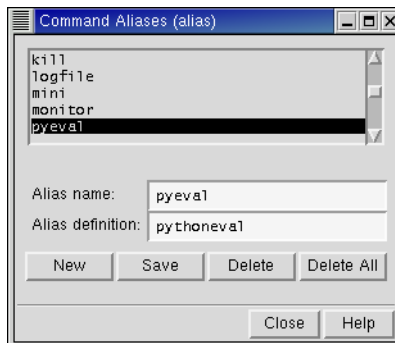


Figure 13.5. The Command Aliases panel.

Saving Aliases

To save your aliases, click **Save**. The aliases are written to the file `alias.cmd` in your preferences directory. Once the aliases have been saved, they can be used in any subsequent Maestro session. If you wish to use an alias only for the current session, it is not necessary to save the alias.

Changing an Existing Alias

Maestro is shipped with a set of command aliases already defined, and these can be edited or deleted. To change an alias, select the alias in the alias list, then type a new alias name or definition in the appropriate box.

To delete an alias from the list, select it and click **Delete**. To delete all currently defined aliases, click **Delete All**.

13.5 Creating Macros

Maestro allows you to create two kinds of macros: named macros and function key macros.

13.5.1 Named Macros

In the Macros panel, you can create a macro, name it, and run the macro from the command input area. To open the Macros panel, choose **Macros** from the **Edit** menu.

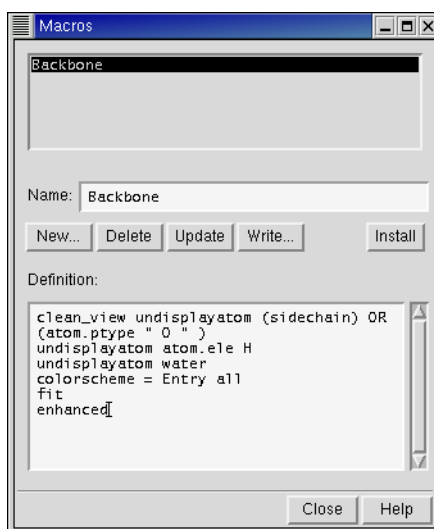


Figure 13.6. The Macros panel.

To create a named macro, click **New**, enter a name for the macro and click **OK**. In the **Definition** text box, you can either type the commands for the macro or copy and paste commands from the **Command Script Editor** panel. When you have finished entering commands, click **Update** to update the macro definition.

Note: Macros are automatically saved when you exit Maestro. However, it is advisable to update the macro whenever you change the definition. This makes the updated macro available during the Maestro session.

You can also add the macro to the **Scripts** menu by clicking **Add to Scripts Menu**. For more information on editing the **Scripts** menu, see [Section 13.1 on page 167](#).

To run the macro, enter the following in the command input area in the main window:

```
macrorun macro_name
```

If the command input area is not visible, choose **Command Input Area** from the **Display** menu.

13.5.2 Function Key Macros

The **Function Key Macros** feature allows you to assign a script command, or series of commands, to a particular keyboard function key (F1-F12).

To open the **Function Key Macros** panel, choose **Function Key Macros** from the **Edit** menu.

To assign a macro to a function key, select the key from the **Macro Key** option menu, and type the command or command sequence into the text area. You can copy and paste commands from the **Command Script Editor** panel into the text area of the **Function Key Macros** panel.

To test your macro, click **Run**.

To save a new macro for use in subsequent Maestro sessions, click **Save**. Saved macros are placed in the `maestro.mac` file in your preferences directory. To delete a macro, edit this file.

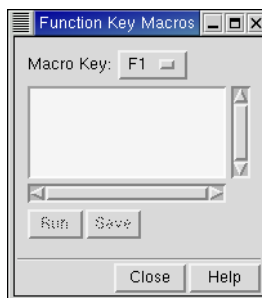


Figure 13.7. The **Function Key Macros** panel.

Printing and Saving Workspace Images

In addition to creating Maestro, MacroModel, PDB, Sybyl Mol2, and MDL SD structure files, Maestro can produce high-quality print and electronic representations of structures depicted in the Workspace. You can use the tools in the Print panel to generate PostScript and Encapsulated PostScript files, and use the Save Image panel to create TIFF and JPEG images. For information on writing (exporting) structure files, see [Chapter 3](#).

14.1 Printing an Image (File or Printer)

From the Print panel, you can print or create PostScript or Encapsulated PostScript images of structures currently rendered in the Workspace. To open the Print panel, choose Print from the Maestro menu. The settings available in the Print panel are described below.

Print To

To print an image, select Printer. To write a PostScript (PS) or encapsulated PostScript (EPS) file, select File.

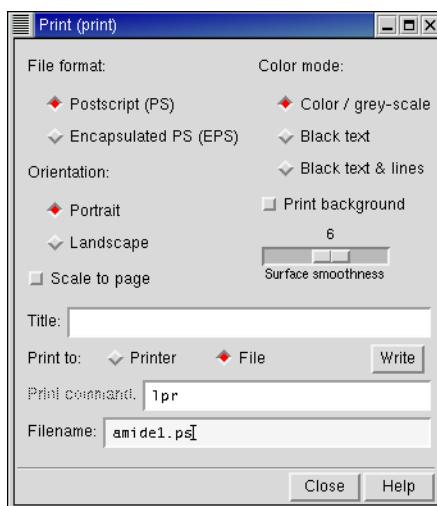


Figure 14.1. The Print panel.

File Format

If you selected Printer in the Print To section, only the Postscript (PS) option is available, because you cannot directly print encapsulated postscript. If you selected File, select Postscript (PS) or Encapsulated PS (EPS).

Orientation

Select Portrait or Landscape to determine the orientation of the image.

Scale To Page

Select Scale to Page to reduce or enlarge the image to occupy the available space on a standard page. When Scale to Page is not selected, the image is printed roughly the same size as it appears on screen.

Color Mode

The Color Mode selections are Color/Grey-scale, Black Text, and Black Text & Lines. Selecting the Color/Grey-scale option allows you to print or save images exactly as they appear on screen. If atom labels and other text appears in color, they are printed in color. It can be difficult to read atom labels and other text in color. To force atom labels and any other text to print in black, select Black Text. To force text plus any lines, such as markers and wire frame bonds, to print in black, select Black Text & Lines. If you print to a black-and-white printer, the images are printed in grey scale instead of color.

Print Background

This button controls the printing of the background color of the Workspace. It is not selected by default because printing a background requires more ink (or toner) and time, and the default background color of the Workspace is black.

Surface Smoothness

Increasing the surface smoothness for an image helps to remove “blocky” edges on curved surfaces in your structure. However, there is a trade-off between degree of smoothing and amount of print time required. Particularly with small curved surfaces, the minute improvement in image clarity obtained from increasing the smoothness factor may not justify the increased time it takes to print the image.

To change the surface smoothing, drag the Surface Smoothness slider.

Title

You can enter a title in the Title text box. The title is printed in black text in the lower left corner of the page.

Print Command

If the UNIX shell you use requires a different print command from the default given in the Print Command text box, change the command in the text box to the required command.

When printing an image, Maestro creates a temporary PostScript file then prints this file using the specified command. If your printer requires the `lp` command instead of `lpr`, use the `-c` option to ensure that the temporary copy of the file to be printed is preserved until printing is complete.

Filename

If you are saving an image to a file, enter the file name in the Filename text box, including either the `.ps` or the `.eps` file extension.

14.2 Creating TIFF and JPEG Image Files

You can save electronic images of structures rendered in the Workspace in TIFF or JPEG format using the Save Image panel. To open the Save Image panel, choose Save Image from the Maestro menu.

When are you are ready to save an image, make sure the Workspace is unobstructed:

- Undisplay the clipping planes window by right-clicking it and choosing Hide or choose Clipping Planes from the Display menu.
- Make sure there are no panels or dialog boxes overlapping the Workspace. This includes the Save Image panel.
- Make sure the cursor is not in the Workspace.

The Save Image panel allows you to set the following options:

- Format—choose Tiff or JPEG:

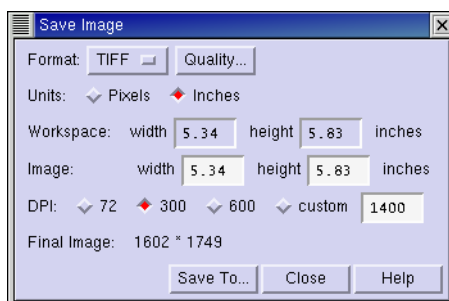


Figure 14.2. The Save Image panel.

- TIFF is an acronym for “Tag(ged) Image File Format”, a file structure that allows pixel data to be stored in a variety of different internal formats. Format parameters such as byte order, image dimensions, and compression method are specified in labeled variables, or “Tags”, which can either be placed within the file header or distributed throughout the file. When Maestro saves an image as a TIFF file, it uses no compression. This means that the image is “lossless” (all of the information in the original image is retained), but the file can be quite large. The TIFF format is recommended for saving high quality Maestro images.
- JPEG stands for “Joint Photographic Experts Group”. JPEG is actually a compression method, not an image format. The JPEG file format specifies a means of compressing natural images, that does not severely impair the perceived image “quality”. The JPEG method is almost always “lossy”, meaning that the reconstituted pixel data differs from that of the original.

For natural images, these differences are not generally noticeable to a human observer. However, in general, the JPEG method is not recommended for saving synthetic images like those Maestro produces, because there may be noticeable distortions in the saved image, such as subtle color changes. This option is, nevertheless, provided because it allows Maestro users to create extremely compact screen shots for sending with an e-mail message or inserting in web pages. Maestro uses the JPEG-based JPEG-File Interchange Format (JFIF), the most commonly used JPEG file format.

Both TIFF and JPEG files can be imported into image manipulation programs.

- **Quality**—open the Image Quality dialog box, which allows you to adjust the quality of the image:
 - **Smooth**—Select this option to turn on scene anti-aliasing. This creates smoother edges, but takes longer to save the image.
 - **Quality**—(for JPEG only) Use the slider to select a quality level from 1 to 100. The default is 75. The higher the quality level, the longer it takes to save the image. For large resolution images, this can take awhile.
- **Units**—Choose pixels or inches. This choice determines the Workspace and Image units (see below).
- **Workspace**—Displays the measurement of the Workspace, not the whole main window. These values can not be edited directly. However, if you resize the Workspace, the values are updated automatically. To change the units, select either Pixels or Inches under Units.

- **Image**—Enter the desired width and height of the image. These proportions are constrained, so if you change one parameter, the other adjusts automatically. To change the units, select either **Pixels** or **Inches** under **Units**.
- **DPI**—Choose a resolution for the image or enter a specific value. If you choose **Pixels** for the units, this option is unavailable and the specified pixels are used.
- **Final Image**—Displays the final resolution of the image (in pixels). This value is not editable, but you can change it by changing the options above.

To save an image:

1. Choose **TIFF** or **JPEG** from the **Format** option menu.
2. If you want to increase the quality of the image:
 - a. Click **Quality** to choose **Smooth**.
 - b. For **JPEG** files, use the slider to select the quality level.
3. Select the units you want to use: **Pixels** or **Inches**.
4. Adjust the image size (width and height).

The Workspace size can only be adjusted by actually resizing the main Maestro window.

5. Select a resolution under **DPI**. (If you chose **Pixels** above, this option is unavailable.)

The final image resolution is displayed in pixels near the bottom of the panel.

6. Click **Save To**, navigate to a directory, and select or enter a file name.
7. Click **Save**.

Note: When saving high resolution images in wire or Ball & Stick molecular representation, you may need to adjust the wire width so the bonds remain visible at the higher resolution. This may take some experimentation, depending on the purpose of the image. See [Section 6.2.3 on page 74](#) for information on changing the molecular representation attributes.

Help

This chapter describes the various ways in which you can obtain information on Schrödinger's products. Maestro has an online help facility, automatic context-sensitive help (Auto-Help), and Balloon help. Manuals in PDF format are included with the product distribution, and each product has an FAQ page on the Schrödinger web site.

15.1 Maestro Online Help

If you have any difficulty performing a particular operation within the Maestro interface, look for answers first in the online help, which is presented in the Help panel.

The Help panel contains two tabs for locating help topics, the Search tab and the Categories tab. Below these tabs, the topic name appears in the Topic text box, and the help topic is displayed in the text area below the topic title.

Most topics contains a list of contents at the top that you can click to go to the relevant section, and a list of related topics at the end that you can click to select a new topic.

For general information, open the Help panel and start a search. To open the Help panel, choose Help from the Help menu, or press CTRL+H in the Workspace.

For help on a particular panel's features, click the Help button in the lower right corner of the panel. The Help panel opens with the relevant topic displayed.

15.1.1 Locating a Topic by Searching

The Search tab allows you to search for a topic using a search string. Type the string in the Search For text box and press the RETURN key or click the Search button. If there are any topics in the help file that contain the search string you entered, they are listed in the Topics which match search list. The search results are sorted in alphabetical order. By default, the last entry in the list is highlighted. This does not mean that the highlighted topic most closely matches the query word. To view a topic, double-click on its name in the list, or click on the name once and then press the RETURN key.

You can restrict or enhance your search by selecting the following options:

- **Title Only:** The search is restricted to the topic titles.
- **Case Sensitive:** The case is matched exactly in the search.
- **Regular Expression:** Allows the use of UNIX regular expressions in the search.

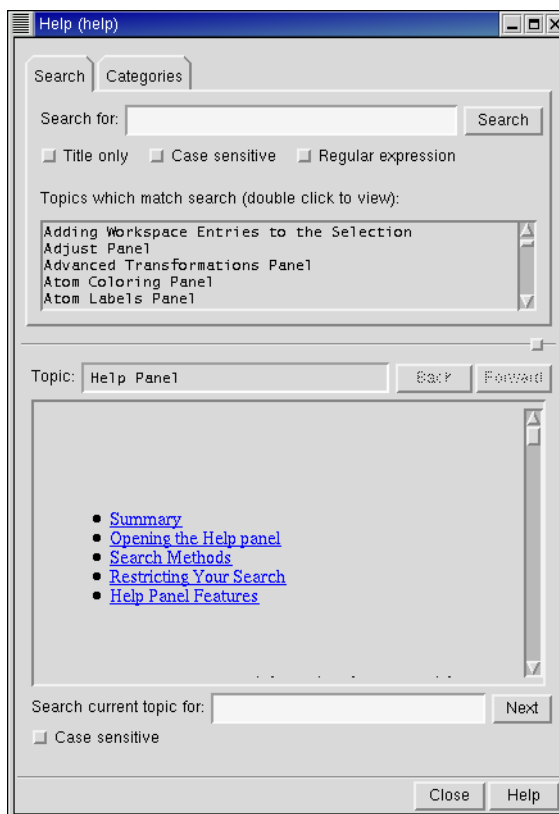


Figure 15.1. The Help panel showing the Search tab.

15.1.2 Locating a Topic by Category

Searching for help by category is helpful when you do not know which keywords may produce the desired topic. Select the option from the Help Categories menu that corresponds to the group of topics you would like to browse. A list of relevant topics is displayed in the Topics in this category list. The topics in the list are arranged in alphabetical order and the last entry in the list is highlighted. To select a topic, double-click its name, or click its name and press ENTER.

15.1.3 Searching Within a Topic

If you want to locate a string in the topic text, type the string in the Search current topic for text box and press ENTER. To locate the next occurrence of the string, press ENTER again or click Next. You can match the exact case of the search string by selecting Case sensitive.

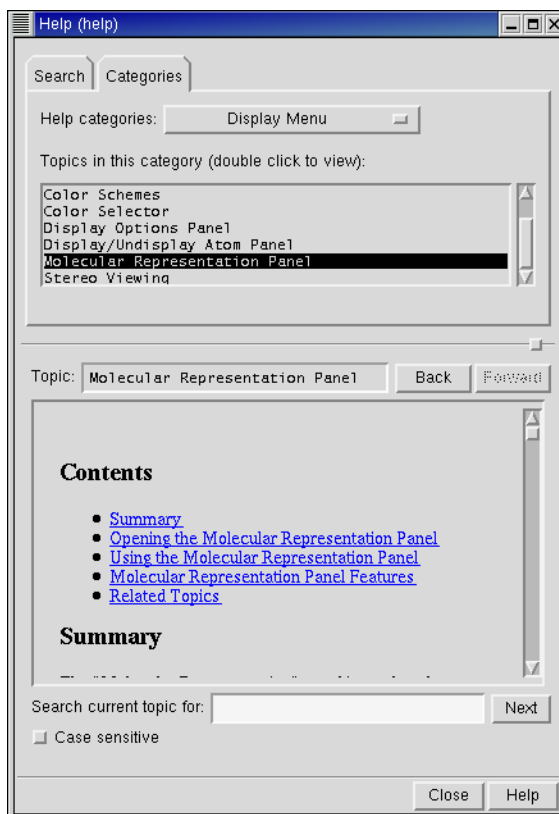


Figure 15.2. The Categories *tab* of the Help *panel*.

15.2 Context-Sensitive Help

Maestro has two other help features, which are described below.

Auto-Help

The Auto-Help text area is located at the foot of the main window. The Auto-Help facility provides automatic context-sensitive help. For example, when you open the Atom Labels panel, Auto-Help prompts you to pick an atom to label. Auto-Help is not available for all operations.

Balloon Help

Balloon Help (tooltips) consists of a floating text box containing context-sensitive help that appears when the pointer is allowed to rest on an object. The balloon disappears when you

move the mouse. You can temporarily disable Balloon Help using the Show Balloon Help item on the Help menu. Balloon Help is not available for all interface components.

15.3 Manuals

If you have difficulty installing Maestro, see the *Installation Guide*. If you have questions about a computational program that you have purchased, see the *User Manual* or *Reference Manual* for that program. Manuals are available in PDF format on the program CD and the Schrödinger web site, and are installed in `$SCHRODINGER/docs`. You can browse to the manuals by opening `$SCHRODINGER/docs/index.html` in a web browser or by opening `$SCHRODINGER/schrodinger_documentation.pdf` in a PDF reader. You can also open the manuals by choosing Manuals from the Help menu in Maestro.

15.4 The Schrödinger Web Site

The Schrödinger web site, <http://www.schrodinger.com>, contains information about all our products. In the [Support Center](#), you can find the manuals in PDF format and a list of frequently asked questions for each of the products, including known issues. The manuals on the web site are the most up-to-date. The manuals page contains links to corrections and additions to the documentation. The web site also contains information about product releases and software updates. Consult the web site for the most recent information on a product.

15.5 Technical Support

If you have questions that have not been answered from any of the above sources, contact Schrödinger using the information below. Generally, e-mail correspondence is best because you can send relevant machine output.

Schrödinger

E-mail: help@schrodinger.com

USPS: 101 SW Main Street, Suite 1300, Portland, OR 97204

Phone: (503) 299-1150

Fax: (503) 299-4532

WWW: <http://www.schrodinger.com>

FTP: <ftp://ftp.schrodinger.com>

When sending e-mail messages, please include the following:

- All relevant user input and machine output
- Software purchaser (company, research institution, or individual)
- Primary software user

- Computer platform type
- Operating system with version number
- Software version number
- Maestro version number
- mmshare version number

Much of the machine and system information listed above can be produced by entering the following command at a shell prompt:

```
$SCHRODINGER/machid
```


The Maestro File Format

All Schrödinger products use the Maestro file format as their primary method of storing molecular structure information. The Maestro file format is extensible. That is, users and third-party programs can add additional fields that will be accepted and retained when the files are read by Schrödinger products. Maestro files are also self-describing; data items are named rather than simply identified by position, and the data type (string, real number, integer and Boolean) is included as a single-digit code letter at the beginning of each data name.

Since Maestro files can accept additional user-defined fields, it is not possible to give an all-inclusive example of a Maestro file. Instead, this section will give a general overview of the Maestro file format and describe the components used for most Schrödinger applications.

A.1 Basic File Description

Maestro format files are free-format ASCII text. All white space (spaces, tab characters, end of lines etc.) is ignored. Schrödinger applications will write files in a consistent layout in order to enhance human readability but the files are not, in fact, required to follow this layout. Except for the requirement that string type data items which contain blanks be surrounded by double quotes, there are no restrictions, such as width or precision of the numeric fields, on the data items.

A.2 Data Blocks

The basic unit of a Maestro file is a block, a series of data items contained within a pair of curly braces { }. For instance, the following is the first block of the example file in shown in this appendix:

```
{  
    s_m_m2io_version  
    :::  
    1.0.0  
}
```

Most blocks in the Maestro file are preceded by a name, for example, `f_m_ct`, where ‘f’ stands for ‘full’ (see below), ‘m’ generally stands for ‘Maestro’ because the data originated from Maestro, and ‘ct’ stands for ‘CT’—short for ‘Connection Table’, which can generally be thought of as a collection of atoms.

Blocks can be nested within other blocks, and if a block contains a list of data values, the block may be indexed. That is, the number of data values contained in the block can be appended to the name. Indexed blocks are identified by having names of the form: 'name[<number of items>]'. For example, the file shown in [Section A.5](#) contains an indexed block of the name `m_atom[18]`. In this block name, 'm' means that the data originated in Maestro, 'atom' means that the data values included in the block belong to atoms, and '18' indicates that there are 18 sets of data in the indexed block (representing 18 atoms).

Within each block there are two sections: a list of the names corresponding to the data fields included in each data item, and the actual data items. In the example file, the first three data field names in the `m_atom[18]` block are:

```
i_m_mmod_type
r_m_x_coord
r_m_y_coord
```

The end of the first section is denoted by a separator ':::', and followed by the actual data items. In the example file, the first three data items are:

```
1 3 0.547623 1.262401 -0.990300 1 " " X " " 2
0.00000 0.00000 CHEX " " " " 6 0 0 1 ""

2 3 -0.930177 1.296701 -1.411700 1 " " X " " 2
0.00000 0.00000 CHEX " " " " 6 0 0 1 ""

3 3 -1.821477 1.772701 -0.253300 1 " " X " " 2
0.00000 0.00000 CHEX " " " c1" 6 0 0 1 ""
```

There must be the same number of data values as there are data names, and the fields within the data items must appear in the same order as the data names are listed. In the above example, for instance, the first field in the first line of data ('1') is the index number. The second number ('3') is the MacroModel atom type, indicated by the data name `i_m_mmod_type`. Index numbers are required in an indexed block and do not count against data names.

The end of the data item section is indicated by a second ':::' separator. A file may contain any number of structures (CT blocks), and each structure block may contain any number of atom and bond data items.

A.3 Compressed Format

Files may use a compressed format, most commonly when storing information for a set of conformers. In a compressed file, the first CT block must be a “full” block containing all the information described above. The name of this block is 'f_m_ct'. This full block can be followed by any number of “partial” CT blocks with the name 'p_m_ct'. The partial blocks contain only information which is different from that in the last preceding full CT block. For example, a file of conformers will have only one copy of the bond table (the m_bond[36] block in the example file), and it will be stored in the full CT. Subsequent partial CT blocks will read bond information from the preceding full CT block.

A.4 Data Item Names

The names of data items in a Maestro format file follow a convention which allows the type of the data item and the “owner” of the data to be identified. Names have the form 't_o_d', where t is the type descriptor, o is the “owner” and d is the actual data name. For instance, the first data name in the m_bond[36] block is 'i_m_mmod_type'. The first character of the name represents the data type. The acceptable types are:

i – integer	s – string value
r – real number	b – Boolean value

Owner values indicate the application from which the data item most likely originated. Including this field in the data name allows multiple applications to store identically titled quantities. For instance, two applications could store data fields named “energy.” Currently, basic geometrical and connectivity information, which is shared by most Schrödinger applications, have 'm' (for Maestro) owner fields. Data introduced by MacroModel have 'mmod' owner fields, data introduced by QikProp have 'qp' owner fields.

In the 'i_m_mmod_type' example, 'i' indicates the value in the corresponding data field is an integer. The letter 'm' in the owner field means that the data originated from Maestro, and 'mmod_type' is an abbreviation for “MacroModel atom type,” which is what the data field represents.

A.5 Example Maestro File

Below is extract from a Maestro format file. The complete file is not shown. Where data has been exempted, it has been noted with an ellipse (...). The text in *italic type* is provided as annotation and is not a part of the file.

The first block in the file is unnamed. This contains required information that is relevant to the whole file.

```
{  
  s_m2io_version  
  :::  
  1.0.0  
}
```

f_m_ct {	<i>The “CT” block. Each structure in the file is contained in such a block. This is a full block.</i>
s_m_title	<i>The only CT-level data name. A string value representing the title of this structure.</i>
:::	<i>The separator between the data names and data values.</i>
"Cyclohexane"	<i>The value of the “title” data item.</i>
m_atom[18] {	<i>The start of the atom block for this CT block. There are 18 atoms in this block.</i>
i_m_mmod_type	<i>The MacroModel atom type.</i>
r_m_x_coord	<i>The X-coordinate.</i>
r_m_y_coord	<i>The Y-coordinate.</i>
r_m_z_coord	<i>The Z-coordinate.</i>
i_m_residue_number	<i>The residue number.</i>
s_m_insertion_code	<i>The PDB insertion code.</i>
s_m_mmod_res	<i>The one-letter MacroModel residue code.</i>
s_m_chain_name	<i>The PDB chain name.</i>
i_m_color	<i>The color for this atom.</i>
r_m_charge1	<i>The partial atomic charge.</i>
r_m_charge2	<i>The partial atomic charge.</i>
s_m_pdb_residue_name	<i>The PDB residue name.</i>
s_m_pdb_atom_name	<i>The PDB atom name.</i>

s_m_grow_name	<i>The name used by the Maestro structure builder.</i>
i_m_atomic_number	<i>The atomic number.</i>
i_m_formal_charg	<i>The formal charge.</i>
i_m_representation	<i>The representation used to draw this atom.</i>
i_m_visibility	<i>A flag to indicate whether this atom is displayed in Maestro or not.</i>
s_m_atom_name	<i>The user-specified atom name.</i>
:::	<i>The separator for the end of the data items.</i>

Next follow the data values. The first column contains an index number that is assigned automatically. Remaining columns represent the data values in the same order as data names given above.

```

1  3  0.547623  1.262401  -0.990300  1 " " X " "  2  0.00000
0.00000 CHEX " " " "  6  0  0  1  ""

2  3  -0.930177  1.296701  -1.411700  1 " " X " "  2  0.00000
0.00000 CHEX " " " "  6  0  0  1  ""

3  3  -1.821477  1.772701  -0.253300  1 " " X " "  2  0.00000
0.00000 CHEX " " " c1"  6  0  0  1  ""

4  3  -1.624777  0.890901  0.990300  1 " " X " "  2  0.00000
0.00000 CHEX " " " "  6  0  0  1  ""

...

17 41 0.493223  -0.677999  0.001900  1 " " X " "  21  0.00000  0.00000
CHEX " " " n3"  1  0  0  1  ""

18 41 1.817023  0.395501  0.566000  1 " " X " "  21  0.00000
0.00000 CHEX " " " n2"  1  0  0  1  ""

:::
}

```

The separator for the end of the data values.

The end of the atoms block.

m_bond[36] {	<i>The bond block. There are 36 bonds.</i>
i_m_from	<i>The atom the bond is from.</i>
i_m_to	<i>The atom the bond is to.</i>
i_m_order	<i>The bond order.</i>
i_m_from_rep	<i>The graphical representation for the “from” half-bond.</i>

i_m_to_rep *The graphical representation for the “to” half-bond.*
::: *The end of the data names for bond data.*

What follows are the data items for the bond block. The first column contains an automatically assigned index number for each item.

```
1  1  2  1 1 1
2  1  6  1 1 1
...
35 17  6  1 1 1
36 18  6  1 1 1
:::
}
}
```

The end of the bond data items.
The end of the bond block.
The end of the CT block.

Atom Types

Maestro atom types are not “hard-wired” into the program, but are read in at run-time from a file called `atom.typ`. We continue to use the original MacroModel atom types for most applications, but now have a mechanism for easily adding new types. The list of explicitly supported atom types has been expanded to include:

Li^+ , Na^+ , K^+ , Rb^+ , Cs^+ , Ca^{2+} , Ba^{2+} , Mg^{2+} , Mn^{2+} , Mn^{3+} , Mn^{4+} , Mn^{5+} , Mn^{6+} , Mn^{7+} , Fe^{2+} , Fe^{3+} , Co^{2+} , Co^{3+} , Ni^{2+} , Ni^{3+} , Cu^+ , Cu^{2+} , Zn^{2+} , Mo^{3+} , Mo^{4+} , Mo^{5+} , Mo^{6+} , and Cl^- .

You can add your own atom types and supply parameters for these types in the force-field file. The default `atom.typ` file is the one located in `$SCHRODINGER/mmshare-vversion/data`. A local file of this name overrides the default, as does a local file whose name is *filename.typ*, where *filename* is the stem of the *filename.mae* file. This system of default version with local overrides is analogous to that used for force field and solvation files.

Currently Maestro accepts atom types up through number 300. We restrict the type numbers that we supply to the range 0 through 199. When adding your own types, you should begin with the number 200. Negative atom types are not accepted.

Inside Maestro you can change atoms in existing structures to new types, including user-defined types. However, no provision is made in the current version for the builder to obtain the geometric information it needs from the `atom.typ` file, which would require extensions to the format of this file. Also, the original atom types remain “hard-wired” in Maestro. Therefore, you should not attempt to redefine types which we supply. Instead, add new types. The `atom.typ` file itself contains comments describing its format and usage.

The atom types and their equivalents in other force fields are listed in [Table B.1](#) and [Table B.2](#). These atom types are defined in the `atom.typ` file. Note that 00 atoms are stored as type 0 in the atom connection table, and that all current force fields require explicit hydrogens on heteroatoms.

Table B.1. Atom types and equivalents.

No.	Symbol	Description	Equivalencies		
			MM2/ MM3	Charmm	Amber
1	C1	Carbon - sp			
2	C2	Carbon - sp ²	2,3	CT,C	C, C*, CA, CB, CC, CM, CF, CG, CN
3	C3	Carbon - sp ³	1	CT	CT
4	CA	United atom CH - sp ³		CH1E	CH
5	CB	United atom CH ₂ - sp ³		CH2E	C2
6	CC	United atom CH ₃ - sp ³		CH3E	C3
7	CD	United atom CH - sp ²		CR1E	CD, CE, CJ, CP
8	CE	United atom CH ₂ - sp ²			
9	CF	United atom CH - sp			
10	CM	Carbanion (C ⁻)			
11	CP	Carbocation (C ⁺)	30		
12	CR	Carbon free radical	29		
...					
14	C0	Any carbon			
15	O2	Oxygen - double bond	7	O	O, O2
16	O3	Oxygen - single bond	6	OH1, OH2	OS, OH
17	OA	United atom OH		OH1E	
18	OM	O ⁻ (alkoxide, carboxylate)	47		OC
19	OW	United atom H ₂ O		OH2E	
20	OP	Oxonium (sp ²) =[O ⁺]-			
21	OQ	Oxonium (sp ³) R ₃ O ⁺			
...					
23	O0	Any oxygen			
24	N1	Nitrogen - sp	10		

Table B.1. Atom types and equivalents. (Continued)

No.	Symbol	Description	Equivalencies		
			MM2/ MM3	Charmm	Amber
25	N2	Nitrogen - sp^2	9	N, NR, NP, NH1, NH2	N, NB, NC, N*, N2
26	N3	Nitrogen - sp^3	8	NH3	NT
27	NA	United atom NH - sp^3			
28	NB	United atom NH_2 - sp^3		NH1E	
29	NC	United atom NH - sp^2			
30	ND	United atom NH_2 - sp^2		NH2E	
31	N4	N^+ - sp^2		NC2	NA
32	N5	N^+ - sp^3	39		N3
33	NE	United atom NH^+ - sp^3			
34	NF	United atom NH_2^+ - sp^3			
35	NG	United atom NH_3^+ - sp^3		NH3E	
36	NH	United atom NH^+ - sp^2			
37	NI	United atom NH_2^+ - sp^2		NC2E	
38	NM	N^- - sp^3			
39	NP	N^- - sp^2			
40	N0	Any nitrogen			
41	H1	H-Electroneut(e.g., C,S)	5, 44	HA	HC, HS
42	H2	H-O (Neutral)	21, 24	H,HC	HO
43	H3	H-N (Neutral)	23, 28	H,HC	H, H2
44	H4	H-Cation	48 MM3 (none in MM2)		H3
45	H5	H-Anion			
...					
48	H0	Any hydrogen			
49	S1	Sulfur	15,17,18	S	S, SH
50	SA	United atom SH		SH1E	

Table B.1. Atom types and equivalents. (Continued)

No.	Symbol	Description	Equivalencies		
			MM2/ MM3	Charmm	Amber
51	SM	S ⁻ (thiolate anion)			
52	S0	Any Sulfur		S	
53	P0	Phosphorus	25		P
54	B2	Boron - sp ²	26		
55	B3	Boron - sp ³	27		
56	F0	Fluorine	11		
57	Cl	Chlorine	12		
58	Br	Bromine	13		
59	I0	Iodine	14		
60	Si	Silicon	19		
61	Du	Dummy atom for FEP			
62	Z0	Special atom to be defined			
63	Lp	Lone electron pair			LP
64	00	Any atom	0	*	
65	Li	Li ⁺			Li
66	Na	Na ⁺		SOD	Na
67	K0	K ⁺		POT	K0
68	Rb	Rb ⁺			Rb
69	Cs	Cs ⁺		CES	
70	Ca	Ca ⁺²		CAL	
71	Ba	Ba ⁺²			
72	Mg	Mg ⁺²		MG	
73	M2	Mn ⁺²			
74	M3	Mn ⁺³			
75	M4	Mn ⁺⁴			
76	M5	Mn ⁺⁵			
77	M6	Mn ⁺⁶			

Table B.1. Atom types and equivalents. (Continued)

No.	Symbol	Description	Equivalencies		
			MM2/ MM3	Charmm	Amber
78	M7	Mn ⁺⁷			
79	f2	Fe ⁺²			
80	f3	Fe ⁺³			
81	o2	Co ⁺²			
82	o3	Co ⁺³			
83	n2	Ni ⁺²			
84	n3	Ni ⁺³			
85	c1	Cu ⁺			
86	c2	Cu ⁺²			
87	Zn	Zn ⁺²			
88	m3	Mo ⁺³			
89	m4	Mo ⁺⁴			
90	m5	Mo ⁺⁵			
91	m6	Mo ⁺⁶			
...					
100	SP	S ⁺			
101	S2	Sulfur - sp ²			
102	Cm	Cl ⁻		CLA	
...					
150	PI	Ligand dummy atom			

Table B.2. Generalized atom types.

No.	Symbol	Description
151	GA	Isolated atom
152	GB	Linear-single coordinate
153	GC	Linear-two coordinate
154	GD	Trigonal-two coordinate
155	GE	Trigonal-three coordinate
156	GF	Tetrahedral-three coordinate
157	GG	Tetrahedral-four coordinate
158	GH	Trigonal bipyramid-three coordinate
159	GI	Trigonal bipyramid-four coordinate
160	GJ	Trigonal bipyramid-five coordinate
161	GK	Octahedral-four coordinate
162	GL	Octahedral-five coordinate
163	GM	Octahedral-six coordinate
164	GN	Pentagonal bipyramid-seven coordinate
165	GO	Twisted cube-eight coordinate
166	GP	Nine coordinate
167	GQ	Ten coordinate
168	GR	Eleven coordinate
169	GS	Icosahedron-twelve coordinate
170	GT	Thirteen coordinate
171	GU	Fourteen coordinate
172	GV	Fifteen coordinate
173	GW	Sixteen coordinate

Entry Suffixes

The entry names within a project must be unique. Maestro uses various suffixes to ensure that the names are unique. Numbers are added after the suffix to make the name unique. Maestro adds a suffix when certain operations are performed on entries. These suffixes are listed in [Table C.1](#).

Table C.1. Entry name suffixes.

Operation	Suffix	When Added
merge projects	-mrg	name conflict
import	-imp	always
job incorporation	-inc	always
create from Workspace	-wse	name conflict
rename entry	-re	name conflict, or multiple entries selected
combine entries	-com	name conflict
disassociate	-dis	always
duplicate	-dup	always

Utilities

Schrödinger software is distributed with a number of utility programs and scripts that can be used in conjunction with any of the computational programs. These UNIX utility programs are installed in `$SCHRODINGER/utilities`. UNIX conventions are used for the command syntax—see [page xiii](#) for more information.

D.1 Structure Conversion

This section describes the utilities available for conversion between various file formats for molecular structures. These utilities are used by Maestro for structure conversion when importing or exporting structures.

D.1.1 Conversions to and From MacroModel Format

maemmod

Convert files from Maestro format (`.mae`) to MacroModel format (`.dat`).

Syntax:

```
maemmod maestro-file mmod-file
```

mmodmae

Convert files from MacroModel format to Maestro format.

Syntax:

```
mmodmae mmod-file maestro-file
```

mmio_convert

Compress and uncompress MacroModel format files.

Syntax:

```
mmio_convert [infile [outfile]]
```

D.1.2 Conversions to and From Mol2 Format

mol2convert

Convert structure files between Maestro and Sybyl Mol2 format.

Syntax:

```
mol2convert [-noarom] {-imae|-imol2} inputfile
             {-omae|-omol2} outputfile
```

-imae *inputfile* Input file in Maestro format
-imol2 *inputfile* Input file in Mol2 format
-omae *outputfile* Output file in Maestro format
-omol2 *outputfile* Output file in Mol2 format

The input and output files must be specified. The `-noarom` option suppresses conversion of aromatic type 4 bonds to single and double bonds according to Maestro and OPLS conventions.

Example:

```
mol2convert -imae myfile.mae -omol2 myfile.mol2
```

mmodmol

Convert files from MacroModel format to Sybyl Mol2 format.

Syntax:

```
mmodmol [input-file] [output-file]
```

If the input file is not specified, input is taken from standard input. Likewise, if the output file is not specified, output is sent to standard output.

Example:

```
mmodmol myfile.dat > myfile.mol2
```

molmmod

Convert files from Sybyl Mol2 format to MacroModel format.

Syntax:

```
molmmod [input-file] [output-file]
```

If the input file is not specified, input is taken from standard input. Likewise, if the output file is not specified, output is sent to standard output.

Example:

```
cat myfile.mol2 | molmmmod > myfile.dat
```

D.1.3 Conversions to and From PDB Format: pdbconvert

This section describes the command-line utility version of the program that converts files between PDB, MacroModel, and Maestro formats. For information on PDB conversion within the Maestro GUI, see [Section 3.1.3 on page 27](#).

Syntax:

```
pdbconvert [options] {-ipdb|-imm|-imae} inputfile
               {-opdb|-omm|-omae} outputfile
```

-ipdb <i>inputfile</i>	Input file in PDB format
-imm <i>inputfile</i>	Input file in MacroModel format
-imae <i>inputfile</i>	Input file in Maestro format
-opdb <i>outputfile</i>	Output file in PDB format
-omm <i>outputfile</i>	Output file in MacroModel format
-omae <i>outputfile</i>	Output file in Maestro format

Options:

-data <i>datafile</i>	Read <i>datafile</i> for more default data
-model <i>number</i>	Specify an input structure other than the first
-no_color	Color atoms by type, not “error code”
-no_dup_connect	Don’t take duplicate CONECT records to mean multiple bonds
-no_renum	Retain atom indices from the MacroModel file
-warn_h	Include hydrogens in missing atom and unrecognized atom name warnings

Conversion Behavior for Importing PDB

Where multiple atomic coordinates exist for a single PDB entry, Maestro chooses the atoms with the highest occupancy ratio.

A few amino acids require choices to be made for placement of double bonds and formal charges. Schrödinger's conventions for these choices are given in [Section 3.1.3 on page 27](#).

Note: Unlike Maestro and SD files, PDB files do not support multiple structures, only multiple sets of coordinates for the same molecule, e.g. different NMR refinements. Certain PDB structures do not follow the PDB standard and have varying number of atoms in different MODELS. Do not import these PDB structures, as they will cause Maestro to fail.

Conversion Error Codes for `pdbconvert`

When running `pdbconvert` from UNIX, you will receive numerical messages (2, 1, 0) indicating the status of an attempted PDB conversion. These messages are equivalent to Maestro's dialog box warnings. The message numbers are defined as follows:

- 2 (ERROR): A fatal error occurred, and no Maestro file was generated. Check disk permissions and disk space.
- 1 (WARNING): A Maestro file was created, but an error at or above the base error level was returned. At the default base error level, this value is returned when red or blue atoms are present, and for orange atoms if the duplicate CONECT record convention is not followed. See [Table 3.2 on page 28](#) for a description of the atom colors.
- 0 (OK): A Maestro file was generated without errors at or above the base error level.

Using Templates for Non-Standard Residues

The `pdbconvert` program has been supplied with connectivity and bond order information for standard residues. For non-standard residues, the program predicts connectivity, but not bond order. However, if the connectivity and bond order information for a non-standard residue is known, it can be used to read in the PDB file. Create a new file that contains a template with this information (see example below).

```
TEMPLATE{
"ALA "
" N " " CA " 1
" N " " H " 1
" CA " " CB " 1
" CA " " C " 1
" C " " O " 2
" HA " " CA " 1
"1HB " " CB " 1
"2HB " " CB " 1
"3HB " " CB " 1
}
```

A connectivity template must begin with `TEMPLATE{` and end with a closing brace `}`. Each line in the template should specify two atom names, followed by a bond order. Templates for ligand molecules can also be created. To specify a template file to be used by `pdbconvert`, use the `-data` option.

The standard residue template file for `pdbconvert` can be viewed at

```
$SCHRODINGER/mmshare-vversion/data/mmpdb.ini.
```

D.1.4 Conversions to and From SD Format: `sdconvert`

Convert between MDL SD, Maestro, and MacroModel format files.

Syntax:

```
sdconvert [options] {-isd|-imm|-imae} inputfile
           {-osd|-omm|-omae} outputfile
```

<code>-isd</code>	Input file in SDfile format
<code>-imm</code>	Input file in MacroModel format
<code>-imae</code>	Input file in Maestro format
<code>-osd</code>	Output file in SDfile format
<code>-omm</code>	Output file in MacroModel format
<code>-omae</code>	Output file in Maestro format

Options:

<code>-a</code>	Append structures to the output file
<code>-all</code>	Convert all structures (synonym for <code>-n 1:</code>). This is the default action.
<code>-h</code>	Show usage summary
<code>-n <i>structs</i></code>	Specify a comma delimited series of colon-separated ranges and single values to convert. For example: <code>1:10,14</code> (structures 1 through 10 and 14) <code>2:</code> (structure 2 through end of file) <code>:5,13:18</code> (structure 1 through 5, and 13 through 18)
<code>-noarom</code>	Do not convert aromatic type 4 bonds to single and double bonds, according to Maestro and OPLS conventions.
<code>-nolewis</code>	Do not fix bond orders and formal charges of nitro and carboxylate groups as appropriate for OPLS use.

- st When files are converted from SD to Maestro format, convert parity and bond-direction information into chirality properties, which are recorded in the structures produced.
- title *prop* For SD to Maestro conversion, define SD property *prop* as the source for assigning the Maestro title.

For multiple -n options, the right-most takes precedence.

Example:

```
$SCHRODINGER/utilities/sdconvert -n 1: -isd lig.sdf -omae     lig.mae
```

D.2 Structure Extraction

This section describes utilities for extracting a subset of structures from a structure file. The first two utilities, `maesubset` and `sdssubset`, extract structures by index. The third utility, `propfilter`, extracts structures from a Maestro-formatted file based on values of properties.

maesubset

Extract a subset of structures from a Maestro format input file, by index or by title.

Syntax:

```
maesubset -n range full.mae > subset.mae
maesubset -t filename full.mae > subset.mae
maesubset [-h] [-v]
```

Options:

- h Show usage summary
- n *range* Extract structures by order in the input file. *range* is a comma separated list of indices or index ranges, in ascending order. Examples:
 - 1:200 (structures 1 through 200)
 - 50: (structures 50 through end of file)
 - "1:5, 14" (structures 1 through 5 and 14)
- t *filename* Extract structures whose titles are listed in the text file *filename*.
- v Show the version number of this utility.

sdssubset

Extract a subset of structures from an SD format input file.

Syntax:

```
sdssubset -n range full.sdf > subset.sdf  
sdssubset [-h] [-v]
```

<i>range</i>	1:200	(structures 1 through 200)
	50:	(structures 50 through end of file)
	"1:5, 14"	(structures 1 through 5 and 14)

Options:

-h	Show usage summary.
-v	Show the version number of this utility.

propfilter

Filter Maestro-formatted input file based on properties. Can filter on QikProp property names or internal Maestro property names.

Syntax:

```
propfilter [options] input-file  
propfilter [-h] [-n] [-v]
```

Options:

-e <i>filter-condition</i>	A single filter condition. Multiple such conditions may be specified.
-f <i>filter-file</i>	Input file listing filter conditions to be used.
-h	Show usage summary.
-n	List the recognized property names. (The input file is not processed.)
-o <i>output-file</i>	Filename for output Maestro- formatted file. If no output file is specified, structures are written to the file <code>propfilter.mae</code> .
-v	Show the version number of this utility.

Filter conditions must be specified in the form “*name op value*”, or just *name*. Here, *name* is the name of a property or descriptor.

name

- For QikProp properties, the standard CSV-file header name may be used, e.g., *MW* or *#nstars* or *dip^2/V*
- The internal Maestro property names may also be used, e.g., *r_qp_mol_wt*, or *i_qp_n_stars* or *r_rp_d2ov*.

op

The conditional operator *op* must be one of the following:

```
>    >=   <    <=
==   !=    ~    !~
```

and it must be surrounded by white space. The *~* and *!~* operators perform pattern matching using regular expressions.

If a condition is simply *name*, then the named property is required to exist, but it may have any value.

If filter conditions are supplied in an input filter file, there must be one condition per line. Lines that start with # are treated as comments. Blank lines are ignored.

D.3 Display of Properties

proplist

List properties in Maestro-formatted files.

Syntax:

```
proplist [options] input-files
```

Options:

- | | |
|----|--|
| -a | Report <i>all</i> the properties found in the input files. |
| -c | Write properties to a CSV file. |
| -h | Show usage summary. |
| -l | List the names of all the properties found in the input files. |

-
- n List the property names for which user-friendly names are known.
(Input files are ignored).
 - p *property-names* Names of properties to report. If more than one property is given, the names must be quoted. e.g., -p "title MW". This is the same as
-p title -p MW.
 - v Show the version number of this utility.

D.4 Structure Preparation

applyhtreat

Apply a hydrogen treatment (add or delete hydrogens) to one or more structures in a Maestro format file.

Syntax:

```
applyhtreat input-file output-file [ "treatment-name" ] [ "ASL-expression" ]
```

Options:

Allowed treatment names:

```
All-atom with Osp3/Nsp3-Lp
All-atom with No-Lp
Csp3 United-atom with S-Lp
Csp3 United-atom with No-Lp
All-atom with S-Lp
C sp2/sp3 United-atom with No-Lp
C sp2/sp3, N,O,S United-atom with No-Lp
```

Default treatment: All-atom with No-Lp

The treatment name must be enclosed in quotes. The ASL expression can be used to specify a range of atoms to which the treatment is to be applied. It must also be enclosed in quotes.

The `applyhtreat` utility can be used to delete or add hydrogens to one or more structures in a Maestro format file, consistent with the hydrogen treatment option that you provide. A *hydrogen treatment* is a protocol that determines which atoms are to have hydrogens and/or lone pairs attached. Several different treatments are supplied. Each of these treatments is associated with a particular molecular mechanics force field. However, some treatments are suitable for several force fields. See [Table 4.1 on page 52](#) for information on the correspondence of treatments to force fields.

The `applyhtreat` utility is used by Maestro when a hydrogen treatment is applied in the Hydrogen Treatment tab of the Build panel. See [Section 4.8 on page 52](#) for more information.

Note:

- Hydrogen addition or deletion takes place only as necessary to make the structure consistent with the selected hydrogen treatment. For example, if the specified treatment option calls for all atoms to have hydrogens and the structures in the input file already have hydrogens on all atoms, no changes will be performed.
- Maestro changes the default hydrogen treatment automatically if the force field is changed. The new hydrogen treatment will be one of those acceptable for the newly selected force field.

D.5 Unique Names

unique_names

This program takes files in Maestro format and adjusts the numeric suffixes for structure titles and entry names to make them unique. This is useful for post-processing output from programs, such as LigPrep and MacroModel, that produce structures with identical titles and entry names.

Note: If the structures in the input file already have unique names, these names might not be preserved if they only differ from a non-unique name by the numeric suffix.

Syntax:

```
unique_names input-file output-file
```

Example:

The file `ligprep.mae` contains the following structures:

Entry Name	Title
entry	title
entry	title2
entry	title2
entry	title3

To create unique titles and entry names for these four structures, enter the following command:

```
unique_names ligprep.mae ligprep-out.mae
```

The structures in the resulting file `ligrpep-out.mae` now have unique names:

Entry Name	Title
entry	title
entry1	title2
entry2	title3
entry3	title4

Customizing Colors

Colors, color schemes, and color ramps are defined in resource files that are installed in the resources subdirectory of the Schrödinger software installation,

```
maestro-vversion/data/res
```

where *version* is the 5-digit Maestro version number. This directory will be referred to as the *installation resources subdirectory*.

You can customize these resources by copying them to your *Maestro user resources directory*:

UNIX: `$HOME/.schrodinger/maestrovv`

Windows: `%USERPROFILE%\schrodinger\maestrovv`

where *vv* is the 2-digit Maestro version, and then editing them. Maestro reads the files in this directory when it starts. Instructions for editing these files are given in the subsections below.

You can also change the standard color resources by editing these files in the distribution. If you do, the colors are changed for all users.

E.1 Color Definitions

The colors available to Maestro are defined in the file `colors.res` in the installation resources subdirectory. To define your own colors, copy this file to your Maestro user resources directory and edit it. You can also change the resource file in the installation, but this will change the colors for all users.

Lines beginning with a `#` are comment lines. Blank lines are ignored. Each line that defines a color has an index number, a triplet of RGB color values in decimal format, and a color name. There must be no spaces after the name.

Maestro supports 256 colors. Colors 1–128 are used by various parts of Maestro. Colors 129–256 are user-defined colors: these colors can be used for new color schemes and color ramps, for example.

E.2 Color Schemes

Each of the standard color schemes is defined in a separate file in the `schemes` subdirectory of the installation resources subdirectory. By convention, color scheme files end with `.sch`. The schemes are called by Maestro using the file `scheme.res`, which is in the installation resources subdirectory

To define your own color schemes:

1. Create a `scheme.res` file in your Maestro user resources directory, or copy `scheme.res` from the installation resources subdirectory and edit it.
2. Create a `schemes` subdirectory in your Maestro user resources directory.
3. Create color scheme files `scheme.sch` in the `schemes` subdirectory, or copy the `.sch` files from the `schemes` subdirectory of the installation resources subdirectory and edit them.

The `scheme.res` file

The `scheme.res` file contains a three-line description of each scheme. The three lines contain the following information:

1. The short name of the color scheme
2. The long name of the color scheme. This is the name used in the Atom Coloring panel.
3. The name of the file that contains the scheme definition

Lines beginning with a `#` are comment lines. Either the short name or the long name can be used to refer to the color scheme when using the `colorscheme scheme=scheme-name` command in Maestro.

Color scheme files

The color scheme files contain a number of rules, one rule per line, which are applied in the order they appear in the file. Each rule consists of the name of a color as it appears in the `colors.res` file and an ASL specification that describes which atoms are to have that color. Color names must be a single token (there must be no spaces). If the color names do have an embedded space, they must be quoted (for example, "spring green").

The following is an example of a possible color scheme:

```
white atom.ele *
blue atom.ele N
red atom.ele O
```

This scheme first sets all atoms to white, then colors all the nitrogen atoms blue and the oxygen atoms red. Any number of such rules can be included in color scheme files.

E.3 Color Ramps

Color ramps are used in Maestro to map volume data onto a surface. In addition to the default color ramps provided with Maestro, you can also create your own color ramps. Like color schemes, each of the standard color ramps is defined in a separate file. These ramp files are stored in the `ramps` subdirectory of the installation resources subdirectory.

By convention, color ramp files end with `.rmp`. The color ramps are called by Maestro using the file `ramp.res`, which is stored in the installation resources subdirectory.

To define your own color ramps:

1. Create a `ramp.res` file in your Maestro user resources directory, or copy `ramp.res` from the installation resources subdirectory and edit it.
2. Create a `ramps` subdirectory in your Maestro user resources directory.
3. Create color ramp files `ramp.sch` in the `ramps` subdirectory, or copy the `.sch` files from the `ramps` subdirectory of the installation resources subdirectory and edit them.

The `ramp.res` file

The `ramp.res` file contains a three-line description of each ramp. The three lines contain the following information:

1. The short name of the color ramp
2. The long name of the color ramp. This is the name used in the Display Options dialog box of the Surface Table panel
3. The name of the file that contains the ramp definition

Lines beginning with a `#` are comment lines. Either the short name or the long name can be used to refer to the color ramp when using the `surfacescheme colorramp=ramp-name` command in Maestro.

Color ramp files

The color ramp files contain a list of color names, in the order that they will be used in the color ramp. The names of the colors must be taken from the `colors.res` file. Color names must be a single token (there must be no spaces). If the color names do have an embedded space, they must be quoted (for example, "spring green"). Lines beginning with a `#` are comment lines. The following example defines a simple 5-color ramp:

```
# Number of colors
# 5
# Color names for all colors of this color ramp
red
orange
yellow
green
blue
```

When a color ramp is applied, the range of values that define the ramp is divided evenly between the colors. The first color on the list is applied to the most negative value, the last color to the most positive. If you want to extend the range of values mapped to a single color, you can duplicate the color in the list. For example, to extend the range for the endpoints in the above example, you could duplicate red and blue, as follows:

```
# Number of colors
# 7
# Color names for all colors of this color ramp
red
red
orange
yellow
green
blue
blue
```

Setting Environment Variables on Windows

To create a new environment variable:

1. Right-click on My Computer, and choose Properties from the shortcut menu.

The System Properties dialog box opens.

2. In the Advanced tab, click Environment Variables.

The Environment Variables dialog box opens.

3. In the User variables section, click New.

The New User Variable dialog box opens.

4. Enter the name of the variable and its value, and click OK.

The New User Variable dialog box closes, and the variable is added to the User variables section of the Environment Variables dialog box.

5. Click OK in the Environment Variables dialog box, and again in the System Properties dialog box.

To modify an existing environment variable:

1. Right-click on My Computer, and choose Properties from the shortcut menu.

The System Properties dialog box opens.

2. In the Advanced tab, click Environment Variables.

The Environment Variables dialog box opens.

3. In the User variables section, select the environment variable you want to modify.

4. Click Edit.

The Edit User Variable dialog box opens.

5. Change the value of the variable and click OK.

The Edit User Variable dialog box closes, and the variable is updated in the User variables section of the Environment Variables dialog box.

6. Click OK in the Environment Variables dialog box, and again in the System Properties dialog box.

Copyright Notices

NCSA HDF5 Software Library and Utilities

Copyright Notice and Statement for NCSA HDF5 (Hierarchical Data Format 5) Software Library and Utilities

NCSA HDF5 (Hierarchical Data Format 5) Software Library and Utilities

Copyright 1998, 1999, 2000, 2001, 2002, 2003, 2004 by the Board of Trustees of the University of Illinois

All rights reserved.

Contributors: National Center for Supercomputing Applications (NCSA) at the University of Illinois at Urbana-Champaign (UIUC), Lawrence Livermore National Laboratory (LLNL), Sandia National Laboratories (SNL), Los Alamos National Laboratory (LANL), Jean-loup Gailly and Mark Adler (gzip library).

Redistribution and use in source and binary forms, with or without modification, are permitted for any purpose (including commercial purposes) provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions, and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions, and the following disclaimer in the documentation and/or materials provided with the distribution.
3. In addition, redistributions of modified forms of the source or binary code must carry prominent notices stating that the original code was changed and the date of the change.
4. All publications or advertising materials mentioning features or use of this software are asked, but not required, to acknowledge that it was developed by the National Center for Supercomputing Applications at the University of Illinois at Urbana-Champaign and to credit the contributors.
5. Neither the name of the University nor the names of the Contributors may be used to endorse or promote products derived from this software without specific prior written permission from the University or the Contributors, as appropriate for the name(s) to be used.

6. THIS SOFTWARE IS PROVIDED BY THE UNIVERSITY AND THE CONTRIBUTORS “AS IS” WITH NO WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED. In no event shall the University or the Contributors be liable for any damages suffered by the users arising out of the use of this software, even if advised of the possibility of such damage.

Portions of HDF5 were developed with support from the University of California, Lawrence Livermore National Laboratory (UC LLNL). The following statement applies to those portions of the product and must be retained in any redistribution of source code, binaries, documentation, and/or accompanying materials:

This work was partially produced at the University of California, Lawrence Livermore National Laboratory (UC LLNL) under contract no. W-7405-ENG-48 (Contract 48) between the U.S. Department of Energy (DOE) and The Regents of the University of California (University) for the operation of UC LLNL.

DISCLAIMER: This work was prepared as an account of work sponsored by an agency of the United States Government. Neither the United States Government nor the University of California nor any of their employees, makes any warranty, express or implied, or assumes any liability or responsibility for the accuracy, completeness, or usefulness of any information, apparatus, product, or process disclosed, or represents that its use would not infringe privately-owned rights. Reference herein to any specific commercial products, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or favoring by the United States Government or the University of California. The views and opinions of authors expressed herein do not necessarily state or reflect those of the United States Government or the University of California, and shall not be used for advertising or product endorsement purposes.

C and C++ Libraries for Parsing PDB Records

The C and C++ libraries for parsing PDB records are Copyright (C) 1989 The Regents of the University of California. All rights reserved.

Redistribution and use in source and binary forms are permitted provided that the above copyright notice and this paragraph are duplicated in all such forms and that any documentation, advertising materials, and other materials related to such distribution and use acknowledge that the software was developed by the University of California, San Francisco. The name of the University may not be used to endorse or promote products derived from this software without specific prior written permission. THIS SOFTWARE IS PROVIDED “AS IS” AND WITHOUT ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE.

Glossary

ASD—Atom Selection dialog box.

ASL—Atom Specification Language.

button menu—The menu available from a toolbar menu button, which you open by holding down the left mouse button.

entry—A structure or set of structures and associated properties. Entries are represented as a row in the Project Table, and can be used as input for jobs.

ESL—Entry Specification Language.

exporting—Writing structures and their associated data from a Maestro project to a structure file, or writing data from a Maestro project to a spreadsheet file.

filter—An ESL expression that is used to select, or filter, entries in the Project Table.

grow bond—The bond that is replaced when you add a fragment to a structure in Grow mode. The grow bond is displayed as a green arrow in the Workspace.

importing—Reading structures and their associated data from a structure file into a Maestro project, or reading data from a spreadsheet file into a Maestro project.

inclusion—Display of entries in the Workspace. Multiple entries can be included in the Workspace. The inclusion state (included, excluded, or fixed) is indicated in the In column of the Project Table.

isovalue—the value of a property that defines a surface. The property is a function of the three spatial coordinates, and is represented as a set of points on a three-dimensional grid.

menu button—A toolbar button that has a menu, which you open by holding down the left mouse button. The button has a black triangle in the lower right corner.

option menu—A drop-down menu from which you can select an option. Clicking the button on the menu displays the menu without having to hold down the left mouse button.

picking—Selecting atoms in the Workspace.

project—A collection of related data, such as structures with their associated properties. The project may include data that does not appear in the Project Table.

Project Table—The table that displays the contents of a project, with rows of entries and columns of properties. The table is located in the Project Table panel.

selection—(1) A group of atoms in the Workspace that has been picked or otherwise chosen. The selection is highlighted by markers when the appropriate option is selected. (2) A group of entries in the Project Table. The entries are highlighted in the table.

shortcut menu—The menu that is displayed when you right-click and hold on a component of a panel, such as the sequence viewer in the main window.

synchronization—updating the project with changes made in the Workspace. The default is to perform synchronizations automatically.

volume—a set of data points on a three-dimensional grid for a property that is a function of the three spatial coordinates.

Workspace—The open area in the center of the Maestro main window in which structures are displayed.

A

active atom, changing	43
Add Property panel	111
Adjust panel	50
Advanced Transformations panel	85
aliases	174
angles	
adjusting	50
measuring	117
ASD	223
ASL	223
Atom Coloring panel	71
atom labels	
adding user text to	78
applying	79
applying while building	38
coloring	155
composing	78
hetero atoms	78
including descriptions in	156
placement	155
removing	79
updating	79
Atom Labels panel	79
atom selection	
by property	65
by proximity	66
examples	67
Atom Selection dialog box	64
atom types (MacroModel)	
setting	54
table	196
united, option for use	158
atoms	
active, changing	43
changing element	44
coloring all	69
coloring with schemes	71
deleting	59
dummy	59, 121
hydrogen, adding	52
labeling	78
molecular representation	72
moving	38, 44
names	56
retyping	44
selecting for measurements	118

valency violation	43
Auto-Help facility	185
axis display in plots	134

B

Babel structure file format conversion	
program	25, 31
Ball & Stick representation	
applying	72
modifying	74
Balloon Help	185
bond order	
automatic assignment	126
changing	46
increasing in Draw mode	44
representation	73
template for pdbconvert	206
bonds	
adjusting	50
changing order of	45
displaying order of	73
hydrogen, displaying	119
representation	73
smoothing of lines	73
Build panel	37
Atom Properties tab	45, 53–57
Fragments tab	40, 42
Residue Properties tab	58
button menu	11
button menu, item types	12

C

Calculate Properties panel	112
centroid, defining	59
Centroids panel	59
chain names	58
charges, assigning	54
chirality information from SD files	208
Choose Entry dialog box	48
Choose Import Key dialog box	36
citing Maestro in publications	6
clipping planes window	9
hiding and redisplaying	174
color schemes	
atom labels	78
atoms	69–71

- customizing..... 216
 - PDB import..... 28
 - resource files..... 216
 - ribbons..... 77
 - surfaces..... 149
 - Workspace background..... 160
 - columns, Project Table
 - aligning data..... 99
 - moving..... 115
 - resizing..... 99
 - Command Aliases panel..... 174
 - command aliases, creating..... 174
 - command history..... 171
 - command input area..... 10
 - hiding..... 173
 - command script..... 171
 - Command Script Editor panel..... 170
 - command scripts..... 170–174
 - customizing Maestro with..... 172
 - running with ePlayer..... 110
 - compressing MacroModel files..... 203
 - Connect & Fuse panel..... 47
 - connectivity template for pdbconvert..... 206
 - Connolly surface..... 139
 - contacts
 - criteria for..... 120
 - displaying..... 119
 - Continuous mode, operating ePlayer in..... 108
 - conventions
 - document..... xiii
 - entry names..... 28
 - entry renaming..... 201
 - converting file formats
 - MacroModel..... 203
 - Mol2..... 204
 - PDB..... 205
 - SD..... 207
 - using Babel..... 25
 - coordinates, updating..... 83
 - Copy Property panel..... 113
 - copy-and-paste..... 113
 - CPK representation
 - applying..... 72
 - as van der Waals surface..... 140
 - modifying..... 74
 - current working directory..... 5
 - customization
 - colors..... 215
 - perspective..... 162
 - stereo viewing..... 161
 - with command scripts..... 172
 - Workspace appearance..... 160
 - See also* preferences
- ## D
- decimal places displayed, setting for
 - properties..... 114
 - Delete Property dialog box..... 114
 - dihedrals
 - adjusting..... 50
 - measuring..... 117
 - directories
 - changing..... 5
 - creating with Make Dir..... 90
 - current working..... 5, 7
 - file import..... 25
 - installation resources..... 215
 - launch..... 7
 - Maestro user resources..... 215
 - preferences..... 151
 - Display Options panel
 - Background tab..... 161
 - Fog tab..... 164
 - Perspective tab..... 163
 - Stereo tab..... 162
 - Display Options panel, surfaces..... 148
 - Display/Undisplay Atoms panel..... 80
 - distances, measuring..... 117
 - Draw mode..... 43
 - dummy atoms..... 59
 - defining for measurements..... 121
- ## E
- electrostatic potential, displaying on a
 - surface..... 149
 - element, changing..... 44
 - entries..... 3
 - adding properties to..... 111
 - combining..... 106
 - copying properties..... 113
 - creating from Workspace..... 91
 - deleting..... 106

-
- displaying in Workspace..... 104
 - duplicating 105
 - excluding 104
 - excluding, using plots 131
 - filtering 101–103
 - groups 97, 99
 - importing 91
 - including 104
 - including, using plots 131
 - incorporating from job output 92
 - merging 92
 - merging on connection or fusion 48
 - moving 106
 - name conflict resolution 92
 - properties 110–115
 - random selection 101
 - renaming 105
 - renaming conventions for 201
 - selecting 100–104
 - selecting from plots 131, 134
 - sorting 94, 107
 - splitting by molecule 106
 - suffixes 201
 - superimposing 122
 - tiling 86
 - See also* properties of entries
 - Entry Selection dialog box 102
 - Entry Selection Language (ESL) 101
 - environment variables
 - BABEL_DIR 25
 - DISPLAY 8
 - SCHRODINGER 8
 - TEMP_PROJECT 7, 87
 - USERPROFILE (Windows) 7
 - ePlayer 108–110
 - reference structure 110
 - speed 109
 - toolbar controls 95
 - ePlayer Options panel 109
 - error codes, pdbconvert utility 206
 - error reporting, in Maestro PDB conversion 27
 - ESL—*see* Entry Selection Language (ESL)
 - excluded entries 94
 - Export panel 31
 - Export Spreadsheet panel 33
 - Extended Radius Surface panel 140
 - extended radius surfaces—*see* surfaces
 - extracting subset of structures
 - Maestro input file 208
 - SD input file 209
 - F**
 - file extensions
 - imported files 24
 - preferences 157
 - file format
 - export, selecting 30
 - import, selecting 24
 - MacroModel, conversion to/from 203
 - Maestro, description 189
 - Maestro, example 191
 - Mol2, conversion to/from 204
 - PDB, conversion to/from 205
 - recognized 25
 - SD, conversion to/from 207
 - files
 - exporting data to 33
 - importing data from 34
 - importing structures from 91
 - JPEG, writing 179
 - overwriting 31
 - selecting for import 24
 - TIFF, writing 179
 - writing structures to 30
 - fill volume, CPK representation 75
 - filters, entry 101
 - Find panel 100
 - fixed entries 94
 - fog 163
 - font size, selecting 152
 - force fields, hydrogen treatment for 52
 - format conversion—*see* converting file formats
 - formats for plot output 136
 - fragment libraries 39–41
 - Fragment mode, building structures in 39
 - fragments
 - adding 39–41
 - building with 39
 - connecting and fusing 46
 - specifying growth settings for 42
 - frequency plots 127
 - full screen mode 10, 20
 - Function Key Macros panel 176

G

- global transformations 83
 - switching to 85
- groups, entry 5, 97, 99
- grow bond 42, 223
- Grow mode, building in 41
- growth settings, specifying 42

H

- help
 - Auto-Help 185
 - Balloon Help 185
 - online 183
- Help panel
 - Categories tab 185
 - Search tab 184
- hetero atoms, labeling 78
- hydrogen bonds
 - criteria for 120
 - displaying 119
- hydrogen treatment
 - applying 52, 211
 - force field compatibility 52
- Hydrogen Treatment panel 53
- hydrophobic and hydrophilic regions,
 - displaying 141
- Hydrophobic/philic map panel 141
- hydrophobic/philic surfaces 141–144

I

- image output
 - PostScript and EPS files 177
 - TIFF and JPEG files 179
- image quality
 - JPEG files 180
 - surfaces 138
- Import panel 23
- Import Spreadsheet panel 35
- importing structures 23
- included entries 94
- inclusion state 93
- incorporation of job results 4
- Install Scripts panel 168
- isosurface 137
- isovalue 137, 147
 - hydrophobic/philic map defaults 143

- restrictions on changing 145

J

- Jaguar files, importing 29
- Jaguar plot files, importing 145
- jobs
 - launching, introduction to 4
 - output incorporation 92
- joining geometry, specifying 42
- JPEG files 136, 179

L

- Label Points dialog box 131
- labels
 - adding user text to 78
 - applying 79
 - applying while building 38
 - coloring 155
 - displaying for hidden atoms 156
 - field descriptions in 156
 - hetero atoms 78
 - placement 155
 - removing 79
 - updating 79
- launch directory 7
- Limit panel, surfaces 148
- line plots 127
- local transformations 84–86
 - markers for 86
 - rotations 86
 - switching to 85
 - translations 86

M

- MacroModel file format, converting to/from .. 203
- Macros panel 175
- macros, creating 175
- Maestro
 - citing use of 6
 - file format 189
 - file format example 191
 - help 183
 - main window 9, 10
 - menus 11
 - starting 8
 - workflow paradigm 3

-
- Maestro command language 170
 - main window 10
 - components, hiding and displaying 9
 - Make Dir, creating directories with 90
 - Manage Scripts panel 167
 - markers
 - Advanced Transformation 86
 - displaying 117
 - H-bond and contact measurement 120
 - measurement, geometric 118
 - MDL SD file format, converting to/from 207
 - measurements
 - angles 117
 - creating properties from 121
 - dihedrals 117
 - distances 117
 - dummy atoms 121
 - H-bond and contact 119
 - writing to a file 121
 - Measurements panel
 - Contacts tab 120
 - Distances tab 118
 - H-Bonds tab 119
 - menu button 11
 - Mol2 file format, converting to/from 204
 - molecular representation
 - atoms 72
 - bonds 73
 - default 75
 - saving 31
 - Molecular Representation panel
 - Atoms tab 72
 - Bonds tab 74
 - opening 72
 - Ribbons tab 76
 - Settings tab 75
 - Molecular Surface panel 139
 - molecular surfaces—*see* surfaces
 - molecules
 - connecting and fusing 46
 - drawing 43–44
 - See also* structures
 - Mopac file format 25
 - mouse functions
 - general 2, 19
 - Project Table 96
 - Workspace 84
 - mouse, two-button 2
 - multiple bonds, creating in Draw mode 44
 - multi-structure files, exporting to 31
 - mutating residues 41
- ## N
- New Plot dialog box 128
 - NMR spin-spin coupling constants,
 - displaying 121
- ## O
- object types 12, 61, 62
 - online help 183
 - opening state of project, reverting to 116
- ## P
- panel layout, saving 151
 - panning plots 130
 - parities 208
 - partial charge, assigning 54
 - pasting values into the Project Table 113
 - PDB atom names, assigning 55
 - PDB database, importing from 26
 - PDB file format conversion 205
 - PDB files, color scheme used on import 28
 - PDB residue name, assigning 58
 - pdbconvert utility 205
 - perspective, enabling 162, 173
 - Pick option menu 62
 - pick states—*see* object types
 - Place mode, building in 41
 - Plot Settings panel
 - Plot tab 132
 - Series tab 133
 - X Axis tab 135
 - Y Axis tab 135
 - Plot XY panel 129
 - plots
 - adding data series to 133
 - axis display 134
 - creating 127
 - deselecting entries from 131
 - editing data 133
 - excluding entries from 131
 - file formats 136
 - including entries from 131

- labeling data points 130
 - line, frequency 127
 - modifying attributes 132
 - panning and zooming 130
 - saving 136
 - selecting 129
 - selecting entries from 131, 134
 - updating from project table 135
 - viewing 128
 - preferences
 - atom labels 155
 - builder 157
 - feedback 160
 - file suffix 157
 - font size 152
 - mouse 156
 - setting 152
 - storage location 7
 - synchronization 110, 153
 - working directory 155
 - preferences directory 151
 - Preferences panel
 - Builder tab 158
 - Directory tab 155
 - Feedback tab 159
 - File Suffix tab 158
 - Font Size tab 152
 - Label tab 156
 - Mouse tab 157
 - Project tab 153
 - Print panel 177
 - printing
 - command for 179
 - enhancing surface smoothing for 178
 - PostScript and EPS files 177
 - titles 178
 - Workspace background color 178
 - Workspace image 177
 - probe radius for surfaces 138
 - project entries—*see* entries
 - project selectors 89
 - Project Table 224
 - Project Table panel 4, 93
 - projects 3, 223
 - annotation 90
 - availability indicators for 89
 - closing 88
 - creating snapshot 116
 - deleting 89
 - locked 89
 - merging 92
 - opening 87
 - opening state, reverting to 116
 - renaming 88
 - saving scratch 88
 - scratch 3, 87
 - properties of entries 110–115
 - adding 111
 - copying 113
 - creating from measurements 121
 - creating from superposition 124
 - decimal places 114
 - pasting values 113
 - plotting 127
 - selecting entries by 101
 - selecting for display 114
 - sorting by 107
 - properties, listing in file 210
 - propfilter utility 209
 - proplister utility 210
 - proteins
 - active site mapping 141
 - aligning 126
 - assigning secondary structure 126
 - atom selection examples 67
 - capping groups 39
 - displaying as ribbons 76
 - mutating residues 41
 - viewing sequence 9
 - Proximity dialog box 66
 - Python commands
 - pyrun 169
 - pythonimport 169
 - pythonrun 169
 - Python scripting tutorial 169
- ## Q
- Quit panel 22
- ## R
- random selection of entries 101
 - Read panel 23
 - reference structure, choosing for ePlayer 110

-
- Rename Property panel 114
 - residues
 - assigning name 58
 - mutating 41
 - non-standard, templates for 206
 - one-letter names 156
 - renumbering 57
 - resource files
 - color 215
 - location 7
 - Ribbon options dialog box 77
 - ribbons
 - color scheme 77
 - display options 78
 - displaying 76
 - style 77
 - rotations
 - from superposition 124
 - local 86
- S**
- Save Image panel 179
 - Save PlotXY Image panel 136
 - Save Project As project selector 89
 - Schrödinger contact information 186
 - Schrödinger web site 186
 - scratch entry 4
 - scratch projects 3, 87
 - location 7
 - Scripting with Python 169
 - scripts, *see* command scripts, utility scripts
 - SD (MDL) file format, converting to/from 207
 - SD (MDL) file, extracting structures from 209
 - SD-formatted files
 - parities and bond directions 208
 - searching for text 100
 - Select menu 101
 - selecting atoms 61
 - selection state 93
 - sequence viewer 9
 - series—*see* plots
 - sets 124
 - reading and writing 125
 - selecting atoms for 125
 - Sets panel 125
 - Show Property panel 115
 - Show/Hide Commands panel 171
 - smoothing wire bonds 73
 - snapshot
 - project 116
 - undoing with 116
 - Sort Project Table panel 107
 - Spartan file format 25
 - spin-spin coupling constants, displaying 121
 - spot-centering on an atom 83
 - spreadsheet
 - exporting data to 33
 - importing data from 34
 - properties for export 115
 - status bar 9
 - Step mode, operating ePlayer in 108
 - stereo viewing method, specifying 161, 173
 - stereochemistry, changing 51
 - structures
 - appending to pre-existing files 31
 - building from fragments 39
 - building in Draw mode 43–44
 - building in Grow mode 41
 - building in Place mode 41
 - coloring 71
 - comparing 122
 - displaying sequentially 108
 - ePlayer reference 110
 - extracting by property from file 209
 - extracting subset from Maestro file 208
 - extracting subset from SD file 209
 - importing multiple 26
 - local translations 86
 - symmetrizing 51
 - tiling 86
 - substructure notation 112
 - superimposing structures 122
 - using the ePlayer 110
 - Superposition panel 123
 - Surface Table panel 146
 - surfaces
 - appearance of 148
 - atoms for display of 138
 - changing isovalue 147
 - color schemes 149
 - Connolly 139
 - context 138
 - deleting 147
 - displaying and undisplaying 148

- extended radius 139
 - generating 137–138
 - hydrophilic, hydrophobic 141
 - hydrophobic/philic 141–144
 - image quality 138
 - importing 145, 147
 - limiting 147
 - molecular 139
 - naming 137
 - probe radius 138
 - transparency 148
 - van der Waals 140
 - Sybyl (Mol2) file format, converting to/from 204
- T**
- templates, for pdbconvert 206
 - terminal atoms, deleting 60
 - text, finding in the Project Table 100
 - TIFF files 136, 179
 - tiling structures 86
 - toolbar
 - Build panel 38
 - main 11–14
 - main, hiding and redisplaying 173
 - Plot 129
 - Project Table panel 94–95
 - translations
 - from superposition 124
 - local 86
 - transparency, surface 144, 148
 - tube representation, adjusting width of 75
 - two-button mouse 2
- U**
- uncompressing MacroModel files 203
 - undoing
 - individual operations 21
 - molecular transformations 21
 - using snapshots 116
 - using the Undo menu option 21
 - unique names, creating 212
 - utility scripts
 - applyhtreat 211
 - maemmod 203
 - maesubset 208
 - mmio_convert 203
 - mmodmae 203
 - mmodmol 204
 - mol2convert 204
 - molmmmod 204
 - pdbconvert 205
 - propfilter 209
 - proplister 210
 - sdconvert 207
 - sdssubset 209
 - unique_names 212
- V**
- valency, violation of 43, 46, 47
 - van der Waals surfaces 140
 - use for molecular surface generation 139
 - VDW Surface panel 140
 - viewing transform, saving and restoring 21
 - volumes 137
 - importing 145
- W**
- warning messages, in Maestro PDB
 - conversion 27
 - window components (main), hiding and displaying 9
 - wire representation
 - changing appearance of 74
 - smoothing 73
 - Workspace 9, 224
 - adding to the selection 104
 - background color 160, 173
 - creating entries from 91
 - full screen mode 10, 20
 - image output 177, 179
 - scratch entry 4
 - toolbar 11
- X**
- XYZ file format 25
- Z**
- zooming
 - in plots 130
 - in Workspace 84

120 West 45th Street
32nd Floor
New York, NY 10036

101 SW Main Street
Suite 1300
Portland, OR 97204

3655 Nobel Drive
Suite 430
San Diego, CA 92122

Dynamostraße 13
68165 Mannheim
Germany

QuatroHouse, Frimley Road
Camberley GU16 7ER
United Kingdom

SCHRÖDINGER.